


lifelines neuro v2.8.4.7 Trackit Plus File Manager Software User Guide

Home » lifelines neuro » lifelines neuro v2.8.4.7 Trackit Plus File Manager Software User Guide 🖺

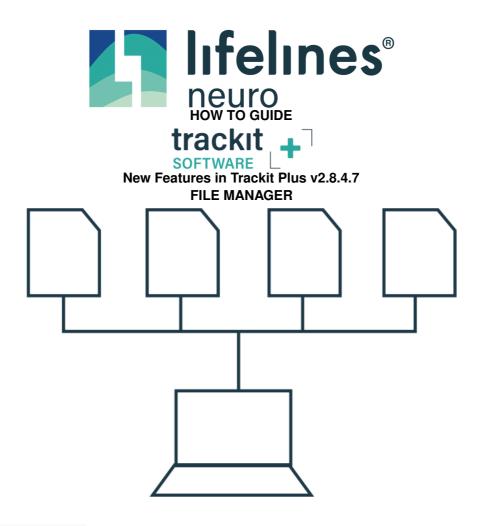

Contents


- 1 File Manager
- 2 Documents /
- **Resources**
- 2.1 References
- **3 Related Posts**

File Manager

At the end of a recording, the File Manager window and dialogue is displayed.

The File Manager automatically collects all the relevant associated files for a study and assembles them into appropriate folders on the PC.

This makes it much easier for local review and/or the upload of files for remote review to take place.

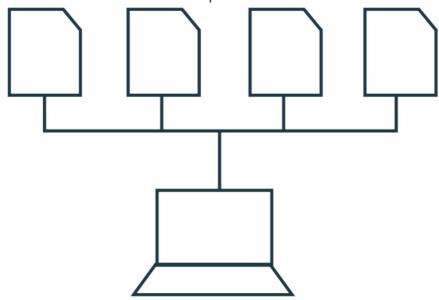

File Manager works for the following recording types:

· Ambulatory (amplifier) recordings with video



• Dual recordings with or without video




File Manager automatically assembles files on the PC as follows:


• Ambulatory folder, e.g. SmithJoh_220311134555

The File Manager assembles the ambulatory files, Events files and all video and synchronization files in this


folder.

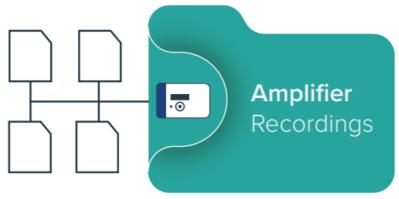

• PC folder with the [PCFiles] identifier, e.g. SmithJoh_220311134555[PCFiles] The File Manager assembles the PC files into this folder.

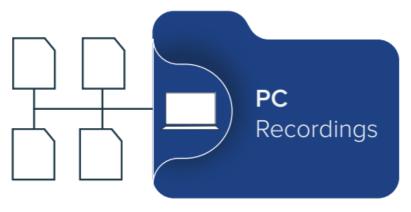

Following Stop Recording, the File Manager steps are as follows:

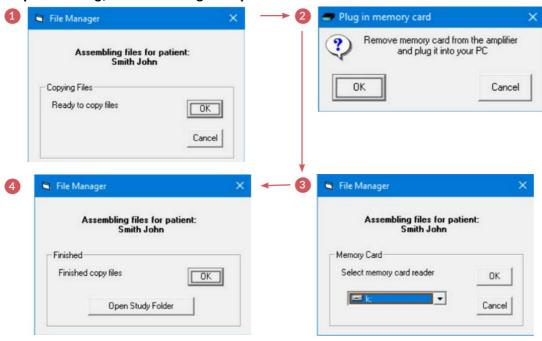

Note: The File Manager can be launched from the main File tab as well as automatically at Stop Recording. The options available here are as follows:


Either

Proceed as above for the last patient recording

Or via the drop down list:

· Select folder and files for File Manager

· Manually manage files

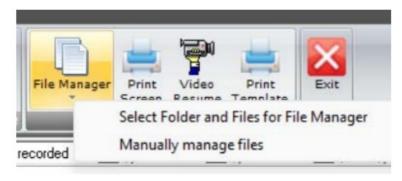



Tel: 866-889-6505
Email: support@lifelinesneuro.com
Customer Experience Department May 2023

Documents / Resources



<u>lifelines neuro v2.8.4.7 Trackit Plus File Manager Software</u> [pdf] User Guide v2.8.4.7 Trackit Plus File Manager Software, Trackit Plus File Manager Software, File Manager Software, Software, Manager Software

References

• User Manual

Manuals+, Privacy Policy

This website is an independent publication and is neither affiliated with nor endorsed by any of the trademark owners. The "Bluetooth®" word mark and logos are registered trademarks owned by Bluetooth SIG, Inc. The "Wi-Fi®" word mark and logos are registered trademarks owned by the Wi-Fi Alliance. Any use of these marks on this website does not imply any affiliation with or endorsement.