

Wireless accessories pairing guide



Phonak ComPilot Air II

Solution for:

- Connecting wireless Brio 3 hearing aids to cell phones or Bluetooth-enabled music sources

Features:

- No neckloop
- Built-in demonstration capability
- Clip for easy wearing
- Binaural streaming of audio signals in stereo
- Built-in remote control for hearing aids
- Spoken messages and caller identification for easy interaction
- Compatible with Phonak RemoteControl App



Phonak ComPilot II

Solution for:

- Connecting wireless Brio 3 hearing aids to cell phones or Bluetooth-enabled music sources

Features:

- Built-in demonstration capability
- Bluetooth and wired audio inputs
- 24 hours of continuous streaming time
- Binaural streaming of audio signals in stereo
- Built-in remote control for hearing aids
- Spoken messages and caller identification for easy interaction
- Exchangeable neckloop available in two colors and lengths
- Compatible with Phonak RemoteControl App



TVLink II

Solution for:

- Listening to TV or other audio devices

Features:

- Supports analog and digital (Optical, Coaxial) audio outputs
- Quick and easy installation
- Test mode for a quick audio check
- Complete set including power plug and cables
- Charging slot for Phonak ComPilot Air II or Phonak ComPilot II
- Compatible with all Phonak streamers



Phonak RemoteMic

Solution for:

- One-on-one conversations over distance

Features:

- Lightweight and small
- Omnidirectional microphone
- Rechargeable battery
- Up to 8 hours of usage time
- Compatible with all Phonak streamers



Phonak DECT II

Solution for:

- Home phone with direct binaural streaming to both hearing aids for maximum speech understanding

Features:

- Unpack and use: no fitting required
- Simultaneous audio transmission to both hearing aids
- 1-click volume boost (up to +15 dB)
- Direct dial button for preferred number
- Suitable for people without hearing loss; may be used like a regular cordless phone



Phonak EasyCall II

Solution for:

- Binaural streaming of a cell phone call to the hearing aids for maximum speech understanding

Features:

- No fitting required – simply pair to cell phone
- Works with any Bluetooth®-enabled cell phone, even non-smartphones
- Attaches onto the cell phone so it's never lost or forgotten
- 10 hour talk time
- Standby time > 1 week



Wireless accessories pairing guide



Phonak ComPilot Air II

Place Phonak ComPilot Air II in pairing mode

- Switch Phonak ComPilot Air II off
- Press and hold the volume up (+) button **A**
- Turn Phonak ComPilot Air II on while pressing the volume up (+) button **B**
- Release the volume up (+) button when the indicator light starts rapidly blinking blue

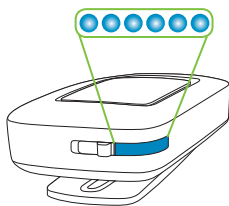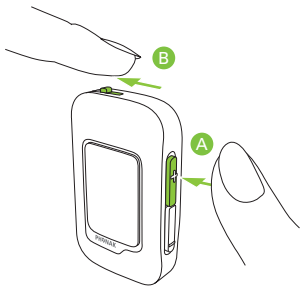

Bluetooth pairing: + on

Place Bluetooth device in pairing mode

- Select "Phonak ComPilot Air II" from the list of found Bluetooth devices
- Enter code "0000" if prompted
- Phonak ComPilot Air II requires access to your phone book. When prompted, accept the request
- After successful pairing, the indicator light will stop rapidly blinking blue



A quick pairing instruction is also shown on the back of Phonak ComPilot Air II. It is only necessary to perform the pairing procedure once with each Bluetooth device. Phonak ComPilot Air II will connect automatically to the phone after the initial pairing.



Phonak ComPilot II

Place Phonak ComPilot II in pairing mode

- Switch on the Phonak ComPilot II
- Press and hold the connect (<=>) and volume up (+) buttons simultaneously (2 sec) until the audio indicator starts to rapidly blink blue

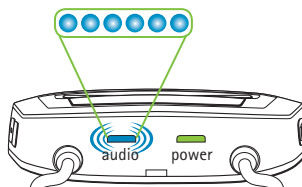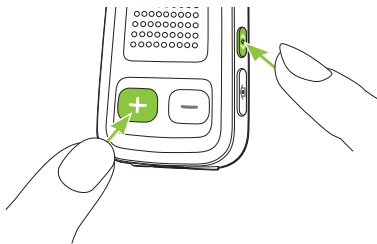

Place Bluetooth device in pairing mode

- Select "Phonak ComPilot II" from the list of found Bluetooth devices
- Enter code "0000" if prompted
- Phonak ComPilot II requires access to your phone book. When prompted, accept the request
- After successful pairing, the indicator light will stop rapidly blinking blue

A quick pairing instruction is also shown on the back of Phonak ComPilot II. It is only necessary to perform the pairing procedure once with each Bluetooth device. Phonak ComPilot II will connect automatically to the phone after the initial pairing.




TVLink II

Place Phonak ComPilot Air II or Phonak ComPilot II in pairing mode

Place TVLink II in pairing mode:

- For the first pairing, TVLink II will automatically be in pairing mode after switching on. This is indicated by a rapidly blinking blue main indicator light
- To manually start pairing mode, switch TVLink II on
- Press the pairing (P) button on the back and the main indicator light will start rapidly blinking blue
- Make sure TVLink II and the streamer are within range (3 ft)
- Lights will turn solid blue once pairing is complete

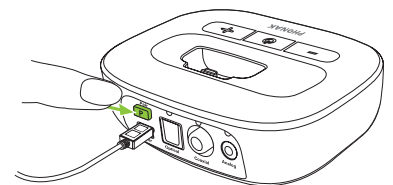

To test the connection (troubleshooting):

- Have the client wear the streamer in the on position
- With TVLink II off, press and hold the volume down (-) and power buttons simultaneously until the light blinks in multiple colors
- Test sound will be heard through the hearing aids if connection between TVLink II and streamer is successful




Phonak RemoteMic

Place Phonak ComPilot Air II or Phonak ComPilot II in pairing mode

Place Phonak RemoteMic in pairing mode:

- For the first pairing, Phonak RemoteMic will automatically be in pairing mode for 30 minutes after turning on. This is indicated by a rapidly blinking blue indicator light
- To manually start pairing mode, switch Phonak RemoteMic off
- Press and hold the volume up (A) button
- Turn the RemoteMic on while pushing the (B) button at the same time

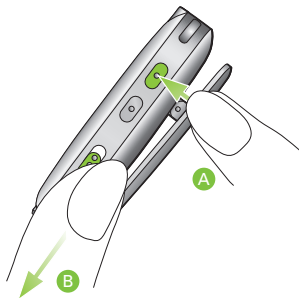

- Release the (A) button when the indicator light starts rapidly blinking blue
- Make sure the streamer and Phonak RemoteMic are within range (3ft)
- Lights will turn solid blue once pairing is complete



Phonak DECT II

No pairing needed

Setup:

- Unpack and use: no setup required
- Plug into a power outlet and home phone line with the included cables

Connecting to your hearing aids:

- Hold the handset close to the hearing aids (within 10")
- The hearing aids will automatically switch to the phone program



max. 25 cm / 10"

Testing the connection to your hearing aids:

- Press and hold the green "Boost" key



- Hold the phone to the ear and the test message will be wirelessly transmitted to the hearing aids
- Use the volume buttons to set the volume to a comfortable level
- Press the "End" key to end the test message

Making a call:

- Enter the phone number on the keypad
- Press the "Talk" key (phone icon) to place the call

Ending a call:

- Press the "End Call" key (phone icon with red X)

Accepting a call:

- Press the "Talk" key (phone icon)




Phonak EasyCall II

Charge for 3 hours before use (a full charge will allow 10 hours of talk time)

Setup:

- Attach EasyCall II to the back of a Bluetooth-enabled cell phone using the enclosed adhesive or Phonak hard case available for popular phone models
- Press and hold the call button until the indicator light begins to flash blue (A)
- Ensure that Bluetooth is enabled on the cell phone and begin searching for devices
- Select EasyCall II from the phone's list of Bluetooth devices found
- Enter "0000" if prompted for a PIN on the device

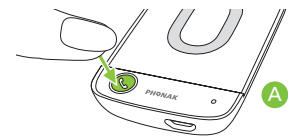

Making a call:

- Enter the phone number on the phone
- Hold the phone (with EasyCall II attached) as usual within 4" of the hearing instruments
- The hearing instruments will automatically switch to the phone program

Ending a call:

- Press the EasyCall II button (B) or use the "End Call" button on the phone
- The hearing instruments will automatically switch back to the last program



Accepting a call:

- Press the EasyCall II button (B) to accept the call