

Medtronic Smart MDI System App User Guide

[Home](#) » [Medtronic](#) » Medtronic Smart MDI System App User Guide 

Medtronic

Engineering the extraordinary



Contents

- 1 Smart MDI System App
- 2 Setting up the Smart MDI system
- 3 Connecting the Guardian™ app with the InPen™ app
- 4 Reminders and notifications
- 5 Documents / Resources
 - 5.1 References

Smart MDI System App

Smart MDI system*

Because CGM alone is not enough



A smart step-by-step user guide

Guardian™ 4 Smart CGM system InPen™ Smart insulin pen system

*Components sold separately




Setting up the Smart MDI system

Ensure you have the following ready:

- A. Guardian™ 4 transmitter and charger
- B. Guardian™ 4 sensor
- C. One Press serter
- D. Guardian™ 4 oval tape
- E. App installed & bluetooth enabled



Setting up the Guardian™4 app

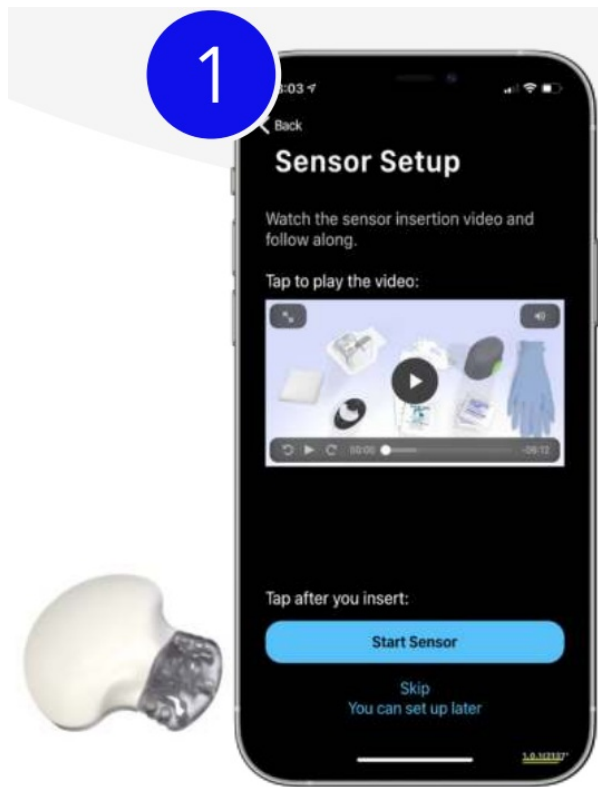

Download the Guardian™ app

Setting up the InPen™ app



Automatically helps track insulin data, administer personalised insulin doses and lets you know if a dose is missed.

- People download the InPen create an InPen™ account
- Pair the InPen †™ app and ™insulin pen to the app
- Follow the app tutorial in the app



Download the InPen™ app

Connecting the Guardian™ app with the InPen™ app


InPen™ system and Guardian™ 4 system are designed to work together.

- Connect the two apps together by following the final InPen™ app instructions
- You will be prompted to enter your CareLink™ Personal username and password




- The Guardian™ app sends the glucose data to the InPen system app so you have a one-stop, clear view™
- Click on Settings in bottom right corner and then on Connections and add all available connections Medtronic CGM and InPen™ Reports in CareLink™
- A Care Partner account can be set up for up to 5 care partners to be invited and be involved



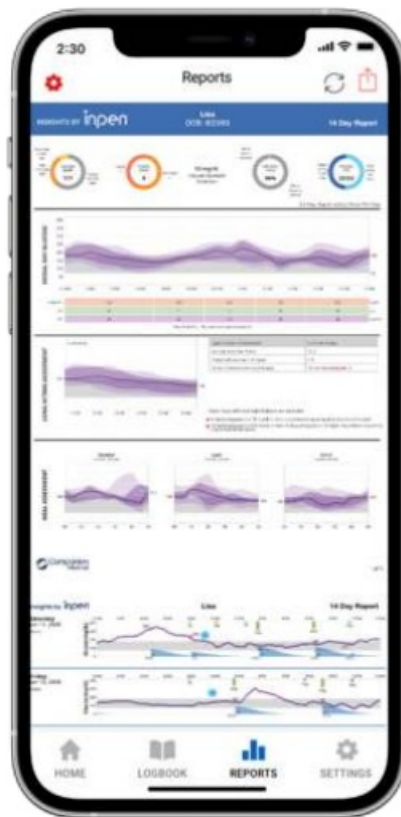



4. Start sharing reports

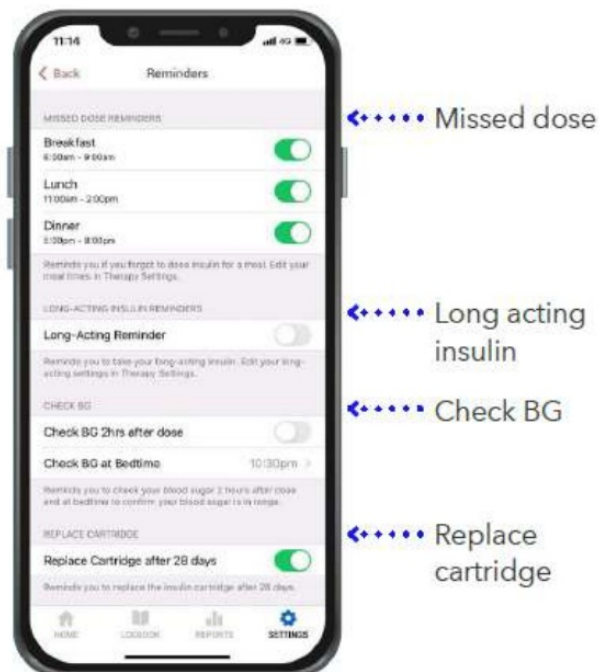

You can send you their InPen™ insight reports by email via the app.
 You can also share your CareLink™ Personal user name with your healthcare team so they can receive automatic data uploads.

Reminders and notifications



	Meal Dose Reminder: Dose Reminder is enabled, and no dose was taken during the time window.
	Long-acting Reminder: Long-acting Reminder is enabled, and no Long-acting dose was logged at the reminder time.
	Check BG Reminder: If enabled you will receive an alarm 2 hours after dose and at bedtime.
	Insulin Age: Cartridge Replace Reminder is enabled, and the InPen™ app has not detected a new cartridge being installed within the past 28 days.
	Low Battery: InPen™ is reaching the end of its 1-year life and needs to be replaced.
	Temperature: InPen™ app detects a very high {above 37°C} or very low (below 0°C) temperature and user may want to consider replacing the insulin cartridge.

Reminders can help you keep on track with your diabetes management plan.

† Some user interaction required.

References:

1. Abraham SB, et al. Improved Real-World Glycaemic Control With Continuous Glucose Monitoring System Predictive Alerts. Journal of Diabetes Science and Technology 2021; 15(1):91–97.
2. Smith M, et al. E5, SIPs Improve Time Below Range in MDI Therapy, AMCP Congress 2020.

This material does not replace or supersede the instructions for use. It should not be considered the exclusive source of information, and should be used in conjunction with the User Guide. See the User Guide for a detailed list of information regarding the instructions for use, indications, contraindications, warnings, precautions, and potential adverse events.

For further information, contact your local Medtronic representative and/or consult the Medtronic website at medtronic.diabetes.com.au

12987-082023 © Medtronic 2023. All rights reserved. Medtronic and the Medtronic logo are trademarks of Medtronic™

* Third party brands are trademarks of their respective owners. All other brands are trademarks of a Medtronic company.

Medtronic

Medtronic Australasia Pty Ltd

2 Alma Road Macquarie Park, NSW 2113 Australia.

Website: medtronic-diabetes.com.au

Tel: +61 2 9857 9000

Fax: +61 2 9889 5167


Toll Free: 1 800 777 808

Email: australia.diabetes@medtronic.com

Facebook: [www.facebook.MedtronicDiabetesAUS](https://www.facebook.com/MedtronicDiabetesAUS)

YouTube: Medtronic Diabetes ANZ

Documents / Resources

 <p>Smart MDI system® InPen™, InPen™, InPen™ A smart step-by-step user guide Medtronic Diabetes Australia</p>	<p>Medtronic Smart MDI System App [pdf] User Guide Guardian™ 4, InPen™, Smart MDI System App, MDI System App, System App, App</p>
--	---

References

- [M Medtronic Diabetes Australia | Insulin Pumps & CGM](#)
- [User Manual](#)

Manuals+, Privacy Policy

This website is an independent publication and is neither affiliated with nor endorsed by any of the trademark owners. The "Bluetooth®" word mark and logos are registered trademarks owned by Bluetooth SIG, Inc. The "Wi-Fi®" word mark and logos are registered trademarks owned by the Wi-Fi Alliance. Any use of these marks on this website does not imply any affiliation with or endorsement.