



addaday Bio6 Glow Massager Pillow with Bluetooth User Manual

[Home](#) » [addaday](#) » addaday Bio6 Glow Massager Pillow with Bluetooth User Manual

addaday Bio6 Glow Massager Pillow with Bluetooth User Manual



For consumers looking for effortless, on-the-go massage, look no further than the GLOW from Addaday. Its cordless design is versatile for use. almost anywhere, including in the car (it comes with a car charger). Six nodules rotate in two directions to knead knotted muscles, with two heated nodules that warm up to soothe stiffness.

Contents [[hide](#)]

- [1 IMPORTANT SAFETY INSTRUCTION](#)
- [2 PRODUCT DIAGRAM](#)
- [3 OPERATING INSTRUCTIONS](#)
- [4 Cleaning & Maintenance](#)
- [5 PRODUCT SPECIFICATIONS](#)
- [6 Battery Disposal](#)
- [7 TROUBLESHOOTING](#)
- [8 Documents / Resources](#)
- [9 Related Posts](#)

IMPORTANT SAFETY INSTRUCTION

When using this device, all instructions and basic precautions outlined in this document should be followed.

Read all instructions before using (this appliance).

DANGER – To reduce the risk of electric shock

1. Always unplug this device from the electrical outlet immediately after use and before cleaning.
2. Do not reach for an appliance that has fallen into water. Unplug immediately.
3. Do not use while bathing or in a shower.
4. Do not place or store appliance where it can fall or be pulled into a tub or sink. Do not place in or drop into water or other liquid.
5. Do not use massager in close proximity to loose clothing or jewelry.
6. Keep long hair away from massager while in use.

WARNING – To reduce the risk of burns, fire, electric shock, or injury to persons:

1. Unplug from outlet when not in use, and before putting on or taking off parts.
2. The massager should never be left unattended when plugged in. Unplug from outlet when it is not in use.
3. Do not operate under blanket or pillow. Excessive heating can occur that might lead to fire, electric shock, or injury.
4. Avoid operating the massager with the air openings blocked. Keep the air openings free of lint, hair, and other fine particles.
5. Avoid inserting objects into any opening.
6. Avoid using outdoors.
7. Avoid handling the electrical cord when your hands are wet.
8. Avoid carrying this device by its electrical cord or using the cord as a handle.
9. Keep the cord away from heated surfaces. If a fault is found, stop using the massager immediately, disconnect it from the power supply and contact the seller or manufacturer.
10. Remove the power plug from the socket immediately after power failure. Otherwise, it may cause electric shock, short circuit or fire.
11. Do not operate this device if it has a damaged cord or plug, if it is not working properly, if it has been dropped or damaged, or submerged in water. Contact the seller or manufacturer for examination.
12. Use this massager only for its intended use as described in this manual. Do not use attachments not recommended by the manufacturer.
13. The massager must only be supplied by a safe low voltage corresponding to the marking on the product.
14. The massager is only to be used with the power supply unit provided.
15. Close supervision is necessary when this appliance is used by, on, or near children, invalids, or disabled persons.
16. As the surface of this product may get hot, those who are sensitive to heat should pay extra attention while using it to avoid minor burns.
17. This massager contains batteries that are non-replaceable.
18. For the purpose of recharging the battery, only use the detachable supply unit provided with this massager.

If you fall under any of the following categories, it is recommended that you consult a physician prior to use:

1. If you are susceptible to electromagnetic malfunctions such as cardiac pacemakers and implantable medical electronic devices.
 2. If you have an inhibited sense of temperature.
 3. If you have malignant tumor(s).
 4. If you have heart disease.
 5. If you are diabetic.
 6. If you have scarring or tissue damage or have recently had surgery.
 7. Do not operate where aerosol (spray) products are being used or where oxygen is being administered.
 8. To disconnect, turn all controls to the off position, then remove plug from outlet.
 9. Prevent unintentional Ensure the switch is in the off-position before connect- ing to battery pack, picking up or carrying the appliance. Carrying the appliance with your finger on the switch or energizing appliance that have the switch on invites accidents.
 10. Recharge only with the charger specified by the A charger that is suitable for one type of battery pack may create a risk of fire when used with another battery pack.
 11. Under abusive conditions, liquid may be ejected from the battery; avoid If contact accidentally occurs, flush with water. If liquid contacts eyes, additionally seek medical help. Liquid ejected from the battery may cause irritation or burns.
 12. Do not use an appliance that is damaged or modified. Damaged or modified appliances may exhibit unpredictable behavior resulting in fire, explosion or risk of injury.
 13. 12. Do not expose a battery pack or appliance to fire or excessive temperature. Exposure to fire or temperature above 130°C may cause explosion.
 13. Follow all charging instructions and do not charge the battery pack or appliance outside of the temperature range specified in the Charging improperly or at temperatures outside of the specified range may damage the battery and increase the risk of fire.
 14. Have servicing performed by a qualified repair person using only identical replace- ment This will ensure that the safety of the product is maintained.
 15. Do not modify or attempt to repair the appliance except as indicated in the instruc- tions for use and care.
 16. Do not charge the massager in an environment where temperature is below 0 °C/32°F or above 40°C/104°F Charging improperly or at temperatures outside of the specified range may damage the battery and increase the risk of fire.
- Do not use the massager in an environment where the temperature is higher than 50°C/122°F

SAVE THESE INSTRUCTIONS

Certificate chargers list

Xiamen Innov Electronics Tech Co.,Ltd IVP1300-1800

Xiamen Eahunt Electronics Co., Ltd E27-1300-1800PU

Guangdong Shanmusen Technology Co.,Ltd SMS-01130180-S20



The massager is not to be used for commercial use or medical purposes.



If you have health concerns, consult your doctor before using this product.



Do not use the massager in damp rooms (e.g. when bathing or showering)



- Do not put in the washing machine.
- Do not dry clean!

USER MANUAL INSERTS FOR PART 15

Warning: Changes or modifications to this unit not expressly approved by the party responsible for compliance could void the user's authority to operate the equipment.

NOTE: This equipment has been tested and found to comply with the limits for a Class B digital device, pursuant to Part 15 of the FCC Rules. These limits are designed to provide reasonable protection against harmful interference in a residential installation. This equipment generates, uses and can radiate radio frequency energy and, if not installed and used in accordance with the instructions, may cause harmful interference to radio communications.

However, there is no guarantee that interference will not occur in a particular installation. If this equipment does cause harmful interference to radio or television reception, which can be determined by turning the equipment off and on, the user is encouraged to try to correct the interference by one or more of the following measures:

- Reorient or relocate the receiving antenna.
- Increase the separation between the equipment and receiver.
- Connect the equipment into an outlet on a circuit different from that to which the receiver is connected.
- Consult the dealer or an experienced radio/TV technician for help.

This device complies with Part 15 of the FCC Rules. Operation is subject to the following two conditions:


1. this device may not cause harmful interference, and
2. this device must accept any interference received, including interference that may cause undesired operation.

USER MANUAL FOR IC

This device contains licence-exempt transmitter(s)/receiver(s) that comply with Innovation, Science and Economic Development Canada's licence exempt RSS(s). Operation is subject to the following two conditions:

1. This device may not cause interference.
2. This device must accept any interference, including interference that may cause undesired operation of the device.

PRODUCT DIAGRAM

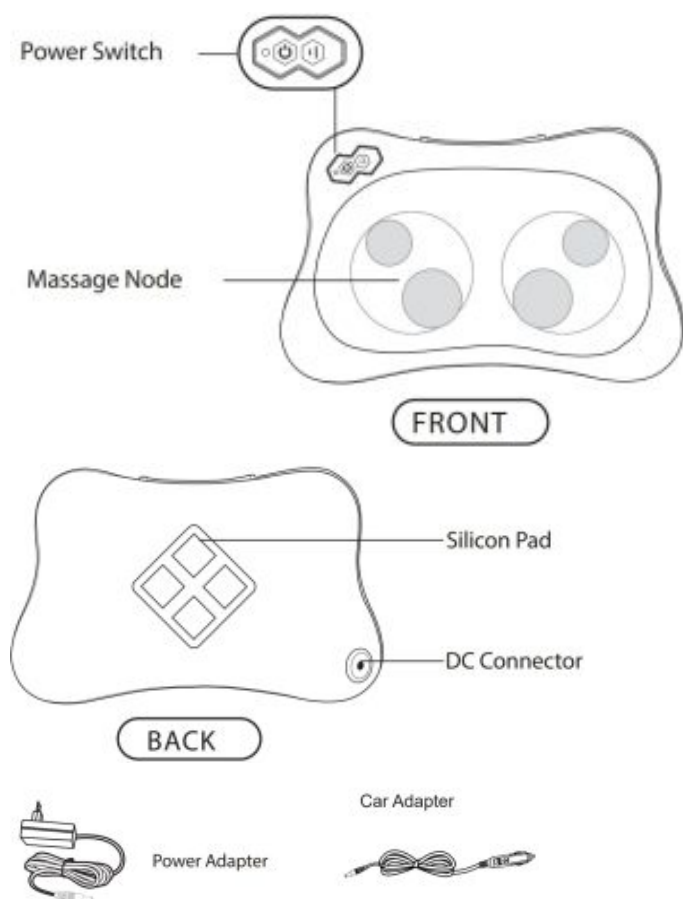

OPERATING INSTRUCTIONS

Ensure the device is charged and ready for use.



Power Switch Button:

Press it once, the product will start its kneading massage and heat therapy for 20 minutes. Press it again, the product will stop its heat therapy while the kneading massage continues. Press for 3 seconds to stop massage session and all functions.



Speed Button:

When the product is turned on, the speed will be on the highest level. Press once to adjust to a lower speed and continue pressing until the desired speed is reached.

Bluetooth :

The GLOW massage pillow has a Bluetooth function that can connect to the Addaday App. When the Bluetooth is successfully connected, the indicator light will be ON after blinking for 10 seconds.

Using Bluetooth connectivity, the Addaday App enables you to:

1. Turn the heating function ON/OFF.
2. Select the Low/High/Auto massage options.
3. Check your massage device usage history.

Charging:

Only use the approved charger that came with your device. The battery light will glow red while charging and

green once the battery is fully charged. When the product is not use, it should be unplugged.

MASSAGE POSITION



Cleaning & Maintenance

To clean and maintain the product

Cleaning

Before cleaning, be sure that the power switch is in the off position and the device is unplugged from its power source. Use a dry cloth to clean the product.

Maintenance

Keep the product clean and away from damp or wet environments.

Attention

Please do not use benzene or thinner and other aggressive detergents to clean the product or it may fade.

PRODUCT SPECIFICATIONS

Item Code:	BIO68095
Product Name	GLOW
Rated Input:	13V 1.8A
Power Consumption	23.4W
Default Time	20 mins
Battery Specifications	11.1V 1800mAh


Battery Disposal

REMOVING THE BUILT-IN RECHARGEABLE BATTERY

Remove the built – in battery when disposing of the massager.

- Please make sure that the battery is disposed of at a qualified disposal location. The instructions and graphic below are intended only for guidance on massager disposal; it is not intended to be used for repairs. Dismantling the appliance yourself is not recommended as it may cause damage or malfunctions and will void the warranty.
Unplug the AC adapter from the appliance before disassembling it.
- Turn on the power and then keep the power on until the battery is completely discharged.

- Perform steps ① to ② and lift the battery, and then remove it.

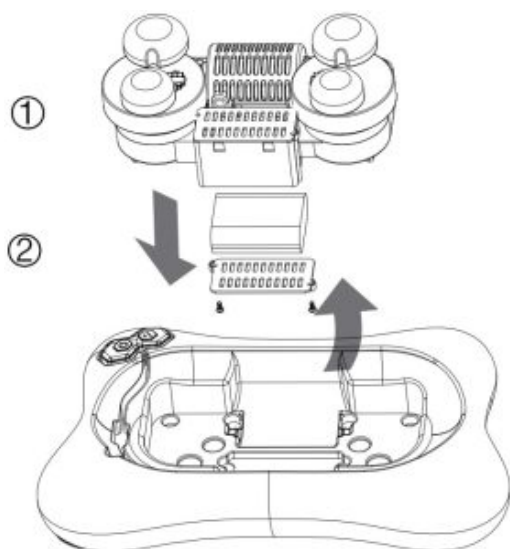

TROUBLESHOOTING

Problem	Cause	Solution
Unable to start the device.	The cord is not connected properly.	Push in the cord harder until you heard a “click”.
	The device is plugged in and charging.	Unplug the device from power source.
The appliance stops during use.	The battery has run down.	Recharge with the AC adaptor.
	The appliance automatically stops after approximately 20 minutes.	This does not indicate a problem.
The main unit does not heat up.	The appliance may not work beyond the recommended temperature.	Use under the recommended temperature of 59 °F – 95 °F

If you still have problems after trying the solutions above, please contact addaday@addaday.com

Correct Disposal of this product



This marking indicates that this product should not be disposed with other household wastes throughout the EU. To prevent possible harm to the environment or human health from uncontrolled waste disposal, recycle it responsibly to promote the sustainable reuse of material resources. To return your used device, please use the return and collection systems or contact the retailer where the product was purchased. They can take this product for environmental safe recycling.

Documents / Resources

	addaday Bio6 Glow Massager Pillow with Bluetooth [pdf] User Manual WSL80233, 2A2WO-WSL80233, 2A2WOWSL80233, Bio6, Glow Massager Pillow with Bluetooth, Bio6 Glow Massager Pillow with Bluetooth
-----------------------------------------------------------------------------------	--------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------