

ADLER AD 8151 Analog Personal Scale User Manual

Home » ADLER » ADLER AD 8151 Analog Personal Scale User Manual



Contents

- 1 ADLER AD 8151 Analog Personal Scale
- **2 Product Usage Instructions**
- **3 DESCRIPTION OF THE DEVICE**
- **4 OPERATION**
- **5 CLEANING AND MAINTENANCE**
- **6 TECHNICAL DATA**
- 7 FAQs
- 8 Documents / Resources
 - 8.1 References



ADLER AD 8151 Analog Personal Scale



Specifications

Capacity: 130 kgAccuracy: 1000g

Product Information

The AD 8151 is a digital scale with the following components:

- · A. Display
- B. Scale platform
- C. Calibration knob

Product Usage Instructions

Calibration:

Before first use, calibrate the scale using knob C. Make sure the red line indicates 0.

Operation:

- 1. Place the scale on a hard and flat surface, avoiding carpets and soft surfaces.
- 2. Step softly on the scale, stand evenly without moving.
- 3. To read the measurement, look at the point marked with a red line when you are on the scale.
- 4. To recalibrate, use knob C to ensure the red line indicates 0.

Cleaning and Maintenance

1. Clean the scale using a moist cloth and wipe it dry.


- 2. Do not use chemical or corrosive cleaning agents.
- 3. Store the scale in a cool and dry place in a horizontal position.

SAFETY CONDITIONS. IMPORTANT INSTRUCTIONS ON SAFETY OF USE PLEASE READ CAREFULLY AND KEEP FOR FUTURE REFERENCE

- 1. Before using this appliance read the user's manual and follow the instructions for use. The manufacturer shall not be held reliable for damage caused by misuse or improper operation of the appliance.
- 2. This appliance is intended for household use only. Do not use it for purposes other than the intended one.
- 3. Take special care when using the appliance near children. Do not allow children to use the appliance as a toy. Do not allow children or other person not familiarized with the appliance to operate it.
- 4. WARNING: This appliance is not intended for use by children under the age of 8 or persons with reduced physical or mental capabilities, or lack of experience and knowledge of the equipment, unless they have been given supervision or instruction concerning the use of the equipment and the related hazards. Children should not play with the equipment. The equipment must not be cleaned or maintained by the children if not supervised.
- 5. To avoid danger of suffocation keep the bag away from babies and children. Do not use the bag in cribs, beds, carriages or play pens. It is not a toy. Knot before throwing away.
- 6. Do not immerse the appliance in water or any other liquid. Do not immerse the appliance in water or other liquid. Do not expose the appliance to weather conditions (rain, sunlight etc.).
- 7. Do not operate any appliance if damaged or malfunctioning in any manner. Do not try to repair the appliance yourself. Return the appliance to a service centre for examination or repair. Any repairs can be carried out by an authorized service facility only. If repaired improperly the appliance can pose serious danger to the user.
- 8. Put the appliance on a cold, stable and even surface away from kitchen appliances such an electric cooker, gas burner etc
- 9. On not operate the appliance near inflammable materials.
- 10. Do not use the appliance near water e.g. under the shower, in a bathtub or above a washbasin with water inside.
- 11. Use the appliance on a stable and flat floor only.
- 12. Keep the appliance away from strong electromagnetic fields.
- 13. Make sure the surface the appliance is used on is not wet or slippery.
- 14. Take care not to drop heavy objects on the display or the glass part of the scale.
- 15. DO NOT hit, DO NOT shake and DO NOT drop the scale.
- 16. Stand still on the scale during the whole weighing time.
- 17. Do not use for commercial or medical purposes.

DESCRIPTION OF THE DEVICE

- · A. display
- · B. scale platform
- · C. calibration knob

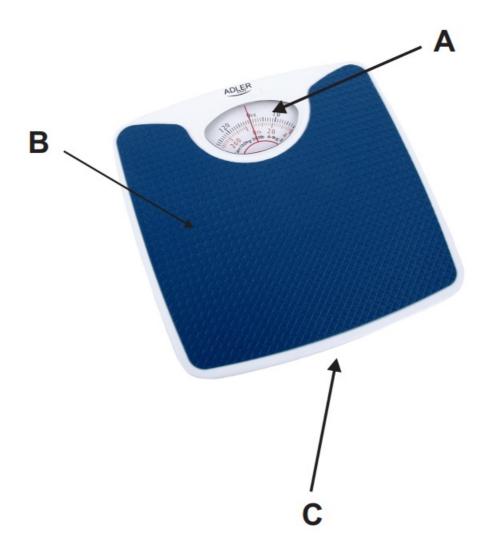

BEFORE FIRST USE

Calibrate the scale using knob C. The red line should indicate "0".

OPERATION

- STEP 1. Place the scale on hard and flat surface (avoid carpets and soft surfaces).
- STEP 2. Softly step on the scale. Stand evenly on the scale, and do not move.
- STEP 3. To read the measurement, look at the point marked with a red line when you are on the scales. NOTE: To calibrate the scales, use knob C. (The red line should indicate "0")

CLEANING AND MAINTENANCE

- 1. Clean the scales using moist cloth. Wipe the scale dry. DO NOT use chemical/corrosive cleaning agents.
- 2. Store the scale in cool and dry place.
- 3. The scale should always be in horizontal position.

TECHNICAL DATA

- · Capacity 130 kg
- · Accuracy: 1000g

To protect your environment: please separate carton boxes and plastic bags and dispose them in corresponding waste bins. Used appliance should be delivered to the dedicated collecting points due to hazarsous components, which may effect the environment. Do not dispose this appliance in the common waste bin.

AFTER SALE SERVICE NOTICE

If you want to buy spare parts or make any complaints, please contact the seller who issued the receipt directly.

FAQs

Q: How do I recalibrate the scale?

A: Use knob C to adjust the calibration so that the red line indicates 0.

Q: What should I do if the scale is not showing accurate measurements?

A: Ensure the scale is on a flat surface, recalibrate using knob C, and stand evenly without moving for accurate readings.

Q: Can I use chemical cleaning agents on the scale?

A: No, only clean the scale using a moist cloth and avoid chemical or corrosive cleaning agents.

Documents / Resources



ADLER AD 8151 Analog Personal Scale [pdf] User Manual

AD 8151, AD 8151 Analog Personal Scale, AD 8151, Analog Personal Scale, Personal Scale, S cale

References

User Manual

Manuals+, Privacy Policy

This website is an independent publication and is neither affiliated with nor endorsed by any of the trademark owners. The "Bluetooth®" word mark and logos are registered trademarks owned by Bluetooth SIG, Inc. The "Wi-Fi®" word mark and logos are registered trademarks owned by the Wi-Fi Alliance. Any use of these marks on this website does not imply any affiliation with or endorsement.