

[Manuals.plus](#) /

› [Waterpik](#) /

› Waterpik FTW-01 Whitening Flosser Tips Instruction Manual

Waterpik FTW-01

Waterpik FTW-01 Whitening Flosser Tips Instruction Manual

INTRODUCTION

Welcome to the instruction manual for your Waterpik FTW-01 Whitening Flosser Tips. These tips are designed to enhance your oral hygiene routine by providing effective flossing and teeth whitening benefits when used with compatible Waterpik Power Flossers. Please read this manual carefully to ensure proper use and maintenance of your flosser tips.

PRODUCT OVERVIEW

The Waterpik FTW-01 Whitening Flosser Tips are disposable replacement tips for compatible Waterpik Power Flossers. Each tip is mint-flavored and coated with silica, an agent that helps brighten teeth while effectively removing plaque and food debris from between teeth and below the gumline.

- **Contents:** This package contains 30 disposable flosser tips.
- **Whitening Feature:** Special mint-flavored tips are coated with a whitening agent containing silica to brighten between teeth.
- **Compatibility:** These tips are compatible with Waterpik flosser models FL-110, FL-310, FLW-110, and FLW-310.

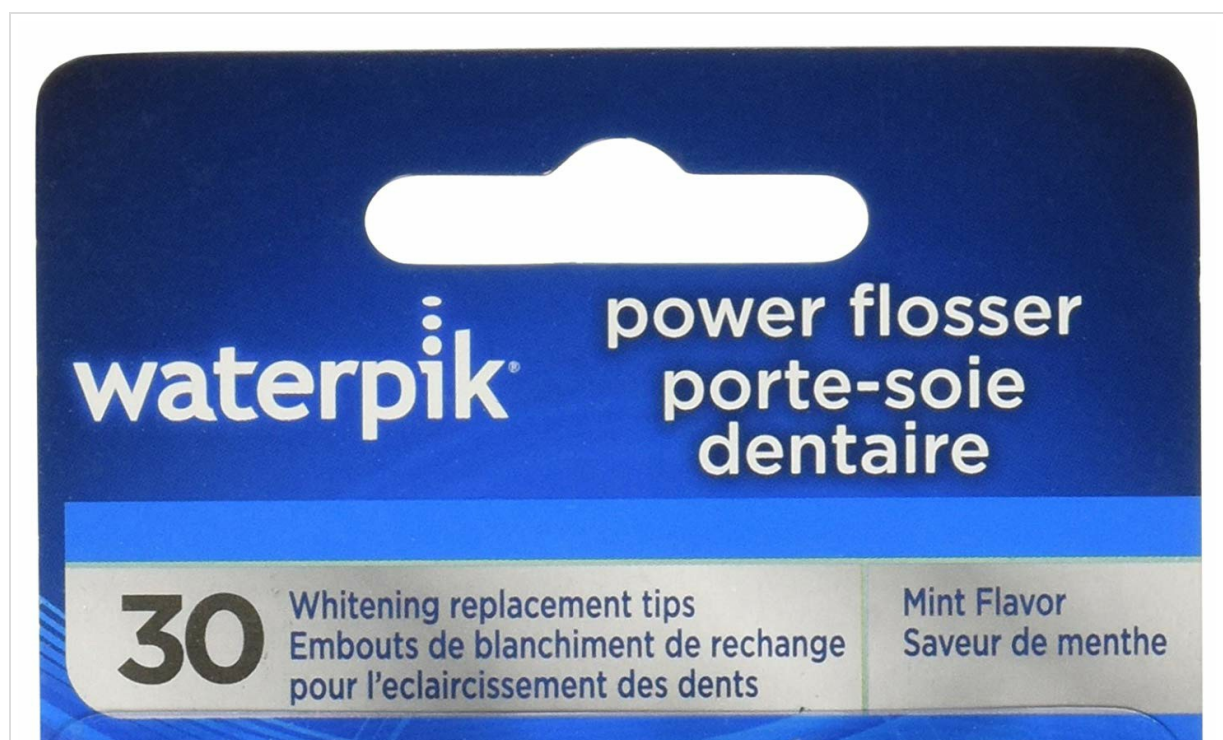




Image showing the retail packaging for Waterpik FTW-01 Whitening Flosser Tips, highlighting the product name, contents (30 disposable tips), and key features like whitening and mint flavor.

SETUP AND INSTALLATION

Replacing Flosser Tips

Changing the Power Flosser tip daily is a simple process to maintain hygiene and effectiveness.

1. Ensure your Waterpik Power Flosser device is turned off before replacing the tip.
2. Gently pull the used flosser tip straight out from the device handle to remove it.
3. Take a new Waterpik FTW-01 Whitening Flosser Tip from its packaging.
4. Align the new tip with the opening on the flosser handle and push it firmly until it clicks securely into place.
5. Dispose of the used tip responsibly in a waste bin.



Changing the Power Flosser tip daily is simple and quick with these convenient cartridges of replacement tips.

Le remplacement quotidien de l'embout de l'appareil est simple et rapide grâce à ces cartouches pratiques d'embouts de remplacement.

Contents: 30 disposable tips.

Contenu: 30 embouts jetables.

Whitening: Special mint-flavored tips are coated with a whitening agent containing silica to brighten between teeth.

Éclaircissement: Les embouts spéciaux à saveur de menthe comportent un agent d'éclaircissement contenant de la silice, qui agit entre les dents.

For Best Results: Depending on your tooth spacing, using more than one Whitening Tip per flossing may be necessary.

Pour de meilleurs résultats: Selon l'espace interdentaire, l'emploi de plusieurs embouts d'éclaircissement peut être nécessaire.

Water Pik, Inc.
1730 East Prospect Road
Fort Collins, CO 80553-0001
www.waterpik.com

Waterpik International, Inc.
Castle Court
4 London Road, Reigate
Surrey, RH2 9RJ, UK

Designed in U.S.A. Made in China. Printed in China.
Conçu aux É.-U. Fabriqué en Chine. Imprimé en Chine.

The Green Dot® is a registered trademark of
Duales System Deutschland and PRO-Europe.
The Green Dot® est une marque déposée de
Duales System Deutschland et PRO-Europe.



©2010 Water Pik, Inc.

Waterpik® (stylized) is a trademark of Water Pik, Inc. registered in Canada, EU, Japan, Mexico, Switzerland and the United States.

©2010 Water Pik, Inc.

Waterpik (stylisé) est une
marque de commerce de Water Pik, Inc.
enregistrée aux États-Unis, au Canada
dans l'Union Européenne, au Japon,
au Mexique et en Suisse.

Model FTW-01
MODÈLE FTW-01
FN: 20014225-F AA



OPERATING INSTRUCTIONS

Daily Flossing Procedure

For optimal results and to maximize the whitening benefits, follow these steps when using your Waterpik Power Flosser with FTW-01 tips:

1. After installing a fresh tip, turn on your Waterpik Power Flosser device.
2. Guide the flosser tip gently between your teeth, ensuring it reaches the gum line.
3. Move the tip slowly along each tooth surface, allowing the gentle vibrations to dislodge food particles and plaque.
4. For best results, especially with varying tooth spacing, it may be necessary to use more than one Whitening Tip per flossing session to ensure thorough cleaning and whitening.
5. After completing your flossing routine, turn off the device.

MAINTENANCE

Tip Disposal

Waterpik FTW-01 Whitening Flosser Tips are designed for single or limited use. For optimal hygiene and effectiveness, it is recommended to change the tip daily or as needed, especially if the tip becomes bent, worn, or loses its mint flavor. Dispose of used tips in a standard waste bin.

TROUBLESHOOTING

Tip Not Fitting Securely

If a tip does not fit securely onto your Waterpik Power Flosser, ensure it is correctly aligned with the flosser handle and pushed in firmly until it clicks. If the issue persists, inspect both the tip and the flosser handle for any visible damage. Replace the tip if damaged.

Reduced Whitening Effect

The whitening effect provided by these tips is gradual and cumulative. Consistent daily use is recommended for noticeable results. Ensure the tips are fresh, as the whitening agent's effectiveness may diminish with prolonged exposure or multiple uses. For significant whitening concerns, consult your dental professional.

SPECIFICATIONS

Model Number	FTW-01
Contents	30 Disposable Tips
Feature	Whitening, Mint-flavored
Compatibility	Waterpik Power Flosser models FL-110, FL-310, FLW-110, FLW-310
Packaging Dimensions	1.5 x 1.5 x 3.3 inches
Packaging Weight	1.1 Pounds
Manufacturer	Water Pik, Inc.
Origin	Designed in U.S.A., Made in China

For more information and detailed product specifications, visit the official [Waterpik website](#).

WARRANTY AND SUPPORT

Waterpik FTW-01 Whitening Flosser Tips are consumable items. For information regarding product warranty, returns, or technical support related to your Waterpik Power Flosser device, please refer to the documentation provided with your specific flosser model or visit the official Waterpik website.

Contact Information:

- **Water Pik, Inc. (USA)**
 - 1730 East Prospect Road, Fort Collins, CO 80553-0001, USA
- **Waterpik International, Inc. (UK)**
 - Castle Court, 4 London Road, Reigate, Surrey, RH2 9RJ, UK
- **Website:** www.waterpik.com

The Green Dot® is a registered trademark of Duales System Deutschland and PRO-Europe. Waterpik® (stylized) is a trademark of Water Pik, Inc. registered in Canada, EU, Japan, Mexico, Switzerland and the United States.