

GIGI 953

GIGI Digital Paraffin Wax Warmer Instruction Manual

Model: 953

INTRODUCTION

Thank you for choosing the GIGI Digital Paraffin Wax Warmer. This professional-grade paraffin bath is designed to provide a luxurious and therapeutic experience for hands, feet, and elbows. This manual provides essential information for the safe and effective operation, maintenance, and care of your device. Please read all instructions carefully before initial use.

IMPORTANT SAFETY INFORMATION

WARNING: To reduce the risk of burns, electrocution, fire, or injury to persons:

- Always unplug the unit from the electrical outlet immediately after use and before cleaning.
- Do not reach for an appliance that has fallen into water. Unplug it immediately.
- Do not place or store the appliance where it can fall or be pulled into a tub or sink. Do not place in or drop into water or other liquid.
- Close supervision is necessary when this appliance is used by, on, or near children or individuals with certain disabilities.
- Use this appliance only for its intended use as described in this manual. Do not use attachments not recommended by the manufacturer.
- Never operate this appliance if it has a damaged cord or plug, if it is not working properly, or if it has been dropped or damaged.
- Keep the cord away from heated surfaces.
- Never use while sleeping or drowsy.
- Do not use outdoors or operate where aerosol (spray) products are being used or where oxygen is being administered.
- Do not use the unit without wax.
- Always test the wax temperature on a small area of skin before full application to ensure it is comfortable and safe.

PACKAGE CONTENTS

Your GIGI Digital Paraffin Wax Warmer package typically includes:

- 1 x GIGI Digital Paraffin Wax Warmer Unit
- Paraffin Wax (quantity may vary, e.g., 6 lbs of Peach Paraffin Wax)
- Removable Slanted Grid

Please note that accessories such as plastic hand/foot protectors or mitts are often sold separately and may be required for optimal use.



Image: GIGI Digital Paraffin Wax Warmer shown with bags of peach paraffin wax.

PRODUCT FEATURES

The GIGI Digital Paraffin Wax Warmer is designed for convenience and effectiveness, offering:

- **Digital Temperature Control:** Allows for precise adjustment of wax temperature for comfort and safety.
- **Oversize Capacity:** Commercial-grade plastic bath can hold up to 8 lbs of paraffin, suitable for full wrist and ankle coverage.

- **See-Through, Spill-Resistant Locking Lid:** Keeps wax contained and allows monitoring of melting progress.
- **Removable Slanted Grid:** Provides ergonomic foot positioning and filters out sediment.
- **Quick Melt Feature:** For fast initial melting of paraffin wax.
- **Cool-to-Touch Exterior:** Unlike traditional metal tubs, enhancing user safety.



A professional-grade paraffin wax warmer designed with all the features you need to give your clients an ultra-luxurious paraffin wax treatment or to enjoy indulgent self-care at home like a beauty pro.

Image: Diagram illustrating key features of the GIGI Digital Paraffin Wax Warmer.

SETUP

1. **Unpack:** Carefully remove all components from the packaging.
2. **Placement:** Place the warmer on a stable, heat-resistant surface. Ensure it is away from water sources and out of reach of children.
3. **Insert Grid:** Place the removable slanted grid inside the warmer.
4. **Add Wax:** Place the desired amount of paraffin wax into the warmer. The unit can hold up to 8 lbs of wax.
5. **Connect Power:** Plug the power cord into a suitable electrical outlet.

OPERATING INSTRUCTIONS

Follow these steps for a relaxing paraffin treatment:

1. **Melt the Wax:** Turn on the warmer and set the desired temperature using the digital controls. The quick melt feature can be used for faster initial melting. Allow sufficient time for the wax to fully melt to a smooth, liquid consistency.
2. **Test Temperature:** Before applying to skin, always test a small amount of wax on your wrist to ensure the temperature is comfortable and safe.
3. **Prepare Skin:** Clean and dry the area to be treated (hands, feet, or elbows). Apply a light, non-greasy moisturizer for easier wax removal.
4. **Dip & Coat:** Slowly dip the hand, foot, or elbow into the melted wax. Remove and allow the wax to set for a few seconds. Repeat this process 2-3 times to build a thick, even coating.
5. **Wrap & Relax:** After the final dip, immediately cover the waxed area with a plastic protector (if using) and then wrap it with a warm towel or insulated mitt to retain heat. Relax for 10-15 minutes.
6. **Peel Away:** Gently peel the hardened wax off the skin. The wax should come off in one or a few large pieces. Discard the used wax.



Image: Visual guide for the paraffin wax treatment process.

The GiGi Digital Paraffin Wax Warmer is ideal for treating hands, wrists, feet, ankles, elbows, and other areas requiring moisturizing and therapeutic heat.

HOW TO USE:



1. Prep Your Skin – Clean and dry the area. Apply a light moisturizer for easier wax removal.

2. Melt the Wax – Plug in the warmer, add wax, and wait for it to melt to a silky, warm consistency.

3. Dip & Coat – Dip your hand or foot into the wax, remove, and repeat 2-3 times for a thicker layer.

4. Wrap & Relax – After the last dip, wrap the area with plastic wrap or a warm towel. Relax for 10-15 minutes.

5. Peel Away – Gently remove the wax to reveal soft, moisturized skin.



Image: Examples of areas suitable for paraffin treatment.

Demonstration Video

Your browser does not support the video tag.

Video: A detailed demonstration of how to use the GiGi Paraffin Bath, including wax melting, dipping, and removal, highlighting the benefits of the treatment.

MAINTENANCE

- **Cleaning:** Unplug the unit and allow it to cool completely before cleaning. Use a soft, damp cloth to wipe the exterior. Do not immerse the unit in water.
- **Wax Replacement:** Over time, the paraffin wax may become contaminated or its therapeutic properties may diminish. Replace the wax as needed.
- **Storage:** Store the unit in a cool, dry place with the lid securely closed when not in use.

TROUBLESHOOTING

Problem	Possible Cause	Solution
---------	----------------	----------

Problem	Possible Cause	Solution
Unit does not turn on.	Not plugged in; power outlet malfunction; unit malfunction.	Ensure unit is securely plugged in. Test outlet with another appliance. If still not working, contact customer support.
Wax not melting or melting slowly.	Insufficient heating time; temperature setting too low; too much wax.	Allow more time for melting. Increase temperature setting. Ensure wax quantity does not exceed capacity. Use the quick melt feature for initial melting.
Wax is too hot.	Temperature setting too high.	Lower the temperature setting. Allow wax to cool slightly before use. Always test temperature on skin before full application.
Wax does not peel off easily.	Wax layer too thin; skin not properly moisturized.	Apply 2-3 thick layers of wax. Ensure skin is lightly moisturized before dipping.

SPECIFICATIONS

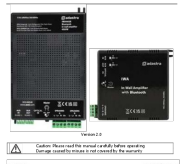
- **Brand:** GIGI
- **Model Number:** 953
- **Product Dimensions:** 19.2 x 12 x 9.3 inches
- **Weight:** Approximately 8 Pounds (empty unit)
- **Bowl Material:** Plastic
- **Capacity:** Up to 8 lbs of paraffin wax
- **UPC:** 073930095309, 884867385650
- **Manufacturer:** American International Industries

WARRANTY AND SUPPORT

For warranty information, technical support, or to inquire about replacement parts, please contact GIGI customer service. Keep your purchase receipt as proof of purchase for any warranty claims. For additional product information and support, visit the official GIGI website or contact your retailer.

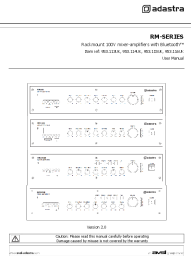


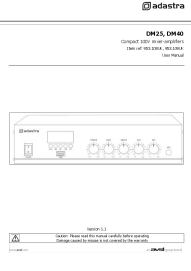
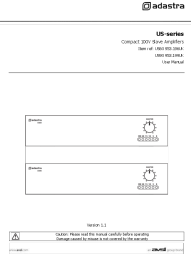
© 2025 GIGI. All rights reserved.

Related Documents - 953



[Adastra IWA-series In-wall Amplifier with Bluetooth User Manual](#)

User manual for the Adastra IWA-series in-wall amplifier with Bluetooth. This guide covers installation, connections, operation, specifications, and safety information for models IWA215B, IWA230B, IWA415B, IWA2100B, and IWA460B.

	<p>Adastra RM-Series Rackmount 100V Mixer-Amplifiers with Bluetooth™ User Manual</p> <p>This user manual provides comprehensive instructions for the Adastra RM-Series rackmount 100V mixer-amplifiers, featuring Bluetooth connectivity. Learn about setup, operation, safety, and specifications for models RM60B, RM120B, RM240SB, and RM360SB.</p>
	<p>ZANE ARTS Inner Tent Footprint & Living Sheet Footprint Manual</p> <p>User manual for ZANE ARTS Inner Tent Footprint and Living Sheet Footprint, covering setup, usage precautions, maintenance, and warranty information.</p>
	<p>Hygger HG-953 Aquarium Magnet Cleaner User Manual</p> <p>User manual for the Hygger HG-953 Aquarium Magnet Cleaner, detailing product features, specifications, application methods, and safety precautions for cleaning freshwater and saltwater tanks.</p>
	<p>Adastra DM25 & DM40 Compact 100V Mixer-Amplifiers User Manual</p> <p>User manual for the Adastra DM25 and DM40 compact 100V mixer-amplifiers. Covers features, setup, operation, media player, Bluetooth, FM tuner, specifications, and troubleshooting.</p>
	<p>Adastra US-series Compact 100V Slave Amplifiers User Manual</p> <p>User manual for Adastra US60 and US90 US-series compact 100V slave amplifiers. This document details installation, connections, operation, safety guidelines, and technical specifications for public address systems.</p>