

Manuals+

[Q & A](#) | [Deep Search](#) | [Upload](#)

manuals.plus /

- › [Veridian Healthcare](#) /
- › [Veridian Healthcare Digital Thermometer Instruction Manual](#)

Veridian Healthcare 08-350

Veridian Healthcare Digital Thermometer Instruction Manual

Model: 08-350

Brand: Veridian Healthcare

INTRODUCTION

Thank you for choosing the Veridian Healthcare 60-Second Digital Thermometer. This manual provides essential information for the safe and effective use of your new thermometer. Please read these instructions thoroughly before use and keep them for future reference.

This digital thermometer is designed for fast and accurate temperature readings, featuring a rigid tip, fever alert, and memory recall. It is suitable for oral, rectal, and underarm use.

IMPORTANT SAFETY INFORMATION

- This product is intended for household use only. It is not intended to diagnose or treat any health problem or disease and should not be used as a substitute for the consultation of a physician or other medical professional.
- Thermometers should only be used under the supervision of an adult.
- Always observe the instructions for use before you begin.
- Do not store the thermometer in direct sunlight, dusty or humid environments, or extreme temperatures.
- This unit is sealed for safety and cannot be replaced. Do not attempt to open the thermometer casing.

PRODUCT OVERVIEW AND COMPONENTS

The Veridian Healthcare Digital Thermometer is designed for ease of use and reliability. Familiarize yourself with its components:

- **Digital Display:** Shows temperature readings.
- **Power Button:** Turns the thermometer on/off and recalls last reading.
- **Rigid Tip:** For accurate and comfortable measurement.
- **Battery Compartment:** Contains a sealed LR41 battery (included).

IN THE BOX

- THERMOMETER
- STORAGE CASE
- BATTERY
- INSTRUCTIONS



Image: The Veridian Healthcare Digital Thermometer and its retail packaging, illustrating the contents included in the box: the thermometer, a storage case, battery (pre-installed), and instructions.





Image: A clear view of the Veridian Healthcare Digital Thermometer, both in its retail packaging highlighting features like 60-second readout and fever alert, and a separate close-up of the thermometer itself.

SETUP

Your Veridian Healthcare Digital Thermometer comes with a sealed LR41 battery pre-installed. No initial battery installation is required.

1. **Clean the Thermometer:** Before first use, and after each subsequent use, clean the thermometer probe as described in the "Maintenance and Cleaning" section.
2. **Prepare for Measurement:** Ensure the thermometer is at room temperature for accurate readings.

OPERATING INSTRUCTIONS

Follow these steps to take a temperature reading:

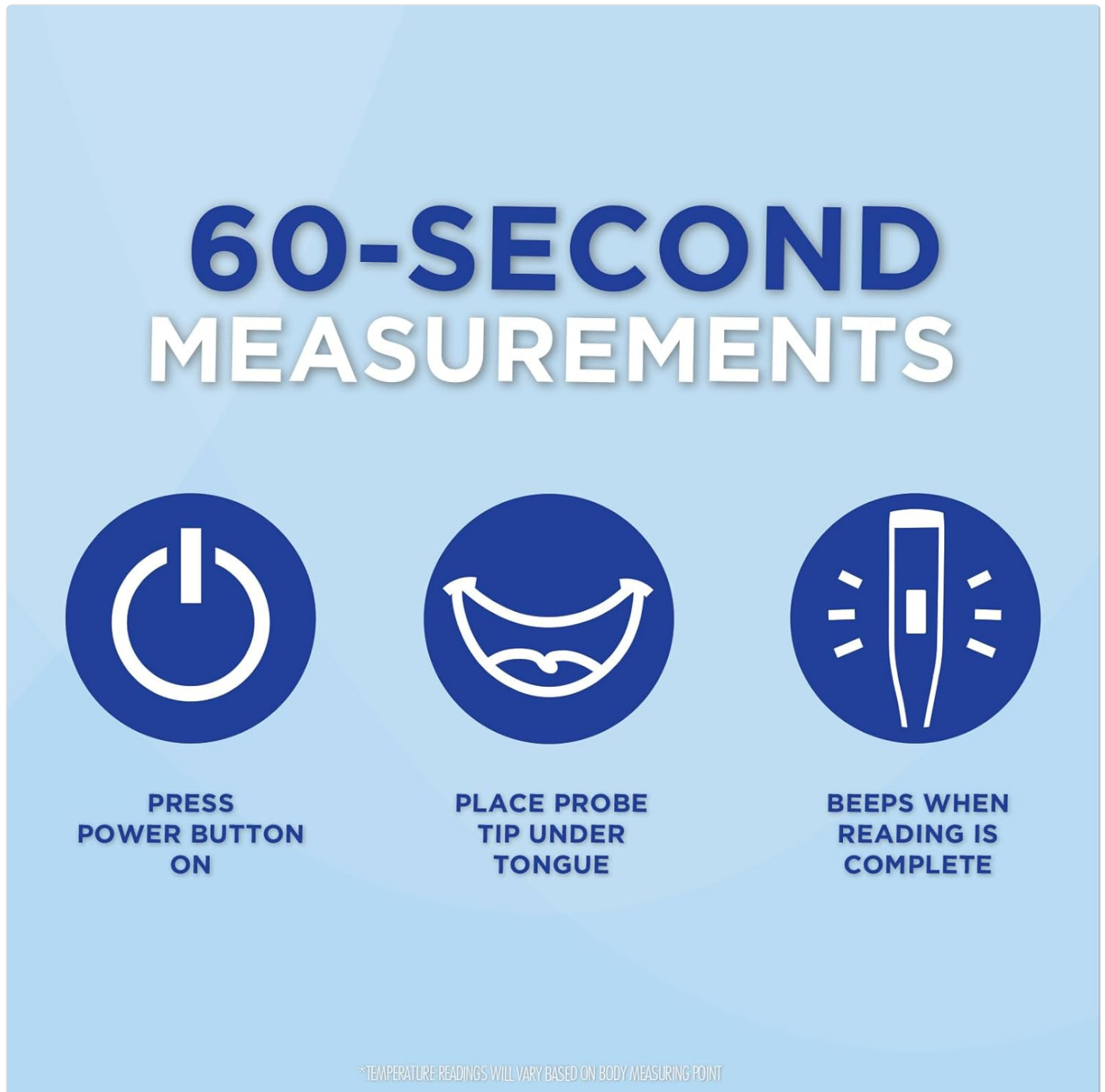


Image: A visual guide demonstrating the three simple steps for taking a temperature: pressing the power button, placing the probe tip under the tongue, and waiting for the beep indicating the reading is complete.

1. **Turn On:** Press the **POWER** button. The last reading will briefly display for 2 seconds, followed by a flashing °F or °C icon, indicating it's ready for a new measurement.
2. **Position the Thermometer:**
 - **Oral Use:** Place the probe tip well under the tongue. Close the mouth and hold the thermometer in place.

- **Rectal Use:** Apply a water-based lubricant to the probe tip. Gently insert the tip no more than ½ inch into the rectum. This method is typically used for infants and young children.
 - **Underarm (Axillary) Use:** Place the probe tip in the armpit so the tip touches the skin. Hold the arm firmly against the body. This method is less accurate and may take longer.
3. **Wait for Reading:** The measurement will take approximately 60 seconds. The °F or °C icon will stop flashing, and a tone will sound for approximately 10 seconds when the peak temperature has been reached.
 4. **Read Temperature:** The temperature will be displayed on the digital screen.
 5. **Fever Alert:** If the reading is 100.0°F (37.7°C) or higher, a fever alert alarm will sound.
 6. **Turn Off:** The unit will automatically shut off after 10 minutes of inactivity to conserve battery life. You can also press the **POWER** button to turn it off manually.

Memory Recall

To recall the last measured temperature, simply press the **POWER** button. The last reading will briefly display for 2 seconds before the thermometer enters measurement mode.

Switching Between Fahrenheit (°F) and Celsius (°C)

While the unit is still on (after pressing the power button and before taking a new reading), hold the **POWER** button down. The unit will switch between °F and °C modes. Release the button when your desired mode is displayed.

MAINTENANCE AND CLEANING

Proper cleaning and maintenance ensure the longevity and accuracy of your thermometer.

- **Cleaning:** The thermometer is water resistant for easy cleaning. Before and after each use, clean the probe tip with 70% isopropyl alcohol diluted with water, or with mild soap and water. Rinse thoroughly with clean water.
- **Disinfection:** For thorough disinfection, wipe the probe with an alcohol swab (70% isopropyl alcohol) or a cotton ball moistened with alcohol.
- **Storage:** Store the thermometer in its hygienic storage half-case in a cool, dry place, away from direct sunlight and extreme temperatures.
- **Battery:** The unit contains a sealed LR41 battery with an approximate life of 600 hours of use. The battery is not user-replaceable.

TROUBLESHOOTING

Problem	Possible Cause	Solution
Thermometer does not turn on.	Battery depleted or unit malfunction.	Since the battery is sealed and not user-replaceable, if the unit does not turn on, it may need to be replaced. Contact customer support.
Inaccurate readings.	Improper placement, movement during measurement, or extreme ambient temperature.	Ensure proper placement (oral, rectal, or underarm) and remain still during measurement. Allow thermometer to acclimate to room temperature before use.
"Lo" or "Hi" displayed.	Temperature is outside the measurable range.	Ensure the thermometer is correctly placed and the body temperature is within the normal range for measurement.

Problem	Possible Cause	Solution
Display shows "Err".	System error.	Turn the thermometer off and then on again. If the error persists, contact customer support.

SPECIFICATIONS

- **Model Number:** 08-350
- **Display Type:** Digital
- **Measurement Time:** Approximately 60 seconds
- **Measurement Range:** Varies by model (typically 90.0°F - 109.9°F / 32.0°C - 42.9°C)
- **Accuracy:** Clinical accuracy (ASTM compliant)
- **Battery:** LR41 (sealed, included)
- **Battery Life:** Approximately 600 hours of use
- **Dimensions:** 5 x 1 x 0.5 inches
- **Weight:** 0.32 ounces
- **Material:** Plastic
- **Special Features:** Fever Alert, Last Reading Memory Recall, Water Resistant, Automatic Shut-Off
- **Included Components:** Digital Thermometer, Storage Case

OUTSTANDING FEATURES



FEATURES

- 60-SECOND READOUT
- RIGID TIP
- PEAK TEMPERATURE TONE
- FEVER ALERT
- LAST READING MEMORY RECALL
- WATERPROOF
- AUTOMATIC SHUT-OFF
- LR41 BATTERY INCLUDED
- ORAL, RECTAL OR UNDERARM USE

*READOUT TIME AND RESULTS WILL VARY BY INDIVIDUAL AND METHOD OF MEASUREMENT. SEE INSTRUCTIONS.

CONTENTS

THERMOMETER, SEALED BATTERY AND STORAGE CASE

CUSTOMER CARE HELP LINE

TOLL-FREE 1-866-326-1313

MON-FRI 8:30 AM-4:30 PM CST



PLEASE READ AND SAVE THESE INSTRUCTIONS

This product is intended for household use; this product is not intended to diagnose or treat any health problem or disease and should not be used as a substitute for the consultation of a physician or other medical professional. Thermometers should only be used under the supervision of an adult. This thermometer should be held during use. Remove the protective case before you begin.

1. The thermometer should be cleaned before and after each use. It can be cleaned with 70% isopropyl alcohol diluted with water or mild soap and water.
 2. Press green POWER button to turn thermometer on. The last reading in memory will briefly display.
 3. "Lo" will appear on the display, along with a flashing °F(°C). You can now take a reading. Measurement should take approximately 60 seconds.
 4. When complete and the peak temperature has been reached, the °F(°C) will cease flashing.
 5. To turn thermometer off, press Power button. Unit will automatically shut off after 10 minutes. To switch between °F or °C mode, while the unit is off press and hold the POWER button. Continue to hold the POWER button, and unit will change the mode and display the new temperature reading mode. Release the POWER button and begin your reading.
- When taking a rectal temperature, use a probe cover and apply water-based lubricant to the probe tip. Do not insert further than 1/2" into the rectum.
 - Do not store in direct sunlight, dusty or humid environments, or extreme temperatures.
 - This unit is sealed for safety; battery cannot be replaced. Battery life ~ 600 hours of use.

MODEL 08-350

MANUFACTURED FOR: VERIDIAN HEALTHCARE
1175 LAKESIDE DRIVE, GURNEE, IL 60031 USA
MADE IN CHINA

92-08001 06/22 VCS10004

©2022 VERIDIAN HEALTHCARE, LLC



Image: The reverse side of the product packaging, detailing key features such as 60-second readout, rigid tip, fever alert, and waterproof design, along with contents, customer care information, and important usage instructions.

WARRANTY

The Veridian Healthcare Digital Thermometer comes with a **1-year limited warranty** from the date of purchase. This warranty covers defects in materials and workmanship under normal use. It does not cover damage resulting from misuse, abuse, accidents, or unauthorized alterations.

For warranty service, please contact Veridian Healthcare Customer Support with your proof of purchase.

CUSTOMER SUPPORT

If you have any questions, concerns, or require assistance with your Veridian Healthcare Digital Thermometer, please contact our customer support team:

Toll-Free: 1-866-326-1313

Hours: Monday - Friday, 8:30 AM - 4:30 PM CST

Manufacturer: Veridian Healthcare LLC

Address: 1175 Larkin Drive, Gurnee, IL 60031 USA

You can also visit the [Veridian Healthcare Store on Amazon](#) for more product information.

DISPOSAL

Please dispose of this product and its battery responsibly according to local regulations. Do not dispose of electronic devices or batteries with household waste. Contact your local waste management authority for information on recycling and disposal options.

© 2024 Veridian Healthcare LLC. All rights reserved.
Product Model: 08-350