

Oral-B 3000

Oral-B PRO 3000 Electric Rechargeable Power Toothbrush User Manual

Model: 3000

INTRODUCTION

The Oral-B PRO 3000 electric rechargeable power toothbrush is designed to provide a superior clean compared to a regular manual toothbrush. Its professionally inspired Cross Action brush head features bristles angled at 16 degrees, and the 3D cleaning action oscillates, rotates, and pulsates to effectively remove plaque. This model includes a pressure sensor to prevent excessive brushing force and offers three distinct brushing modes. An integrated timer assists in achieving the dentist-recommended two-minute brushing duration.

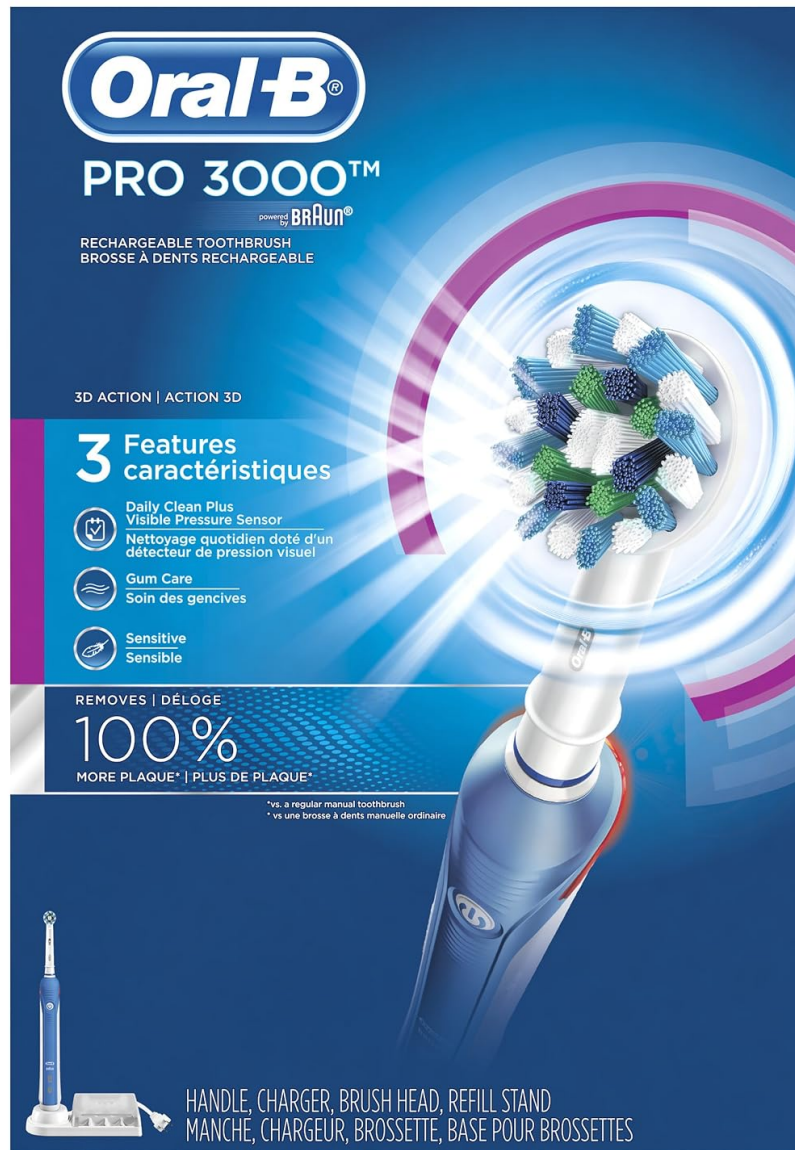


Image: Packaging for the Oral-B PRO 3000 Electric Rechargeable Power Toothbrush, highlighting its 3D Action and three features: Daily Clean Plus Visible Pressure Sensor, Gum Care, and Sensitive mode.

SETUP

- Unpack the Contents:** Carefully remove the toothbrush handle, brush head, charging station, and brush head stand from the packaging.
- Attach the Brush Head:** Align the brush head with the metal shaft on the toothbrush handle and push it down until it clicks into place.
- Initial Charge:** Place the toothbrush handle on the charging station. For optimal battery life, charge the toothbrush for at least 24 hours before its first use. The charging indicator light will show the charging status.



Image: The Oral-B PRO 3000 electric toothbrush handle placed on its charging station, with a brush head attached.

OPERATING INSTRUCTIONS

To use your Oral-B PRO 3000 electric toothbrush, follow these steps:

1. **Wet the Brush Head:** Moisten the brush head and apply a small amount of toothpaste.
2. **Power On:** Press the power button to turn on the toothbrush.
3. **Select Mode:** Press the mode button (if available, or cycle through modes with the power button on some models) to choose your desired brushing mode:
 - **Daily Clean:** Standard mode for everyday cleaning.
 - **Gum Care:** Gentle massage for gums.
 - **Sensitive:** Lower speed for sensitive areas.
4. **Brush Your Teeth:** Guide the brush head slowly from tooth to tooth. Hold the brush head in place for a few seconds on each tooth surface before moving to the next.
5. **Pressure Sensor:** The visible pressure sensor will illuminate if you are brushing too hard. Reduce pressure to protect your gums and teeth.
6. **Built-in Timer:** The in-handle timer pulses every 30 seconds to indicate when to move to a different

quadrant of your mouth. A longer pulse indicates the completion of the dentist-recommended 2-minute brushing time.

7. **Power Off:** Press the power button again to turn off the toothbrush after brushing.

Your browser does not support the video tag.

Video: An animation demonstrating the Oral-B electric toothbrush's cleaning action compared to a manual toothbrush, showing its effectiveness in plaque removal.

Your browser does not support the video tag.

Video: A short explanation detailing the benefits of Oral-B's round toothbrush head design for effective cleaning.

MAINTENANCE

- **After Each Use:** Rinse the brush head thoroughly under running water with the handle switched on. Remove the brush head and clean both the handle and the brush head separately. Wipe the handle dry.
- **Cleaning the Charging Station:** Unplug the charging station before cleaning. Wipe it with a damp cloth. Do not immerse the charging station in water.
- **Brush Head Replacement:** Dentists recommend replacing your brush head every 3 months, or sooner if the bristles become splayed or worn. Indicator bristles will fade halfway to signal when to replace your brush head.



Image: Illustration of the in-handle timer feature on the Oral-B PRO 3000 electric toothbrush, indicating brushing time.

TROUBLESHOOTING

Problem	Possible Cause	Solution
Toothbrush does not turn on.	Battery is discharged.	Place the toothbrush on the charging station for at least 24 hours.
Toothbrush does not charge.	Charging station not plugged in or faulty.	Ensure the charging station is securely plugged into a working outlet. Check for damage to the charger.
Pressure sensor light illuminates frequently.	Brushing too hard.	Reduce the pressure applied while brushing. The toothbrush is designed to clean effectively with minimal pressure.
Reduced brushing performance.	Worn brush head or low battery.	Replace the brush head if bristles are splayed. Ensure the toothbrush is fully charged.

SPECIFICATIONS

- **Model Number:** 3000
- **Brand:** Oral-B (Powered by Braun)
- **Power Source:** Rechargeable Battery
- **Brushing Modes:** 3 (Daily Clean, Gum Care, Sensitive)
- **Special Features:** Visible Pressure Sensor, In-handle Timer
- **Dimensions:** Approximately 4 x 6.88 x 10 inches
- **Weight:** Approximately 1.04 Pounds
- **Bristle Material:** Nylon
- **Handle Material:** Plastic
- **Color:** White and Blue

SAFETY INFORMATION

Please read all instructions before using this appliance. When using electrical products, especially when children are present, basic safety precautions should always be followed, including the following:

- This product is not intended for children 3 years and under.
- Do not use on infants or persons with limited physical, sensory, or mental capabilities unless supervised.
- Do not use if the product appears damaged.
- Do not modify or repair the product. This may result in fire, electric shock, or injury.
- Keep the charging unit away from water. Do not place or store the charger where it can fall or be pulled into a tub or sink.
- Unplug the charger during storms or when not in use for extended periods.
- Only use this product for its intended household use as described in this manual.

WARRANTY AND SUPPORT

For warranty information or product support, please contact Oral-B customer service. You can find contact details on the official Oral-B website or refer to the documentation included with your product.

Manufacturer Contact Information: P&G Plaza, Cincinnati, OH 45202

Please retain your proof of purchase for warranty claims.