

[Manuals.plus](#) /

› [Beurer](#) /

› Beurer BG 17 Glass Body Analysis Scale User Manual

## Beurer 769.25

# Beurer BG 17 Glass Body Analysis Scale User Manual

MODEL: 769.25



## 1. Introduction

---

Thank you for choosing the Beurer BG 17 Glass Body Analysis Scale. This device is designed to accurately measure your body weight, body fat percentage, body water percentage, and muscle ratio. It features an easy-to-read XXL display and stores data for up to 10 users.



Image 1.1: The Beurer BG 17 Glass Body Analysis Scale, showcasing its sleek design and large display.

## Important Safety Information

- This scale must not be used by individuals with medical implants (e.g., pacemakers). Its function may be impaired.
- Do not use during pregnancy.
- Always place the scale on a hard, flat surface to ensure accurate measurements. Avoid carpets or uneven floors.

## 2. Setup

### 2.1 Unpacking and Battery Installation

1. Carefully remove the scale from its packaging.
2. Locate the battery compartment on the underside of the scale.
3. Insert the included CR 2032 3V lithium metal battery, ensuring correct polarity (+/-).
4. Close the battery compartment securely.

### 2.2 Initial Placement

For accurate readings, place the Beurer BG 17 scale on a firm, level surface. Soft surfaces like carpets can lead to inconsistent measurements.

## 2.3 User Profile Configuration

The scale allows for 10 user memories to store personal data and track progress. To set up a user profile:

1. Press the 'SET' button to enter setup mode.
2. Use the arrow buttons to select a user memory slot (P1-P10). Press 'SET' to confirm.
3. Enter your gender (male/female). Press 'SET'.
4. Enter your height. Press 'SET'.
5. Select your activity level (1-5, where 1 is low activity and 5 is high activity). Press 'SET'.
6. The scale will save your profile and return to measurement mode.

You can also select the unit of measurement (kg/lb/st) by pressing the unit button on the underside of the scale.

## 3. Operating the Scale

### 3.1 Taking a Measurement

The Beurer BG 17 scale features a convenient Tap-On activation.

1. Gently tap the scale surface with your foot to turn it on. Wait for '0.0' to appear on the display.
2. If you have a user profile, select your profile number using the arrow buttons before stepping on.
3. Step onto the scale barefoot, ensuring your feet are positioned on the brushed stainless steel electrodes. Stand still.
4. The scale will first display your weight, then cycle through your body fat percentage, body water percentage, and muscle ratio.
5. Step off the scale once all measurements have been displayed. The scale will automatically shut off after a few seconds of inactivity.



Image 3.1: Close-up of the XXL display showing a weight reading of 62.4 kg.

## 4. Maintenance

### 4.1 Cleaning

Clean the scale regularly with a damp cloth and a mild detergent if necessary. Do not use abrasive cleaning agents or immerse the scale in water.

### 4.2 Battery Replacement

When the battery indicator appears on the display, replace the CR 2032 3V battery as described in Section 2.1.

### 4.3 Storage

Store the scale in a dry place, away from extreme temperatures and direct sunlight.

## 5. Troubleshooting

- **Inaccurate Readings:** Ensure the scale is placed on a hard, flat, and stable surface. Avoid carpets or uneven flooring. Re-calibrate by stepping on and off the scale once before taking a measurement.
- **Display Error:** If the display shows an error message, try removing and re-inserting the battery. If the issue persists, the battery may need replacement.
- **No Display:** Check if the battery is correctly inserted and has sufficient charge. Replace the battery if necessary.
- **Body Analysis Not Working:** Ensure you are stepping on the scale barefoot and your feet are making good contact with both stainless steel electrodes. Your user profile must be correctly configured.

## 6. Specifications

| Feature                | Specification                                   |
|------------------------|---|
| Model Number           | 769.25  |
| Display                | LCD, XXL digits (40 mm)                         |
| Capacity               | 150 kg (330 lb / 23 st)                         |
| Accuracy               | 100 Grams                                       |
| Electrodes             | Brushed Stainless Steel                         |
| User Memories          | 10  |
| Activity Levels        | 5   |
| Power Source           | 1 x CR 2032 3V Lithium Metal Battery (included) |
| Dimensions (L x W x H) | 30 x 30 x 1 cm                                  |
| Weight                 | 1 Kilogram                                      |

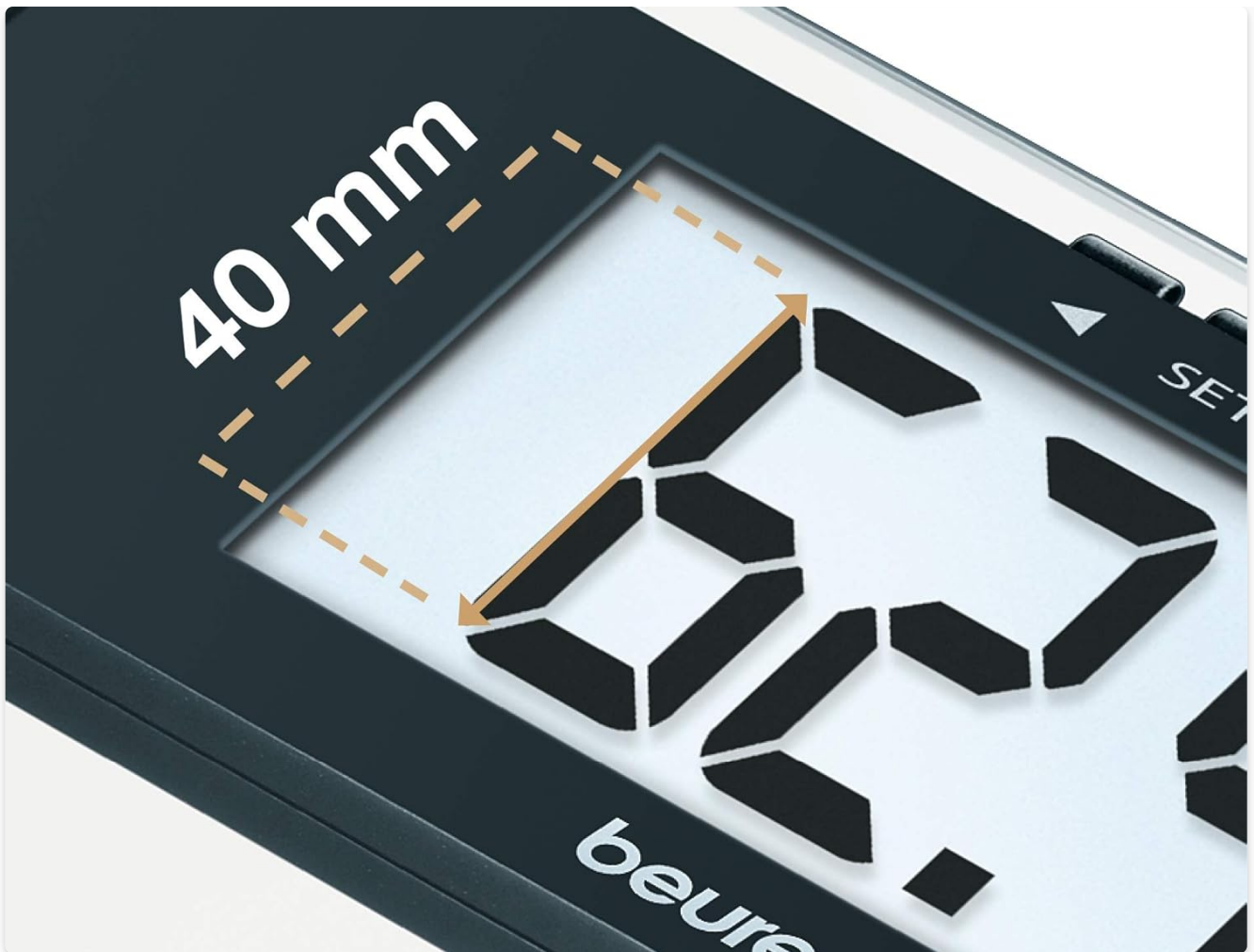


Image 6.1: Detail of the 40mm digit size on the Beurer BG 17 display.

## 7. Warranty and Support

### 7.1 Warranty

The Beurer BG 17 Body Analysis Scale comes with a 5-year warranty. Please retain your proof of purchase for warranty claims.

### 7.2 Customer Support

For further assistance, technical support, or warranty inquiries, please contact Beurer customer service. Refer to the official Beurer website or product packaging for contact details specific to your region.