

Taylor 7407

Taylor Model 7407 Digital Body Weight Scale Instruction Manual

Brushed Stainless Steel Thin Glass Platform with Blue LCD

INTRODUCTION

This manual provides instructions for the proper use and care of your Taylor Model 7407 Digital Body Weight Scale. Please read all instructions carefully before use to ensure accurate measurements and product longevity. This scale is designed to provide precise body weight measurements for personal use.



Image: The Taylor Model 7407 Digital Body Weight Scale with a brushed stainless steel platform and a blue backlit LCD displaying '134.6 lb'.

SAFETY INFORMATION

WARNING: This product contains chemicals known to the State of California to cause cancer and birth defects or other reproductive harm.

Always use the scale on a hard, flat surface. Avoid using on carpet or uneven flooring, as this can affect accuracy. Do not use in excessively wet environments. Keep out of reach of children.

PACKAGE CONTENTS

- Taylor Model 7407 Digital Body Weight Scale
- 3 AAA batteries (included)

SETUP

1. **Battery Installation:** Locate the battery compartment on the underside of the scale. Open the compartment cover and insert the 3 AAA batteries, ensuring correct polarity (+/-). Close the cover securely.
2. **Initial Placement:** Place the scale on a hard, flat, and stable surface. Avoid soft surfaces like carpets or rugs, as they can interfere with accurate readings.
3. **First Use Calibration:** For the first use, or if the scale has been moved, gently press down on the center of the scale with your foot until digits appear on the display. Remove your foot and wait for the display to show '0.0' or '0.00'. The scale is now calibrated and ready for use.



Image: The bottom view of the scale, highlighting the battery compartment and the anti-skid foot pads.

OPERATING INSTRUCTIONS

1. **Step On:** Gently step onto the center of the scale with bare feet. Stand still and distribute your weight evenly.

2. **Read Measurement:** Your weight will appear on the blue backlit LCD display. The display will hold the reading for a few seconds before automatically turning off.
3. **Automatic Shut-off:** The scale features an automatic shut-off function to conserve battery life. It will turn off after a few seconds of inactivity.



Image: A user stepping onto the digital scale to take a weight measurement.

MAINTENANCE

- **Cleaning:** Clean the scale platform with a damp cloth and mild detergent. Do not use abrasive cleaners or immerse the scale in water.
- **Storage:** Store the scale in a dry place at room temperature. Avoid extreme temperatures or high humidity.
- **Battery Replacement:** If the display shows 'Lo' or becomes dim, replace the batteries with new AAA batteries.

TROUBLESHOOTING

Problem	Solution
Scale does not turn on.	Check battery installation and ensure batteries are new and correctly inserted.
Inaccurate readings.	Ensure the scale is on a hard, flat surface. Recalibrate by stepping on and off the scale, waiting for '0.0' before weighing. Avoid moving the scale frequently.
Display shows 'Err'.	This indicates an overload. The scale's maximum capacity is 400 lbs (180 kg). Remove weight immediately to prevent damage.
Display shows 'Lo'.	Batteries are low. Replace with new AAA batteries.

SPECIFICATIONS

- **Model Number:** 7407
- **Capacity:** 400 lbs (180 kg)
- **Increments:** 0.2 lbs
- **Platform Material:** Brushed Stainless Steel, Thin Glass
- **Display:** Blue Backlit LCD (3" x 1.5")
- **Power:** 3 x AAA batteries (included)
- **Product Dimensions:** 13.25 x 1.75 x 13.75 inches
- **Item Weight:** 4 Pounds

LARGE SIZE

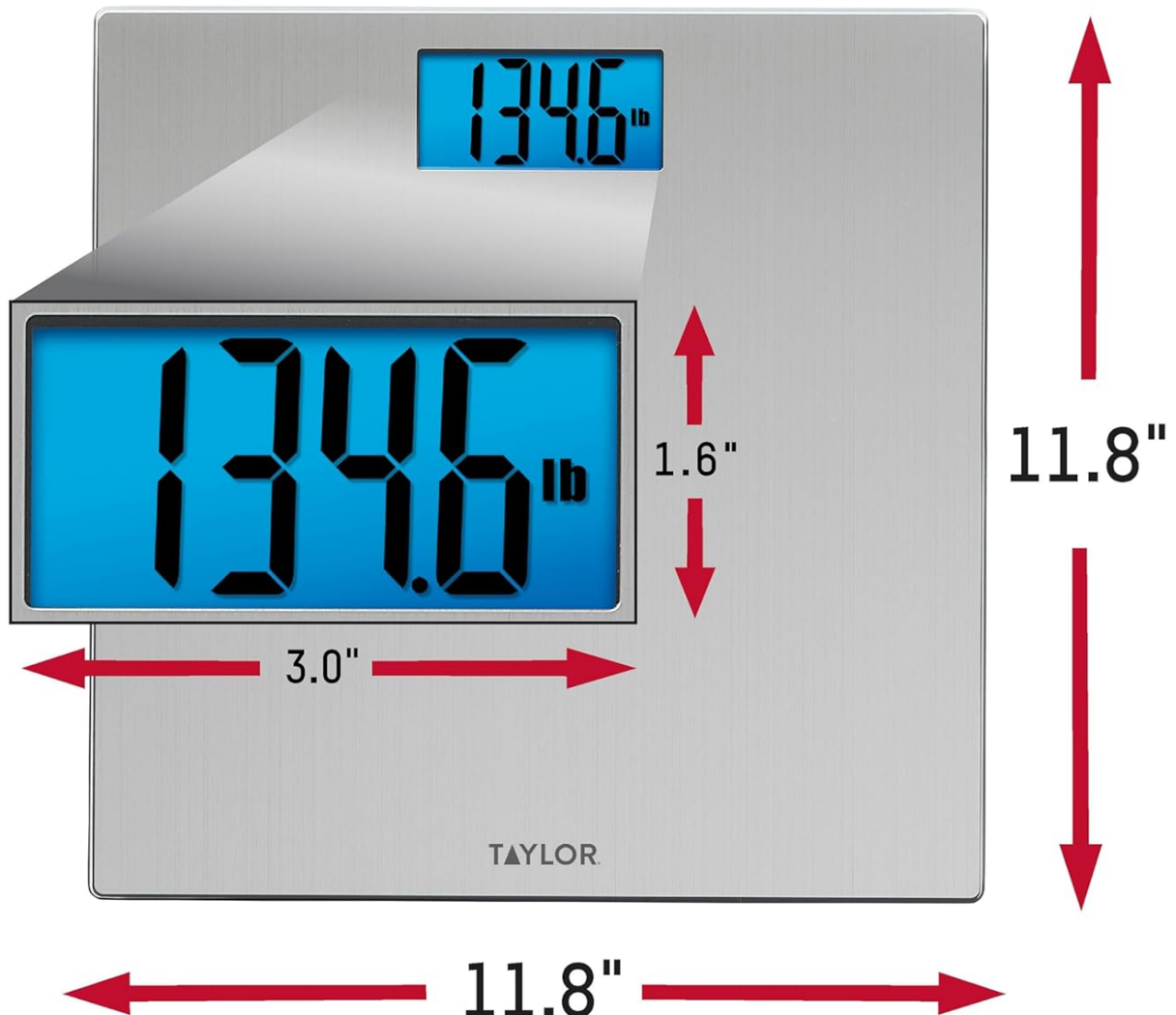


Image: A diagram illustrating the scale's dimensions and display size.

WARRANTY AND SUPPORT

Taylor Precision Products are manufactured to high-quality standards. This product is backed by a limited warranty against defects in material and workmanship. For specific warranty details or technical support, please refer to the warranty card included with your purchase or visit the official Taylor Precision Products website. You may also contact customer service for assistance.

For further information, visit the [Taylor Store](#).



Related Documents - 7407



A quick guide for setting up and using the Taylor 5721F Body Composition Scale, including user profile setup and measurement instructions.



Comprehensive user manual for the Taylor Body Fat Analyzer and Scale (Model 5743). Learn how to use Bioelectrical Impedance Analysis (BIA) technology, set up your scale, take readings, understand results, and find troubleshooting tips and warranty information.



Comprehensive guide for the Taylor 5324506KHL Body Composition Smart Scale. Learn how to set up, use, and troubleshoot your scale, connect with the AFit Health app, and understand its features and warranty.



Quick use guide and frequently asked questions for the Taylor Body Fat Scale (models 5593, 5594). Covers setup, weight-only operation, body fat/water analysis, and memory functions.



Comprehensive guide for the Taylor 5721F Body Composition Scale, covering setup, usage, features, maintenance, troubleshooting, and product specifications. Learn how to measure body fat, body water, and BMI.



This guide provides step-by-step instructions for using the Taylor Model # 5771F Body Composition Scale to measure body fat, body water, and muscle mass. Learn how to set up, enter personal data, and obtain accurate readings.