Manuals+

[Q & A](#) | [Deep Search](#) | [Upload](#)

manuals.plus /

- › [Solac](#) /
- › [Solac ME7711 Sculptural - Professional Lipo-modeling Firming System User Manual](#)

Solac ME7711

Solac ME7711 Sculptural Professional Lipo-modeling Firming System

USER MANUAL

Product Overview

The Solac ME7711 Sculptural is a professional lipo-modeling firming system designed to combat sagging skin and improve circulation. Inspired by professional massage techniques, this device offers an effective therapeutic massage at home, specifically recommended for cellulite, water retention, and tissue relaxation. It features a motor with eccentric rotation, pneumatic pressure, and infrared heat for deep massage.

Figure 1: Solac ME7711 Sculptural device with its interchangeable treatment heads.

The system includes 6 interchangeable heads for different skin and muscle actions across various body zones, enhancing the penetration of anti-cellulite creams.

Important Safety Warnings

Please read all safety instructions carefully before using the device.

- Do not use the device on injured or broken skin.
- Consult a medical professional before use if you have any particular medical conditions.
- Ensure the usage environment is safe and free from obstacles to prevent accidents.
- Keep the device away from water and other liquids.
- Do not use the device for prolonged periods without breaks, as it may overheat and shut off automatically. Allow it to cool down before resuming use.
- Keep out of reach of children.

Setup

1. **Unpack the Device:** Carefully remove the Solac ME7711 Sculptural and all its accessories from the packaging.
2. **Inspect for Damage:** Check the device and all attachments for any signs of damage. Do not use if damaged.
3. **Select an Attachment:** Choose the appropriate interchangeable head for the desired treatment area and purpose.
4. **Attach the Head:** Securely attach the selected head to the device. Ensure it clicks into place or is firmly seated.
5. **Connect to Power:** Plug the power cable into a suitable electrical outlet (220-240V).

Operating Instructions

Follow these steps for effective use of your Solac ME7711 Sculptural device:

1. **Prepare the Skin:** For best results, apply an anti-cellulite cream or gel to the area you intend to treat. This can aid in absorption and improve the massage effect.
2. **Turn On and Adjust Power:** Turn on the device using the power switch. Use the adjustable power dial to select your desired intensity level. Start with a lower setting and gradually increase as comfortable.

Figure 2: Adjusting the power setting on the device.

3. **Apply to Skin:** Place the chosen attachment gently onto the skin of the area to be treated.

4. **Perform Massage:** Move the device in a circular motion over the skin. The eccentric rotation, pneumatic pressure, and infrared heat work together for a deep massage.

Figure 3: Demonstrating the circular eccentric massage technique.

5. **Target Areas:** The device is suitable for various body parts including arms, abdomen, buttocks, thighs, and back.

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Figure 4: Examples of target areas for treatment.

6. **Changing Attachments:** To change an attachment, turn off the device and carefully pull the current head off. Then, firmly push the new head into place until it is secure.

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Figure 5: Easy interchangeability of massage heads.

7. **Session Duration:** Use for recommended session durations (e.g., 10-15 minutes per area). If the device overheats and shuts off, allow it to cool down completely before continuing.

Regular use can lead to visible reduction of cellulite and improved blood circulation.

Figure 6: Potential benefits of consistent use.

Maintenance

- **Cleaning:** After each use, unplug the device and wipe the main unit and attachments with a soft, damp cloth. Do not use abrasive cleaners or immerse the device in water.
- **Storage:** Store the device and its attachments in a cool, dry place, away from direct sunlight and extreme temperatures. Use the provided storage pouch if available.
- **Cord Care:** Do not wrap the power cord tightly around the device. Store it loosely to prevent damage.

Troubleshooting

Problem	Possible Cause	Solution


Problem	Possible Cause	Solution
Device does not turn on.	Not plugged in; power outlet issue; device malfunction.	Ensure the power cord is securely plugged into a working outlet. Try a different outlet. If the problem persists, contact customer support.
Device overheats and shuts off.	Prolonged continuous use; motor protection activated.	This is a safety feature. Turn off the device and allow it to cool down completely (at least 15-20 minutes) before resuming use. Avoid continuous operation for more than 10-15 minutes.
Attachments are difficult to change.	Attachment not aligned correctly; residue build-up.	Ensure the device is off. Apply firm, even pressure when attaching or detaching heads. Clean any residue from the attachment mechanism.
Reduced massage intensity.	Incorrect power setting; worn attachment head.	Adjust the power dial to a higher setting. Ensure the attachment head is clean and not worn out. Replace if necessary.

Specifications

Model Number: ME7711

Brand: Solac

Dimensions (L x W x H): 31 x 21 x 13.2 cm (12.2 x 8.3 x 5.2 inches)

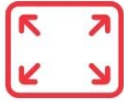**18 x 18 x 14.7
cm**

Figure 7: Device dimensions.

Weight: 1.3 Kilograms (2.87 lbs)

Power Source: Corded Electric (220-240V)

Material: Silicone (for some attachments)

Color: White/Pink

Special Feature: Infrared Heat

Massage Form: Suction cup massage (referring to the action, not literally a suction cup)

Specific Uses: Cellulite treatment, firming, circulation improvement

Interchangeable Heads: 6

Warranty and Support

The Solac ME7711 Sculptural device comes with a **2-year warranty** from the date of purchase.

For technical support, warranty claims, or any questions regarding your product, please contact Solac customer service through their official website or the retailer where you purchased the device. Please have your proof of purchase and

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Solac

[Solac Taste Slim Pro Inox CE4522 Safety Instructions](#)

Comprehensive safety instructions and warnings for the Solac Taste Slim Pro Inox espresso coffee maker (Model CE4522), covering safe usage, maintenance, and environmental considerations.