

## Manuals+

[Q & A](#) | [Deep Search](#) | [Upload](#)

manuals.plus /

- › [Waterpik](#) /
- › [Waterpik Waterflosser NANO WP-250 Oral Irrigator Instruction Manual](#)

## Waterpik WP-250UK

# Waterpik Waterflosser NANO WP-250 Oral Irrigator Instruction Manual

Model: WP-250UK

## 1. INTRODUCTION

The Waterpik Waterflosser NANO WP-250 is an advanced oral irrigator designed to enhance daily oral hygiene. It utilizes a combination of water pressure and pulsations to clean deep between teeth and below the gumline, areas traditional brushing and string floss may miss. Regular use can significantly improve gum health and remove plaque from treated areas.

## 2. PRODUCT OVERVIEW



The Waterpik Waterflosser NANO WP-250 is a compact oral irrigator featuring a main unit with a water reservoir, a handle with a water flossing tip, and a pressure control dial. The reservoir is transparent blue with

measurement markings, and the unit is primarily white with blue accents. The handle connects to the main unit via a coiled hose.

### Included Components:

- Waterpik Waterflosser NANO WP-250 Main Unit
- 1 Classic Jet Tip
- 1 Plaque Seeker Tip



The unit comes with two specialized tips: a Classic Jet Tip for general cleaning and a Plaque Seeker Tip designed for enhanced plaque removal around dental work.

## 3. KEY FEATURES



- **Effective Plaque Removal:** Eliminates up to 99.9% of plaque from treated areas.
- **Improved Gum Health:** Up to 2 times more effective than string floss for improving gum health.
- **Compact Design:** Approximately half the size of larger Waterpik models (e.g., WP100 or WP70), ideal for small bathrooms and travel.
- **Quiet Operation:** Engineered for reduced noise during use.
- **Pressure Control System:** Features 3 pressure settings for customized cleaning.
- **Reversible Reservoir:** The water reservoir can be inverted for compact storage and travel, also helping to keep out dust and debris when not in use.
- **Water Capacity:** Provides 60 seconds of water capacity for a full cleaning session.
- **Power Source:** Operates via a cable, compatible with shaver sockets.

## 4. SAFETY WARNINGS

- **Important:** Do not use this device for massaging or stimulating gums. This model is specifically designed

for plaque removal and gum health improvement through water flossing.

- Always read all instructions before initial use.
- Keep out of reach of children.
- Do not direct water under the tongue, into the ear, nose, or other delicate areas.
- Use only as directed for oral hygiene.

## 5. SETUP INSTRUCTIONS

---

1. **Unpack the Unit:** Remove all components from the packaging.
2. **Place the Unit:** Position the Waterflosser on a flat, stable surface near a power outlet and a sink.
3. **Fill the Reservoir:** Lift the hinged lid on the reservoir and fill it with lukewarm water. You may also add mouthwash as recommended by your dental professional. Place the reservoir back onto the base.
4. **Insert Tip:** Press the tip firmly into the center of the handle until it clicks into place. To remove a tip, press the eject button on the handle.
5. **Prime the Unit (First Use):** If using for the first time, set the pressure control dial to high. Direct the tip into the sink. Turn the unit ON until water flows from the tip. Turn the unit OFF.

## 6. OPERATING INSTRUCTIONS

---

1. **Adjust Pressure:** Select the desired pressure setting using the dial on the base unit. Start with the lowest setting and gradually increase as comfortable.
2. **Position the Tip:** Lean over the sink and place the tip in your mouth. Close your lips enough to prevent splashing, but allow water to flow from your mouth into the sink.
3. **Start Flossing:** Turn the unit ON. Direct the water stream at a 90-degree angle to your gumline. Begin with the back teeth and work your way to the front. Glide the tip along the gumline, pausing briefly between teeth.
4. **Clean All Areas:** Continue until you have cleaned the inside and outside of both your upper and lower teeth. The 60-second reservoir capacity is designed for a thorough cleaning.
5. **Turn Off:** When finished, turn the unit OFF.
6. **Empty Reservoir:** Empty any remaining water from the reservoir after each use to prevent bacterial growth.

## 7. MAINTENANCE

---

### Cleaning the Unit:

- **Exterior:** Wipe the exterior of the unit with a soft cloth and mild, non-abrasive cleaner. Do not immerse the unit in water.
- **Reservoir:** The reservoir can be removed and cleaned with warm soapy water. For thorough cleaning, it is top-rack dishwasher safe. Ensure it is completely dry before reattaching.

### Hygienic Storage:



For hygienic storage and travel, the reservoir can be inverted and placed over the main unit. This compact configuration helps protect the unit from dust and debris when not in use.

Replacing Tips:

- Waterpik recommends replacing flossing tips every 3-6 months, or sooner if the tip becomes damaged.

8. TROUBLESHOOTING

Problem	Possible Cause	Solution
No water flow	Reservoir empty; unit not primed; tip clogged.	Fill reservoir; prime unit (see Setup); remove and clean tip.
Low pressure	Pressure dial set too low; mineral deposits.	Increase pressure setting; descale unit with vinegar solution (refer to full manual for detailed descaling instructions).
Unit leaking	Reservoir not seated correctly; hose or tip damaged.	Ensure reservoir is properly seated; inspect hose and tip for damage and replace if necessary.

9. SPECIFICATIONS

- **Model Number:** WP-250UK
- **Brand:** Waterpik
- **Dimensions (L x W x H):** 22.9 x 15.2 x 15.2 cm (9 x 6 x 6 inches)
- **Weight:** 710 grams (1.56 lbs)
- **Power Source:** Corded Electric
- **Special Feature:** Compact design
- **Water Capacity:** 60 seconds
- **Color:** White






10. WARRANTY AND SUPPORT


The Waterpik Waterflosser NANO WP-250 comes with a **2-year warranty** from the date of purchase. This warranty covers manufacturing defects under normal use. Please retain your proof of purchase for warranty

claims.


For technical support, troubleshooting assistance, or warranty inquiries, please refer to the official Waterpik website or contact their customer service department. Contact information can typically be found in the product packaging or on the manufacturer's website.

## Related Documents - WP-250UK

	<p><a href="#">Waterpik Waterflosser Quick Start Guide: WP-440-466 Models</a></p> <p>Get started quickly with your Waterpik Waterflosser. This guide covers unit preparation, charging, filling the reservoir, and proper water flossing techniques for models WP-440 through WP-466.</p>
	<p><a href="#">Waterpik Waterflosser WP-440 Series Quick Start Guide</a></p> <p>A quick start guide for the Waterpik Waterflosser WP-440 series, covering unit preparation and water flossing instructions. Includes model numbers WP-440, 450, 455, 460, 462, 463, 464, 465, and 466.</p>
	<p><a href="#">Waterpik Waterflosser WP-860 Series Quick Start Guide</a></p> <p>Concise quick start guide for the Waterpik Waterflosser WP-860 series, covering first use, starting water flossing, technique, and how to watch instructional videos.</p>
	<p><a href="#">Waterpik Waterflosser WP-560 Series Quick Start Guide</a></p> <p>Concise quick start guide for the Waterpik Waterflosser WP-560 series, covering unit preparation, charging, and water flossing instructions. Includes model numbers WP-560, WP-562, WP-563, WP-564, WP-565, WP-566, and WP-569.</p>
	<p><a href="#">Waterpik Waterflosser WP-400 Series Quick Start Guide</a></p> <p>Get started quickly with the Waterpik Waterflosser WP-400 Series. This guide provides essential steps for unit preparation, including charging and filling, and instructions on how to effectively use the water flosser for optimal oral hygiene.</p>



Waterpik® Water Flosser  
Model WP-150 Series  
Hydropulseur Waterpik®  
Modèle WP-150 Série




www.waterpik.com

[Waterpik® Water Flosser WP-150 Series User Manual & Guide](#)

Official user manual and guide for the Waterpik® Water Flosser Model WP-150 Series. Learn about important safeguards, product description, getting started, tip usage, cleaning, troubleshooting, and warranty information.

## Documents - Waterpik – WP-250UK

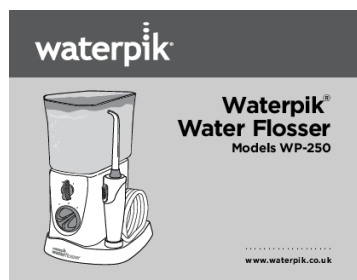

[\[pdf\]](#) User Manual Quick Start Guide Instructions Troubleshooting Guide Warranty  
 Instruction Manual Nano Water Flosser WP 250UK Series Waterpik® Owner s EN Manuals Oral Health  
 Products 20019689 FAA IM rdF2 LR Web waterpik co uk media wysiwyg Waterpik 2 TABLE OF  
 CONTENTS Important Safeguards 3 Product Description 5 Getting Started 7 Cleaning and  
 Troubleshooting Your 9 Limited Two Year Warranty 10 |||  
 Waterpik Water Flosser Models WP-250 O I 1 2 3 ..... www.waterpik.co.uk 2  
 TABLE OF CONTENTS Important Safeguards 3 Product Description 5 Getting  
 Started 7 Cleaning and Troubleshooting Your Waterpik Water Flosser 9 Limited Two-  
 Year Warranty 10 3 IMPORTANT SAFEGUARDS ...  
 lang:en **score:40** filesize: 498.81 K page\_count: 11 document date: 2019-12-12