Omron HBF-214

Omron HBF-214 Body Composition Meter Instruction Manual

Model: HBF-214

1. INTRODUCTION

This manual provides essential instructions for the safe and effective use of your Omron HBF-214 Body Composition Meter. Please read this manual thoroughly before using the device and retain it for future reference. The Omron HBF-214 is designed to measure body weight and estimate various body composition metrics.

2. SETUP

2.1 Unpacking

Carefully remove the body composition meter from its packaging. Ensure all components are present and undamaged. Dispose of packaging materials responsibly.

2.2 Battery Installation

The Omron HBF-214 requires four (4) AAA batteries. Locate the battery compartment on the underside of the device. Insert the batteries according to the polarity indicators (+/-). Close the battery compartment securely.

2.3 Device Placement

For accurate measurements, place the body composition meter on a hard, flat, and stable surface. Avoid placing the device on soft surfaces such as carpets, rugs, or tatami mats, as this can affect measurement accuracy.

2.4 Initial User Registration

Upon first use or after battery replacement, you will need to register user profiles. The device supports up to four user profiles. Follow the on-screen prompts to enter personal data such as age, height, and gender. Please note that the device display may be in Japanese. Refer to the operating instructions in Section 3 for guidance on navigating the interface.

Figure 1: The Omron HBF-214 Body Composition Meter displaying a weight measurement. The device features a clear digital display and control buttons below it.

3. OPERATING INSTRUCTIONS

3.1 Taking a Measurement

1. **Power On:** Gently step onto the device. It will automatically power on.
2. **User Selection:** If multiple users are registered, the device may attempt to auto-recognize you. If not, use the arrow buttons (▲ / ▼) to select your user profile. Press the 'SET' or 'OK' button (often represented by a circle icon) to confirm.
3. **Stand Still:** Stand still with bare feet on the metal electrodes. Ensure your feet are clean and dry for proper contact.
4. **Measurement Process:** The device will begin measuring. Remain still until all measurements are displayed.
5. **Power Off:** The device will automatically power off after a period of inactivity.

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Figure 2: Proper stance for taking a measurement on the Omron HBF-214. Ensure bare feet are placed on the electrodes.

3.2 Understanding Display Readings

The Omron HBF-214 provides various body composition metrics. The display may show icons or abbreviations. Below are common measurements:

- **Weight:** Your current body weight (e.g., 52.3 kg).
- **Body Fat Percentage:** The proportion of fat in your body.
- **Skeletal Muscle Percentage:** The percentage of skeletal muscle in your body.
- **BMI (Body Mass Index):** A measure of body fat based on height and weight.
- **Visceral Fat Level:** An estimation of fat surrounding your internal organs.
- **Body Age:** An estimation of your body's physiological age.
- **Basal Metabolism:** The minimum amount of energy required to maintain basic bodily functions at rest.

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Figure 3: Visual representation of the different body composition metrics measured by the device.

3.3 Comparing Measurements

The device can display previous measurement values, allowing you to track changes over time. Use the 'Previous Value' button (often labeled with an icon like ' ' or a left arrow) to review past data. This feature helps in monitoring progress effectively.

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Figure 4: The device display, emphasizing large text for readability and the ability to compare with previous measurements.

4. MAINTENANCE

4.1 Cleaning

To clean the device, wipe it with a soft, damp cloth. Do not use abrasive cleaners, solvents, or immerse the device in water. Ensure the surface is dry before storage or next use.

4.2 Storage

Store the body composition meter in a cool, dry place away from direct sunlight and extreme temperatures. Its compact design allows for convenient vertical storage in narrow spaces, such as between furniture or on a shelf.

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28mm

Figure 5: The Omron HBF-214's slim profile allows for easy vertical storage, saving space.

5. TROUBLESHOOTING

- **No Display/Device Not Turning On:** Check battery installation and ensure batteries are not depleted. Replace batteries if necessary.
- **Inaccurate Readings:** Ensure the device is placed on a hard, flat surface. Soft surfaces like carpets can cause inaccuracies. Ensure bare feet are clean and dry and make full contact with the electrodes.
- **Display in Japanese:** The device's user interface may be in Japanese. While this manual provides English instructions, you may need to use a translation app for on-screen prompts during initial setup or specific functions.
- **Error Messages:** If an error code appears, refer to the specific error code section in the full product manual (if available) or contact customer support.

6. SPECIFICATIONS

Feature	Detail
Model Number	HBF-214
Product Dimensions	11.22 x 11.02 x 1.1 inches
Weight	3.53 Pounds
Batteries	4 AAA batteries required
Manufacturer	OMRON
Color	White
Recommended Uses	Humans
Special Feature	Backlit Display, High Weight Capacity
Display Type	LED or LCD

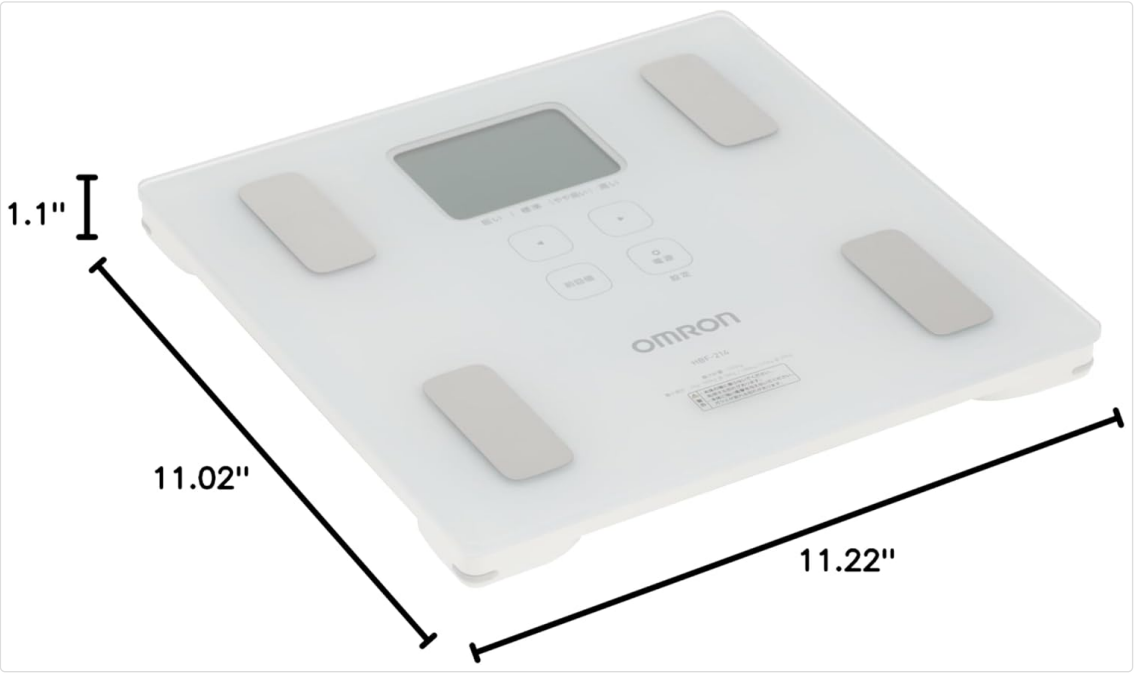
Figure 6: Dimensions of the Omron HBF-214 Body Composition Meter.

7. WARRANTY AND SUPPORT

7.1 Warranty Information

The Omron HBF-214 Body Composition Meter typically comes with a 2-year manufacturer's warranty. Please retain your proof of purchase for warranty claims. For specific warranty terms and conditions, refer to the documentation included with your product or contact Omron customer support.

7.2 Customer Support

For technical assistance, troubleshooting, or warranty inquiries, please contact Omron customer support. When contacting support, it is helpful to have the following information ready:

- Product Model Number (HBF-214)
- Manufacturing Number (if available)
- Date and Place of Purchase
- Detailed description of the issue

Contact Number: 0120-30-6606 (Please verify operating hours and regional availability.)

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
Figure 7: Omron customer support contact details. Prepare product information before calling.

Related Documents - HBF-214

	OMRON HBF-214 Karada Scan Body Composition Monitor Instruction Manual Comprehensive instruction manual for the OMRON HBF-214 Karada Scan Body Composition Monitor, detailing setup, usage, safety precautions, and interpretation of results.
	OMRON HBF-224 Body Composition Monitor Instruction Manual Comprehensive instruction manual for the OMRON HBF-224 Body Composition Monitor. Learn how to set up, use, and interpret results for accurate body composition analysis.


	OMRON HBF-222T Body Composition Monitor User Guide & Instructions Comprehensive instruction manual for the OMRON HBF-222T Body Composition Monitor, covering setup, operation, measurement, troubleshooting, and technical specifications for home use.
	OMRON HBF-212 Karada Scan Body Composition Monitor Instruction Manual This manual provides instructions for using the OMRON HBF-212 Karada Scan Body Composition Monitor. It covers setup, measurement, understanding results, troubleshooting, and specifications.
	OMRON HBF-702T Body Composition Monitor Instruction Manual Comprehensive instruction manual for the OMRON HBF-702T Body Composition Monitor, detailing setup, usage, measurement parameters (weight, body fat, BMI, etc.), and troubleshooting for home health monitoring.
	OMRON BCM-500 Body Composition Monitor and Scale with Bluetooth® Instruction Manual OMRON BCM-500 Instruction Manual: A guide to using the OMRON Body Composition Monitor and Scale with Bluetooth® for accurate body composition analysis, including weight, body fat, visceral fat, skeletal muscle, resting metabolism, and BMI. Features Bluetooth connectivity for smart device integration.

Documents - Omron – HBF-214

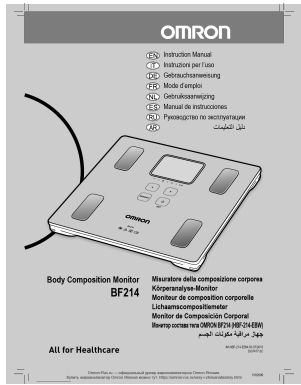[\[pdf\]](#) User Manual


Provided by ManualsAndMore Download User Manual Omron BF214 Instruction omron bf214 instruction manual 2 manualsandmor Instruction Manual Istruzioni per l'uso Gebrauchsanweisung Mode d'emploi Gebruiksaanwijzing Manual .. composition corporelle Lichaamscompositiemeter Monitor de Composicin Corporal OMRON BF214 **HBF-2** EBW IM-**HBF-214**-EBW-03-07/2015 5337477-2C Omron-Rus.ru -- Omron . Omron : <https://www.omron.com> lang:it score:50 filesize: 2.66 M page_count: 13 document date: 2019-07-24

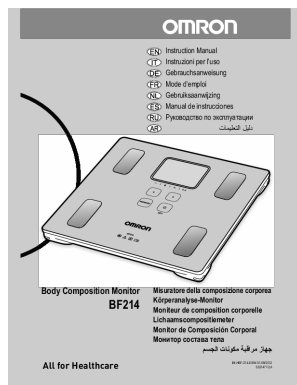[\[pdf\] User Manual](#)

Provided by ManualsAndMore Download User Manual Omron BF214 Instruction omron bf214 instruction manual 1 manualsandmor Instruction Manual Istruzioni per l'uso Gebrauchsanweisung Mode d'emploi Gebruiksaanwijzing Manual .. Moniteur de composition corporelle Lichaamscompositiemeter Monitor de Composicin Corporal IM-**HBF-21** EBW-01-08/2012 5337477-2A RG RPSRVLWLRQ 0RQLWRU 0LVXUDWRUH GHOD FRPSRVL LRQH FRUSR...

lang:it **score:47** filesize: 519.91 K page_count: 12 document date: 2019-07-23

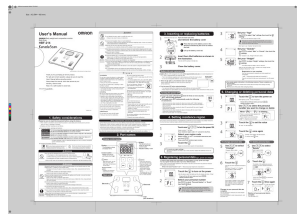
[\[pdf\] User Manual](#)

HBF 214 C m sideA Manual NEW STAR COMPANY LTD OMRON 222T Body Composition Monitor Instruction drive google file d 1BKsaxXDbjWvBaRGgLTe5r4tRO6d6Q0Q view usp sharing |||

N TJEF QEG Size : A2 594 420 mm User s Manual Body weight and composition monitor Body Scan **HB** Caution 3. Inserting or replacing batteries 3 Do not attempt to open, repair or modify the main unit. Disasse Doing so may result in injury or malfunction. prohibited Do not use i...


lang:en **score:44** filesize: 2.1 M page_count: 2 document date: 2012-02-06

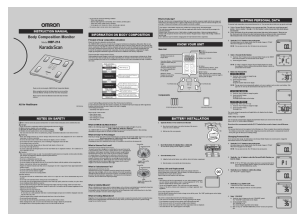[OMRON HBF-214 Karada Scan Body Composition Monitor Instruction Manual](#)

Comprehensive instruction manual for the OMRON HBF-214 Karada Scan Body Composition Monitor, dete setup, usage, safety precautions, and interpretation of results.

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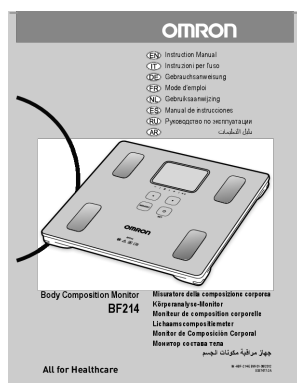
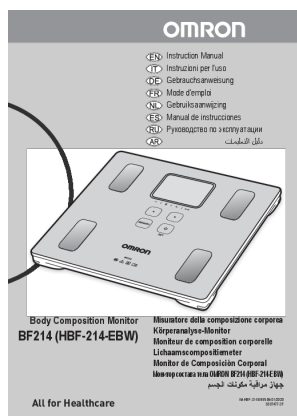
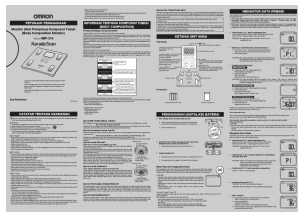
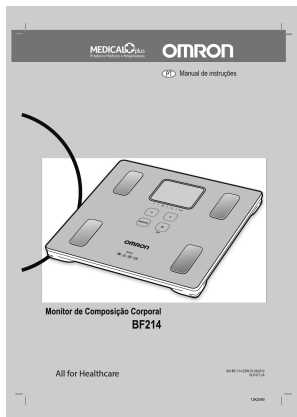
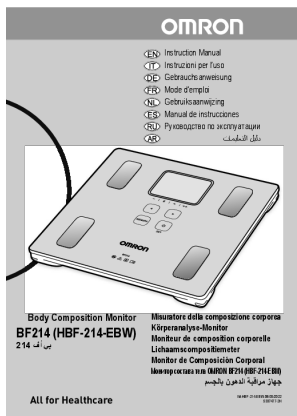
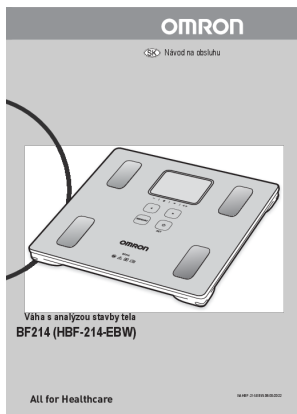
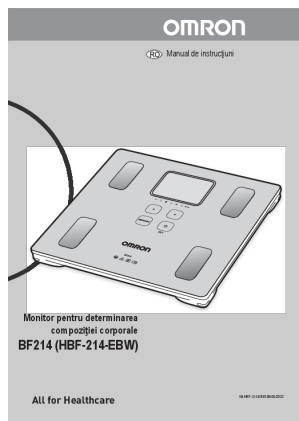
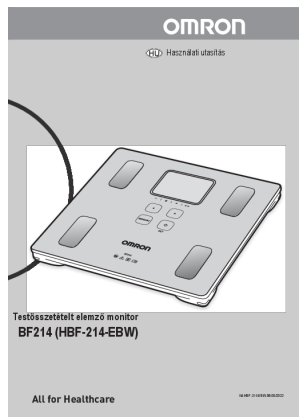
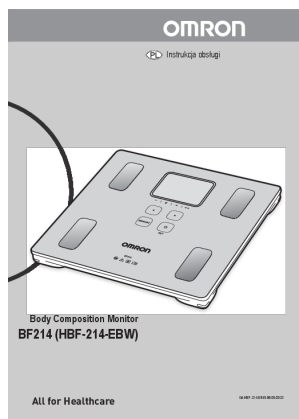
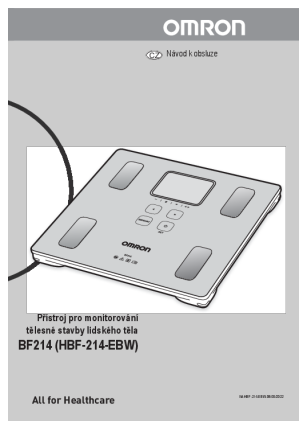
lang:i-klngon **score:30** filesize: 574.27 K page_count: 16 document date: 2022-10-13

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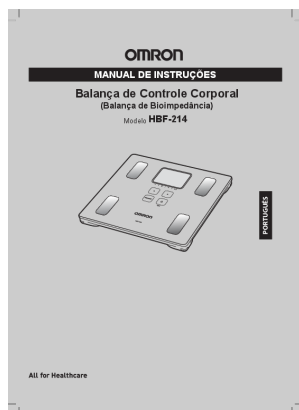
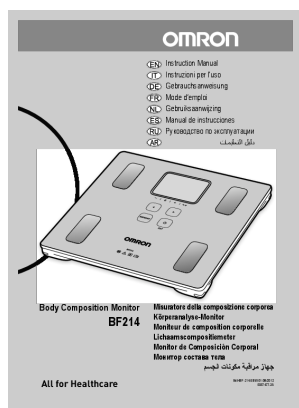
lang:it **score:28** filesize: 7.28 M page_count: 12 document date: 2013-04-11

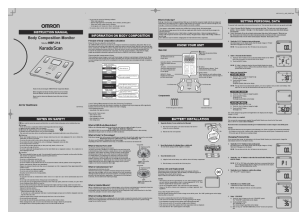
[\[pdf\]](#) User Manual Instructions

HBF 214 omronhealthcare ap 80 sg resources

INSTRUCTION MANUAL Body Composition Monitor Model**HBF-214** Thank you for purchasing the OMRO
Body Composition Monitor. Before using this unit for the first time, please be sure to read this Instruction M
carefully and use the unit safely and properly. Please keep this Instruction Manual at hand...

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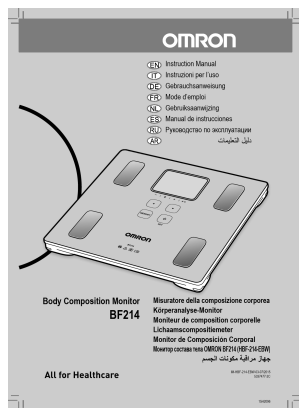
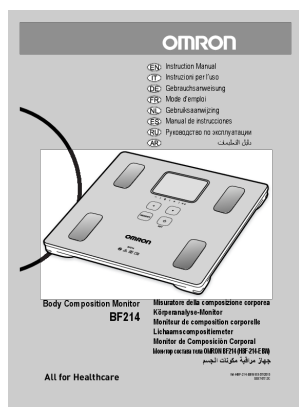
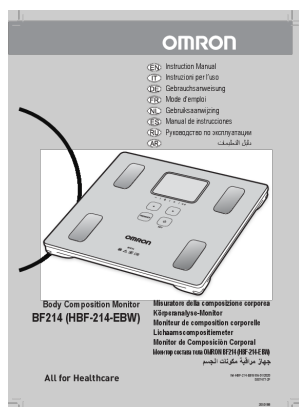
[pdf] User Manual

BF212 214 06 indd enzo Body Composition g ec2 images amazon eu ssl Monitors BF214 Fat Your body fat % is calculated using th Bioelectrical Impedance BI method The monitor uses an extremely weak electrical current passed through your to determine amount tissue 4 AAA batteries R03 You may also use alkaline LR03 Battery life Approx 12 months* when used measurement a day at room temperature 23°C *Supplied battery for trial This can run out within months Display Weight 2 150 kg with increment 0 1 750 2771 Fe 2772 Email information@eu omron BOX CONTAINS Monitor Manganese Batteries Instruction Manual bcm bf212 bf214 V37515710 uk health and beauty Omron

Body Composition Monitors BF214**HBF-214**-EBW Also available BF212 HBF-212-EW Small footprint EAN 15672 10703 8 RRP 44.95 CODE: EAN RRP: TYPE: **HBF-214**-EBW brown 40 15672 10704 5 59.95 BOD' COMPOSITION MONITOR DESCRIPTION BF214 It couldn t be simpler to check your body composition. T pro...

lang:en **score:26** filesize: 693.77 K page_count: 2 document date: 2013-02-27

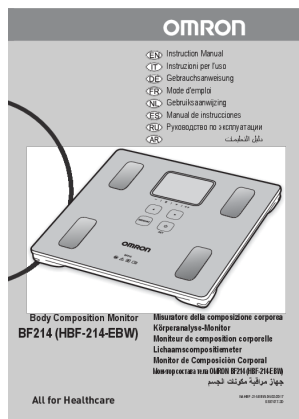
[\[pdf\]](#) User Manual

BF214 HBF 214 EBW IM EN 04 02 2017 omrondoc ru |||

Instruction Manual Istruzioni per l'uso Gebrauchsanweisung Mode d'emploi Gebruiksaanwijzing Manual d'instrucciones Body Composition Monitor BF214 **HBF-214**-EBW All for Healthcare Misuratore della composizione corporea Krperanalyse-Monitor Moniteur de composition corporelle Lichaamscompositiemeter...

lang:it **score:25** filesize: 1.28 M page_count: 28 document date: 2017-08-03

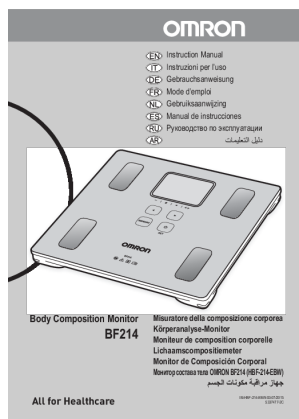[\[pdf\]](#) User Manual

Omron BF214 electronic body analyser scale Manual 879cfe156f812a216cc3e015cdf8e3c4dc43 gzhlis at blob ldb 5 9 2 |||

Instruction Manual Istruzioni per l'uso Gebrauchsanweisung Mode d'emploi Gebruiksaanwijzing Manual .. sition corporelle Lichaamscompositiemeter Monitor de Composicin Corporal 20521 : IM-**HBF-214**-EBW-03-07/2015 5337477-2C Thank you for purchasing this body composition monitor. It s int...

lang:it **score:25** filesize: 652.23 K page_count: 12 document date: 2017-01-05


[\[pdf\]](#) User Manual Declaration of Conformity Borchure

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lang:id **score:24** filesize: 845.85 K page_count: 2 document date: 2015-10-27

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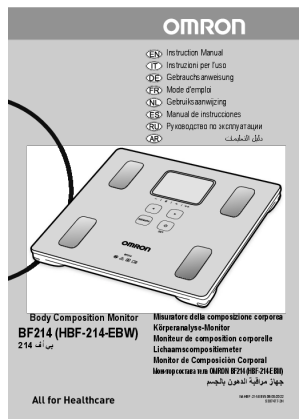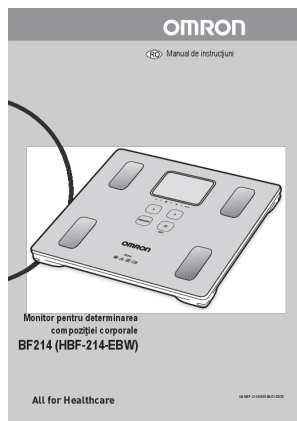
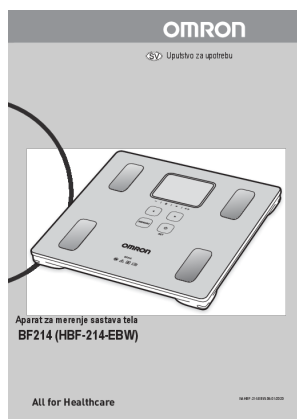
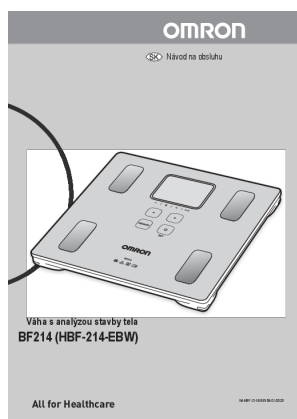
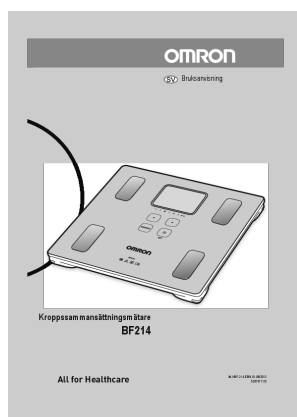
lang:no **score:22** filesize: 891.37 K page_count: 12 document date: 2020-10-07


[\[pdf\]](#) User Manual Catalog

Omron blood pressure monitors wrist BLOOD PRESSURE MONITORING OMRON offers you true quality accuracy and reliability. We are looking for an upper arm or type device we have the omronkatalog 2013 ajmedical.se PRODUCT CATALOGUE HEALTHCARE 637 Manual Pdf Wrist Blood Pressure Monitor with Advanced Positioning |||

OMRON PRODUCT CATALOGUE 2013 HEALTHCARE AJ Medical Uddvgen 32A 181 30 Liding Tel: 08-76 70 70 Fax: ... linically tested to ensure accurate results, helping you to track your body composition over time.

214-EBW HBF-212-EW * only BF 214 BF306 Speed accuracy = total control Regularly used by prof...

lang:en **score:21** filesize: 4.35 M page_count: 64 document date: 2013-08-26

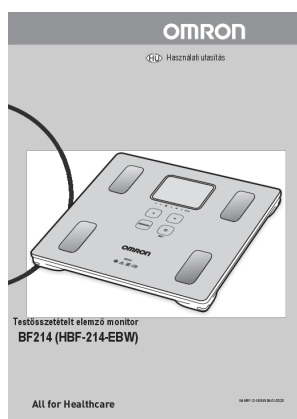
[\[pdf\]](#) Catalog

OMRON Healthcare Weight Management BF214 web IM HBF 214 EBW HU 06 01 2020 omron healthcare on demandware static Si master catalog default dw1260432c s |||

Használati útmutató 7HVW VV HWpWHOWHOHP PRQLWRU BF214 **HBF-214-EBW** All for Healthcare IM-**HBF-214-EBW-06-01/2020** . V QMNKRJ H WDWHVW VV HWpWHOWHOHP PpUNpV OpNHW YiODV WRWWC OpNDWHVW VV HWpWHON YHWNH pUWPnHLWWXGMD PpUQLpVPHJMHQHtWHQLWHVWW PHJWHVW VtUEDQ Yi L PRN W PHJEDQ pV7HVWW PHJLQGH 0, V...

lang:tl **score:21** filesize: 568.56 K page_count: 16 document date: 2020-04-24

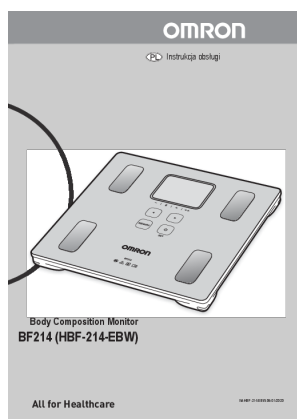[\[pdf\]](#) Catalog

OMRON Healthcare Weight Management BF214 web IM HBF 214 EBW PL 06 01 2020 omron healthcare on demandware static Si master catalog default dw5dc3b7d7 s |||

,QVWUXNFMMDREVXJL Body Composition Monitor BF214 **HBF-214-EBW** All for Healthcare IM-**HBF-214-EBW-06-01/2020** LNXMHP D DNXSDQDOL DWRUDVNDGXFLDD6X RQGRSRPLDUX LZ ZLHWODQLDQDVWSXMF FKSDUDPHWUyZVNDGXFLDDPDV FLDD DZDUWRFLWNDQNLWXV F RZHMZFLHOHZ PLQLV NLHOHWZRZ FK Z RUD ZVNDQLNDPDV FLDDERG PDVVLQGH ...

lang:sw **score:21** filesize: 592.02 K page_count: 16 document date: 2020-04-24


[\[pdf\]](#) Catalog

OMRON Healthcare Weight Management BF214 web IM HBF 214 EBW CS 06 01 2020 omron healthcare on demandware static Si master catalog default dwf97f6eea s |||

CZ Nvod k obsluze 3tVWURMSURPRQLWRURYiQt WcOHVQpVWDYE OLGvNpKRWCOD : All for Health IM-**HBF-214-EBW-06-01/2020** CZ CNXMHPHYiP

DQINXSWRKRWRStVWURMHSURPRQLWRURYiQtWCOHVQp VWDYE OLGvNpKRWCOD6ORXtNPCl REUD HQtQiVOHGXMtFtFK SDUDPHWUWCOHVQpVWDYE OLGvNpKRWCODKPRWQRVWWCOHVC Y NRVWHUQtV...


lang:i-klingon **score:20** filesize: 574.24 K page_count: 16 document date: 2020-04-24

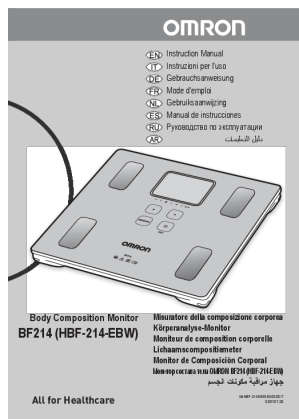[pdf] User Manual

90237 Omron KoerperfettanalysewaageBF214 GA medplus24 de out media |||

Instruction Manual Istruzioni per l'uso Gebrauchsanweisung Mode d'emploi Gebruiksaanwijzing Manual d'instrucciones Body Composition Monitor BF214 **HBF-214**-EBW All for Healthcare Misuratore della composizione corporea Krperanalyse-Monitor Moniteur de composition corporelle Lichaamscompositieme...
lang:it **score:19** filesize: 973.34 K page_count: 28 document date: 2021-09-17

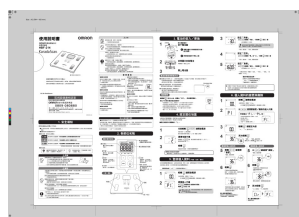
[pdf]

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67a07bb4e81d5omronhealthcare tw storage

67a07bb4e81d5 |||

Size : A2 594 420 mm eef 3. 3 cdc **HBF-214** f f df f df df f f 1 1 2 d 1 d f 2 1080f df df df 4 cdc f df eddf f (-df df ddf f 2 f 3 df 5 cdc 100.0199.5 cmf deed 0.0 kgf f OMRONf ddf df ee ef f f d f f d f d f ...
lang:nl **score:18** filesize: 2.19 M page_count: 2 document date: 2012-04-16

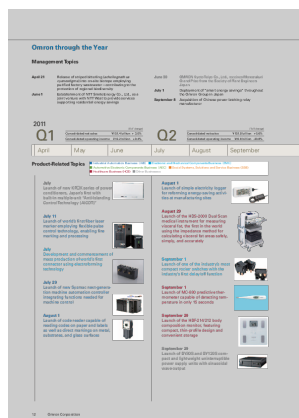[pdf]

Omron through the Year Apr 21 2025 · Launch of HJA 310 Calorie Scan activity monitor for indicating calorie reduction necessary to monthly dieting targets and determine ar12 09e omron about ir irlib s ar12e |||


Omron through the Year Management Topics April 21Release of striped bitterling acheilognathus cyano ... well as direct markings on metal, substrates, and glass surfaces September 20 Launch of the **HBF-214/214** composition monitor, featuring compact, thin-profile design and convenient storage ...
lang:en **score:18** filesize: 627.51 K page_count: 2 document date: 2012-09-19

[pdf]

Microsoft PowerPoint HCB FY12 1Q Presentation ppt omronir Healthcare OMRON Global No 1 blood pressure monitor global shar Brand building through academic activity Sales network covering 110 nations 20120730 presentation healthcare omron en assets fi 80 about s

OMRON HEALTHCARE Business Strategies OMRON HEALTHCARE Co., Ltd. President and Chief Executive Office ... - Electronic meter Pulse Massager 19/21 Tooth Brush Massager Stylish body composition monitor **214** Mass-production Technology to Achieve Cost Competitiveness Full automation and easy-to-mak...
lang:en **score:17** filesize: 1.99 M page_count: 21 document date: 2012-07-30

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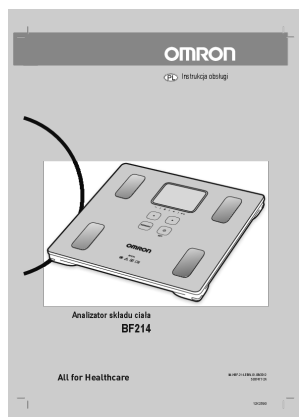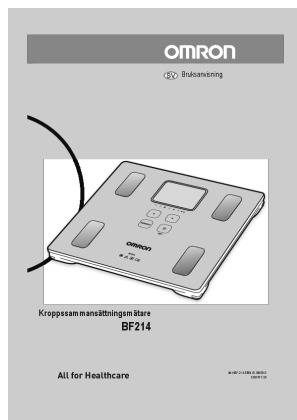
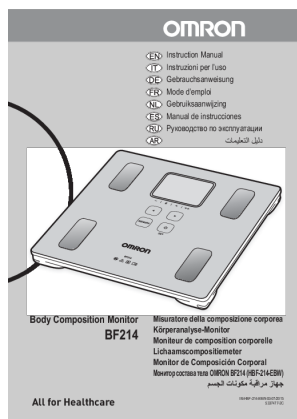
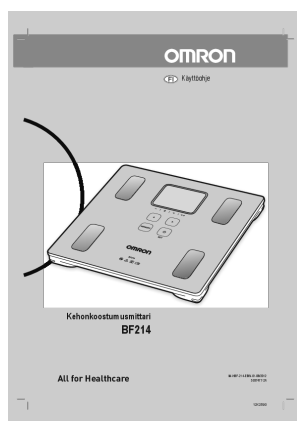
lang:fr **score:14** filesize: 5.66 M page_count: 35 document date: 2020-12-14

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ce0af90516354bff909cfa7b63c4632b omronhealthcare ap Content products |||

MY O TNHH PHN C TH HBF-255T Hng dn s dng 5649687-9A Cm n bn mua my o thnh phn c th OMRON. s

my ng cch v an ton, vui lng c k Hng dn s dng trc khi dng my ln u tin. Hy lu gi cn thn bn Hng dn s dng tham kl ny. Thit b ny c dng o v hin th cc ch s thnh phn c th sau y. - Trng lng c th ...

lang:tl **score:13** filesize: 12.98 M page_count: 10 document date: 2021-11-24

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