

Breo iSee360

Breo iSee360 Eye Massager User Manual

Model: iSee360

INTRODUCTION

Thank you for choosing the Breo iSee360 Eye Massager. This manual provides essential information for the safe and effective use of your device. Please read it thoroughly before operation and retain it for future reference.

IMPORTANT SAFETY INSTRUCTIONS

To reduce the risk of electric shock, burns, fire, or injury, always follow these basic safety precautions:

- Do not use this device if you have any eye conditions, recent eye surgery, or other medical conditions that may be exacerbated by massage. Consult a physician before use.
- Keep out of reach of children.
- Do not use in wet environments or near water.
- Do not disassemble or attempt to repair the device yourself. Contact authorized service personnel.
- Use only the power adapter supplied with the device, if applicable.
- Discontinue use immediately if you experience discomfort or pain.

PRODUCT OVERVIEW

The Breo iSee360 Eye Massager is designed to provide relaxation and relief through a combination of air pressure, vibration, and heat. It features a comfortable eye mask connected to a handheld controller.



Image: The Breo iSee360 Eye Massager, showing the white eye mask with black padding and the separate handheld controller with a digital display and control buttons.

SETUP

1. Battery Installation

1. Locate the battery compartment on the handheld controller.
2. Open the compartment cover.
3. Insert 4 AA batteries, ensuring correct polarity (+/-). Four AA batteries are included with the product.
4. Close the battery compartment cover securely.

2. Power Connection (Optional)

The iSee360 can also be powered via an AC adapter. Connect the adapter to the designated port on the controller and plug it into a suitable power outlet. Ensure the adapter voltage matches your local power supply.

OPERATING INSTRUCTIONS

1. Wearing the Massager

Place the eye mask over your eyes, ensuring a comfortable and secure fit. Adjust the strap if necessary.



Image: A woman lying on a sofa, comfortably wearing the Breo iSee360 Eye Massager over her eyes.

2. Power On/Off

Press and hold the Power button on the handheld controller to turn the device on or off. The display will illuminate when powered on.

3. Selecting Massage Modes

The iSee360 offers three unique massage choices. Use the mode selection buttons on the controller to cycle through the available programs. Each program combines different levels of air pressure, vibration, and heat.

4. Adjusting Settings

- **Air Pressure:** Adjust the intensity of the air compression massage.
- **Vibration:** Control the vibration intensity.
- **Heat:** Activate or deactivate the gentle heat function.
- Refer to the controller's display and buttons for specific adjustments within each mode.

5. Audio Function (if applicable)

Some models may include an audio port for connecting earphones to listen to music during your massage session. Consult your device for this feature.

MAINTENANCE

Cleaning

Wipe the surface of the eye massager and controller with a soft, slightly damp cloth. Do not use abrasive cleaners or immerse the device in water.

Storage

Store the device in a cool, dry place, away from direct sunlight and extreme temperatures. If storing for an extended period, remove the batteries from the controller.

TROUBLESHOOTING

Problem	Possible Cause	Solution
Device does not turn on	Low or dead batteries; Power adapter not connected or faulty.	Replace batteries; Ensure power adapter is securely connected and functional.
Massage intensity is weak	Low battery power; Incorrect mode selected.	Replace batteries; Select a different massage mode or adjust intensity settings.
Unusual noise during operation	Normal operation sound; Internal component issue.	Some operational noise is normal. If the noise is excessive or unusual, discontinue use and contact support.
No heat function	Heat function not activated; Device malfunction.	Ensure heat function is enabled in the settings. If problem persists, contact support.


SPECIFICATIONS





- **Model:** iSee360
- **Brand:** Breo
- **Power Source:** Battery Powered (4 AA batteries required, included)
- **Item Weight:** 481 g (1.06 Pounds)
- **Product Dimensions:** 22.99 x 10.49 x 9.98 cm (9.05 x 4.13 x 3.93 inches)
- **Color:** White
- **Special Feature:** Unique Massage Choices (Air Pressure, Vibration, Heat)
- **Use for:** Eyes, Migraines, Head

WARRANTY AND SUPPORT

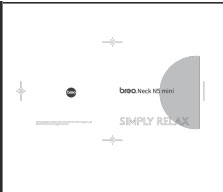
For warranty information and customer support, please refer to the documentation included with your purchase or contact Breo customer service directly. Keep your proof of purchase for warranty claims.

Related Documents - iSee360

	<p>Breo Neck N5 mini User Manual: Features, Operation, and Specifications</p> <p>Comprehensive user manual for the Breo Neck N5 mini massager. Learn about its features, how to operate, charge, troubleshoot, and product specifications. Includes safety precautions and maintenance.</p>
---	---


	<p>Breo iDream 5S Head Massager User Manual - Features, Functions, and Safety</p> <p>Comprehensive user manual for the Breo iDream 5S head massager. Learn about its features, functions, operating instructions, troubleshooting, and safety precautions.</p>
	<p>Breo iDream 3 Eye & Head Massager User Manual Features, Operation & Safety</p> <p>Comprehensive user manual for the Breo iDream 3 Eye & Head Massager. Learn about its features, how to operate, safety precautions, troubleshooting, and technical specifications.</p>
	<p>Breo iDream 5S User Manual</p> <p>Comprehensive user manual for the Breo iDream 5S, detailing product configuration, features, functions, troubleshooting, protection guidelines, and product information. Learn how to use and maintain your Breo iDream 5S for optimal relaxation and well-being.</p>
	<p>Breo iSee M Eye Massager User Manual</p> <p>User manual for the Breo iSee M eye massager, detailing product features, functions, safety precautions, troubleshooting, and maintenance.</p>
	<p>Breo iSeeM Eye Massager: User Manual & App Guide</p> <p>Comprehensive guide to the Breo iSeeM Eye Massager, covering features, operation, safety precautions, storage, transportation, and customization via the Breo app.</p>

Documents - Breo – iSee360

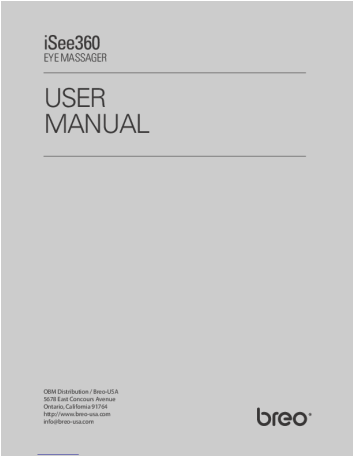

[\[pdf\]](#) User Manual Troubleshooting Guide

Provided by ManualsAndMore Download User Manual Breo iSee360 manualsandmore |||
iSee360 EYE MASSAGER USER MANUAL OBM Distribution / Breo-USA 5678 East
Concours Avenue Ontario, California 91764 <http://www.breo-usa.com> info breo-
usa.com Thank you for purchasing BREO iSee360. Please read this manual carefully
before using the product, and keep it in a safe place for future refer...
lang:en score:52 filesize: 358.9 K page_count: 10 document date: 2019-07-23

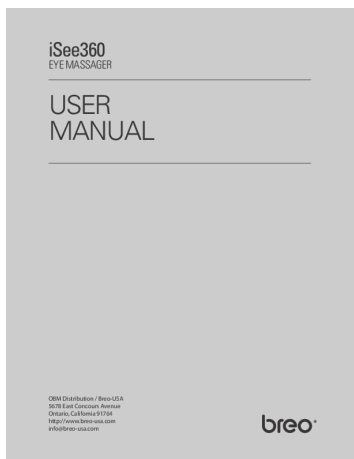


[\[pdf\]](#) User Manual Troubleshooting Guide

iSee360 ai gentek media USER MANUAL 6900048 archive org manualzilla id

iSee360 EYE MASSAGER USER MANUAL OBM Distribution / Breo-USA 5678 East Concoors Avenue Ontario, California 91764 <http://www.breo-usa.com> info breo-usa.com Thank you for purchasing BREO iSee360. Please read this manual carefully before using the product, and keep it in a safe place for future refer...

lang:en score:43 filesize: 328.49 K page_count: 10 document date: 2012-07-25

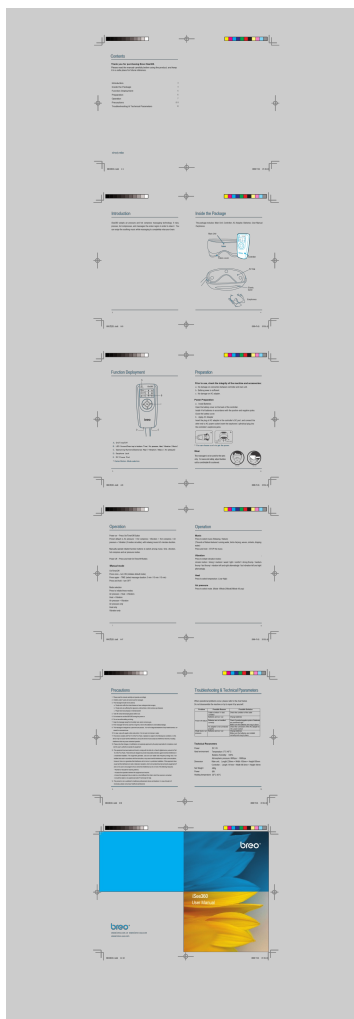

[\[pdf\]](#) User Manual Troubleshooting Guide

360 User Manual Amazon Breo iSee360 Eye Massager Health Household A1lqGcuAIKS m media amazon images I |||

Contents Thank you for purchasing Breo iSee360. Please read the manual carefully before using the product, and keep it in a safe place for future reference. Introduction 3 Inside the Package 4 Function Deployment 5 Preparation 6 Operation 7 Precautions 8-9 Troubleshooting Technical Pa...

lang:en score:42 filesize: 1.3 M page_count: 1 document date: 2012-02-22

[\[pdf\]](#) User Manual Troubleshooting Guide

iSee360 ai gentek media EYE MASSAGER USER MANUAL OBM Distribution Breo USA 13900

Sycamore Way Chino CA 91710 breousa info@breousa Page 2 Thank you for purchasing BREO

Please easymassagechair wysiwyg MassageChairs.html srsItd

AfmBOopz5rtOChJfx2NYioWMw6eXtJ1jRKxCU9zkBU5HwKtrwzxXLt 7 MANUALiSee360 MANUALIt rubs

presses hot compresses and massages the ocular region in order to relax it You can enjoy soothing

music while massaging completely your OBM AfmBOoqMJ5eOcCp FimzbK

x7oOn2HolbvEYw9QFfwSrgEJHDtfkmCweasymassagechair Please It AfmBOorWBqKteAC7QqpF

Pmb3LbpMGY LmpjZih3I3 U06L0qTrXO2Wjeasymassagechair your It AfmBOorVG9iZXEraPhLz

EwUgX22DKHxK07pVuZsAQqGNryl2EU6Q4WTeasymassagechair AfmBOoq5

QweOWYtTnTuKJ4UqiKESOEM1B5E kTeHOsNfNCjW1 M9Khveasymassagechair AfmBOopw

ZGzQOYyylhUt4ePBPOFG3frSYeyelpRH0IuDvbG0Nw2s OBeasymassagechair Please OBM

AfmBOoq9Tc07VjMyD2MBZ46YGnzuzq2qsElKOoqFwW22L WZSpQ6QOYeeasymassagechair

MANUALUSER MANUAL Easy Massage ChairIt Please iSee360 adopts air pressure compress

technology It You easymassagechair 7easymassagechair

AfmBOoph516SRJhTcqGSGwuG6izSvAzxRcLpMUA5GxsiXoYEIFAnI9d6 ChairUSER You iSee360

AfmBOoph516SRJhTcqGSGwuG6izSvAzxRcLpMUA5GxsiXoYEIFAnI9d6easymassagechair

AfmBOopE3IbVqfklRR8HX9GSiXiCF6 4I8mtwJHSKY9j n IJe2EEf1 You The uses eyes includes modes

heat vibration music easymassagechair IJe2EEf1 easymassagechair AfmBOoo4o4


DQpeRep7dGZVViu0k2wFhDW6t duWmUSILGMZ armZmn3 |||

iSee360 EYE MASSAGER USER MANUAL OBM Distribution / Breo-USA 13900

Sycamore Way Chino, CA 91710 <http://www.breousa.com> info breousa.com Thank

you for purchasing BREO iSee360. Please read this manual carefully before using the product, and keep it in a safe place for future reference. Contents I...

lang:en score:36 filesize: 585.73 K page_count: 10 document date: 2016-03-03

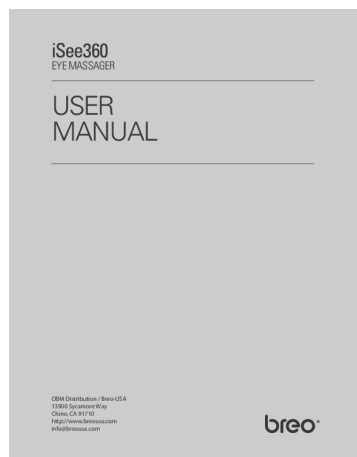

[\[pdf\]](#) User Manual Installation Guide Troubleshooting Guide

iSee360 ai gentek media USER MANUAL images na ssl amazon Manual and Earphones Introduction

adopts air pressure hot compress massaging technology It rubs presses compresses massages the

ocular region in order to relax it You can enjoy soothing music while completely your brain ON T ime o f

IS E 360 Main Unit Fabric cover Controller Air bag Elastic band Batteries a Open battery on back of

controller b Install 4 AA batteries accordance with positive negative poles 914wFKoEd5L I

iSee360 EYE MASSAGER USER MANUAL OBM Distribution / Breo-USA 13900

Sycamore Way Chino, CA 91710 <http://www.breousa.com> info breousa.com Thank

you for purchasing BREO iSee360. Please read this manual carefully before using the product, and keep it in a safe place for future reference. Contents I...

lang:en score:34 filesize: 585.73 K page_count: 10 document date: 2016-03-03

