

Beurer BG13

Beurer BG13 Diagnostic Bathroom Scale Instruction Manual

Brand: Beurer | **Model:** BG13

INTRODUCTION

The Beurer BG13 Diagnostic Bathroom Scale is designed to measure body weight, body fat percentage, body water percentage, muscle density, and Body Mass Index (BMI). This scale features a weighing area made from safety glass and brushed stainless steel electrodes for accurate measurements. It supports up to 10 user memories and offers 5 activity levels for personalized analysis.



Image: The Beurer BG13 Diagnostic Bathroom Scale, showcasing its sleek design in a bathroom environment.

SETUP

1. **Battery Installation:** The scale requires 1x 3V CR 2032 Lithium Metal battery, which is included. Open the battery compartment on the underside of the scale and insert the battery, ensuring correct polarity.
2. **Placement:** Place the scale on a firm, level surface. Avoid carpets or uneven flooring, as this can affect measurement accuracy.
3. **Initial Power On:** The scale features automatic switch-on. Simply step onto the weighing surface to activate it.
4. **Unit Selection:** To change the unit of measurement (kg/lb/st), locate the unit switch on the underside of the scale or use the 'SET' button if available on the display.
5. **User Profile Setup:** The scale supports 10 user memories. To set up a user profile, press the 'SET' button. Use the arrow buttons to input personal data such as age, height, gender, and activity level (1-5). Confirm each entry with the 'SET' button.



Image: Visual representation of key features including Body Mass Index calculation, 10 user memory slots, 150 kg weight capacity, and a slim 18 mm profile.

OPERATING INSTRUCTIONS

- Weighing:** Step onto the scale barefoot. Ensure your feet are positioned on the brushed stainless steel electrodes for diagnostic measurements. The scale will automatically switch on and display your weight.
- Diagnostic Measurement:** After displaying your weight, the scale will proceed to calculate body fat, body water, muscle percentage, and BMI. Remain still until all measurements are displayed.
- User Recognition:** If multiple user profiles are set up, the scale will attempt to recognize the user based on their weight. If recognition is successful, the stored personal data will be used for diagnostic calculations.
- Reading the Display:** The scale features an extra-large LCD display with a digit size of 38 mm for clear readability.



Image: A detailed view of the scale's large digital display, showing a weight measurement in kilograms.

MAINTENANCE

- **Cleaning:** Clean the scale surface with a damp cloth and mild detergent. Do not use abrasive cleaning agents or immerse the scale in water.
- **Battery Replacement:** When the battery indicator appears on the display, replace the 3V CR 2032 Lithium Metal battery.
- **Storage:** Store the scale in a dry place, away from extreme temperatures and direct sunlight.

TROUBLESHOOTING

- **"Err" or Overload Indication:** If the display shows "Err" or an overload symbol, it indicates that the maximum weight capacity of 150 kg has been exceeded. Step off the scale immediately to prevent damage.
- **Inaccurate Readings:** Ensure the scale is placed on a hard, flat surface. Avoid moving the scale between measurements. Ensure feet are clean and dry and positioned correctly on the electrodes for diagnostic measurements.

- **No Display:** Check if the battery is correctly inserted and has sufficient charge. Replace the battery if necessary.

SPECIFICATIONS

Model Number	BG13 (Item no.: 760.30)
Weight Capacity	150 kg
Measurement Increments	100 g for body weight, 0.1% for body fat, body water, and muscle density
User Memories	10
Activity Levels	5
Display Type	LCD, 38 mm digit size
Weighing Surface	Safety glass (6 mm)
Electrodes	Brushed stainless steel
Power Source	1x 3V CR 2032 Lithium Metal battery (included)
Product Dimensions	Approx. 30 x 30 x 1.8 cm
Product Weight	Approx. 1.55 kg
Manufacturer	Beurer GMBH, Soflinger, Ulm, Germany



Image: The Beurer BG13 scale illustrating its compact dimensions of 30 cm length, 30 cm width, and 1.8 cm height.

WARRANTY

The Beurer BG13 Diagnostic Bathroom Scale comes with a **5-year guarantee** from the date of purchase. Please retain your proof of purchase for warranty claims. The warranty covers manufacturing defects and material faults under normal use conditions.



BG 13


beurer
wellbeing

BG 13 Diagnostic bathroom scale

 **german**
engineering



With extra large display

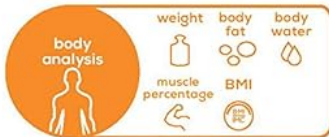

Digit size 38 mm

Weighing area made from safety glass (6 mm)

Electrodes made of brushed stainless steel

Determination of body weight, fat, water, muscle density and BMI

Capacity 150 kg

Introduction:

- 0.1% for body fat, body water and muscle density
- 100 g for body weight

Change over to kg/lb/st

10 user memories

5 activity levels

Automatic switch-on and switch-off

Overload indication

1 x 3V CR 2032 battery included

Product measurements: 30 x 30 x 1,8 cm

Product weight: approx. 1549 g

5 year guarantee

EAN no.: 4211125 76030 4

Item no.: 760.30



760.30_BG13_2018-07-09_01_DS_EN_BEU Subject to errors and changes

Image: Product information sheet for the Beurer BG13, confirming a 5-year warranty period.

SUPPORT

For technical assistance, warranty claims, or further inquiries, please contact Beurer customer support:

- **Manufacturer:** Beurer GMBH, Sofflinger, Ulm, Germany
- **Email:** info@beurer.com
- **Phone:** +49 (731) 3989-0
- **Indian Importer Support:** Beurer India Pvt Ltd, 523-524, JMD Megapolis, Sector - 48, Gurugram, Haryana-122018
- **Indian Support Email:** support@beurerindia.com
- **Indian Support Phone:** 01244477111