

SANITAS 755.18

Sanitas SGS03 Pure White Glass Bathroom Scales User Manual

Model: 755.18

1. INTRODUCTION

The Sanitas SGS03 is a modern digital bathroom scale designed for accurate body weight measurement. It features a pure white safety glass platform and a slim-line profile. The large 40mm LCD display is backlit in blue for enhanced visibility. This scale offers easy switching between metric (kilograms) and imperial (stones, pounds) units. It is equipped with an automatic on/off function via a vibration sensor to conserve battery life. The scale can measure weights up to 150 kilograms with 100-gram increments.



Image 1.1: Sanitas SGS03 Pure White Glass Bathroom Scale showing a weight reading.

2. SETUP

2.1. Battery Installation

The scale typically comes with batteries included. Before first use, locate the battery compartment on the underside of the scale. Remove any protective film or tab from the batteries to activate them. Ensure batteries are inserted with the correct polarity (+/-) as indicated inside the compartment.

2.2. Initial Placement

Place the scale on a firm, flat surface. Uneven or soft surfaces (like carpet) can affect accuracy. Avoid placing the scale near heat sources or in direct sunlight.

2.3. Unit Selection (kg/st/lb)

To switch between measurement units (kilograms, stones, or pounds), locate the unit switch button on the underside of the scale, usually near the battery compartment. Press this button to cycle through the available units until your desired unit is displayed on the screen.

3. OPERATING INSTRUCTIONS

3.1. Taking a Measurement

1. Ensure the scale is on a flat, stable surface.
2. Gently step onto the scale. The scale will automatically turn on due to its vibration sensor.
3. Stand still with your weight evenly distributed until the display shows a stable reading.
4. Step off the scale. The reading will typically hold for a few seconds before the scale automatically shuts off to conserve power.



Image 3.1: Close-up view of the digital LCD display showing a weight in kilograms.

4. MAINTENANCE

4.1. Cleaning

Clean the scale's glass surface with a damp cloth and a mild detergent if necessary. Do not use abrasive cleaners or immerse the scale in water. Dry thoroughly after cleaning.

4.2. Battery Replacement

When the display shows a 'Lo' indicator, it is time to replace the batteries. Open the battery compartment on the underside of the scale, remove the old batteries, and insert new ones of the same type, ensuring correct polarity. Dispose of old batteries responsibly.

4.3. Storage

Store the scale in a dry place at room temperature. Avoid extreme temperatures or high humidity, which can affect electronic components.

5. TROUBLESHOOTING

- **No Display / 'Lo':** Replace batteries. Ensure they are inserted correctly.
- **'Err' / Overload:** The maximum weight capacity of 150 kg has been exceeded. Step off the scale immediately.
- **Inaccurate Readings:**
 - Ensure the scale is on a hard, flat, and level surface.
 - Step onto the scale gently and stand still.
 - Ensure your feet are clean and dry.
 - Avoid moving the scale between measurements.
 - Recalibrate by stepping on and off the scale quickly, then wait for it to turn off before taking a measurement.
- **Display Flashes:** This may indicate an unstable reading. Stand still until the reading stabilizes.

6. SPECIFICATIONS

Brand	SANITAS
Model	755.18 (SGS03)
Color	White
Material Type	Glass
Weight Capacity Maximum	150 Kilograms (330 lbs / 23.5 st)
Readout Accuracy	100 Grams (0.2 lb)
Display Type	LCD with Backlight
Special Features	Auto Shut Off, Backlit Display, Large Display, Unit Conversion
Item Weight	798 Grams
Power Source	Battery (type not specified, typically CR2032 or AAA)

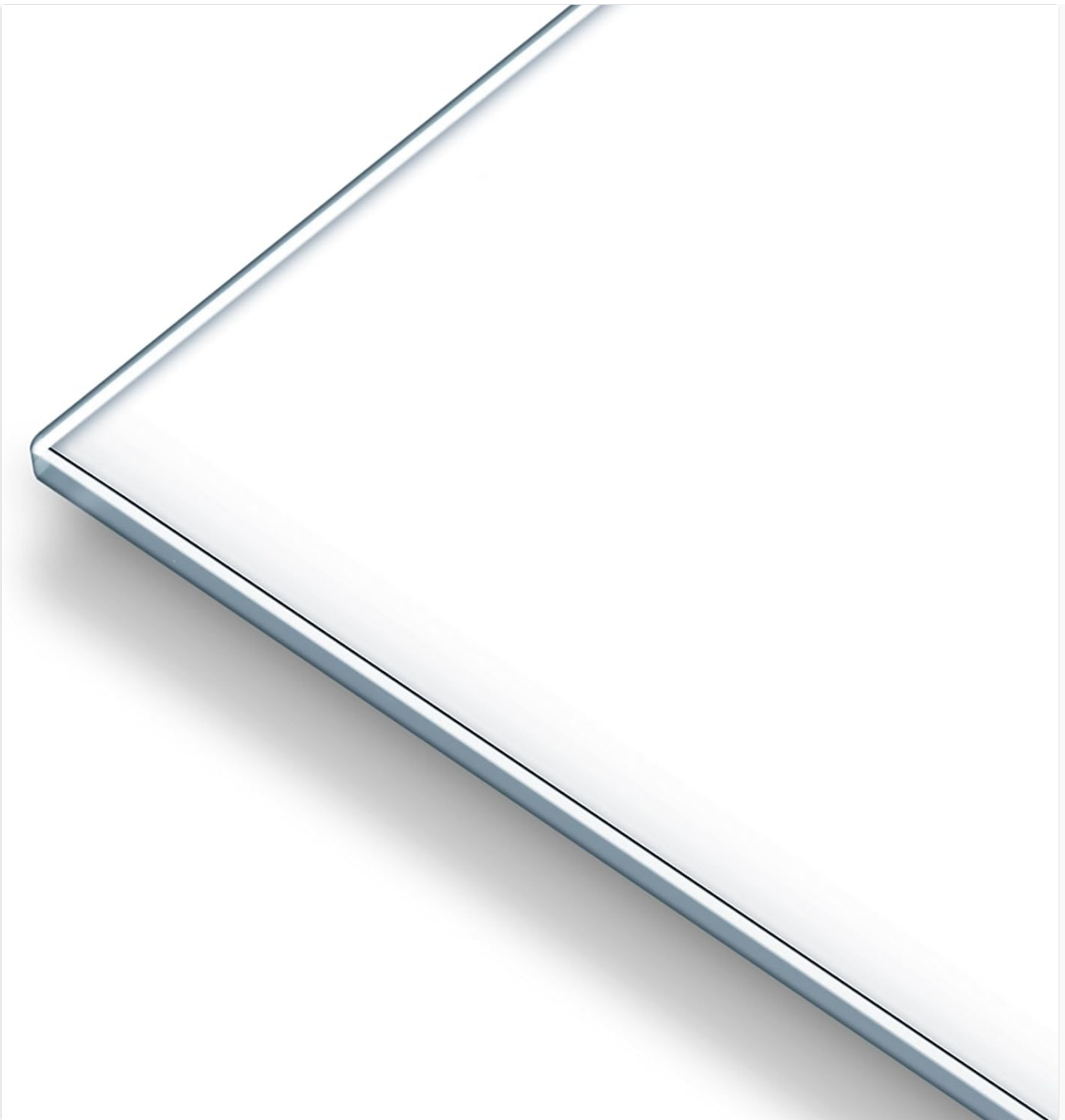


Image 6.1: Close-up view of the safety glass edge of the scale, highlighting its material and slim design.

7. WARRANTY

The Sanitas SGS03 Pure White Glass Bathroom Scales come with a 3-year guarantee from the date of purchase. This warranty covers manufacturing defects under normal use. Please retain your proof of purchase for any warranty claims. For further details or to initiate a claim, please contact Sanitas customer support.

8. CUSTOMER SUPPORT

For technical assistance, troubleshooting not covered in this manual, or warranty inquiries, please refer to the contact information provided with your product packaging or visit the official Sanitas website. When contacting support, please have your model number (755.18) and proof of purchase readily available.

© 2023 SANITAS. All rights reserved.