

Project E Beauty PE094

Project E Beauty Portable Rechargeable Galvanic Device

MODEL: PE094

Introduction

The Project E Beauty Portable Rechargeable Galvanic Device is designed to provide a natural and non-invasive method for skin care. It utilizes galvanic current to deeply cleanse the skin, reduce impurities, and enhance the absorption of nutrients and creams. This device is beneficial for improving skin texture, firming the skin, reducing redness, smoothing out wrinkles, increasing blood circulation, and eliminating acne-causing oil. It offers a safe and effective way to achieve beautiful skin without the need for surgeries, chemical peels, or injections.

Safety Information

Before using the device, ensure all jewelry and metal objects are removed from your person. This is a crucial step to prevent interference with the galvanic current. Always follow the instructions carefully to ensure safe and effective use.

Product Components

The Project E Beauty Galvanic Device comes with two interchangeable probes to cater to different facial areas:

- **Roller Probe:** Ideal for broad areas such as the face and neck.
- **Ball Probe:** Perfect for smaller, more delicate areas like around the nose, eyes, and mouth.



Figure 1: The Project E Beauty Portable Rechargeable Galvanic Device, including the main unit, roller probe, and ball probe.

How Does *Galvanic Current* Work on The Skin?

Scientifically proven skincare

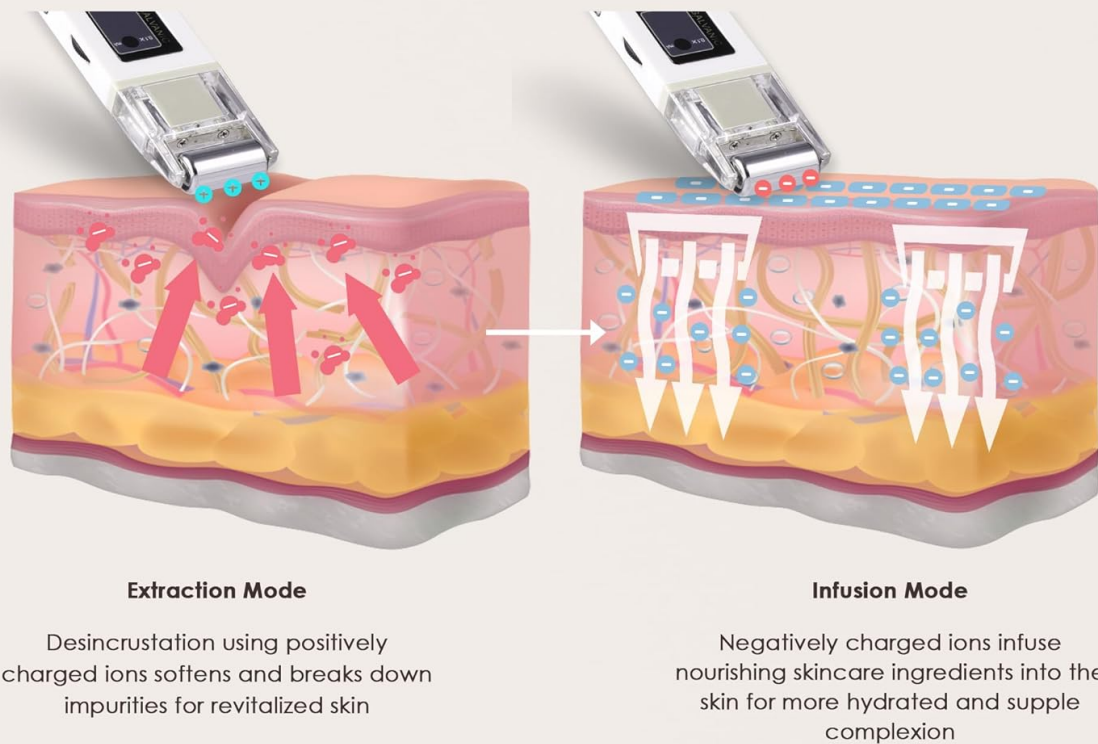


Figure 2: The device with the ball probe attached, demonstrating its compact design and highlighting its dual function and multiple intensity settings.

Setup

1. **Cleanse Skin:** Thoroughly wash and dry your face before beginning any treatment.
2. **Disinfect Probe:** Use 70% alcohol to disinfect the chosen probe (roller or ball) before each use.
3. **Attach Probe:** The probes are attached by aligning them with the unit and pressing down onto the USB-like jack until securely in place.
4. **Power Source:** The device can be connected to an electrical outlet for continuous use or charged for wireless operation.

Operating Instructions

The device features two main modes: Desincrustation (Extraction) and Iontophoresis (Infusion).

Premium Beauty Technology

Powerful extract + infusion therapy

1. EXTRACTION MODE

Negative galvanic current softens skin, extracting impurities, oil, dirt, and dead skin, revitalizing congested pores for better skincare absorption.

2. INFUSION MODE

Positive galvanic current boosts deeper skincare penetration, enhancing hydration and nourishment for a radiant glow.



Figure 3: Visual representation of the Extraction Mode (negative galvanic current) and Infusion Mode (positive galvanic current).

Mode 1: Disincrustation (Deep Cleansing)

This mode uses negative galvanic current to soften skin, emulsify sebum and keratin, and prepare pores for deep cleansing by eliminating dead skin cells, makeup, and dust.

1. Press the **ON/OFF** button to power on the device.
2. Press the **Negative (-)** button to activate the Disincrustation mode. The red LED light will illuminate, indicating the device is operating in this mode.
3. Gradually increase the intensity to a comfortable level. You do not need to use the highest setting to achieve results.
4. Ensure your fingers are touching the side electrodes of the device for proper current flow. Place the probe firmly on the skin.
5. Move the device in upward motions across the treatment areas. Never move the unit downwards.
 - **Forehead:** Move from your eyebrows upward to the hairline.
 - **Cheeks:** Begin at the nose and move across the cheek towards the ear.
 - **Chin & Neck:** Move down the chin to the neck, then from the bottom of the neck up to the face and then the forehead.

Mode 2: Iontophoresis (Nutrient Infusion)

This mode uses positive galvanic current to boost deeper skincare penetration, enhancing hydration and nourishment for a radiant glow.

1. Decrease the intensity to 0.5.
2. Apply an ion gel or your preferred serum/cream to the areas you will be working on.
3. Press the **Positive (+)** button. The blue LED light will illuminate, indicating the device is operating in this mode.
4. Adjust the intensity, starting at a lower level, to a comfortable setting.
5. Ensure your fingers are touching the side electrodes and place the probe firmly on the skin.
6. Move the device in the same upward directions as in the Disincrustation mode:
 - **Forehead:** Move from your eyebrows upward to the hairline.
 - **Cheeks:** Begin at the nose and move across the cheek towards the ear.
 - **Chin & Neck:** Move down the chin to the neck, then from the bottom of the neck up to the face and then the forehead.
 - **Around Eyes/Temples:** Slowly move from the chin to the bottom of the ear, lift off the skin, and move from the chin to the corners of the mouth. Lastly, move from the corners of the eyes to the temples, to the top of the neck. Lift the device and move from the chin to the top of the neck and down to the armpits.

Treatment Frequency & Duration

Each treatment session should last for 10 minutes. For initial results, perform treatments four times a week. Once desired goals are achieved, reduce usage to twice a week for maintenance.

Video 1: An instructional video demonstrating the usage of the Project E Beauty Portable Rechargeable Galvanic Device, including probe attachment, mode selection, and proper movement techniques for facial treatment.

Maintenance

After each use, clean the probes and the device surface with 70% alcohol to ensure hygiene and longevity of the product. Store the device in a cool, dry place.

Troubleshooting

If you experience any discomfort or a sensation of shock during use, ensure that you have applied sufficient conductive gel or serum to the skin and that your fingers maintain proper contact with the side electrodes. The device requires a complete circuit to function correctly. If issues persist, discontinue use and contact customer support.

Specifications

Item Model Number	PE094
Product Dimensions	7.48 x 1.97 x 0.59 inches
Item Weight	0.15 Kilograms (5.29 ounces)
Color	White




Power Source	Rechargeable (can be used plugged in or wirelessly)
Intensity Settings	11 levels




Warranty and Support

The Project E Beauty Portable Rechargeable Galvanic Device typically comes with a 1-year warranty. For any product inquiries, technical support, or warranty claims, please visit the official Project E Beauty store or contact their customer service directly.

[Visit the Project E Beauty Store on Amazon](#)

Related Documents - PE094

	<p>ageLOC Galvanic Spa: Manuale Utente per Cura della Pelle e Capelli</p> <p>Scopri il dispositivo ageLOC Galvanic Spa di Nu Skin per migliorare l'aspetto della pelle e dei capelli. Questa guida fornisce istruzioni dettagliate sull'uso, la manutenzione e la garanzia per ottenere risultati visibili.</p>
	<p>Nu Skin ageLOC Galvanic Spa Owner's Manual: Enhance Skin & Hair Vitality</p> <p>Explore the Nu Skin ageLOC Galvanic Spa, a personal care device designed to visibly improve skin and hair appearance. Learn how to use it for anti-aging benefits, radiance, and healthier-looking hair.</p>
	<p>Photon Skin Rejuvenation Face & Neck Mask (PE706) User Manual</p> <p>User manual for the Project E Beauty Photon Skin Rejuvenation Face & Neck Mask (Model No. PE706). Learn about LED photon therapy, its benefits, safety precautions, and how to use the device for optimal skin rejuvenation.</p>

	<p>Project E Beauty VANA PE711 Blackhead Removal Pore Vacuum User Guide</p> <p>User guide for the Project E Beauty VANA PE711 Blackhead Removal Pore Vacuum, detailing safety precautions, features, usage instructions, probe types, technical specifications, and disposal information.</p>
	<p>FOREO LUNA 2 Professional User Manual: High-Performance Deep Cleansing & Firming Facial Massage</p> <p>Discover the FOREO LUNA 2 Professional with this comprehensive user manual. Learn about its advanced T-Sonic™ technology for deep cleansing and firming facial massage, step-by-step usage, cleaning, charging, troubleshooting, and warranty information.</p>
	<p>PALSAR7 Galvanic and Cleansing Face Iron: User Manual & Guide</p> <p>Comprehensive user manual for the PALSAR7 Galvanic and Cleansing Face Iron. Learn how to use its modes (Cleansing, Nourishing, Warming), LED light therapy, and maintenance for optimal skin care results.</p>

Documents - Project E Beauty – PE094



[\[pdf\]](#) User Manual

20210331 PE094 User Manual All Download Multiple Languages Project E Beauty Galvanic Facial Lifting Roller SKU v 1618275155 cdn shopify s files 1 0732 7299 |||

...

lang: **score:22** filesize: 6.87 M page_count: 96 document date: 2021-04-12