

Manuals+

[Q & A](#) | [Deep Search](#) | [Upload](#)

[manuals.plus](#) /

> [Homedics](#) /

> [Homedics Me My Elos Soft Quartz Lamp Cartridge 120,000 Light Pulses Instruction Manual](#)

Homedics 4051123000909

Homedics Me My Elos Soft Quartz Lamp Cartridge 120,000 Light Pulses Instruction Manual

Model: 4051123000909

INTRODUCTION

This document provides instructions for the Homedics Me My Elos Soft Quartz Lamp Cartridge, designed as a replacement for compatible hair removal devices. This cartridge utilizes patented elōs Quartz technology to deliver 120,000 light pulses.

The cartridge is intended for use with the following devices:

- ME Smooth Hair Removal Device
- Illuminage Touch
- ME Soft
- ME Touch
- Tanda ME Elos



Image: Packaging for the Me My Elos Soft Quartz Lamp Cartridge.

SETUP AND INSTALLATION

The Me My Elos Soft Quartz Lamp Cartridge is designed for easy replacement. Ensure your hair removal device is turned off and unplugged before proceeding.

1. Locate the cartridge slot on your compatible hair removal device. Refer to your device's original instruction manual for specific guidance on cartridge removal.
2. Carefully remove the old cartridge by following the device's instructions.
3. Unpack the new Homedics Me My Elos Soft Quartz Lamp Cartridge.
4. Align the new cartridge with the slot and gently push it into place until it securely snaps. Do not force the cartridge.
5. Verify that the cartridge is properly seated before plugging in and powering on your device.



Image: The Me My Elos Soft Quartz Lamp Cartridge, ready for installation.

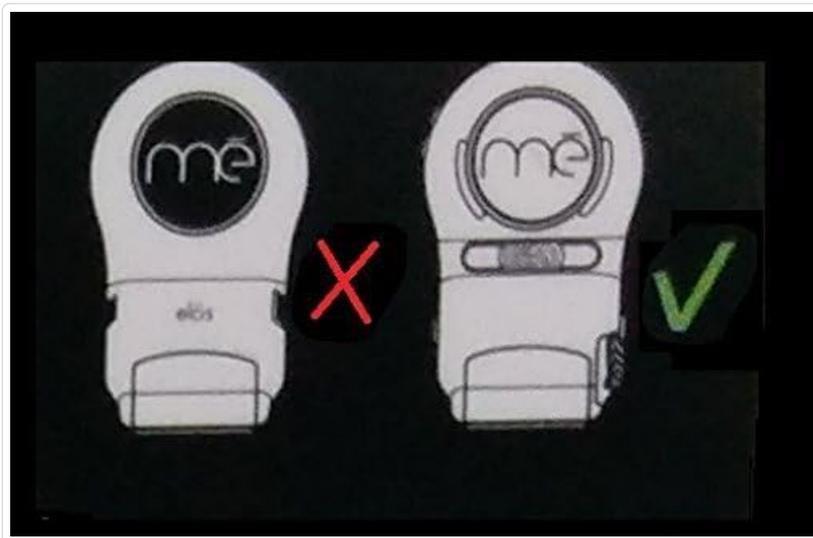


Image: Compatibility diagram illustrating the correct device shape for cartridge installation.

OPERATING INSTRUCTIONS

Once the new cartridge is installed, your hair removal device is ready for use. The cartridge functions as the light source for the elos technology.

- Refer to your specific hair removal device's instruction manual for detailed operating procedures, safety guidelines, and treatment schedules.
- The cartridge provides 120,000 light pulses, designed for long-term use.

- Ensure the treatment window of the cartridge is clean before each use.



Image: A Me My Elos hair removal device with the cartridge installed, shown on a black base.

MAINTENANCE

Proper maintenance of the cartridge ensures optimal performance and longevity.

- **Cleaning:** After each use, ensure the light emission window of the cartridge is clean. Gently wipe it with a soft, dry cloth. Do not use abrasive cleaners or solvents.
- **Storage:** When not in use, store the cartridge within the device or in its original packaging in a cool, dry place, away from direct sunlight and extreme temperatures.
- **Inspection:** Periodically inspect the cartridge for any visible damage or debris. If any damage is observed, discontinue use and contact customer support.

TROUBLESHOOTING

If you encounter issues with your Me My Elos Soft Quartz Lamp Cartridge, consider the following troubleshooting steps:

- **Cartridge Not Recognized:**
 - Ensure the cartridge is fully and correctly inserted into the device. Remove and re-insert it firmly.
 - Verify that your device is compatible with this specific cartridge model (refer to the "Setup and Installation" section).

- **No Light Pulses:**

- Check if the device is properly powered on and connected to an electrical outlet.
- Confirm that the cartridge has not exceeded its 120,000 pulse limit.
- Inspect the light emission window for any obstructions or visible damage (e.g., a "defective bubble" or "blown fuse" appearance in the lamp). If damage is present, the cartridge may need replacement.
- Consult your hair removal device's manual for device-specific error indicators or troubleshooting steps.

- **Reduced Effectiveness:**

- Ensure the treatment window is clean and free of debris.
- Follow the recommended treatment schedule and technique outlined in your device's manual.

If problems persist after following these steps, please contact customer support for further assistance.

SPECIFICATIONS

Model Number	4051123000909
Pulses	120,000
Product Dimensions	5.12 x 3.54 x 1.18 inches
Item Weight	1.76 ounces
Manufacturer	Syneron
ASIN	B00CSKCDNM

WARRANTY AND SUPPORT

Specific warranty information for the Homedics Me My Elos Soft Quartz Lamp Cartridge may be provided with your purchase documentation or by the manufacturer. Please retain your proof of purchase for warranty claims.

For technical support, product inquiries, or warranty service, please contact the manufacturer, Syneron, or the seller, Illuminage Me. Contact details can typically be found on the product packaging or their official websites.

Manufacturer: Syneron

Seller: Illuminage Me