

GIMA 28321

GIMA MT BASE PLUS Low-Frequency Magnetotherapy Device User Manual

Model: 28321

1. INTRODUCTION

The GIMA MT BASE PLUS is a low-frequency magnetotherapy device specifically designed and manufactured for convenient home applications. This device utilizes pulsed magnetic fields to support various therapeutic needs.

The primary effects of low-frequency magnetic fields include:

- **Anti-edema:** Helps reduce swelling.
- **Analgesic:** Provides pain relief.
- **Anti-phlogistic:** Reduces inflammation.
- **Regenerative:** Supports tissue regeneration.

These effects are primarily applied in fields such as traumatology, rheumatology, angiology, neurology, dermatology, and endocrinology.

The device incorporates an electric current generator with a square wave output, offering a variable frequency from 10 to 99 Hz and adjustable pulse duration. This allows for a pre-determinable field intensity ranging from 10 to 80 Gauss. It also features an audible magnetic field detector and a safety system that halts operation in case of a short circuit.

Both included plates (solenoids) can be connected simultaneously; in this configuration, the magnetic field produced by each plate will be half of the total intensity.



Image 1.1: The GIMA MT BASE PLUS device kit, highlighting its suitability for home magnetotherapy applications.

2. SAFETY INFORMATION

Please read all safety instructions carefully before using the device. Failure to do so may result in injury or device malfunction.

- **Consult a Physician:** Always consult your doctor or a qualified healthcare professional before starting any magnetotherapy treatment, especially if you have pre-existing medical conditions.

- **Contraindications:** Magnetotherapy generally has no specific contraindications; however, its use is **NOT RECOMMENDED** for individuals with pacemakers or other implanted electronic devices, and for pregnant women.
- **Children:** Keep the device out of reach of children.
- **Electrical Safety:** Do not use the device near water or in damp environments. Ensure the power cord is not damaged.
- **Device Integrity:** Do not attempt to open or repair the device yourself. Contact qualified service personnel for any maintenance or repair needs.
- **Short Circuit Protection:** The device is equipped with a safety system that blocks distribution in case of a short circuit. If this occurs, immediately unplug the device and check all connections before attempting to restart. If the issue persists, contact customer support.

3. PRODUCT OVERVIEW

The GIMA MT BASE PLUS system includes the main control unit, two solenoids (applicator plates), and a carrying case for portability and storage.



Image 3.1: Complete GIMA MT BASE PLUS kit, including the main unit, solenoids, and carrying case.

Le applicazioni con la Magnetoterapia si effettuano con i solenoidi (trasduttori) in dotazione che trasformano l'impulso elettrico in campo magnetico



Image 3.2: The main control unit and the two solenoids (transducers) that generate the magnetic field.

Components:

- **Main Control Unit:** Features a digital display for frequency and intensity settings, and control buttons.
- **Solenoids (Applicator Plates):** Two transducers that convert electrical impulses into a magnetic field for application.
- **Power Adapter:** For connecting the device to a power source.
- **Carrying Case:** For easy transport and protection of the device.

4. SETUP

1. **Unpack:** Carefully remove all components from the carrying case.
2. **Placement:** Place the main control unit on a stable, flat surface.
3. **Connect Solenoids:** Plug the cables from the solenoids into the designated ports on the main control unit. Ensure a secure connection.
4. **Power Connection:** Connect the power adapter to the main control unit, then plug it into a suitable electrical outlet (230 V / 50-60 Hz).
5. **Initial Check:** Before starting treatment, ensure all connections are firm and the device is powered on.

5. OPERATING INSTRUCTIONS

The GIMA MT BASE PLUS is designed for simple, fast, and lightweight operation, offering professional magnetotherapy at home.



Semplice, veloce e leggero

La potenza di una magnetoterapia professionale, dove vuoi tu

Image 5.1: Applying magnetotherapy to the knee using the GIMA MT BASE PLUS.

1. **Power On:** Turn on the device using the power button. The digital displays will illuminate.
2. **Set Frequency:** Use the control buttons to adjust the frequency (Hz) to the desired value, ranging from 10 to 99 Hz.
3. **Set Intensity:** Use the control buttons to adjust the magnetic field intensity (Gauss) to the desired value, ranging from 10 to 80 Gauss.
4. **Position Solenoids:** Place the solenoids directly on the area of the body requiring treatment. They can be secured with elastic bands or wraps if needed.
5. **Audible Field Detector:** The device has an integrated audible magnetic field detector. A sound indicates that the magnetic field is being generated. If no sound is heard, check connections and settings.
6. **Simultaneous Use:** If using both solenoids simultaneously, be aware that the intensity delivered by

each will be half of the total set intensity.

7. **Treatment Duration:** The duration of a session varies depending on the individual, the condition being treated, and its severity. Always follow your healthcare professional's recommendations.
8. **Power Off:** After the treatment session, turn off the device and unplug it from the power outlet.



Image 5.2: Applying magnetotherapy to the arm/wrist area.

6. APPLICATIONS

Magnetotherapy is a form of physiotherapy that uses electromagnetic waves for therapeutic purposes. The GIMA MT BASE PLUS is suitable for a wide range of applications, including:



Image 6.1: Overview of the effects and primary applications of low-frequency magnetic fields.



Image 6.2: Specific applications of magnetotherapy, including sports injuries and various pathologies.

Common Applications:

- **Traumatology:** For injuries such as sprains, strains, and fractures.
- **Rheumatology:** For conditions like arthritis and other rheumatic pains.
- **Angiology:** To improve circulation.
- **Neurology:** For certain neurological conditions.
- **Dermatology:** For skin-related issues.
- **Endocrinology:** For endocrine system support.
- **Sports Injuries:** Recovery and pain management for athletes.
- **Foot Pathologies:** Addressing various foot-related conditions.
- **Cervical Spine:** Relief for neck pain and related issues.
- **Knee Pathology:** Support for knee conditions and pain.
- **General Pain Relief:** Effective for various types of localized pain.
- **Low Back Pain:** Alleviating discomfort in the lumbar region.

7. MAINTENANCE AND CARE

- **Cleaning:** Disconnect the device from the power supply before cleaning. Wipe the main unit and solenoids with a soft, slightly damp cloth. Do not use abrasive cleaners or solvents.
- **Storage:** Store the device and its components in the provided carrying case in a cool, dry place, away from direct sunlight and extreme temperatures.
- **Inspection:** Periodically inspect the power cord and solenoid cables for any signs of damage. If damage is found, discontinue use and contact customer support.

8. TROUBLESHOOTING

- **Device Does Not Power On:**
 - Ensure the power adapter is securely connected to both the device and the electrical outlet.
 - Check if the electrical outlet is functional by plugging in another appliance.
- **No Magnetic Field Detected (No Audible Sound):**
 - Verify that the solenoids are correctly plugged into the main unit.
 - Check that the frequency and intensity settings are not set to zero.
 - Ensure the device is powered on and functioning.
- **Device Stops During Operation (Short Circuit Protection):**
 - Immediately unplug the device from the power source.
 - Inspect all cables and connections for any visible damage or loose connections.
 - If no obvious issue is found, wait a few minutes, then reconnect and restart. If the problem persists, contact customer support.

For any issues not covered here, please contact GIMA customer support.

9. TECHNICAL SPECIFICATIONS



Magnetoterapia

- ✓ Magnetoterapia e sport
- ✓ Magnetoterapia e patologie del piede
- ✓ Magnetoterapia e fratture
- ✓ Magnetoterapia applicata alla cervicale
- ✓ Magnetoterapia e patologia del ginocchio
- ✓ Magnetoterapia e dolore
- ✓ Magnetoterapia e lombalgie

Image 9.1: Detailed technical specifications of the GIMA MT BASE PLUS device.

Feature	Specification
Model Number	28321
Waveform	Square
Max Emission/Pause Ratio	50%
Adjustable Magnetic Field	10 to 80 Gauss
Adjustable Frequency	10 to 99 Hz
Magnetic Field Intensity	80 Gauss (8 mT)
Operating Voltage	230 V / 50-60 Hz
Dimensions (L x W x H)	193 x 138 x 68 mm
Weight	1 kg
Safety Class	II BF - Norms: EN60601-1 (CEI)

Manufacturer	LED SPA (Brand: GIMA)
Country of Origin	Italy

10. WARRANTY AND SUPPORT

For warranty information, service, or technical support, please contact GIMA SPA, the manufacturer of this device. Please refer to the product packaging or the official GIMA website for the most up-to-date contact details and warranty terms.

Manufacturer Contact: Gima SPA

It is recommended to keep your purchase receipt as proof of purchase for warranty claims.

© 2024 GIMA. All rights reserved. This manual is for informational purposes only.

Related Documents - 28321

	GIMA KURA TENS AD-2126 User Manual and Operation Guide Comprehensive user and operation guide for the GIMA KURA TENS AD-2126 device, detailing intended use, specifications, setup, operation, troubleshooting, and maintenance for pain relief.
	GIMA Safety Belt - Patient Restraint Device GIMA offers a range of safety belts designed to secure patients and prevent them from sliding or falling from stretchers. Available in various models with different lengths and materials, these belts are essential for patient safety during transport.
	GIMA Pressure Activated Safety Lancets (PA) - Medical Device for Capillary Blood Sampling Information on GIMA Pressure Activated Safety Lancets (Model PA), a sterile, single-use device for capillary blood sampling. Includes intended use, instructions, warnings, storage, and disposal.

[illegible]

[\[pdf\]](#)

printgimaitaly print sku 28321 lcid 1040 |||

MAGNETOTERAPIA BASE Codice: **28321** Categoria: Magnetoterapia Unit di vendita:

1 pz. Quantitativi minimi: 1 Dispositivo: Classe: Dispositivo medico II A NSIS: 11709

CND: Z120606 EAN13: 8033638326299 Descrizione: Caratteristiche tecniche:

Dotazione standard: GIMA MT BASE PLUS - magneto...

lang:it score:24 filesize: 97.8 K page_count: 1 document date: 2025-09-11

[illegible]

[\[pdf\]](#) User Manual Accessories

011 Electromedical Dev 2015 indd macpro MT BASE MAGNETOTHERAPY Manual GB IT 0051 Intensity

of magnetic field 80 Gauss 8mT Operating voltage 230V 50 60 Hz Size 193x138xh 68 mm Weight 1 kg

Safety class II BF Norms EN60601 CEI Multilanguage manual and box FR ES PT DE HU SK LV CZ EE

PL RU BG Adjustable angle Box size 185x140x225 150 eng 2018 19 gimaitaly DocumentiGIMA Pagine

EN

PHYSIOTHERAPY - NEBULIZERS - SUCTION PUMPS MT BASE:

MAGNETOTHERAPY **28321** GIMA MT BASE PLUS - magnetotherapy Magneto Base

Plus is an ELF low frequency magnetotherapeutic unit designed for domiciliary

applications. The most important effects of the magnetic fields are: Antiedemigen /

Analgesic / ...

lang:en score:22 filesize: 1.12 M page count: 1 document date: 2017-06-29



GIMA

CATALOGUE INTERNATIONAL

PRODUITS PROFESSIONNELS MEDICAUX

1924



TESTS DE DIAGNOSTIC

PROTECTION DES SOIGNANTS

INSTANTANÉ MÉTÉOROLOGIQUE

ET CHIMIQUE DES URINES

TROUSSES, SACS ET RESULTATS

MOSCAUX

BALANCES ET

APPAREILS DE MESURE

SOIN, AMPLIFIET ET

AGE A LA TENDANCE

MOBILIER MEDICAL

MATERNITES

NEBULISATEURS, ASPIRATEURS

OPHTHALMOLOGIE

PROCTOLOGIE

ELECTROCARDIOGRAPHES

CATHETERS

LIQUIDES MEDICAUX

ET ENDOSCOPIE

ECLAIRAGE MEDICAL

DIAGNOSTIC

AUDIMETRIE

ET PROTHETIQUE

ECLIMATRIE

ET THERAPIES

ULTRASONES

ET PREMIERE SECOURS

STIMULATION

APPAREILS ET INSTRUMENTS

ET TENDANCE

MOBILIER MEDICAL

ET TENDANCE



STETHOSCOPES ENVOIE



CHARGES MULTIFONCTIONS




page 158

[GIMA International Catalogue 2024 - Professional Medical Products](#)

The GIMA International Catalogue 2024 presents over 9,500 professional medical products, including diagnostic tests, medical instruments, patient care equipment, and protective gear. GIMA emphasizes fast delivery, quality, innovation, and comprehensive customer service across 145 countries.

lang:fr score:21 filesize: 51.43 M page_count: 468 document date: 2023-06-13

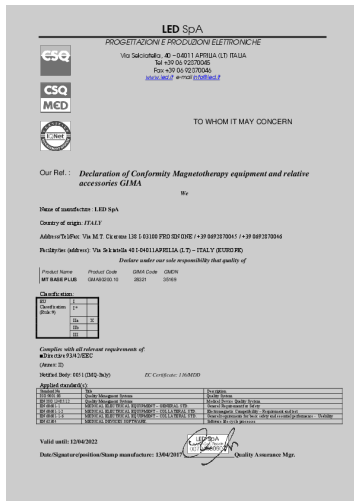

[pdf] Decleration of Conformity Accessories

Microsoft Word DeclConf GIMA MAGNETOTHERAPY EQUIPMENT 04 2017 Claudio Mauti Declaration of conformity MT Base Plus magnetterapi Domiciliary Swemed EU d28321 swemed eu docs gima ||| LED SpA PROGETTAZIONI E PRODUZIONI ELETTRONICHE Via Selciatella, 40 04011 APRILIA LT ITALIA Tel ... esponsibility that quality of Product Name MT BASE PLUS Product Code GMA80200.10 GIMA Code GMDN **28321** 35169 Classification: EU I Classification I* Rule:9 IIa X IIb III Complies with all ... lang:en **score:21** filesize: 84.97 K page_count: 1 document date: 2018-01-31



[pdf] Decleration of Conformity Accessories

Microsoft Word DeclConf GIMA MAGNETOTHERAPY EQUIPMENT 04 2017 Claudio Mauti d28321gimaitaly DocumentiGIMA DC d28321 ||| LED SpA PROGETTAZIONI E PRODUZIONI ELETTRONICHE Via Selciatella, 40 04011 APRILIA LT ITALIA Tel ... esponsibility that quality of Product Name MT BASE PLUS Product Code GMA80200.10 GIMA Code GMDN **28321** 35169 Classification: EU I Classification I* Rule:9 IIa X IIb III Complies with all ... lang:en **score:20** filesize: 84.97 K page_count: 1 document date: 2018-01-31



[pdf] User Manual Accessories

MAGNETOTHERAPY UNIVERSAL ELASTIC BELTS Gima 28 giu 2023 — 28298 ORTHOMAG low frequency magnetotherapy kit A professional PEMF device that can be used everywhere Lightweight 191 ENG gimaitaly DocumentiGIMA Pagine EN ||| PHYSIOTHERAPY - NEBULIZERS - SUCTION PUMPS ORTHOMAG LOW FREQUENCY MAGNETOTHERAPY - 2 PROGRAMMES O ... n GB, FR, IT, ES . On request DE. 28296 Muscle pain Anti-inflammatory MT BASE: MAGNETOTHERAPY **28321** GIMA MT BASE PLUS - magnetotherapy Magneto Base Plus is an ELF low frequency magnetotherapeutic ... lang:en **score:20** filesize: 2.22 M page_count: 1 document date: 2023-06-13



[\[pdf\]](#) User Manual Specifications Catalog

011 Electromedical ING macpro suction pumps FOOT WARMER HEATING BLANKETS 14 giu 2023 — Electrode position manual and user GB FR IT ES On request DE PT Size 150x90x35 mm 0068 6 language software Page 7 08 PHYSIOTHERAPY NEBULIZERS SUCTION PUMPS ENG 2024 gimaitaly en assets cataloghi |||

PHYSIOTHERAPY - NEBULIZERS - SUCTION PUMPS HEATING BLANKETS FOOT WARMER 28660 28660 HEATING UNDE ... n GB, FR, IT, ES . On request DE. 28296 Muscle pain Anti-inflammatory MT BASE: MAGNETOTHERAPY **28321** GIMA MT BASE PLUS - magnetotherapy Magneto Base Plus is an ELF low frequency magnetotherapeutic ...

lang:en **score:17** filesize: 5 M page_count: 24 document date: 2023-07-28



[\[pdf\]](#) Catalog

011 Electromedical ITA andrea Catalogo ELETTRODO GOMMA 50 x mm ET2 Elettrodi in gomma per IONO BASE PLUS IONOFORESIS GALVANICA cod 28306 124 ita gimaitaly it assets prodotti ||| ELETTROMEDICALI IONO BASE PLUS - IONOFORESIS GALVANICA 28306 IONO BASE - PLUS - 2 canali Generator ... uscita: 60 mA Classe: II BF - 93/42/EEC. Produzione Italiana. 0051 28308 MAGNETOTERAPIA BASE **28321** GIMA MT BASE PLUS - magnetoterapia Il Magneto Base un apparato per magnetoterapia a bassa fr...

lang:it **score:16** filesize: 1.73 M page_count: 1 document date: 2014-05-15



[\[pdf\]](#) Specifications Accessories

011 Electromedical ING mbpro IONO BASE PLUS GALVANIC IONOPHORESIS ELECTROMEDICAL DEVICES 124 28308 STIM 12 programs 2 channels Bipolar pulsing square wave current generator for muscle toning and stimulation as well pain relief eng gimaitaly en assets prodotti ||| ELECTROMEDICAL DEVICES IONO BASE PLUS - GALVANIC IONOPHORESIS 28306 IONO BASE - PLUS - 2 channels ... t Current: 60 mA Safety class: II BF - 93/42/EEC. Made in Italy. 28308 MT BASE: MAGNETOTHERAPY **28321** GIMA MT BASE PLUS - magnetotherapy Magneto Base Plus is an ELF low frequency magnetotherapeutic...

lang:en **score:16** filesize: 1.73 M page_count: 1 document date: 2014-05-15



[\[pdf\]](#) User Manual Catalog

untitled AIR COMPRESSION THERAPY SOLUTIONS Manual GB FR IT ES This system can operate with 3 compression modes A B C 28435 LEG 93x32 cm 28434 EXTENSION ZIPPER 10 extra extension 28440 DOUBLE HOSE to use two leg cuffs at the same time spare 28436 ARM CUFF unique size 86x29 28437 WAIST 130x38 28442 L 6 CHAMBERS 09 PHYSIOTHERAPY NEBULIZERS SUCTION PUMPS 2018 19 ENG gimaitaly en assets cataloghi 80 PHYSIOTHERAPY - NEBULIZERS - SUCTION PUMPS AIR COMPRESSION THERAPY SOLUTIONS All human body is sup ... sectors without solenoids PHYSIOTHERAPY - NEBULIZERS - SUCTION PUMPS MT BASE: MAGNETOTHERAPY **28321** GIMA MT BASE PLUS - magnetotherapy Magneto Base Plus is an ELF low frequency magnetotherapi...

lang:en **score:16** filesize: 8.18 M page_count: 22 document date: 2017-06-30



[\[pdf\]](#) Dimension Guide

011 Electromedical FRA andrea ÉLECTROSTIMULATEUR NUMÉRIQUE 12 PROGRAMMES 28308 STIM BASE PLUS programmes 2 canaux Générateur de courant à forme d onde carrée et impulsion bipolaire pour la stimulation tonification 124 fra gimaitaly fr assets prodotti ||| APPAREILS ELECTROMEDICAUX IONO BASE PLUS - IONOPHORESE GALVANIQUE 28306 IONO BASE - plus - 2 canaux ... ax.: 60 mA Classe de scurit II BF - 93/42/CEE. Fabriqué en Italie. 28308 MT BASE: MAGNTOTHRAPIE **28321** GIMA MT BASE PLUS - magntothrapie Appareil de magntothrapie ELF basse frquence destin aux app...

lang:fr **score:16** filesize: 1.73 M page_count: 1 document date: 2014-05-15



[\[pdf\]](#) Catalog

GIMA catalogue part 1 IV GIMAcatalogue part1 gimaitaly 80 en assets cataloghi HEALTH-CARE - PHARMACY Pill Splitter, Pill Pulverizer These easy tools are ideal for elderly peopl ... 60601-1 CE162-5 - 93/42 EEC Plates and Cylinder are not included in standard configuration. Code **28321** GIMA MT BASE PLUS - magnetotherapy Magneto Base is an ELF low frequency magnetotherapeutic unit...

lang:en **score:15** filesize: 4.01 M page_count: 88 document date: 2007-07-10



[\[pdf\]](#) Catalog

Catalogo Generale GIMA 2017 ITA gimaitaly it assets cataloghi

CATALOGO GENERALE ARTICOLI APPARECCHI PER MEDICINA iHEALTH PRODOTTI WIRELESS DIAGNOSTICa rapida PR ... tale e 2 settori passivi senza solenoidi ELEtTROMEDICALi - fisioterapia magnetoterapia base **28321** GIMA MT BASE PLUS - magnetoterapia Il Magneto Base un apparato per magnetoterapia a bassa fr...

lang:it **score:14** filesize: 94.98 M page_count: 387 document date: 2017-05-10



[\[pdf\]](#) Catalog

untitled echipamente medicale gima Index of upload prezentare med shop ro |||

2017 2018 international CATALOGue PROFESSIONAL MEDICAL PRODUCTS iHEALTH WIRELESS PRODUCTS DIAGNOSTIC ... h sector, totally 6 solenoids and 2 passive sectors without solenoids MT BASE: MAGNETOTHERAPY **28321** GIMA MT BASE PLUS - magnetotherapy Magneto Base Plus is an ELF low frequency magnetotherapi...

lang:en **score:14** filesize: 38.76 M page_count: 353 document date: 2016-12-28



[\[pdf\]](#) Catalog

untitled GIMA International Catalogue 2020 ENG gimaitaly en assets cataloghi |||

INTERNATIONAL CATALOGUE PROFESSIONAL MEDICAL PRODUCTS 2020 NEWS 2020 Cautapldoagtuue n3e6wpapgrees ... sectors without solenoids PHYSIOTHERAPY - NEBULIZERS - SUCTION PUMPS MT BASE: MAGNETOTHERAPY **28321** GIMA MT BASE PLUS - magnetotherapy Magneto Base Plus is an ELF low frequency magnetotherapi...

lang:en **score:14** filesize: 67.64 M page_count: 422 document date: 2019-11-11



[\[pdf\]](#) Catalog

Gima articoli medicali apparecchi per medicina elettromedicali endoscopi misuratori arredo ambulatorio elettrobisturi il pronto soccorso la sterilizzazione GIMA International Catalogue 2023 LR ENG gimaitaly en assets cataloghi |||

INTERNATIONAL CATALOGUE PROFESSIONAL MEDICAL PRODUCTS 2023
PREVIEW DIAGNOSTIC TESTS LABORATORY OPE ... tis Fractures Articular pain
Muscle pain 10 Hz 20 Hz Anti-inflammatory MT BASE: MAGNETOTHERAPY **28321**
GIMA MT BASE PLUS - magnetotherapy Magneto Base Plus is an ELF low
frequency magnetotherapeutic ...

lang:en **score:14** filesize: 76.27 M page_count: 468 document date: 2022-11-08



[\[pdf\]](#) SAR Rating Catalog

untitled Gima articoli medicali apparecchi per medicina elettromedicali endoscopi misuratori arredo ambulatorio elettrobisturi il pronto soccorso la sterilizzazione Catalogo Generale GIMA 2024 LR ITA gimaitaly it assets cataloghi |||

CATALOGO GENERALE ARTICOLI APPARECCHI PER MEDICINA TEST
DIAGNOSTICI LABORATORIO PROTEZIONE OPERATO ... lori articolari
Autoscansione Dolori muscolari 28296 Anti-infiammatorio MAGNETOTERAPIA MT
BASE **28321** GIMA MT BASE PLUS - magnetoterapia Apparecchio per
magnetoterapia a bassa frequenza, progettato...

lang:it **score:14** filesize: 51.44 M page_count: 468 document date: 2023-06-13



[\[pdf\]](#) Dimension Guide Catalog

011 Electromedical FRA macpro CHAUFFE PIEDS LITS Gima ITALYChauffe lits électriques en tissu non tissé polyester estampillé 41706 BIBERON 150 ml avec bouchon set de 3 rechange 41705 CAPUCHON POUR 08 PHYSIOTHERAPIE NEBULISEURS ASPIRATEURS 2024gimaitaly fr assets cataloghi 08 2024 |||

PHYSIOTHRAPIE - NBULISEURS - ASPIRATEURS CHAUFFE-LITS CHAUFFE-PIEDS 28660 28660 CHAUFFE-LIT - in ... es Auto-balayage Douleurs musculaires Anti-inflammatoire Sur demande DE. MT BASE: MAGNTOTHRAPIE **28321** GIMA MT BASE PLUS - magntothrapie Appareil de magntothrapie ELF basse frquence destin aux a...

lang:fr **score:14** filesize: 16.11 M page_count: 24 document date: 2023-06-13



[\[pdf\]](#) Catalog

untitled Gima articoli medicali apparecchi per medicina elettromedicali endoscopi misuratori arredo ambulatorio elettrobisturi il pronto soccorso la sterilizzazione GIMA International Catalogue 2024 LR ENG gimaitaly en assets cataloghi |||

INTERNATIONAL CATALOGUE PROFESSIONAL MEDICAL PRODUCTS 2024
DIAGNOSTIC TESTS LABORATORY OPERATOR S ... n GB, FR, IT, ES . On request DE. 28296 Muscle pain Anti-inflammatory MT BASE: MAGNETOTHERAPY **28321** GIMA MT BASE PLUS - magnetotherapy Magneto Base Plus is an ELF low frequency magnetotherapic ...

lang:en **score:14** filesize: 51.45 M page_count: 468 document date: 2023-06-13



[\[pdf\]](#) Catalog

untitled GIMA International Catalogue 2018 19 ENG gimaitaly en assets cataloghi
INTERNATIONAL CATALOGUE PROFESSIONAL MEDICAL PRODUCTS 2018
2019 iHEALTH WIRELESS HEALTH-CARE PRODU ... sectors without solenoids
PHYSIOTHERAPY - NEBULIZERS - SUCTION PUMPS MT BASE:
MAGNETOTHERAPY **28321** GIMA MT BASE PLUS - magnetotherapy Magneto Base Plus is an ELF low frequency magnetotherapi...

lang:en **score:14** filesize: 63.95 M page_count: 388 document date: 2017-07-06



[\[pdf\]](#) Catalog

untitled GIMA Complete catalogue 2017 ENG LD gimaitaly en assets cataloghi
2018 2019 INTERNATIONAL CATALOGUE PROFESSIONAL MEDICAL
PRODUCTS iHEALTH WIRELESS HEALTH-CARE PRODUCT ... sectors without
solenoids PHYSIOTHERAPY - NEBULIZERS - SUCTION PUMPS MT BASE:
MAGNETOTHERAPY **28321** GIMA MT BASE PLUS - magnetotherapy Magneto Base Plus is an ELF low frequency magnetotherapi...

lang:en **score:14** filesize: 40.25 M page_count: 388 document date: 2017-06-14

lang:en score:13 filesize: 1.43 M page count: 152 document date: 2025-06-09

GIMA PRICE LIST 2023					
CODE	DESCRIPTION	UNIT	PRICE	REMARKS	
13000	13000 GIMA PRICE LIST 2023	1	13000		
13001	13001 GIMA PRICE LIST 2023	1	13001		
13002	13002 GIMA PRICE LIST 2023	1	13002		
13003	13003 GIMA PRICE LIST 2023	1	13003		
13004	13004 GIMA PRICE LIST 2023	1	13004		
13005	13005 GIMA PRICE LIST 2023	1	13005		
13006	13006 GIMA PRICE LIST 2023	1	13006		
13007	13007 GIMA PRICE LIST 2023	1	13007		
13008	13008 GIMA PRICE LIST 2023	1	13008		
13009	13009 GIMA PRICE LIST 2023	1	13009		
13010	13010 GIMA PRICE LIST 2023	1	13010		
13011	13011 GIMA PRICE LIST 2023	1	13011		
13012	13012 GIMA PRICE LIST 2023	1	13012		
13013	13013 GIMA PRICE LIST 2023	1	13013		
13014	13014 GIMA PRICE LIST 2023	1	13014		
13015	13015 GIMA PRICE LIST 2023	1	13015		
13016	13016 GIMA PRICE LIST 2023	1	13016		
13017	13017 GIMA PRICE LIST 2023	1	13017		
13018	13018 GIMA PRICE LIST 2023	1	13018		
13019	13019 GIMA PRICE LIST 2023	1	13019		
13020	13020 GIMA PRICE LIST 2023	1	13020		
13021	13021 GIMA PRICE LIST 2023	1	13021		
13022	13022 GIMA PRICE LIST 2023	1	13022		
13023	13023 GIMA PRICE LIST 2023	1	13023		
13024	13024 GIMA PRICE LIST 2023	1	13024		
13025	13025 GIMA PRICE LIST 2023	1	13025		
13026	13026 GIMA PRICE LIST 2023	1	13026		
13027	13027 GIMA PRICE LIST 2023	1	13027		
13028	13028 GIMA PRICE LIST 2023	1	13028		
13029	13029 GIMA PRICE LIST 2023	1	13029		
13030	13030 GIMA PRICE LIST 2023	1	13030		
13031	13031 GIMA PRICE LIST 2023	1	13031		
13032	13032 GIMA PRICE LIST 2023	1	13032		
13033	13033 GIMA PRICE LIST 2023	1	13033		
13034	13034 GIMA PRICE LIST 2023	1	13034		
13035	13035 GIMA PRICE LIST 2023	1	13035		
13036	13036 GIMA PRICE LIST 2023	1	13036		
13037	13037 GIMA PRICE LIST 2023	1	13037		
13038	13038 GIMA PRICE LIST 2023	1	13038		
13039	13039 GIMA PRICE LIST 2023	1	13039		
13040	13040 GIMA PRICE LIST 2023	1	13040		
13041	13041 GIMA PRICE LIST 2023	1	13041		
13042	13042 GIMA PRICE LIST 2023	1	13042		
13043	13043 GIMA PRICE LIST 2023	1	13043		
13044	13044 GIMA PRICE LIST 2023	1	13044		
13045	13045 GIMA PRICE LIST 2023	1	13045		
13046	13046 GIMA PRICE LIST 2023	1	13046		
13047	13047 GIMA PRICE LIST 2023	1	13047		
13048	13048 GIMA PRICE LIST 2023	1	13048		
13049	13049 GIMA PRICE LIST 2023	1	13049		
13050	13050 GIMA PRICE LIST 2023	1	13050		
13051	13051 GIMA PRICE LIST 2023	1	13051		
13052	13052 GIMA PRICE LIST 2023	1	13052		
13053	13053 GIMA PRICE LIST 2023	1	13053		
13054	13054 GIMA PRICE LIST 2023	1	13054		
13055	13055 GIMA PRICE LIST 2023	1	13055		
13056	13056 GIMA PRICE LIST 2023	1	13056		
13057	13057 GIMA PRICE LIST 2023	1	13057		
13058	13058 GIMA PRICE LIST 2023	1	13058		
13059	13059 GIMA PRICE LIST 2023	1	13059		
13060	13060 GIMA PRICE LIST 2023	1	13060		
13061	13061 GIMA PRICE LIST 2023	1	13061		
13062	13062 GIMA PRICE LIST 2023	1	13062		
13063	13063 GIMA PRICE LIST 2023	1	13063		
13064	13064 GIMA PRICE LIST 2023	1	13064		
13065	13065 GIMA PRICE LIST 2023	1	13065		
13066	13066 GIMA PRICE LIST 2023	1	13066		
13067	13067 GIMA PRICE LIST 2023	1	13067		
13068	13068 GIMA PRICE LIST 2023	1	13068		
13069	13069 GIMA PRICE LIST 2023	1	13069		
13070	13070 GIMA PRICE LIST 2023	1	13070		
13071	13071 GIMA PRICE LIST 2023	1	13071		
13072	13072 GIMA PRICE LIST 2023	1	13072		
13073	13073 GIMA PRICE LIST 2023	1	13073		
13074	13074 GIMA PRICE LIST 2023	1	13074		
13075	13075 GIMA PRICE LIST 2023	1	13075		
13076	13076 GIMA PRICE LIST 2023	1	13076		
13077	13077 GIMA PRICE LIST 2023	1	13077		
13078	13078 GIMA PRICE LIST 2023	1	13078		
13079	13079 GIMA PRICE LIST 2023	1	13079		
13080	13080 GIMA PRICE LIST 2023	1	13080		
13081	13081 GIMA PRICE LIST 2023	1	13081		
13082	13082 GIMA PRICE LIST 2023	1	13082		
13083	13083 GIMA PRICE LIST 2023	1	13083		
13084	13084 GIMA PRICE LIST 2023	1	13084		
13085	13085 GIMA PRICE LIST 2023	1	13085		
13086	13086 GIMA PRICE LIST 2023	1	13086		
13087	13087 GIMA PRICE LIST 2023	1	13087		
13088	13088 GIMA PRICE LIST 2023	1	13088		
13089	13089 GIMA PRICE LIST 2023	1	13089		
13090	13090 GIMA PRICE LIST 2023	1	13090		
13091	13091 GIMA PRICE LIST 2023	1	13091		
13092	13092 GIMA PRICE LIST 2023	1	13092		
13093	13093 GIMA PRICE LIST 2023	1	13093		
13094	13094 GIMA PRICE LIST 2023	1	13094		
13095	13095 GIMA PRICE LIST 2023	1	13095		
13096	13096 GIMA PRICE LIST 2023	1	13096		
13097	13097 GIMA PRICE LIST 2023	1	13097		
13098	13098 GIMA PRICE LIST 2023	1	13098		
13099	13099 GIMA PRICE LIST 2023	1	13099		
13100	13100 GIMA PRICE LIST 2023	1	13100		
13101	13101 GIMA PRICE LIST 2023	1	13101		
13102	13102 GIMA PRICE LIST 2023	1	13102		
13103	13103 GIMA PRICE LIST 2023	1	13103		
13104	13104 GIMA PRICE LIST 2023	1	13104		
13105	13105 GIMA PRICE LIST 2023	1	13105		
13106	13106 GIMA PRICE LIST 2023	1	13106		
13107	13107 GIMA PRICE LIST 2023	1	13107		
13108	13108 GIMA PRICE LIST 2023	1	13108		
13109	13109 GIMA PRICE LIST 2023	1	13109		
13110	13110 GIMA PRICE LIST 2023	1	13110		
13111	13111 GIMA PRICE LIST 2023	1	13111		
13112	13112 GIMA PRICE LIST 2023	1	13112		
13113	13113 GIMA PRICE LIST 2023	1	13113		
13114	13114 GIMA PRICE LIST 2023	1	13114		
13115	13115 GIMA PRICE LIST 2023	1	13115		
13116	13116 GIMA PRICE LIST 2023	1	13116		
13117	13117 GIMA PRICE LIST 2023	1	13117		
13118	13118 GIMA PRICE LIST 2023	1	13118		
13119	13119 GIMA PRICE LIST 2023	1	13119		
13120	13120 GIMA PRICE LIST 2023	1	13120		
13121	13121 GIMA PRICE LIST 2023	1	13121		
13122	13122 GIMA PRICE LIST 2023	1	13122		
13123	13123 GIMA PRICE LIST 2023	1	13123		
13124	13124 GIMA PRICE LIST 2023	1	13124		
13125	13125 GIMA PRICE LIST 2023	1	13125		
13126	13126 GIMA PRICE LIST 2023	1	13126		
13127	13127 GIMA PRICE LIST 2023	1	13127		
13128	13128 GIMA PRICE LIST 2023	1	13128		
13129	13129 GIMA PRICE LIST 2023	1	13129		
13130	13130 GIMA PRICE LIST 2023	1	13130		
13131	13131 GIMA PRICE LIST 2023	1	13131		
13132	13132 GIMA PRICE LIST 2023	1	13132		
13133	13133 GIMA PRICE LIST 2023	1	13133		
13134	13134 GIMA PRICE LIST 2023	1	13134		
13135	13135 GIMA PRICE LIST 2023	1	13135		
13136	13136 GIMA PRICE LIST 2023	1	13136		
13137	13137 GIMA PRICE LIST 2023	1	13137		
13138	13138 GIMA PRICE LIST 2023	1	13138		
13139	13139 GIMA PRICE LIST 2023	1	13139		
13140	13140 GIMA PRICE LIST 2023	1	13140		
13141	13141 GIMA PRICE LIST 2023	1	13141		
13142	13142 GIMA PRICE LIST 2023	1	13142		
13143	13143 GIMA PRICE LIST 2023	1	13143		
13144	13144 GIMA PRICE LIST 2023	1	13144		
13145	13145 GIMA PRICE LIST 2023	1	13145		
13146	13146 GIMA PRICE LIST 2023	1	13146		
13147	13147 GIMA PRICE LIST 2023	1	13147		
13148	13148 GIMA PRICE LIST 2023	1	13148		
13149	13149 GIMA PRICE LIST 2023	1	13149		
13150	13150 GIMA PRICE LIST 2023	1	13150		
13151	13151 GIMA PRICE LIST 2023	1	13151		
13152	13152 GIMA PRICE LIST 2023	1	13152		
13153	13153 GIMA PRICE LIST 2023	1	13153		
13154	13154 GIMA PRICE LIST 2023	1	13154		
13155	13155 GIMA PRICE LIST 2023	1	13155		
13156	13156 GIMA PRICE LIST 2023	1	13156		
13157	13157 GIMA PRICE LIST 2023	1	13157		
13158	13158 GIMA PRICE LIST 2023	1	13158		
13159	13159 GIMA PRICE LIST 2023	1	13159		
13160	13160 GIMA PRICE LIST 2023	1	13160		
13161	13161 GIMA PRICE LIST 2023	1	13161		
13162	13162 GIMA PRICE LIST 2023	1	13162		
13163	13163 GIMA PRICE LIST 2023	1	13163		
13164	13164 GIMA PRICE LIST 2023	1	13164		
13165	13165 GIMA PRICE LIST 2023	1	13165		
13166	13166 GIMA PRICE LIST 2023	1	13166		
13167	13167 GIMA PRICE LIST 2023	1	13167		
13168	13168 GIMA PRICE LIST 2023	1	13168		
13169	13169 GIMA PRICE LIST 2023	1	13169		
13170	13170 GIMA PRICE LIST 2023	1	13170		
13171	13171 GIMA PRICE LIST 2023	1	13171		
13172	13172 GIMA PRICE LIST 2023	1	13172		
13173	13173 GIMA PRICE LIST 2023	1	13173		
13174	13174 GIMA PRICE LIST 2023	1	13174		
13175	13175 GIMA PRICE LIST 2023	1	13175		
13176	13176 GIMA PRICE LIST 2023	1	13176		
13177	13177 GIMA PRICE LIST 2023	1	13177		
13178	13178 GIMA PRICE LIST 2023	1	13178		
13179	13179 GIMA PRICE LIST 2023	1	13179		
13180	13180 GIMA PRICE LIST 2023	1	13180		
13181	13181 GIMA PRICE LIST 2023	1	13181		
13182	13182 GIMA PRICE LIST 2023	1	13182		
13183	13183 GIMA PRICE LIST 2023	1	13183		
13184	13184 GIMA PRICE LIST 2023	1	13184		
13185	13185 GIMA PRICE LIST 2023	1	13185		
13186	13186 GIMA PRICE LIST 2023	1	13186		
13187	13187 GIMA PRICE LIST 2023	1	13187		
13188	13188 GIMA PRICE LIST 2023	1	13188		
13189	13189 GIMA PRICE LIST 2023	1	13189		
13190	13190 GIMA PRICE LIST 2023	1	13190		
13191	13191 GIMA PRICE LIST 2023	1	13191		
13192	13192 GIMA PRICE LIST 2023	1	13192		
13193	13193 GIMA PRICE LIST 2023	1	13193		
13194	13194 GIMA PRICE LIST 2023	1	13194		
13195	13195 GIMA PRICE LIST 2023	1	13195		
13196	13196 GIMA PRICE LIST 2023	1	13196		
13197	13197 GIMA PRICE LIST 2023	1	13197		
13198	13198 GIMA PRICE LIST 2023	1	13198		
13199	13199 GIMA PRICE LIST 2023	1	13199		
13200	13200 GIMA PRICE LIST 2023	1	13200		

[pdf] Catalog Price List

Cristiano Calendino GIMA PRICE LIST 2023 SPA REVISION from May 1st up to August 31st Unit of Minimum Code Description products sale EURO order 13000 CATALOGUE anonymous 1 ask for price 10 13001 your address 100 13002 4 cover pages personalized 500 13010 WEB with 3 EN Price list gimaitaly en newsletter Gima lists euro |||

GIMA Code 13000 13001 13002 13010 13100 13102 13103 13104 13105 13106 13107 13108 13109 13110 13118 ... EOMAT MATTRESS for magnetotherapy 28319 NEW ASPEED 3 SUCTION ASPIRATOR - 230V single pump - 2 | **28321** GIMA MT BASE PLUS MAGNETOTHERAPY - domiciliary 283

[illegible][illegible]

[\[pdf\]](#) Catalog

011 Electromedical ING mbpro AIR COMPRESSION THERAPY SOLUTIONS 146 PHYSIOTHERAPY
NEBULIZERS SUCTION PUMPS 28410 I TECH LA 500 LASER Laser therapy device equipped with a
new probe an emitting diode that allows beam laser application directly to 09 ELECTROMEDICAL

DEVICES 2017 ENG gimaitaly en assets cataloghi |||

PHYSIOTHERAPY - NEBULIZERS - SUCTION PUMPS AIR COMPRESSION

THERAPY SOLUTIONS All human body is sup ... sectors without solenoids
PHYSIOTHERAPY - NEBULIZERS - SUCTION PUMPS MT BASE:
MAGNETOTHERAPY **28321** GIMA MT BASE PLUS - magnetotherapy Magneto Base
Plus is an ELF low frequency magnetotherapi...

lang:en score:11 filesize: 4.82 M page_count: 22 document date: 2017-06-14

[\[pdf\]](#) Catalog

untitled 09 ELETTRONICALI 2018 19 ITA gimaitaly it assets cataloghi |||

ELETTROMEDICALI - FISIOTERAPIA SISTEMI DI PRESSOMASSAGGIO Il corpo umano costituito da circolazio ... tale e 2 settori passivi senza solenoidi
ELETTROMEDICALI - FISIOTERAPIA MAGNETOTERAPIA BASE **28321** GIMA MT
BASE PLUS - magnetoterapia Il Magneto Base un apparato per magnetoterapia a bassa fr...

lang:it score:11 filesize: 8.4 M page count: 22 document date: 2017-06-30