



Manuals.plus /

- › Oral-B /
- › Oral-B Pulsar Vibrating Bristles Toothbrush User Manual

Oral-B Pulsar

Oral-B Pulsar Vibrating Bristles Toothbrush User Manual

Model: Pulsar

INTRODUCTION

The Oral-B Pulsar toothbrush is designed to provide an enhanced cleaning experience. It features soft vibrating bristles that assist in dislodging plaque between teeth. The pressure-sensitive split head adapts to the contours of your teeth, facilitating cleaning in hard-to-reach areas and regulating pressure on teeth and gums.

This manual provides instructions for the proper setup, operation, and maintenance of your Oral-B Pulsar toothbrush.



Image: Packaging for the Oral-B Pulsar 3D White Toothbrush, showing a 4-pack with battery-powered feature and stain removal claim.

SETUP

The Oral-B Pulsar toothbrush is battery-operated and comes with a pre-installed battery. No initial battery installation is required for first use.

1. Remove the toothbrush from its packaging.
2. Ensure the bristles are clean before first use.




Image: Four Oral-B Pulsar 3D White Toothbrushes, showing their individual design and clear handles.

OPERATING INSTRUCTIONS

To use your Oral-B Pulsar toothbrush:

1. Apply a small amount of toothpaste to the bristles.
2. Press the **On Button** located on the handle to activate the vibrating bristles.
3. Brush your teeth as you would with a standard manual toothbrush, ensuring all surfaces are cleaned. The vibrating bristles will assist in plaque removal.
4. After brushing, press the **Off Button** to deactivate the vibrations.
5. Rinse your mouth and the toothbrush thoroughly with water.

Oral-B®
PULSAR
3D WHITE™
TOOTHBRUSHES

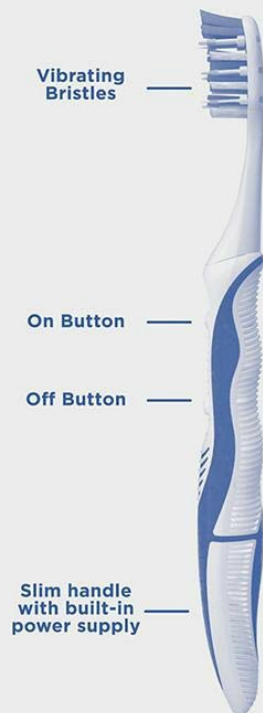

Vibrating
Bristles

On Button

Off Button

Slim handle
with built-in
power supply

WHY IS PULSAR 3D WHITE RIGHT FOR YOU?

It pulsates away plaque to help improve gum health and helps whiten your smile. Pulsar 3D White removes surface stains ON and BETWEEN teeth.

HOW DOES PULSAR 3D WHITE WORK?

It's a unique manual brush with an extra boost of vibrating power. It reaches surface stains even between teeth. Just turn it on and brush as you would with a regular manual brush.

AND IT'S GUARANTEED...

If you're not satisfied, call 1-877-593-4199 within 60 days of purchase with UPC and receipt to receive a pre-paid card in the amount of the purchase.

Recommended to use with:

Crest 3D White™ Toothpaste - Whitens by removing surface stains.

Bristles contain an antimicrobial to keep them cleaner. Does not kill bacteria in the mouth OR protect against disease.

Battery included. No replacement necessary.

Each Oral-B® Pulsar 3D White toothbrush comes with a battery that does not require any replacing. When the battery dies, simply change out your Pulsar brush for a new one. The battery is designed to last the lifetime of the brush head, approximately 3 months. Actual battery life may vary.

**NOT INTENDED FOR CHILDREN 3 YEARS AND UNDER.
REPLACE BRUSH EVERY 3 MONTHS.**

For more oral care information, visit www.OralB.com or call 1-800-56-ORALB (1-800-566-7252).

DISTR. BY ORAL-B, A DIVISION OF PROCTER & GAMBLE, CINCINNATI, OH 45202

Patents: www.pg.com/patents
© 2019 P&G

OM065

90749474



P&G
www.pg.com

Image: The back of the Oral-B Pulsar 3D White Toothbrush packaging, illustrating the vibrating bristles, On/Off buttons, and slim handle with built-in power supply.

MAINTENANCE

- **Cleaning:** After each use, rinse the toothbrush head thoroughly under running water to remove any toothpaste residue and debris.
- **Storage:** Store the toothbrush in an upright position to allow it to air dry.
- **Replacement:** For optimal oral hygiene and performance, replace your Oral-B Pulsar toothbrush every 3 months, or sooner if the bristles appear frayed or worn. The battery is designed to last the lifetime of the brush head.

SAFETY INFORMATION

- This product is not intended for children 3 years of age and under.
- Replace the brush every 3 months for effective cleaning and hygiene.

SPECIFICATIONS

Brand	Oral-B
Model	Pulsar
Power Source	Battery Powered (1 Lithium Ion battery required, pre-installed)
Special Feature	Vibrating Bristles, Pressure-Sensitive Head
Bristle Firmness	Medium
Age Range	Adult
Product Dimensions	9.3 x 2 x 0.4 inches
Item Weight	5.6 ounces
UPC	300416693668

WARRANTY AND SUPPORT

For additional oral care information, product support, or inquiries, please visit the official Oral-B website or contact customer service.

- **Website:** www.OralB.com
- **Customer Service:** 1-800-56-ORALB (1-800-566-7252)

Oral-B offers a satisfaction guarantee. If you are not satisfied, you may call 1-877-593-4199 within 60 days of purchase with your UPC and receipt to receive a pre-paid card in the amount of the purchase.