

Safety 1st TH068

Safety 1st Quick Read Forehead Thermometer User Manual

Model: TH068

[Introduction](#)

[Setup](#)

[Operation](#)

[Care & Maintenance
Information](#)

[Troubleshooting](#)

[Specifications](#)

[Safety](#)

1. INTRODUCTION

The Safety 1st Quick Read Forehead Thermometer provides fast and accurate temperature readings. Designed for ease of use, it allows for quick temperature checks without disturbing a sleeping child. The thermometer displays readings in large numerals and includes fever icons for immediate interpretation of results.



Image 1: The Safety 1st Quick Read Forehead Thermometer.

Key Features:

- Accurate temperature reading in 1 second.
- No-touch method for minimal disturbance, ideal for sleeping children.
- Large digital display with easy-to-understand fever icons.
- Durable construction using high-quality materials.

2. SETUP

Battery Installation:

1. Locate the battery compartment on the back of the thermometer.
2. Use a small screwdriver to open the battery compartment cover.
3. Insert 2 AAA batteries, ensuring correct polarity (+/-).

4. Replace the battery compartment cover and secure it with the screw.

Your browser does not support the video tag.

Video 1: Demonstration of battery installation and basic use of the Safety 1st Quick Read Forehead Thermometer.

3. OPERATING INSTRUCTIONS

Taking a Temperature Reading:

1. Press the power button to turn on the thermometer.
2. Position the thermometer's lens approximately 1/2 inch (1.27 cm) from the center of the forehead. Ensure the lens is aimed directly at the forehead, not downward toward the nose.
3. Press the measurement button. The thermometer will beep once the reading is complete (approximately 1 second).
4. Read the temperature displayed on the screen. The display will show the temperature in large numerals, along with a fever icon if a high temperature is detected.

**Accurate
readings**

in 1 second,
displayed
in °F or °C

**Des lectures
précises**

en 1 seconde
affichées en
°F ou en °C



Image 2: Proper positioning for a no-touch forehead temperature reading.

Understanding Fever Icons:

The thermometer features two easy-to-understand icons to indicate temperature status:

- **Smiling Face Icon:** Indicates a normal temperature.
- **Frowning Face Icon:** Indicates a high temperature (fever).

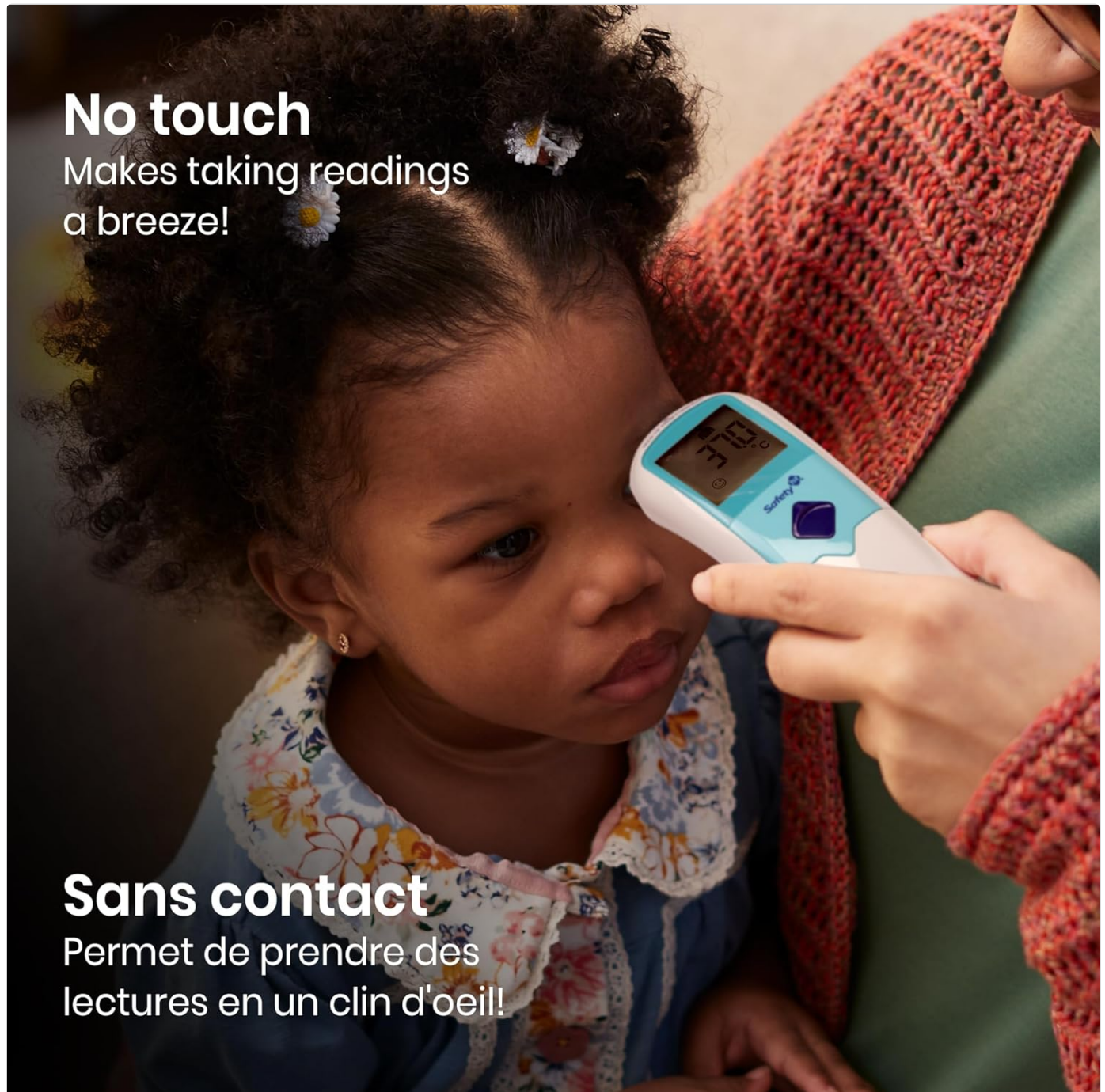


Image 3: Thermometer display showing temperature and fever icons.

Switching Between Fahrenheit and Celsius:

Refer to the detailed instructions provided in the packaging for specific steps on how to switch between Fahrenheit (°F) and Celsius (°C) modes, as this may involve accessing a small switch within the battery compartment.

Your browser does not support the video tag.

Video 2: Overview of the Safety 1st Digital Thermometer, highlighting its ease of use for quick temperature readings.

4. CARE AND MAINTENANCE

Cleaning:

- Wipe the thermometer clean with a soft, dry cloth after each use.
- Do not use abrasive cleaners or immerse the thermometer in water.

Storage:

- Store the thermometer in a dry place, away from direct sunlight and extreme temperatures.
- Remove batteries if the thermometer will not be used for an extended period.

5. TROUBLESHOOTING

Problem	Possible Cause	Solution
No display or weak display	Low or dead batteries, incorrect battery installation.	Replace batteries with new AAA batteries, ensuring correct polarity.
Inaccurate readings	Improper positioning, dirty sensor lens, extreme ambient temperature.	Ensure correct positioning (1/2 inch from forehead, not angled). Clean the sensor lens gently with a soft cloth. Allow thermometer to acclimate to room temperature for 30 minutes before use.
Thermometer does not turn on	Batteries are not installed or are dead.	Install new AAA batteries.

6. SPECIFICATIONS

Feature	Detail
Brand	Safety 1st
Model Name	Safety 1st Easy Read Dual Mode Thermometer (TH068)
Special Feature	High Accuracy
Display Type	Digital
Response Time	1 second
Power Source	Battery Powered (2 AAA batteries included)
Battery Life	Approximately 100 hours
Product Dimensions	7.2 x 2.8 x 1.5 inches
Item Weight	4.2 ounces
Outer Material	Plastic
Product Care	Wipe Clean

7. SAFETY INFORMATION

Legal Disclaimer: This product is not a toy. Do not allow children to play with it. When not in use, please keep out of reach of children. For home use only.

Always consult a healthcare professional for any health concerns or if you receive a high temperature reading.

