

[Manuals.plus](#) /

› [reVive Light Therapy](#) /

› reVive dpl IIa Blue & Red Light Therapy Device User Manual

reVive Light Therapy RVPNLSYS

reVive dpl IIa Blue & Red Light Therapy Device for Face

Model: RVPNLSYS

Brand: reVive Light Therapy

INTRODUCTION

The reVive dpl IIa Blue & Red Light Therapy Device is an OTC Class II medical device designed for comprehensive facial skin treatment. It utilizes advanced LED light therapy to address both anti-aging concerns, such as wrinkles and fine lines, and acne-related issues. This non-invasive device offers a clinical-strength treatment for home use, promoting healthier, clearer, and more youthful-looking skin.

The device features multiple treatment settings to cater to various skin needs, including specific modes for wrinkle reduction, acne treatment, and a combination of both. Its design focuses on ease of use and effectiveness, making it a convenient addition to your skincare routine.

WHAT'S INCLUDED

Upon unboxing your reVive dpl IIa device, you should find the following components:

- dpl IIa Device (Light Panel)
- Instruction Manual
- AC/DC Power Supply
- Goggles for Eye Comfort



Image: The reVive dpl IIa device, shown open to display the LED panels, alongside its retail packaging.

SETUP

Setting up your reVive dpl IIa device is straightforward:

1. **Prepare Skin:** Start with clean, washed, and dried skin.
2. **Connect Power:** Plug the AC/DC Power Supply into the device and then into a suitable electrical outlet.
3. **Wear Goggles:** For eye comfort and protection, always wear the provided goggles during treatment.

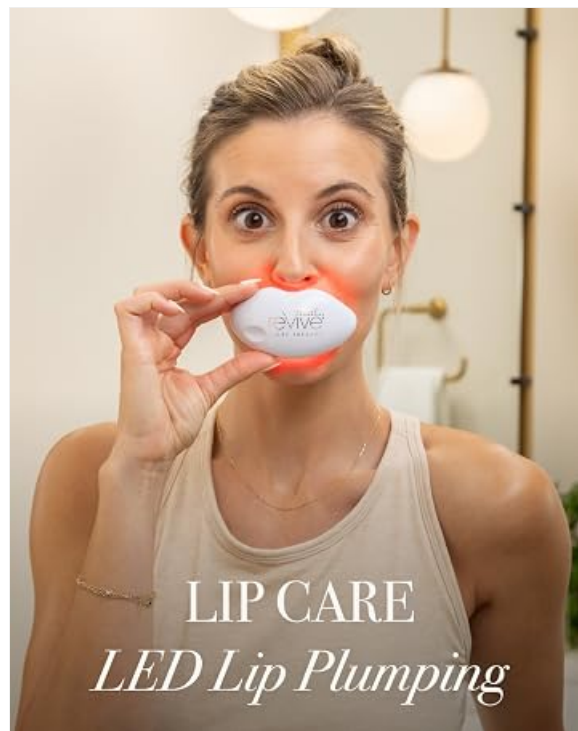


Image: A user wearing the protective goggles, demonstrating the preparation step before using the device.

OPERATING INSTRUCTIONS

The reVive dpl IIa device offers three distinct treatment settings:

Setting One (Anti-Aging): Red, Infrared, and Amber lights. This setting is designed to reduce wrinkles, smooth fine lines, and improve skin texture and tone.

Setting Two (Anti-Acne): Blue and Red lights. This setting targets acne by destroying bacteria under the skin and fighting inflammation.

Setting Three (Anti-Acne & Anti-Aging): Red, Infrared, Amber, and Blue lights. This comprehensive setting targets all skin concerns simultaneously.

To operate the device:

1. **Power On:** Press the On/Off button once to turn on the device.
2. **Select Setting:** Press the On/Off button again to cycle through the available settings (One, Two, Three).
3. **Apply Treatment:** Gently place the device on the desired skin area and hold it in place. Do not move or wave the device during treatment.
4. **Treatment Duration:** Each cycle is approximately 3 minutes. The device will automatically shut off when the treatment is complete.
5. **Repeat (if needed):** Move to the next treatment area and repeat the process.

For optimal results, consistent use is recommended. Anti-aging results can take up to 12 weeks to become noticeable, while acne treatment results may be seen in as little as one day. The device can be used up to 3 times per day.



3 SETTINGS FOR *a Variety of Skin Needs*

SETTING *One*

RED, INFRARED, and AMBER



Anti-Aging

reduce wrinkles, smooth
fine lines, improve skin
texture and tone

SETTING *Two*

BLUE and RED



Anti-Acne

destroys bacteria
under the skin and
fights inflammation

SETTING *Three*

RED, INFRARED, AMBER, and BLUE



Anti-Acne & Anti-Aging

target all skin
concerns at once

Image: The reVive dpl IIa device displaying its various light settings for anti-aging and acne treatment.



Image: A user applying the reVive dpl IIa device to their face, demonstrating proper usage during a light therapy session.

MAINTENANCE

To ensure the longevity and optimal performance of your reVive dpl IIa device, follow these maintenance guidelines:

- **Cleaning:** After each use, gently wipe the LED panels and the device exterior with a soft, dry cloth. Do not use abrasive cleaners or immerse the device in water.
- **Storage:** Store the device in a cool, dry place, away from direct sunlight and extreme temperatures. Use the provided soft carrying bag for protection when not in use or during travel.
- **Handling:** Handle the device with care to avoid dropping it or subjecting it to strong impacts, which could damage the LEDs or internal components.

TROUBLESHOOTING

If you encounter any issues with your reVive dpl IIa device, please refer to the following common troubleshooting tips:

- **Device Not Turning On:**
 - Ensure the AC/DC power supply is securely plugged into both the device and a working electrical outlet.
 - Check if the power outlet is functional by plugging in another appliance.
- **Lights Not Illuminating Properly:**
 - Verify that the correct treatment setting has been selected by pressing the On/Off button to cycle through modes.
 - Confirm the device is receiving power.
- **No Apparent Results:**
 - Ensure consistent daily use as recommended in the operating instructions. Results for anti-aging may take up to 12 weeks.
 - Confirm that the device is being held directly against the skin area being treated for optimal light penetration.

If the issue persists after attempting these steps, please contact customer support for further assistance.

SPECIFICATIONS

Feature	Detail
Product Dimensions	11 x 4.25 x 9.5 inches
Item Weight	2.98 Pounds
Model Number	RVPNLSYS
UPC	857563004958
Manufacturer	LED Technologies, Inc
Country of Origin	China
Device Type	OTC Class II Medical Device

WARRANTY & SUPPORT

Your reVive dpl Ila device comes with a**One Year Warranty** from the date of purchase. This warranty covers manufacturing defects and ensures the quality and performance of your device under normal use. For warranty claims, technical support, or any other inquiries, please contact reVive Light Therapy customer service. It is recommended to register your product upon purchase to facilitate any future support needs. You can find more information and contact details on the official reVive Light Therapy website or through the retailer where you purchased the device.

