

Manuals+

[Q & A](#) | [Deep Search](#) | [Upload](#)

manuals.plus /

› [Waterpik](#) /

› [Waterpik Cordless Freedom Water Flosser WF-03 Instruction Manual](#)

Waterpik WF-03

Waterpik Cordless Freedom Water Flosser WF-03 Instruction Manual

Model: WF-03 | Brand: Waterpik

INTRODUCTION

The Waterpik Cordless Freedom Water Flosser is designed for effective oral hygiene, offering a convenient, cordless solution for daily flossing. It is portable, battery-operated, and suitable for use in smaller bathrooms or while traveling. This device features dual-water pressure control and is accepted by the American Dental Association (ADA) for its proven effectiveness in removing plaque and improving gum health.

IMPORTANT SAFETY INFORMATION

Please read all instructions before using this appliance.

- KEEP OUT OF REACH OF CHILDREN.
- DO NOT USE IF YOU HAVE A HEART CONDITION, EPILEPSY, OR ANY OTHER MEDICAL CONDITION THAT MAY BE AFFECTED BY WATER PRESSURE.
- ALWAYS FOLLOW THE INSTRUCTIONS PROVIDED WITH THIS DEVICE.
- RISK OF ELECTRICAL SHOCK IF USED IN WATER OR OTHER LIQUIDS. DO NOT IMMERSE DEVICE IN WATER OR OTHER LIQUIDS.
- DO NOT USE DEVICE IF IT HAS BEEN DAMAGED OR IF YOU SUSPECT IT HAS BEEN DAMAGED.
- DO NOT USE DEVICE WITH DAMAGED OR WORN-OUT BATTERIES. REPLACE BATTERIES IMMEDIATELY IF THEY LEAK OR SWELL.
- DO NOT DISCHARGE BATTERIES IN FIRE. DISPOSE OF USED BATTERIES ACCORDING TO LOCAL LAWS AND REGULATIONS.

PACKAGE CONTENTS

Your Waterpik Cordless Freedom Water Flosser package includes:

- 1 Waterpik Cordless Freedom Water Flosser Unit
- 3 Water Flossing Tips (Classic Jet Tip, Orthodontic Tip, etc.)
- 1 Microfiber Travel Bag
- 1 Travel Water Plug
- 3 AA Batteries

3 Water Flossing Tips Included



The Cordless Freedom Water Flosser includes a travel bag and a travel water plug for convenience.



Three specialized water flossing tips are included for various oral care needs.

PRODUCT OVERVIEW



The Waterpik Cordless Freedom Water Flosser WF-03, a portable dental oral irrigator.

SETUP AND FIRST USE

1. **Battery Installation**

Remove the reservoir from the main unit. Twist the battery door counter-clockwise and pull to open. Insert three AA batteries, ensuring correct polarity. Replace the battery door by aligning arrows, pushing, and twisting clockwise.

Note: Replace batteries every 1-2 months with daily use, depending on battery quality, or when pressure significantly decreases.

2. **Filling the Reservoir**

The reservoir can be filled in two ways:

- Detach the reservoir from the unit and fill it directly with warm water.
- Alternatively, open the fill door on the reservoir while it's attached to the unit and fill with warm water.

Easy-to-Fill Reservoir



30 Seconds of
Water Capacity



The reservoir is designed for easy refilling and holds approximately 30 seconds of water capacity.

3. **Attaching a Tip**

Insert your chosen tip into the center of the handle until it clicks into place. To remove or change tips, press the eject button on top of the unit.

OPERATING INSTRUCTIONS

1. **Power On/Off and Pressure Settings**

The unit features a dual-water pressure control system with two settings: low and high. Slide the pressure switch to your desired setting.



Select between two pressure settings for personalized comfort.

2. **Flossing Technique**

Lean over a sink. Place the tip in your mouth, aiming it at the gumline at a 90-degree angle. Close your lips enough to prevent splashing but allow water to flow into the sink. Move the tip along the gumline, pausing briefly between teeth. Clean both the front and back surfaces of your upper and lower teeth.

Tip: For optimal results, hold the unit upright and avoid tipping it during use.



Water flossing effectively removes up to 99.9% of plaque from treated areas.

3. **Waterproof Design**

The Waterpik Cordless Freedom is waterproof and can be used in the shower for added convenience.



The waterproof design allows for use in the shower.

MAINTENANCE

- **Cleaning the Reservoir:** After each use, empty any remaining water from the reservoir. It is recommended to leave the fill door open to allow the reservoir to air dry.
- **Cleaning the Unit:** Wipe the exterior of the unit with a soft, damp cloth as needed. Do not immerse the main unit in water.
- **Tip Care:** Rinse tips after each use. Replace tips every 3-6 months, or as recommended by your dental professional.

TROUBLESHOOTING

- **Low Pressure:** Check battery levels. If batteries are old or weak, replace them with new AA batteries. Ensure the tip is securely inserted.
- **Unit Not Turning On:** Verify that batteries are correctly installed and have sufficient charge.
- **Water Leakage:** Ensure the reservoir is properly seated and the travel water plug is securely in place if transporting the unit.

SPECIFICATIONS

Feature	Detail
Power Source	Battery Powered (3 AA batteries included)
Pressure Settings	2 (45 to 75 PSI)
Pulses Per Minute	1200
Water Reservoir Capacity	5 ounces (approx. 30 seconds flossing time)
Special Features	3 Water Flossing Tips, Travel Bag included, Waterproof
Model Number	WF-03
UPC	073950199773
Product Dimensions	12 x 12 x 12 inches; 1 Pound
Manufacturer	Water Pik, Inc.

WARRANTY AND SUPPORT

The Waterpik Cordless Freedom Water Flosser WF-03 comes with a one-year warranty. For further assistance, please refer to the official Waterpik website or contact customer support.

OFFICIAL PRODUCT VIDEOS

How to Use the Waterpik Cordless Freedom Water Flosser

Your browser does not support the video tag.

This video provides a step-by-step guide on how to effectively use the Waterpik Cordless Freedom Water Flosser for optimal oral hygiene.

Waterpik Water Flosser Cordless Dental Oral Irrigator for Teeth with Portable Travel Bag and 3 Jet Tips

Your browser does not support the video tag.

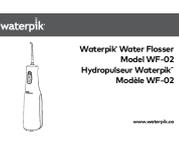
This video highlights the features of the Waterpik Cordless Water Flosser, including its portability and the included travel bag and jet tips.

Related Documents - WF-03



[Waterpik Cordless Express WF-02 Quick Start Guide | Oral Care](#)

Get started quickly with your Waterpik Cordless Express WF-02 water flosser. This guide provides simple instructions for setup, use, and optimal oral hygiene.

 <p>waterpik Waterpik™ Water Flosser WF-20 21 Series Hydropulseur Waterpik™ Modèles WF-20 21 www.waterpik.com</p>	<p>Waterpik WF-20 21 Series Water Flosser User Manual</p> <p>Comprehensive user manual for the Waterpik WF-20 and WF-21 Series Water Flossers, covering setup, operation, cleaning, maintenance, battery care, warranty, and safety information.</p>
 <p>waterpik Waterpik Water Flosser WF-20 21 Series Limpideur dental Waterpik™ WF-20 21 Serie www.waterpik.com</p>	<p>Waterpik Water Flosser WF-20 21 Series User Manual</p> <p>Comprehensive user manual for the Waterpik Water Flosser WF-20 and 21 Series, covering setup, usage, cleaning, troubleshooting, battery care, and warranty information.</p>
 <p>waterpik Waterpik™ Water Flosser Model WF-02 Hydropulseur Waterpik™ Modèle WF-02 www.waterpik.com</p>	<p>Waterpik Water Flosser WF-02 User Manual and Guide</p> <p>Official user manual for the Waterpik Water Flosser Model WF-02. Includes setup, operation, cleaning, troubleshooting, and warranty information for optimal oral hygiene. Learn how to use and maintain your Waterpik water flosser for effective dental care.</p>
 <p>waterpik Waterpik™ Water Flosser Model WF-02 Hydropulseur Waterpik™ Modèle WF-02 www.waterpik.com</p>	<p>Waterpik Water Flosser Model WF-02 User Manual and Guide</p> <p>Comprehensive user manual for the Waterpik Water Flosser Model WF-02, covering setup, usage, cleaning, troubleshooting, and warranty information. Learn how to effectively use your Waterpik for optimal oral hygiene.</p>
 <p>waterpik Waterpik™ Water Flosser Model WF-03 Hydropulseur Waterpik™ Modèle WF-03 www.waterpik.com</p>	<p>Waterpik Water Flosser Model WF-03 Instruction Manual</p> <p>Comprehensive guide for the Waterpik Water Flosser Model WF-03, covering setup, usage, cleaning, troubleshooting, and warranty information.</p>