

Clean + Easy 45026

Clean + Easy Electrolysis Machine Model 45026 User Manual

Brand: Clean + Easy | Model: 45026

1. IMPORTANT SAFETY INFORMATION

WARNING: Before using this device, carefully read and understand all safety instructions. Failure to follow these instructions may result in injury or adverse reactions. Consult a healthcare professional if you have any pre-existing medical conditions or concerns.

Do not use this device if you have any of the following conditions:

- Diabetes, Hepatitis, Blood disorders, Pigmentation problems, Lupus, Scleroderma (Pigmentation).
- Impetigo, Chronic skin disorders, or if you scar easily.
- If you wear a pacemaker or similar device.
- If you are pregnant or taking medication such as Retin-A, Accutane, Glycolic Acid, Alpha Hydroxy treatments, or Steroids that may cause skin sensitivity.
- Do not use on your upper lip if you currently wear braces.
- Do not use to remove moles, warts, and beauty spots, or on nipples, eyelids, eyelashes, tattoos, and sunburned skin or inside nose or ears.

The unit should not be used by anyone under the age of 18 without parental supervision.

2. PRODUCT OVERVIEW

The Clean + Easy Electrolysis Machine Model 45026 is designed for at-home permanent hair removal. This non-laser device utilizes radio frequency transmitted via a stylet tip to weaken and eventually stop hair regrowth by targeting the hair root and shaft.



Image: Clean + Easy Electrolysis Machine packaging and device components.

How Electrolysis Works

Electrolysis is the only non-laser method for permanent hair removal. It works by introducing a mild electrical current into the moist base of the hair follicle. This current stops the blood supply to the root, which is then neutralized by a chemical action, preventing regrowth. This method is effective on a few hairs at a time, making it suitable for various body areas and all hair/skin types.

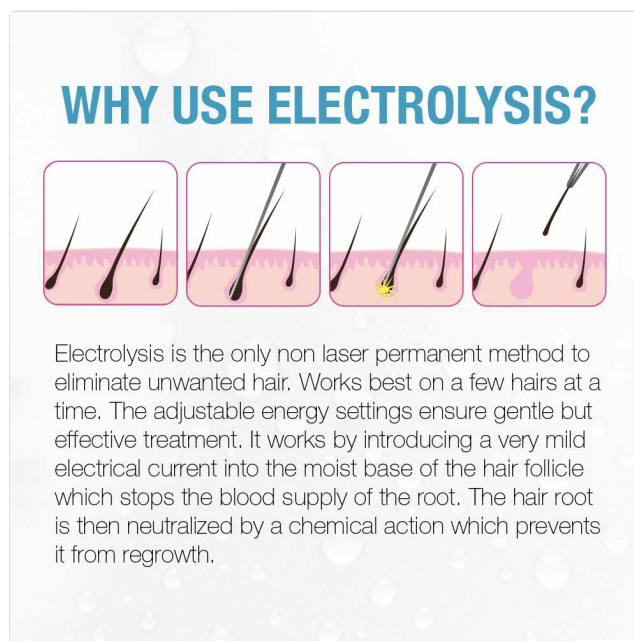


Image: Diagram illustrating the process of electrolysis on hair follicles.

Key Features

- Non-laser permanent hair removal for face and body.

- Weakens hair root and shaft with each use.
- Suitable for coarse and stubborn hairs.
- Economical and convenient for at-home use.
- Retractable stylet tip designed not to pierce the skin.

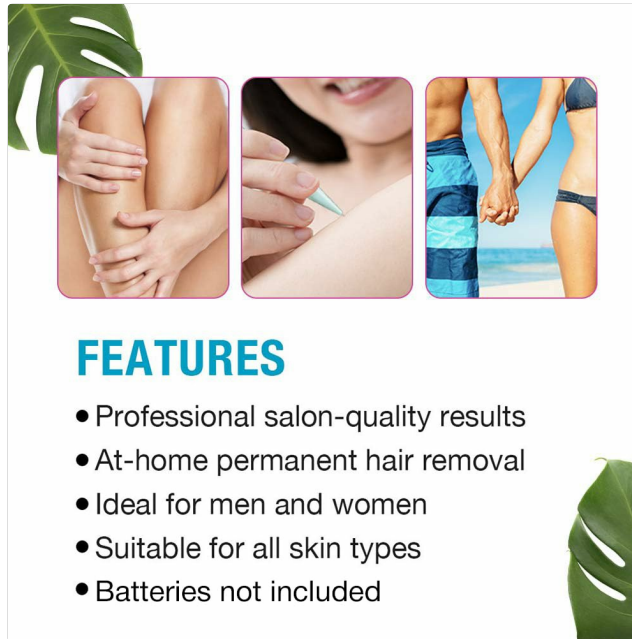


Image: Visual representation of key product features.

3. PACKAGE CONTENTS

Your Clean + Easy Electrolysis Machine package includes:

- Clean + Easy Personal Electrolysis Unit
- Spare Stylet
- Tweezer
- Magnifying Mirror

Note: A 9-volt battery is required for operation and is NOT included in the package. Please purchase separately.

4. SETUP

1. **Insert Battery:** Open the battery compartment and insert a new 9-volt battery. Ensure correct polarity.
2. **Cleanse Skin:** Cleanse the skin area to be treated with isopropyl alcohol. Allow the skin to dry completely.
3. **Prepare Salt Solution:** Mix $\frac{1}{2}$ teaspoon of salt with 1 cup of water to create a salt solution.
4. **Moisten Fingers:** Moisten your fingers with the prepared salt solution. Your fingers must remain moist throughout the treatment to ensure proper conductivity.
5. **Clean Stylet and Tweezers:** Wipe the stylet tip and tweezers with isopropyl alcohol before each session.



Image: Visual guide for preparation steps before using the electrolysis device.

5. OPERATING INSTRUCTIONS

Follow these steps carefully for effective hair removal:

1. **Set Comfort Dial:** Turn the Comfort Control Dial to position one (1). You can increase this setting later if needed.
2. **Expose Stylet Tip:** Push the slide button forward to expose the stylet tip.
3. **Hold Device:** Hold the stylet like a pen with your moistened fingers around the touch band.
4. **Insert Stylet:** Observe the angle of hair growth. Gently insert the stylet tip along the underside of the hair and into the follicle. The tip should slide in easily without resistance.
5. **Confirm Contact:** Hold the stylet tip firmly in place. When proper contact is made, you should feel a slight tingling sensation. If no tingling is felt, increase the Comfort Control Dial setting gradually.
6. **Rotate Gently:** Gently rotate the stylet tip in the follicle a few times to ensure contact with the hair root.
7. **Apply Current:** Once a slight tingling is felt, hold the stylet firmly in place and count to 15 slowly.
8. **Remove Hair:** Remove the stylet and use the provided tweezers to gently remove the treated hair. The hair should slide out easily without pulling or resistance. If resistance is felt, the root may not be sufficiently treated; re-treat the hair.
9. **After Treatment:** After completing the session, cleanse the treated skin area with isopropyl alcohol.



Image: The electrolysis device in use, showing stylet insertion.

Important Tips for Use

- A magnifying mirror and good lighting are essential for precise insertion.
- Do not share stylet tips with others to maintain hygiene.
- Ensure your fingers remain moist throughout the entire treatment session.
- Avoid applying make-up, perfume, or exposing treated areas to direct sunlight for 24 hours after treatment.
- Wait 24 hours before repeating treatment on the same area.

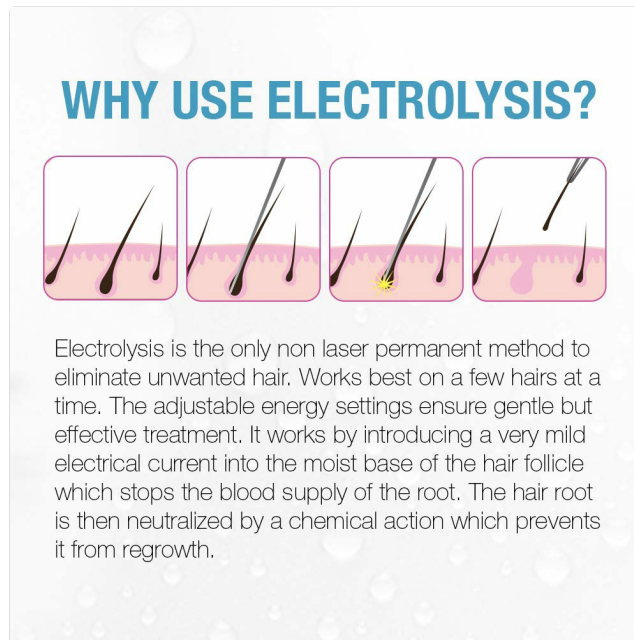


Image: Key tips for optimal use and post-treatment care.

6. MAINTENANCE

- **Cleaning:** After each use, clean the stylet tip and tweezers with isopropyl alcohol.
- **Storage:** Store the device in a cool, dry place, away from direct sunlight and moisture.
- **Battery:** Remove the 9-volt battery if the device will not be used for an extended period to prevent leakage.
- **Stylet Care:** The stylet tip is delicate. Handle with care to prevent bending or breaking.

7. TROUBLESHOOTING

Problem	Possible Cause	Solution
Device not activating / No tingling sensation	<ul style="list-style-type: none">• Battery is dead or incorrectly inserted.• Fingers are not sufficiently moist.• Stylet tip not making proper contact with the follicle or touch band.• Comfort Control Dial set too low.	<ul style="list-style-type: none">• Replace 9-volt battery; check polarity.• Re-moisten fingers with salt solution.• Ensure stylet is correctly inserted into the follicle and fingers are on the touch band.• Gradually increase Comfort Control Dial setting.

Problem	Possible Cause	Solution
Stinging or burning sensation / Skin turning dark	<ul style="list-style-type: none"> • Stylet not inserted correctly into the follicle. • Comfort Control Dial setting is too high. • Holding the stylet in place for too long. 	<ul style="list-style-type: none"> • Remove stylet and reinsert in a different position until a tingling sensation is felt. • Reduce Comfort Control Dial setting. • Limit application time to 15 seconds per hair.
Hair does not remove easily	<ul style="list-style-type: none"> • Hair root not sufficiently treated. • Stylet not making proper contact with the root. 	<ul style="list-style-type: none"> • Re-treat the hair, ensuring proper insertion and sufficient application time. • Gently rotate the stylet to ensure contact.
Brown stains on fingers	<ul style="list-style-type: none"> • Minerals/impurities in tap water reacting with the device. 	<ul style="list-style-type: none"> • Use distilled water for preparing the salt solution.

8. SPECIFICATIONS

- **Model Number:** 45026
- **UPC:** 672153450263
- **Product Dimensions:** 8.82 x 3.31 x 1.97 inches
- **Weight:** 8 ounces
- **Power Source:** 9-volt battery (not included)
- **Manufacturer:** American International Industries

9. WARRANTY AND SUPPORT

Specific warranty information for the Clean + Easy Electrolysis Machine Model 45026 is not provided in this manual. For details regarding warranty coverage, please refer to the product packaging or contact the manufacturer directly.

For technical support or further assistance, please contact Clean + Easy customer service through their official website or the retailer where the product was purchased.