

manuals.plus /

› [Waterpik](#) /

› [Waterpik Ultra Dental Water Jet WP-100W Instruction Manual](#)

Waterpik WP-100W

Waterpik Ultra Dental Water Jet WP-100W Instruction Manual

PRODUCT OVERVIEW

The Waterpik Ultra Dental Water Jet WP-100W is designed to enhance oral hygiene by providing a unique combination of water pressure and pulsations. This device effectively cleans deep between teeth and below the gumline, targeting harmful bacteria and debris that traditional brushing and flossing may miss. It is clinically proven to remove up to 99.9% of plaque from treated areas and is up to 50% more effective for improving gum health compared to string floss. The WP-100W features a contemporary, compact design, a high-volume covered reservoir, and an advanced pressure-control system with 10 settings.



Image 1: Waterpik Ultra Dental Water Jet WP-100W product packaging, showing the device and its key features.

Key Features:

- **Powerful Deep Cleaning:** Effectively removes plaque and debris from between teeth and below the gumline.
- **Clinically Proven:** Up to 50% more effective than string floss for improving gum health.
- **10 Pressure Settings:** Adjustable water pressure for a customized, comfortable clean.
- **Comprehensive Tip Set:** Includes multiple specialized tips for different oral care needs.
- **Compact Design:** 20% smaller than previous models for convenient placement.
- **Quiet Operation:** Engineered for reduced noise during use.

- **90 Seconds Water Capacity:** High-volume reservoir for uninterrupted flossing.

IMPORTANT SAFETY INFORMATION

Please read all instructions before using this product. When using electrical products, especially when children are present, basic safety precautions should always be followed.

- For oral use only.
- Use only as directed in this instruction manual.
- Do not use if the unit, cord, or plug is damaged.
- Keep out of reach of children unless supervised by an adult.
- Do not direct water under the tongue, into the ear, nose, or other delicate areas. This product is capable of producing pressures that may cause serious damage in these areas.
- Do not fill the reservoir with hot water.
- Unplug the product before cleaning.

SETUP GUIDE

1. **Unpack the Unit:** Carefully remove the Waterpik Ultra Dental Water Jet from its packaging.
2. **Clean Removable Parts:** Before first use, clean all removable parts, including the reservoir and tips, with warm, soapy water. Rinse thoroughly.
3. **Fill the Reservoir:** Lift the hinged lid of the reservoir and fill it with warm water. You may also use mouthwash as desired. Ensure the reservoir is securely placed back on the base.
4. **Insert a Tip:** Select the desired tip and insert it into the handle until it clicks into place. The tip rotates 360 degrees for easy access to all areas of the mouth.
5. **Prime the Unit (First Use Only):** If this is the first time using the unit, or if it has not been used for an extended period, prime it. Set the pressure control dial to high (setting 10). Direct the tip into the sink. Turn the unit ON until water flows from the tip. Turn the unit OFF.
6. **Plug In:** Plug the power cord into a grounded electrical outlet.

OPERATING INSTRUCTIONS

1. **Adjust Pressure Setting:** Select your desired pressure setting using the dial on the base of the unit. Start with a low setting (1-3) and gradually increase as you become accustomed to the water jet.
2. **Position the Tip:** Place the tip in your mouth, aiming it at the gumline at a 90-degree angle. Lean over the sink to avoid splashing.
3. **Start Flossing:** Turn the unit ON using the switch on the base. Close your lips enough to prevent splashing, but allow water to flow from your mouth into the sink.
4. **Guide the Water Jet:** Begin with the back teeth and work your way forward. Glide the tip along the gumline and pause briefly between each tooth. Continue until you have cleaned all areas of your mouth.
5. **Pause Function:** The ergonomic handle includes a convenient pause button. Press and hold this button to temporarily stop the water flow without turning off the unit. Release to resume.
6. **Finish:** When finished, turn the unit OFF. Remove the tip by pressing the tip eject button on the handle.

The Waterpik Ultra Dental Water Jet provides approximately 90 seconds of water capacity, ensuring a complete cleaning cycle without the need for refills.

MAINTENANCE

Daily Cleaning:

- After each use, turn off the unit and unplug it.
- Empty any remaining water from the reservoir.
- Remove the tip and rinse it under running water. Store the tip in the unique undercover tip compartment for hygienic storage.
- Wipe the exterior of the unit with a soft, damp cloth.

Deep Cleaning (Monthly or as needed):

- To prevent mineral deposits, add 2 tablespoons of white vinegar to a full reservoir of warm water.
- Run the unit until the reservoir is empty, directing the water into the sink.
- Rinse by filling the reservoir with clean warm water and running the unit until empty.

TROUBLESHOOTING

If you encounter any issues with your Waterpik Ultra Dental Water Jet WP-100W, please refer to the following common solutions:

- **No Water Flow:** Ensure the reservoir is filled and properly seated. Check that the tip is securely inserted. If the unit has not been used recently, it may need priming (refer to Setup Guide).
- **Low Pressure:** Verify the pressure control dial is set to your desired level. Check for any blockages in the tip or hose.
- **Leaking:** Ensure all connections are tight, especially the reservoir to the base. Check for cracks in the hose or reservoir.
- **Unit Not Turning On:** Confirm the power cord is securely plugged into a working electrical outlet.

For further assistance or issues not listed here, please visit the official Waterpik website or contact their customer support.

PRODUCT SPECIFICATIONS

Brand	Waterpik
Model Number	WP-100W
Product Dimensions	12 x 12 x 8 inches
Item Weight	3.8 Pounds
Power Source	Corded Electric
Special Feature	Rotatable Tip, 10 Pressure Settings
UPC	073950278195
Manufacturer	Water Pik, Inc.






WARRANTY AND SUPPORT


The Waterpik Ultra Dental Water Jet WP-100W comes with a **Three Year Warranty** from the date of purchase, covering defects in materials or workmanship. Please retain your proof of purchase for warranty claims.

For warranty service, technical support, or to purchase replacement parts and accessories, please visit the official Waterpik website or contact their customer service department. Contact information can typically be found on the product packaging or the Waterpik website.

Website: www.waterpik.com

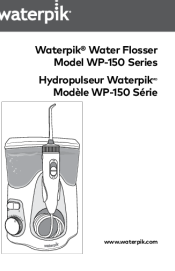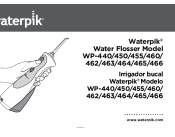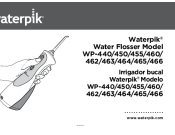Related Documents - WP-100W

 <p>waterpik</p> <p>Waterpik® Water Flosser Model WP-150 Series Hydropulseur Waterpik® Modèle WP-150 Série</p> <p>www.waterpik.com</p>	<p>Waterpik® Water Flosser WP-150 Series User Manual & Guide</p> <p>Official user manual and guide for the Waterpik® Water Flosser Model WP-150 Series. Learn about important safeguards, product description, getting started, tip usage, cleaning, troubleshooting, and warranty information.</p>
 <p>waterpik</p> <p>Waterpik® Water Flosser Model WP-440/450/455/460/ 462/463/464/465/466 Hydropulseur Waterpik® Modèle WP-440/450/455/460/ 462/463/464/465/466</p> <p>www.waterpik.com</p>	<p>Waterpik Water Flosser User Manual: Models WP-440 to WP-466</p> <p>Official user manual for Waterpik Water Flosser models WP-440, WP-450, WP-455, WP-460, WP-462, WP-463, WP-464, WP-465, and WP-466. Covers setup, operation, maintenance, troubleshooting, and warranty information.</p>
 <p>waterpik</p> <p>Waterpik® Water Flosser Model WP-440/450/455/460/ 462/463/464/465/466 Hydropulseur Waterpik® Modèle WP-440/450/455/460/ 462/463/464/465/466</p> <p>www.waterpik.com</p>	<p>Waterpik Water Flosser WP-440 Series User Manual and Guide</p> <p>This user manual provides detailed instructions and safety information for the Waterpik Water Flosser models WP-440, WP-450, WP-455, WP-460, WP-462, WP-463, WP-464, WP-465, and WP-466. It covers product description, getting started, tip usage, troubleshooting, cleaning, and warranty details.</p>
 <p>waterpik</p> <p>Waterpik® Water Flosser Model WP-660/670/ 672/674/676</p> <p>www.waterpik.co.uk</p>	<p>Waterpik Water Flosser WP-660/670/672/674/676 User Manual</p> <p>This manual provides important safeguards, product description, getting started guide, and cleaning/troubleshooting information for the Waterpik Water Flosser models WP-660, WP-670, WP-672, WP-674, and WP-676.</p>
 <p>waterpik waterflosser</p> <p>QUICK START GUIDE / GUIDE DE DÉMARRAGE RAPIDE</p> <p>1. Check the water level in the reservoir. Fill with clean water. Do not use tap water directly from the tap. 2. Turn the power switch on. 3. Select the desired flossing mode. 4. Insert the flossing tip into the mouth and floss. 5. Turn the power switch off. 6. Clean the flossing tip after use. 7. Store the flosser in a clean, dry place.</p>	<p>Waterpik Water Flosser Quick Start Guide - Setup and Usage</p> <p>A concise guide to setting up and using your Waterpik Water Flosser, including priming, tip selection, pressure adjustment, and usage techniques for optimal results. Features instructions for timer and massage modes.</p>

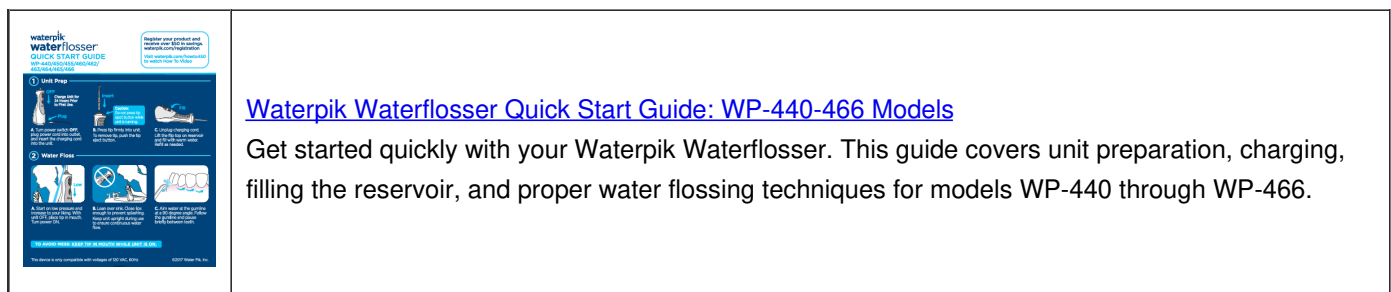

Waterpik Waterflosser Quick Start Guide: WP-440-466 Models

Get started quickly with your Waterpik Waterflosser. This guide covers unit preparation, charging, filling the reservoir, and proper water flossing techniques for models WP-440 through WP-466.

Documents - Waterpik – WP-100W


[pdf] Quick Start Guide Guide

20029786 FAA WP 100W ThankYou QSG F RGB WEB Waterpik Ultra Water Flosser Classic Blue 113
 Best Buy englishopens a new window Quick Start Guide files bbystatic VVm5Ox6r10HeFPo6Qi7UTQ |||
 Register and Save Save over 50 Visit Waterpik.com/WelcomeWF to watch How To
 video. Register your product and receive over 50 in savings.
 Waterpik.com/registration O er includes 50 o your next order of assorted tips tip
 storage case. We d Love to Hear from You Connect with Us Waterpik Ca...
 lang:en score:32 filesize: 818.94 K page_count: 2 document date: 2022-10-24




[pdf]


waterpik seniors and boomers CE professional s education |||
 Seniors and Boomers: Living Longer, Living Healthier. Considerations for Dental
 Professionals Course ... SeekerTM Tip Figure 12 , and string floss.47,49 8 Figure 5:
 Waterpik Ultra Water Flosser, Model **WP-100W**; comes with the Classic Jet Tip,
 Plaque SeekerTM Tip, Pik PocketTM Tip, Orthodontic Tip, Ton...
 lang:en score:18 filesize: 1.38 M page_count: 12 document date: 2010-08-05

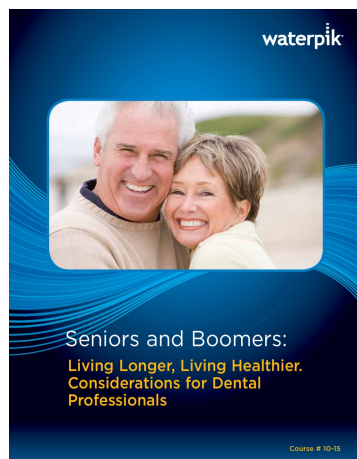

[pdf]

to Optimal Interdental Health waterpik Disclosure Statement The content for this self study course was
 developed and written by Water Pik Inc employee Carol A Jahn RDH MS Waterpik Coaching
 ContinuingEducation oral health pro education Your Patients Dental professionals know that regular
 interdental cleaning is a prerequisite most people in order achieve optimal Achieving patient compliance
 founded on the belief provider expert practitioner educates patients believes if they why how will do as are
 Coaching Your Patients to Optimal Interdental Health Course #12-19 Disclosure
 Statement: The cont ... egardless of type, to toothbrushing was nearly twice Figure 1:
 Waterpik Ultra Water Flosser, Model **WP-100W** 6 Figure 2: Waterpik Complete Care,
 Model WP-900W Other products including toothpicks an...
 lang:en score:18 filesize: 6.06 M page_count: 12 document date: 2017-02-20

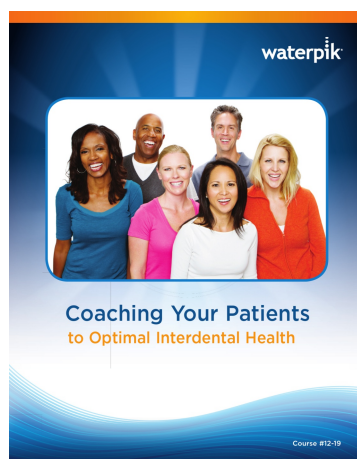

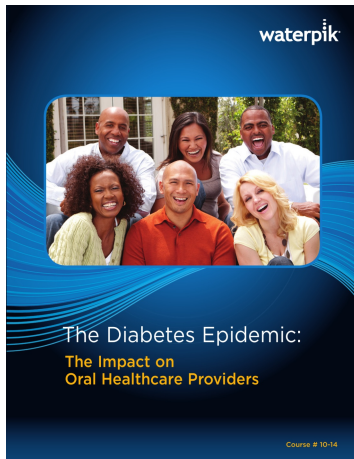

[pdf]

The Diabetes Epidemic Impact on Oral Healthcare Providers Diabetics Water Flosser Waterpik Dentist Hygienist Periodontal Disease Gum Gingivitis Plaque Bleeding Carol Jahn waterpik diabetes epidemic CE oral health pro education

The Diabetes Epidemic: The Impact on Oral Healthcare Providers Course # 10-14 Disclosure Statement ... e to avoid periods when the risk of hypoglycemia 7 Figure 1: Waterpik Ultra Water Flosser, Model **WP-100W**; comes with the Classic Jet Tip, Plaque Seeker™ Tip, Pik Pocket™ Tip, Orthodontic Tip, and...

lang:en score:18 filesize: 2.35 M page_count: 12 document date: 2014-08-14

[pdf] Quick Start Guide Guide

20027611 FAA WP 100W Series QSG F LR web QUICK START GUIDE 1 First Use ON OFF Water Floss original content syndigo asset 7acbaa61 51e6 41a7 b028 68f39e4bc203 FN 20021812 AB EI dispositivo es compatible con voltajes de 120 VAC 60 Hz ©2019 Pik Inc Registre su producto y reciba

AWSAccessKeyId AKIAIIE5CHZ4PRWSLYKQ Expires 1893513942 Signature

2jiVSdrHvH5QGnGQ8Riv8Ddpbjw media3 webcollage net

a6bce136f39538350dc31f839f74cd53af329325 response type application 1893523857

s4VFTZrpj9IzDUGnx 9W43ekNUo 349a12d9 bac5 4960 b2ae e52ed961e7f2 action resource smedia

rvvfp wc cp 1569275430203 33f9cdf9 f9cd 48f2 9919 77ea0dd881f4 module waterpik products

1457037257557 tab 0ce9e8a5 f23c 4491 8d79 2104d3fecc68 d56b91f8 7459 4b4c a4ac 914bd201be2b

1570765112535 20a3e824 966b 4e40 935b e7f20de002b3 1463755513028 34bcb350 182c 4fde 96da

d0f0f8632462

QUICK START GUIDE WP-100 Series Visit waterpik.com/welcomeWF Register your product and receive over 50 in savings. waterpik.com/registration 9 10 9 10 9 10 9 10 1 First Use Fill ON OFF O I Fill reservoir with warm 12345678 water and replace it rmly on the base. Insert ON OFF Insert t...

lang:es score:10 filesize: 272.6 K page_count: 2 document date: 2019-07-03

