

Manuals+

[Q & A](#) | [Deep Search](#) | [Upload](#)

manuals.plus /

› [Biosan](#) /

› [Non-Original Instruction Manual for Biosan Beauty Center 1000 Programs Mp6 Plus](#)

Biosan Beauty Center 1000 Programs Mp6 Plus

Non-Original Instruction Manual for Biosan Beauty Center 1000 Programs Mp6 Plus

Comprehensive guide for operation, maintenance, and troubleshooting.

1. INTRODUCTION

This manual provides detailed instructions for the safe and effective use of the Biosan Beauty Center 1000 Programs Mp6 Plus device. It covers essential information regarding setup, operation, maintenance, and troubleshooting. Please read this manual thoroughly before using the device to ensure proper function and to prevent potential hazards.

Important Safety Information: Always consult a healthcare professional before starting any new electrostimulation program, especially if you have pre-existing medical conditions, are pregnant, or have a pacemaker. Do not use the device on broken skin, open wounds, or sensitive areas. Keep out of reach of children.

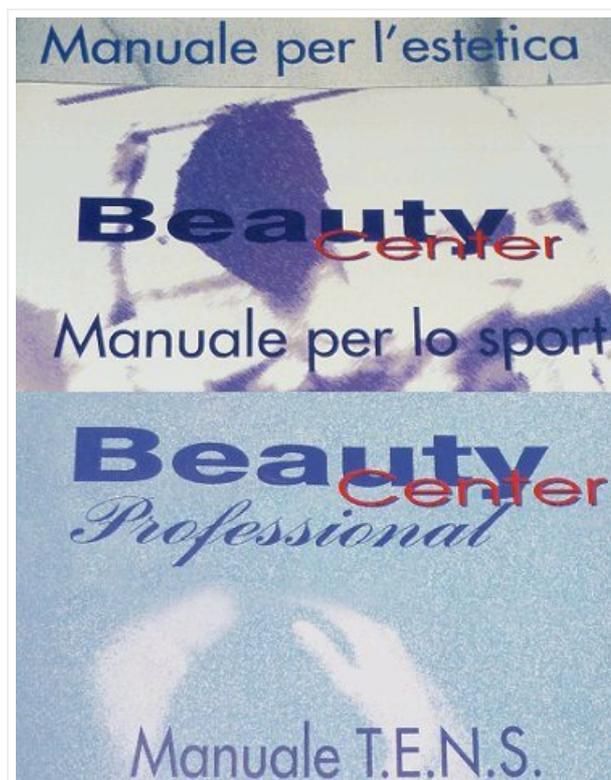


Image 1.1: Representation of the Biosan Beauty Center 1000 Programs Mp6 Plus device and its accompanying manual.
This image illustrates the general appearance of the product for identification purposes.

2. SETUP

2.1 Unpacking and Inspection

- Carefully remove all components from the packaging.
- Inspect the device and accessories for any signs of damage. If any damage is found, do not use the device and contact your supplier.
- Retain the packaging for future storage or transport.

2.2 Component Identification

Familiarize yourself with the main components:

- Main unit (Biosan Beauty Center 1000 Programs Mp6 Plus)
- Electrode cables
- Self-adhesive electrodes
- Power adapter (if applicable)
- User manual (this document)

2.3 Initial Preparation

1. Ensure the device is fully charged or connected to a power source using the provided adapter.
2. Connect the electrode cables to the output ports on the main unit.
3. Attach the self-adhesive electrodes to the ends of the electrode cables.

3. OPERATING INSTRUCTIONS

3.1 Powering On/Off

- To power on: Press and hold the **POWER** button until the display illuminates.
- To power off: Press and hold the **POWER** button until the display turns off.

3.2 Program Selection

The Biosan Beauty Center 1000 Programs Mp6 Plus offers various programs categorized for Sport, Aesthetics, and TENS applications. Use the **MODE** or **PROGRAM** buttons to navigate through the available options.

- **Sport Programs:** Designed for muscle stimulation, recovery, and performance enhancement.
- **Aesthetics Programs:** Focus on toning, firming, and cosmetic applications.
- **TENS Programs:** Transcutaneous Electrical Nerve Stimulation for pain relief.

Refer to the device's display for specific program names and descriptions. Select the desired program by pressing the **SELECT** or **ENTER** button.

3.3 Electrode Placement

Proper electrode placement is crucial for effective and safe treatment. Clean and dry the skin area before applying electrodes. Place electrodes on the target muscle groups or pain areas as indicated by the specific program guidelines or your healthcare professional's advice. Avoid placing electrodes over joints, bony prominences, or directly on the heart area.

3.4 Adjusting Intensity and Duration

- Use the + and - buttons to adjust the intensity of the stimulation. Start with a low intensity and gradually increase until a comfortable, yet effective, sensation is achieved.
- The program duration is typically pre-set. Some programs may allow for manual adjustment of time using dedicated **TIME** buttons.

Do not exceed comfortable intensity levels. If you experience any discomfort or pain, reduce the intensity or stop the session immediately.

4. MAINTENANCE

4.1 Cleaning the Device

- Ensure the device is powered off and disconnected from any power source before cleaning.
- Wipe the main unit with a soft, slightly damp cloth. Do not use abrasive cleaners, solvents, or immerse the device in water.
- Allow the device to air dry completely before storage or next use.

4.2 Electrode Care and Replacement

- After each use, re-apply the electrodes to their protective film to maintain their adhesive quality.
- If the electrodes lose their stickiness or become damaged, replace them with new, compatible electrodes. Using worn-out electrodes can reduce effectiveness and cause skin irritation.
- Store electrodes in a cool, dry place.

4.3 Storage

Store the Biosan Beauty Center 1000 Programs Mp6 Plus and its accessories in a clean, dry environment, away from direct sunlight and extreme temperatures. Keep it in its original packaging or a protective case when not in use.

5. TROUBLESHOOTING

If you encounter issues with your Biosan Beauty Center 1000 Programs Mp6 Plus, refer to the following common problems and solutions:

Problem	Possible Cause	Solution
Device does not power on.	Low battery or power adapter not connected.	Charge the device or ensure the power adapter is securely connected.
No sensation or weak stimulation.	Electrodes not properly applied, worn out electrodes, or low intensity setting.	Ensure electrodes are firmly attached to clean skin. Replace old electrodes. Increase intensity gradually. Check cable connections.
Skin irritation under electrodes.	Skin not clean, electrodes worn out, or allergic reaction.	Clean skin thoroughly before use. Replace electrodes. Discontinue use if irritation persists and consult a doctor.

If the problem persists after attempting these solutions, please refer to the support section.

6. SPECIFICATIONS

Key technical specifications for the Biosan Beauty Center 1000 Programs Mp6 Plus:

- **Model:** Beauty Center 1000 Programs Mp6 Plus
- **Brand:** Biosan
- **Item Weight:** 200 grams
- **ASIN:** B01KVUSROE
- **Manufacturer Reference:** X
- **Country of Origin:** Italy
- **Programs:** 1000 (Sport, Aesthetics, TENS)

Note: Specifications are subject to change without prior notice.

7. WARRANTY AND SUPPORT

7.1 Warranty Information

This non-original manual does not come with a direct product warranty from Biosan. For warranty information regarding the Biosan Beauty Center 1000 Programs Mp6 Plus device itself, please refer to the original manufacturer's documentation or the retailer from whom the device was purchased. Standard consumer rights apply to the purchase of this manual.

7.2 Customer Support

For questions regarding the content of this manual or general inquiries about the Biosan Beauty Center 1000 Programs Mp6 Plus, please contact your retailer or the manual's publisher. For technical support or service related to the device hardware, please contact the official Biosan customer service or an authorized service center.

Please note that this manual is a non-original version and support for the manual's content is provided by its publisher, not the original device manufacturer.