

## Waterpik WP-663

# Waterpik Aquarius Water Flosser Instruction Manual

Model: WP-663

## PRODUCT OVERVIEW

---

The Waterpik Aquarius water flosser is designed for advanced oral hygiene, offering a compact and contemporary design with powerful water flossing technology. It is equipped with various features to enhance your daily dental care routine. This device is clinically proven to remove up to 99.9% of plaque bacteria from treated areas, contributing to improved gum health and fresher breath. It is the #1 water flosser brand recommended by dental professionals and has earned the American Dental Association (ADA) Seal of Acceptance for its safety and effectiveness.



Image: The Waterpik Aquarius Water Flosser in blue, shown alongside its product packaging, highlighting its sleek design and key features.

## KEY FEATURES

---

- **Advanced Water Flossing:** Features enhanced pressure with 10 settings for a custom clean and a massage mode for gum stimulation.
- **Built-in Timer/Pacer:** Pauses briefly at 30 seconds and 1 minute to help track flossing time and ensure thorough cleaning.
- **Large Reservoir:** Provides 90+ seconds of flossing time without needing refills, with a removable 22 oz reservoir that is top-rack dishwasher safe for easy cleaning.
- **Easy to Use:** Small, maneuverable flossing wand with a convenient water on/off switch on the handle.
- **360-Degree Tip Rotation:** Allows the tip to easily reach all areas of the mouth, including back teeth.
- **7 Flossing Tips:** Includes a variety of tips for different dental needs and multiple users.
- **PrecisionPulse Technology:** Delivers a combination of water pressure and pulsations for effective plaque removal.

**7 specialty  
tips included**

**360° tip reaches  
whole mouth**

**flossing wand  
with water control**

**90 seconds of  
water capacity**

**10 safe & effective  
pressure settings**



Image: A visual breakdown of the Waterpik Aquarius's key features, including 7 specialty tips, 360-degree tip rotation, flossing wand with water control, 90 seconds of water capacity, and 10 pressure settings.



## Compact & Customizable

Advanced water flossing technology for **maximum performance.**



Image: An illustration demonstrating the 360-degree tip rotation feature of the Waterpik Aquarius, highlighting its ability to reach all areas of the mouth.

### WHAT'S INCLUDED

---

Your Waterpik Aquarius Water Flosser package includes the following components:

- Waterpik Aquarius Water Flosser Unit
- 7 Flossing Tips:
  - 4 x Precision Tips (JT-100E) - Designed for everyday use.
  - 1 x Plaque Seeker™ Tip (PS-100E) - Perfect for implants, crowns, bridges, and other restorations.
  - 1 x Orthodontic Tip (OD-100E) - Easily cleans around orthodontic brackets and wires.
  - 1 x Pik Pocket™ Tip (PP-100E) - Designed for gentle delivery of therapeutic rinses up to 6mm.

# What's Included in the Box?



Image: A clear view of the Waterpik Aquarius unit and the seven different flossing tips included in the product box.

## smart features



**On/Off**  
**Water Control**  
**On Handle**



Image: An illustration showcasing the four types of specialty flossing tips: Precision, Plaque Seeker, Orthodontic, and Pik Pocket, with descriptions of their uses.

## SAFETY INFORMATION

---

### Important Safety Precautions:

- Do not direct water under the tongue, into the ear, nose, or other delicate areas.
- Only use tips and accessories recommended by Water Pik, Inc.
- Use this product only as indicated in this instruction manual or as recommended by your dental professional.
- This product is designed for use in North America only (120VAC/60Hz).

## SETUP

---

Before first use, ensure the unit is placed on a flat, stable surface near a power outlet and a sink.

1. **Unpack:** Carefully remove all components from the packaging.
2. **Clean Reservoir:** Wash the removable reservoir with warm soapy water. It is top-rack dishwasher safe.
3. **Insert Tip:** Select the desired flossing tip and firmly insert it into the handle until it clicks into place.
4. **Fill Reservoir:** Fill the reservoir with warm water. You may also add a small amount of mouthwash if desired. Place the reservoir back onto the base.
5. **Plug In:** Plug the power cord into a grounded electrical outlet.

## OPERATING INSTRUCTIONS

---

Follow these steps for effective water flossing:

1. **Prepare:** Lean over the sink and place the tip into your mouth.
2. **Turn On:** Turn the unit ON using the power button on the base.
3. **Adjust Pressure:** Start with the lowest pressure setting (dial on the base) and gradually increase it over time to a comfortable level. The Aquarius offers 10 pressure settings.
4. **Aim Tip:** Aim the tip along the gumline at a 90-degree angle.
5. **Floss:** Turn the water ON using the switch on the handle. Pause briefly between teeth to allow water to flow from your mouth into the sink. Floss both the inside and outside surfaces of your upper and lower teeth.
6. **Timer/Pacer:** The built-in timer will briefly pause at 30 seconds and 1 minute to help you track your flossing time and ensure you clean all quadrants of your mouth.
7. **Finish:** When finished, turn the water OFF on the handle, then turn the unit OFF on the base. Empty any remaining water from the reservoir.

## smart features



### Tip Rotation

Reaches your whole mouth, **even back teeth**



Image: A close-up of the Waterpik Aquarius's pressure control dial, illustrating the 10 available settings for customized cleaning.

## smart features



### Large Reservoir

For **90+ seconds** of flossing time

Image: A hand demonstrating the on/off water control switch located directly on the handle of the Waterpik Aquarius flossing wand for convenient operation.

## MAINTENANCE

Proper maintenance ensures the longevity and optimal performance of your Waterpik Aquarius Water Flosser.

- **Daily Cleaning:** After each use, empty any remaining water from the reservoir.
- **Reservoir Cleaning:** The reservoir is top-rack dishwasher safe. For manual cleaning, wash with warm soapy water.
- **Tip Cleaning:** Tips can be cleaned by rinsing them under warm water. For thorough cleaning, soak tips in a solution of 2 parts water and 1 part white vinegar for 5-7 minutes.
- **Unit Exterior:** Wipe the exterior of the unit with a soft cloth and mild, non-abrasive cleaner as needed. Do not immerse the unit in water.
- **Prevent Mineral Deposits:** To prevent mineral deposits, especially if you have hard water, run a solution of 2 tablespoons of white vinegar to 16 ounces of warm water through the unit every 1 to 3 months.

## TROUBLESHOOTING

Problem	Possible Cause	Solution
Unit does not turn on.	Not plugged in; power outlet not working.	Ensure unit is securely plugged into a working electrical outlet.
Low water pressure.	Pressure setting too low; mineral deposits; kinked hose.	Increase pressure setting. Descale unit (see Maintenance). Check hose for kinks.
Water leaks from unit.	Reservoir not seated correctly; damaged hose/tip.	Ensure reservoir is firmly seated. Inspect hose and tips for damage; replace if necessary.
Unit is noisy.	Normal operation; air in pump.	Some noise is normal. Run unit with full reservoir until water flows smoothly to remove air.

If you experience issues not listed here, please refer to the Warranty and Support section for assistance.

## SPECIFICATIONS



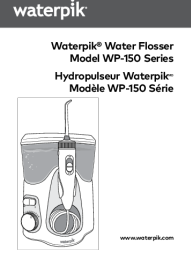

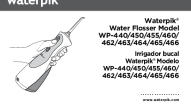

Feature	Detail
Model Number	WP-663
Product Dimensions	3.8 x 4.7 x 10.3 inches
Item Weight	1.46 Pounds
Power Source	Corded Electric (120VAC/60Hz for North America only)
Pressure Settings	10
Reservoir Capacity	22 oz (90+ seconds of flossing time)
Included Tips	7 (4 Precision, 1 Plaque Seeker, 1 Orthodontic, 1 Pik Pocket)
ADA Acceptance	Yes

## WARRANTY AND SUPPORT

The Waterpik Aquarius Water Flosser is backed by a **3-year limited manufacturer's warranty**. For detailed warranty information, please refer to the full product manual. For any product questions or support needs, Waterpik's Colorado-based support team is available to assist you. You can find additional support and contact information on the official Waterpik website: [www.waterpik.com](http://www.waterpik.com)





	<p><a href="#">Waterpik Cordless Water Flosser WP-360: Quick Start Guide and Welcome</a></p> <p>Get started with your Waterpik Cordless Water Flosser WP-360. This guide provides setup instructions, usage tips, and information on registration, support, and special offers from Waterpik.</p>
	<p><a href="#">Waterpik Water Flosser WP-60/65/70 User Manual and Guide</a></p> <p>Discover the benefits of the Waterpik Water Flosser WP-60/65/70. This guide provides essential information on setup, usage, safety precautions, and maintenance for optimal oral hygiene, backed by Waterpik's satisfaction guarantee.</p>
	<p><a href="#">Waterpik® Water Flosser WP-150 Series User Manual &amp; Guide</a></p> <p>Official user manual and guide for the Waterpik® Water Flosser Model WP-150 Series. Learn about important safeguards, product description, getting started, tip usage, cleaning, troubleshooting, and warranty information.</p>
	<p><a href="#">Waterpik Water Flosser WP-900 Quick Start Guide</a></p> <p>A step-by-step guide to setting up and using your Waterpik Water Flosser Model WP-900, covering reservoir preparation, tip insertion, unit priming, pressure adjustment, usage, and tip removal for effective oral hygiene.</p>
	<p><a href="#">Waterpik Water Flosser WP-440/450/455/460/462/463/464/465/466 User Manual</a></p> <p>User manual for Waterpik Water Flosser models WP-440 to WP-466, covering safety, product details, usage, troubleshooting, and warranty information.</p>
	<p><a href="#">Waterpik Waterflosser WP-560 Series Quick Start Guide</a></p> <p>Get started quickly with your Waterpik Waterflosser WP-560 series. This guide provides essential steps for unit preparation, charging, filling, and water flossing for optimal oral hygiene.</p>