

Manuals+

[Q & A](#) | [Deep Search](#) | [Upload](#)

manuals.plus /

- › [Colgate](#) /
- › [Colgate 360 Optic White Battery Powered Sonic Toothbrush User Manual](#)

Colgate CN05776A

Colgate 360 Optic White Battery Powered Sonic Toothbrush User Manual

Model: CN05776A

INTRODUCTION

The Colgate 360 Optic White Battery Powered Sonic Toothbrush combines the familiarity of a manual toothbrush head with the cleaning power of sonic technology. This toothbrush is designed to provide a comprehensive clean for your teeth, cheeks, tongue, and gums, contributing to a brighter, whiter smile.

SAFETY INFORMATION


Battery choke hazard. Do not charge or put in fire. May explode or leak causing injury.

SETUP

1. Unpacking and Battery Installation

Your Colgate 360 Optic White Sonic Toothbrush comes with one AAA battery included. To install or replace the battery:

1. Gently twist and pull off the bottom cap of the toothbrush handle.
2. Insert one fresh AAA battery, ensuring the positive (+) and negative (-) terminals are aligned correctly as indicated inside the compartment.
3. Replace the bottom cap, twisting it securely until it clicks into place.

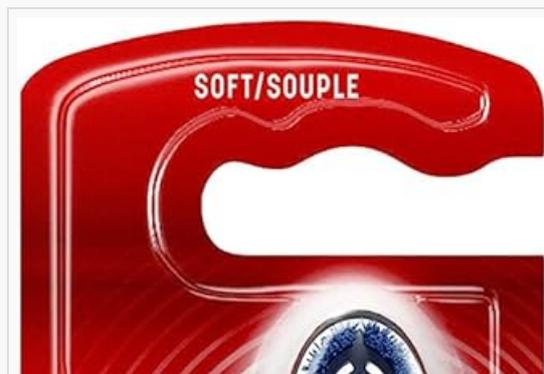

2X

WHITENING
ACTION†
L'ACTION
BLANCHISSANTE†

360°
OPTIC
WHITE®

SONIC POWER
**PUISSANCE
SONIQUE**

Colgate®

Colgate

POWERED TOOTHBRUSH
BROSSE À DENTS ALIMENTÉE



Energizer
MAX

Figure 1: Colgate 360 Optic White Battery Powered Sonic Toothbrush in its original packaging.

2. Initial Use

Before first use, rinse the brush head thoroughly under running water. Apply a small amount of toothpaste to the bristles.

Your browser does not support the video tag.

Video 1: An overview of the Colgate 360 Optic White Sonic Toothbrush, highlighting its features and deep cleaning power.

OPERATING INSTRUCTIONS

1. Turning On/Off

To activate the sonic vibrations, press the larger button located on the handle. To turn off the toothbrush, press the smaller button.

2. Brushing Technique

Dentists and hygienists recommend brushing your teeth for 2 minutes, twice a day. The Colgate 360 Optic White Sonic Toothbrush features soft Floss Tip bristles designed to clean effectively between teeth. The powerful sonic technology delivers 20,000 strokes per minute for a thorough clean.

- Hold the toothbrush at a slight angle towards your gum line.
- Gently move the brush in small circular motions, allowing the sonic vibrations to do the work.
- Brush all surfaces of your teeth: outer, inner, and chewing surfaces.



Figure 2: The soft Floss Tip bristles are designed to reach between teeth for a comprehensive clean.




*vs. brushing teeth alone with an ordinary flat-trim toothbrush

Figure 3: The vibrating bristles are powered by sonic technology for effective cleaning.

3. Tongue and Cheek Cleaning

For a healthier whole mouth clean, the toothbrush also features a textured pad on the back of the brush head for cleaning your tongue and cheeks. This helps remove bacteria that can cause bad breath.

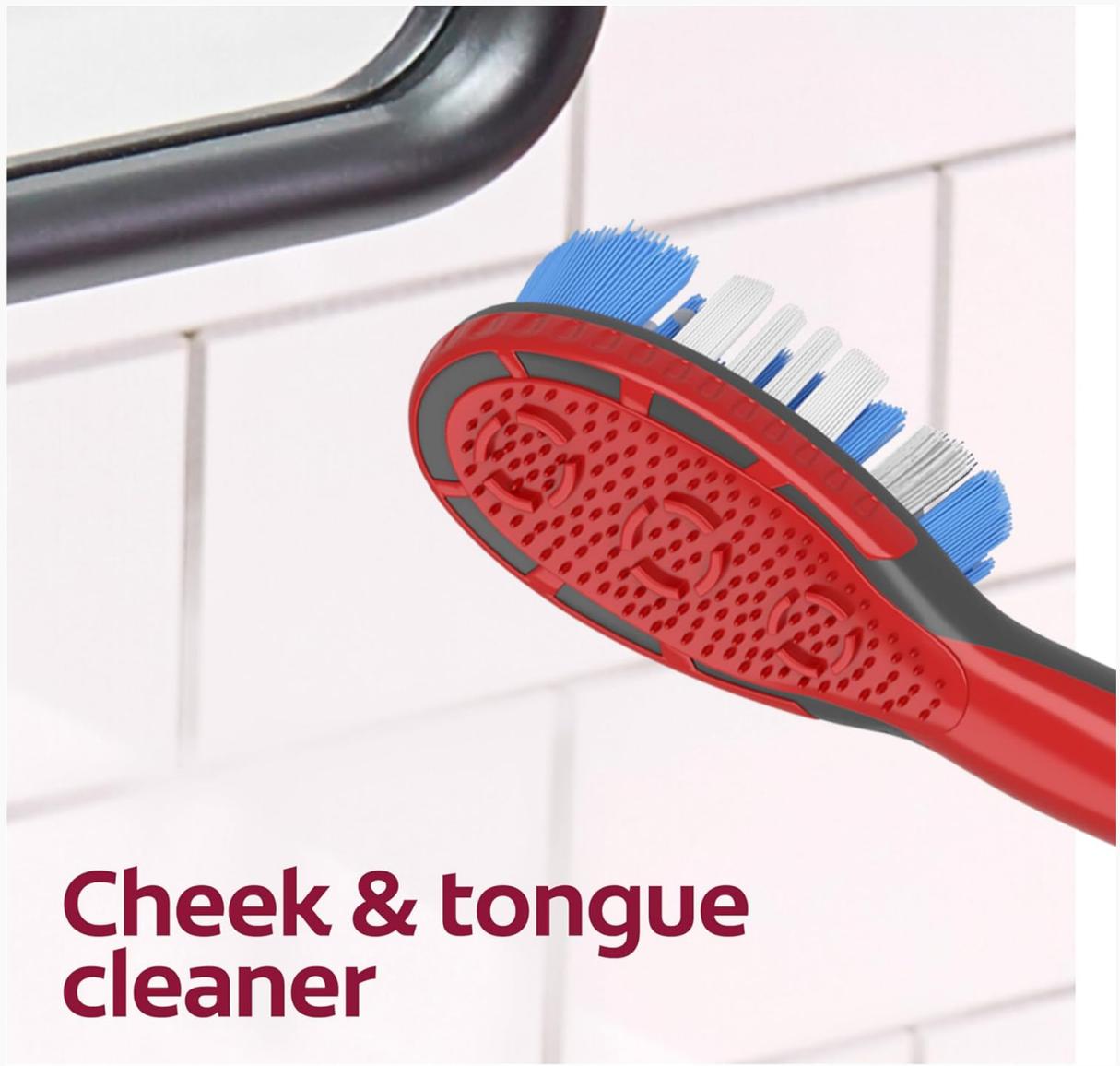

Cheek & tongue cleaner

Figure 4: The textured back of the brush head serves as a cheek and tongue cleaner.



Figure 5: This toothbrush is designed to clean the entire mouth, including teeth, cheeks, tongue, and gums.

Your browser does not support the video tag.

Video 2: A short demonstration of the Colgate 360 Optic White toothbrush in action, highlighting its whitening capabilities and ease of use.

MAINTENANCE

1. Replacing the Brush Head

Dentists and hygienists recommend replacing your toothbrush every 3 months, or sooner if the bristles are frayed. The brush head is designed to be easily replaceable.

2. Replacing the Battery

The toothbrush operates on one AAA battery. When the vibrations weaken, it's time to replace the battery. Follow the battery installation steps in the Setup section.

Vibrating bristles powered by sonic technology



Figure 6: The bottom cap can be twisted off to access the battery compartment for replacement.

3. Cleaning the Toothbrush

After each use, rinse the brush head thoroughly under running water to remove any toothpaste residue and debris. Wipe the handle with a damp cloth as needed.


TROUBLESHOOTING

If your toothbrush is not vibrating or performing as expected:

- **Check the battery:** Ensure the AAA battery is correctly installed and not depleted. Replace with a fresh battery if necessary.
- **Secure the cap:** Make sure the bottom cap is twisted on tightly to ensure proper contact for the battery.
- **Inspect the brush head:** Ensure the brush head is securely attached to the handle.


SPECIFICATIONS

Feature	Detail
---------	--------

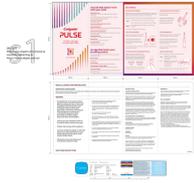

[Colgate Pulse Electric Toothbrush User Manual and Guide](#)

Comprehensive guide for the Colgate Pulse electric toothbrush, covering setup, usage, features, maintenance, warranty, and safety information.

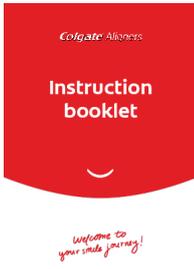

[Colgate Aligners Instruction Booklet: Your Guide to a Straighter Smile](#)

Comprehensive instruction booklet for Colgate Aligners, guiding you through the process of achieving a straighter smile with clear steps, care instructions, and important precautions.



[Colgate My Superhero Moment: A Smile Transformation Case Study](#)

Detailed case study by Colgate Professional showcasing a patient's smile transformation through advanced dental treatments, including prosthodontics and cosmetic dentistry, to achieve desired results quickly.