

## DOMAS SM9065N

# DOMAS Ab Belt Replacement Pads User Manual

Model: SM9065N

## INTRODUCTION

---

Thank you for choosing DOMAS Ab Belt Replacement Pads. These pads are designed to replace worn-out electrode pads on compatible DOMAS Ab toning belts, ensuring continued effective use of your device. This manual provides essential information for proper identification, installation, and maintenance of your replacement pads.

## IMPORTANT COMPATIBILITY INFORMATION


---

It is crucial to verify the compatibility of these replacement pads with your existing DOMAS Ab toning belt before installation. Incorrect pads may not function properly or could damage your device.

- **Correct Compatibility:** These replacement pads are designed for DOMAS Ab toning belts where the connect wires are **both black in color**.
- **Incorrect Compatibility:** If your Ab toning belt has connect wires that are **black and red in color**, these pads are NOT compatible. Please do not purchase or attempt to use these pads with such a belt.

Refer to the image below to confirm your belt's wire configuration.

## Ab Belt with black connect wire, then it's the right replacement pads for your belt.

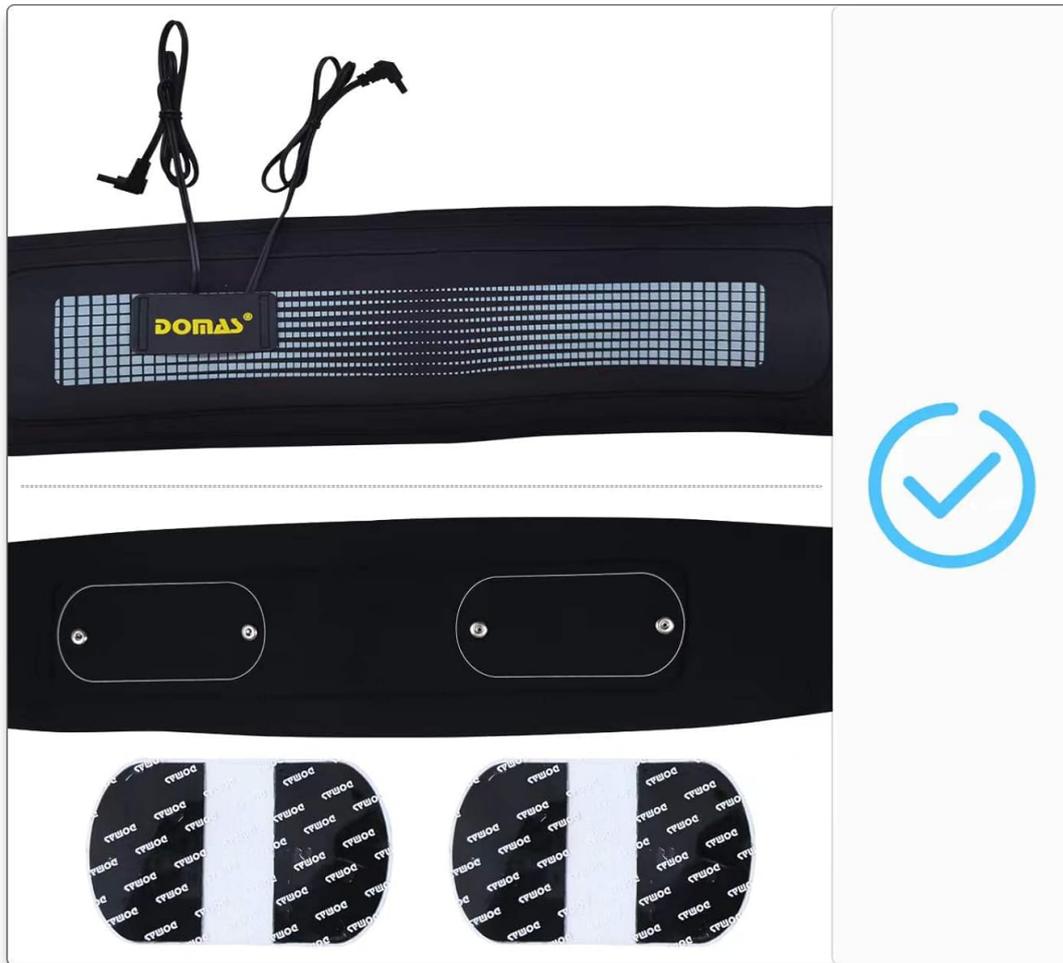

**Image Description:** This image displays a DOMAS Ab toning belt with two black connection wires clearly visible. Below the belt, two compatible replacement pads are shown, indicating the correct type of pads for this belt configuration. A green checkmark signifies compatibility.

### SETUP AND INSTALLATION

---

1. **Prepare the Belt:** Ensure your DOMAS Ab toning belt is turned off and disconnected from any power source.
2. **Remove Old Pads:** Carefully peel off the old, worn-out electrode pads from the belt. Dispose of them responsibly.
3. **Prepare New Pads:** Take a new replacement pad from its packaging. Gently peel off the protective film from the adhesive gel side of the pad.



**Image Description:** A hand is shown carefully peeling back the clear protective film from the black adhesive gel surface of a DOMAS replacement electrode pad. This illustrates the preparation step before applying the pad to the belt.

4. **Attach New Pads:** Align the new pad with the electrode connection points on your Ab toning belt. Press firmly and evenly to ensure good contact and adhesion.
5. **Connect Wires:** If applicable, ensure the black connection wires are securely attached to the metal studs on the new pads.
6. **Repeat:** Repeat the process for all pads requiring replacement.

## OPERATING INSTRUCTIONS

---

These replacement pads are passive components of your DOMAS Ab toning belt system. Their operation is directly tied to the functionality of the belt itself. Once installed, operate your Ab toning belt according to its specific user manual. The pads facilitate the electrical impulses from the belt to your muscles.

- Ensure the pads are clean and properly adhered to the belt for optimal conductivity.
- Apply the belt to clean, dry skin for best results.

## MAINTENANCE

---

Proper maintenance extends the lifespan and effectiveness of your replacement pads.

- **Cleaning:** After each use, gently wipe the gel surface of the pads with a damp, lint-free cloth. Do not use harsh chemicals or abrasive materials.
- **Storage:** After cleaning, re-apply the protective film to the gel surface of the pads. Store the pads in

their original packaging or a sealed bag to prevent drying out and contamination. Keep them in a cool, dry place away from direct sunlight.

- **Reusability:** The gel pads are reusable. Their lifespan depends on skin type, frequency of use, and proper maintenance. Replace pads when they lose their adhesiveness or when you notice a decrease in conductivity.



**Image Description:** This image shows two views of a single DOMAS replacement electrode pad. The left side displays the white cloth backing with two metal studs for connection. The right side shows the black adhesive gel surface, partially covered by a protective film, indicating its reusable nature.

## TROUBLESHOOTING

Problem	Possible Cause	Solution
Pads not sticking well	Gel surface is dirty or dry; pads are old.	Clean the gel surface with a damp cloth. Ensure skin is clean and dry before application. If pads are old, replace them.
Weak or inconsistent stimulation	Poor contact with skin; pads are worn out; incorrect pad type.	Ensure pads are firmly pressed onto the skin. Check pad adhesiveness and replace if necessary. Verify pad compatibility with your belt (black wires only).
Pads do not fit the belt	Incorrect pad model purchased.	Refer to the "Important Compatibility Information" section. If you have a belt with black and red wires, these pads are not suitable. Contact customer service for guidance on finding the correct pads.

If you are unsure about choosing the correct pads or encounter issues not listed here, please contact DOMAS customer service for assistance.

## SPECIFICATIONS

---

- **Model Number:** SM9065N
- **Pad Type:** Reusable Gel Electrode Pads
- **Backing Material:** Cloth
- **Compatibility:** DOMAS Ab toning belts with **black connect wires only**
- **Item Weight:** Approximately 159 g (for a pack of 4 pads)
- **Parcel Dimensions:** 17.53 x 13.72 x 1.78 cm

## CUSTOMER SUPPORT

---

For any questions, concerns, or assistance with your DOMAS Ab Belt Replacement Pads, please contact DOMAS customer service. Refer to your original product packaging or the DOMAS official website for contact details.

**Note:** This manual pertains specifically to the replacement pads. For instructions on operating your DOMAS Ab toning belt, please consult the manual provided with your belt device.