

CallToU CC01US+01Ax2

CallToU Wireless Caregiver Pager System User Manual

Model: CC01US+01Ax2

Brand: CallToU

1. INTRODUCTION

The CallToU Wireless Caregiver Pager System is designed to provide reliable communication between caregivers and individuals requiring assistance. This system is ideal for elderly individuals, patients, or those with disabilities, ensuring help is readily available at the touch of a button. The system includes a plug-in receiver and two portable call buttons.



Figure 1: CallToU Wireless Caregiver Pager System components.

WIRELESS CAREGIVER PAGER



The elderly need help
press the CALL button

Receiver receive
signal and ring

The elderly get help

CallToU



Figure 2: The CallToU system facilitates communication between an elderly individual and their caregiver.

2. SAFETY INFORMATION

- Do not expose the receiver to water or excessive moisture. The call button is waterproof, but the receiver is not.
- Keep the device away from extreme temperatures and direct sunlight.
- Do not attempt to disassemble or repair the device yourself. Contact customer support for assistance.
- Ensure the receiver is plugged into a properly functioning electrical outlet.
- Keep batteries out of reach of children.

3. PACKAGE CONTENTS

- 1 x Plug-in Receiver
- 2 x Call Buttons (with pre-installed 12V batteries)
- 2 x Lanyards
- 1 x Mini Screwdriver (for battery replacement)

- 2 x Double-sided Adhesive Stickers (for mounting call buttons)
- 1 x User Manual

4. SETUP

4.1 Receiver Setup

1. Unpack the receiver from its packaging.
2. Plug the receiver into a standard AC electrical outlet in a central location where the caregiver can easily hear the alerts.
3. The receiver will emit a sound when powered on, indicating it is ready for use.

FULLY EXPANDABLE SYSTEM

500ft distance coverage, Add extra receivers and transmitters for a wider range of security



Figure 3: Simple plug-and-play setup for the CallToU system.

4.2 Call Button Preparation

1. The call buttons come with batteries pre-installed.
2. **Attaching the Lanyard:** If desired, attach the lanyard to the call button. Use the provided mini screwdriver to gently pry open the back cover of the call button. Thread the small loop of the lanyard through the designated slot on the back cover, then pass the larger part of the lanyard through the small loop to secure it. Reattach the back cover.

Video 1: Instructions on how to attach the lanyard to the CallToU Wireless Call Button.

3. **Mounting the Call Button:** The call button can be worn around the neck using the lanyard, placed on a bedside table, or mounted to a wall using the provided double-sided adhesive stickers. Ensure the button is within reach of the user.



Figure 4: The call button features an IP55 waterproof rating, making it suitable for various environments.

4.3 Pairing Call Buttons (Pre-paired from factory)

The call buttons and receiver are pre-paired at the factory. No manual pairing is typically required for initial use. If you need to add new call buttons or re-pair existing ones, follow these steps:

1. Plug the receiver into an AC socket.
2. On the receiver, press the 'Forward' or 'Backward' key to select your desired ringtone (there are 55 ringtones to choose from).
3. Press and hold the 'Volume' key on the receiver for 5 seconds until the LED flashes and you hear a "ding" sound. This indicates the receiver is in pairing mode.
4. Press the call button you wish to pair. The receiver will make a "ding-ding" sound.
5. Press the call button again to confirm pairing. The receiver will emit the selected ringtone.

5. OPERATING INSTRUCTIONS

5.1 Using the Call Button

To request assistance, simply press the large red 'CALL' button on the portable transmitter. The receiver will immediately sound an alert and flash its LED indicator.



Figure 5: The call button provides immediate alert in case of an emergency.


5.2 Adjusting Receiver Settings

The receiver offers customizable settings for ringtones and volume levels.

- **Ringtone Selection:** Use the 'Forward' (▶▶) and 'Backward' (◀◀) buttons on the side of the receiver to cycle through the 55 available ringtones.
- **Volume Adjustment:** Use the 'Volume' () button to select one of the 5 adjustable volume levels, ranging from 0dB to 110dB.

How to Use The Caregiver Pager?

Plug-and-play

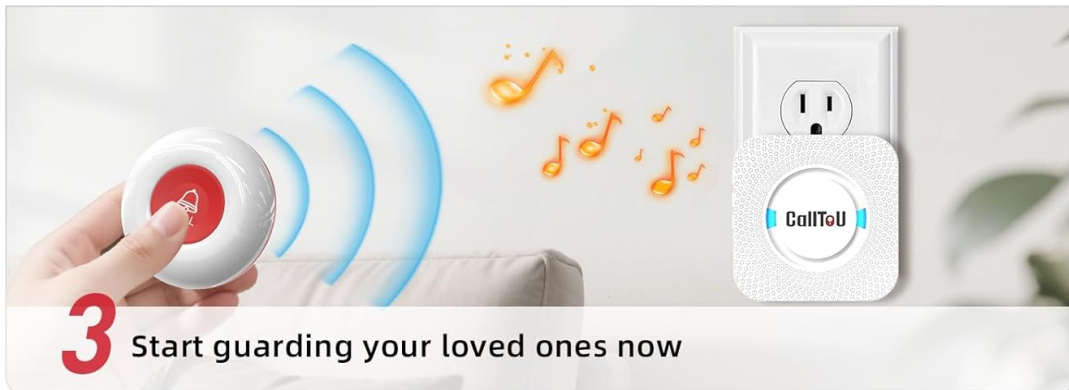

Products are pre-paired at the factory. No need to pair them yourself.

Figure 6: Receiver controls for ringtone and volume adjustments.

6. MAINTENANCE

6.1 Call Button Battery Replacement

The call button uses a 12V battery. When the button's indicator light no longer flashes or the signal range decreases significantly, it's time to replace the battery.

1. Use the provided mini screwdriver to gently pry open the back cover of the call button.
2. Carefully remove the old 12V battery.
3. Insert a new 12V battery, ensuring correct polarity (+/-).
4. Snap the back cover securely into place.

Your browser does not support the video tag.

Video 3: Instructions on how to replace the battery in the CC01 Wireless Call Button.

6.2 Cleaning

Wipe the devices with a soft, dry cloth. Do not use abrasive cleaners or solvents.

7. TROUBLESHOOTING

- **No alert from receiver:**

- Ensure the receiver is properly plugged into a live electrical outlet.
- Check if the receiver's volume is set to an audible level (not 0dB).
- Verify the call button battery is not depleted. Replace if necessary.
- Re-pair the call button with the receiver following the steps in Section 4.3.
- Ensure the call button is within the operating range (up to 500 feet in open areas). Walls and other obstacles can reduce range.

- **Reduced signal range:**

- Replace the call button battery.
- Relocate the receiver to a more central position.
- Avoid placing the receiver near large metal objects or other electronic devices that may cause interference.

- **Call button not responding:**

- Check the battery.
- Ensure the button is not physically damaged.
- Attempt to re-pair the button with the receiver.

8. SPECIFICATIONS

Feature	Detail
Model Number	CC01US+01Ax2
Voltage (Receiver)	110-260 Volts
Control Method	Touch
Noise Level (Max)	110 Decibels
Mounting Type	Wall Mount (Receiver)
Maximum Range	500+ Feet (Open Area)
Sensor Technology	Contact Sensor
Call Button Battery	2 x 12V Alkaline batteries (included)
Item Weight	8.4 ounces
Product Dimensions	3.23 x 3.23 x 0.94 inches
Waterproof Rating (Call Button)	IP55
Ringtone Options	55
Volume Levels	5 (0dB to 110dB)







9. WARRANTY AND SUPPORT

CallToU is committed to providing quality products and customer satisfaction. While specific warranty details are not provided in this manual, a dedicated customer service and professional technology team are available to assist you with any questions or issues you may encounter.

For support, please refer to the contact information provided with your product packaging or visit the official CallToU website.

© 2025 CallToU. All rights reserved.

Related Documents - CC01US+01Ax2

	<p>CallToU Caregiver Pager System CC01 Series Instruction Manual</p> <p>Comprehensive guide to the CallToU CC01 Series Caregiver Pager System, detailing setup, features, troubleshooting, and specifications for enhanced home care and elderly monitoring.</p>
	<p>CallToU Caregiver Pager System User Manual</p> <p>User manual for the CallToU Caregiver Pager System, including receiver model CC28, transmitter model BT009, and wrist call button E-05W. Learn about features, operation, pairing, troubleshooting, and specifications.</p>
	<p>CallToU CC28 Receiver & BT009-WH Transmitter User Manual</p> <p>User manual for the CallToU CC28 receiver and BT009-WH transmitter, a home care and elderly caregiver pager system. Learn about features, operation, installation, and troubleshooting.</p>
	<p>CallToU CC28 Caregiver Pager FAQ</p> <p>Frequently Asked Questions about the CallToU CC28 Caregiver Pager, covering topics like range, battery issues, pairing, and water resistance.</p>
	<p>CALLTOU CC01 Series Caregiver Pager System: Instruction Manual & Guide</p> <p>Comprehensive guide to the CALLTOU CC01 Series Caregiver Pager System. Learn about setup, pairing, features, specifications, troubleshooting, and precautions for reliable elderly care.</p>
	<p>CallToU TY01 Wireless Doorbell Operation Instructions and Specifications</p> <p>Comprehensive guide to operating the CallToU TY01 wireless doorbell system, covering setup, features, Tuya APP integration, and technical specifications for optimal performance.</p>

