

AUVON PRO18AB

AUVON PRO18AB Rechargeable TENS Unit Muscle Stimulator User Manual

Model: PRO18AB

INTRODUCTION

This manual provides comprehensive instructions for the safe and effective use of your AUVON PRO18AB Rechargeable TENS Unit Muscle Stimulator. Please read this manual thoroughly before operating the device and retain it for future reference.



Image: The AUVON TENS Unit, electrode pads, lead wires, USB cable, and storage pouch.



Image: Visual representation of common body areas where the TENS unit can be applied for pain relief, including shoulder, arm, wrist, ankle, neck, elbow, knee, and lower back.

SAFETY INFORMATION AND WARNINGS

WARNING: Failure to follow these safety instructions may result in injury or adverse reactions.

This product is not suitable for individuals with certain medical conditions or circumstances. Consult with your physician before using this device if you have any of the following conditions:

- Implanted defibrillator, cardiac pacemaker, or any other implanted metallic or electronic device.
- Heart rhythm problems, epilepsy, brain diseases, or neurological disorders.
- Undiagnosed chronic pain.
- Pregnancy or labor.
- Cancer.
- Sensitivity to electrotherapy.
- Lack of skin sensation or an abnormal skin condition.

- Any condition requiring the active supervision of a physician.

Additional Precautions:

- Avoid excessive use of electrotherapy.
- Avoid use above the neck, especially on the head or neck area.
- Do not use this product if you are under 22 years of age.
- Individual symptoms, situations, circumstances, and results may vary. The information provided is not intended to be used for medical diagnosis or treatment or as a substitute for professional medical advice. Please consult your physician or qualified health care provider regarding your condition and appropriate medical treatment.



Image: Illustration of a person experiencing neck discomfort, indicating a potential area for TENS application.

PACKAGE CONTENTS

Upon opening your AUVON TENS Unit package, please verify that all the following items are included:

- 1 x AUVON TENS Muscle Relaxer (Main Unit)
- 8 x Square Electrode Pads (2"x2" Premium)
- 2 x Lead Wires (1/2 Wire)
- 1 x Pads Holder
- 1 x USB Cable
- 1 x Pouch
- 1 x Professional User Manual

What You Receive



Image: All components included in the AUVON TENS Machine package, laid out for inspection.

PRODUCT FEATURES

- **Professional TENS Machine:** Designed to stimulate sensory nerves and muscles for pain management.
- **Upgraded Gel Pads:** Low impedance solid gel for improved self-stick performance and extended lifespan (up to 45 uses). Compatible with 2mm pigtail TENS unit pads.
- **24 Pre-programmed Modes:** Offers 24 user-friendly massage modes and adjustable timing for varied therapy.
- **Rechargeable Lithium Battery:** Provides up to 10 hours of continuous use. Features a power indicator for battery life monitoring. Fully recharges in approximately 1 hour.

Up to 24 Mode Choices

24 mode options and controls to give you a relaxing massage with little fuss.



24 Modes



20 Intensity



Pocket Size




10h Usage



Steady Current



USB Rechargeable

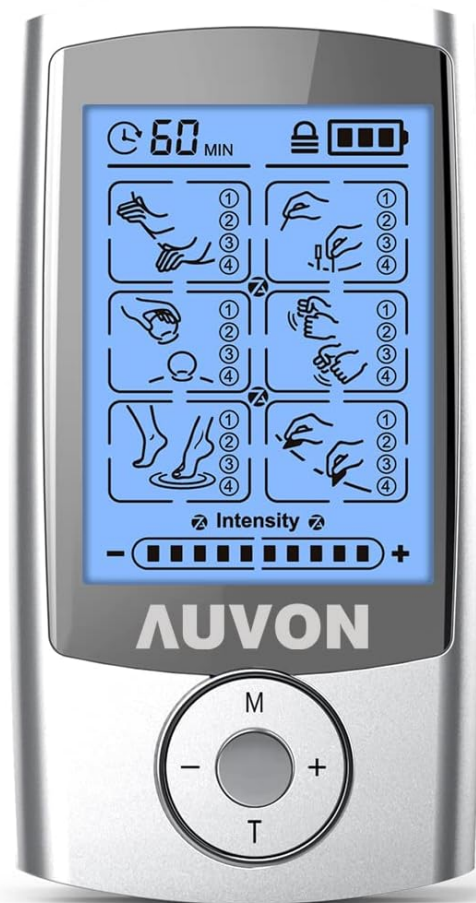

*** The 24 sub-mode choices selection is under the 6 massage modes.**

Image: AUVON TENS Unit display highlighting 24 modes, 20 intensity levels, 10-hour usage, steady current, and USB rechargeable capability.



Good Stickiness

American Gel, More Reuses

Lower Impedance, Better Conductivity.


	AUVON GEL	OTHERS GEL
Country of origin	America 	/
Number of reuse	45 times	30 times
Resistance	45 Ω	50 Ω
Biocompatibility	√	×
Production Tool	Automated Machine	Hands
Certificate Report	510K Biocompatibility Report	×

Image: Table comparing AUVON gel pads with other brands, detailing country of origin, reuses, resistance, biocompatibility, production tool, and certification.

Current Stability

Classic pin-style electrodes are better than other button-style ones on current stability and anti-corrosion.

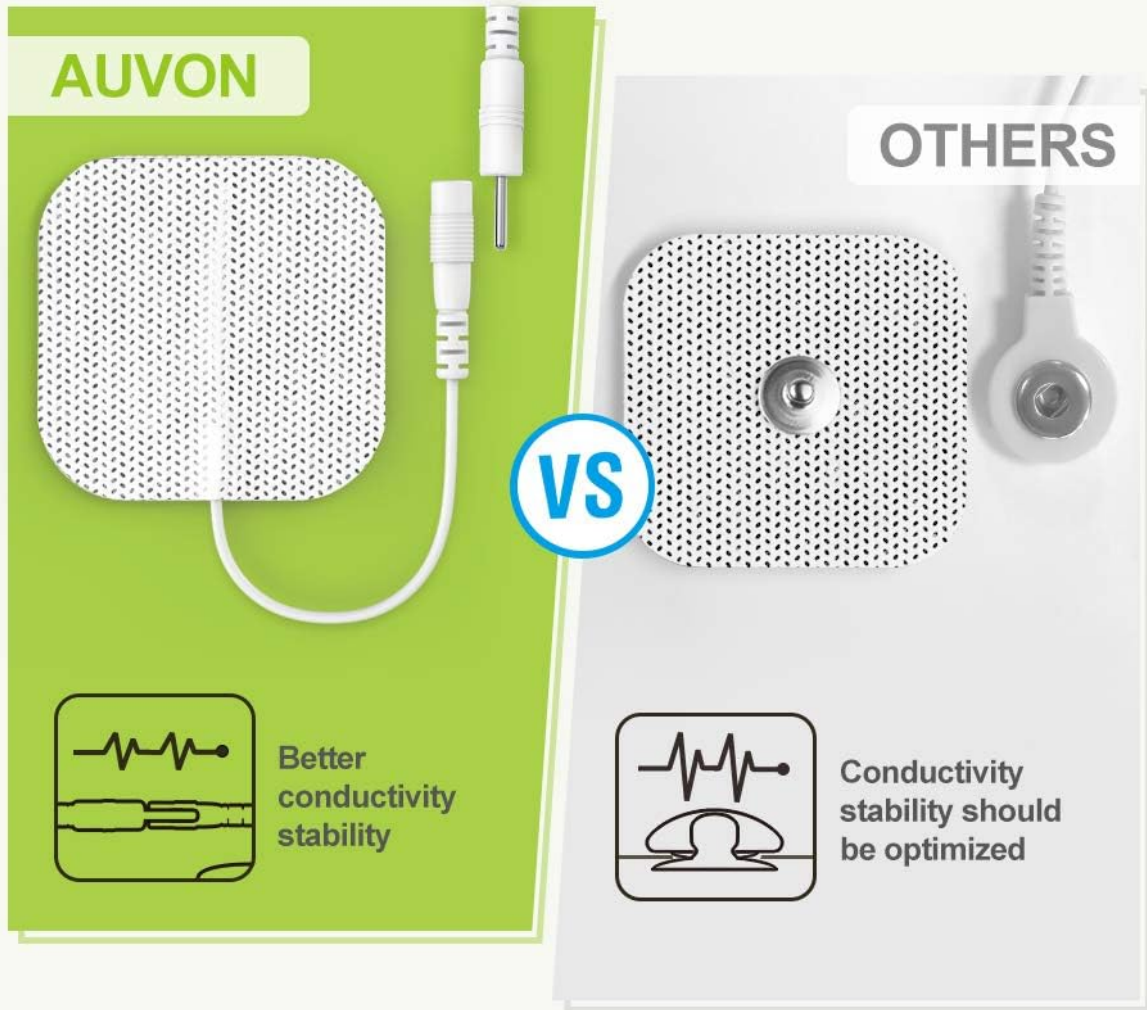


Image: Visual comparison between AUVON pin-style electrodes and button-style electrodes, emphasizing AUVON's better conductivity stability.

Superior Quality

Precise electronic circuitry highly ensures for optimal safety (CE / ROHS Cleared).

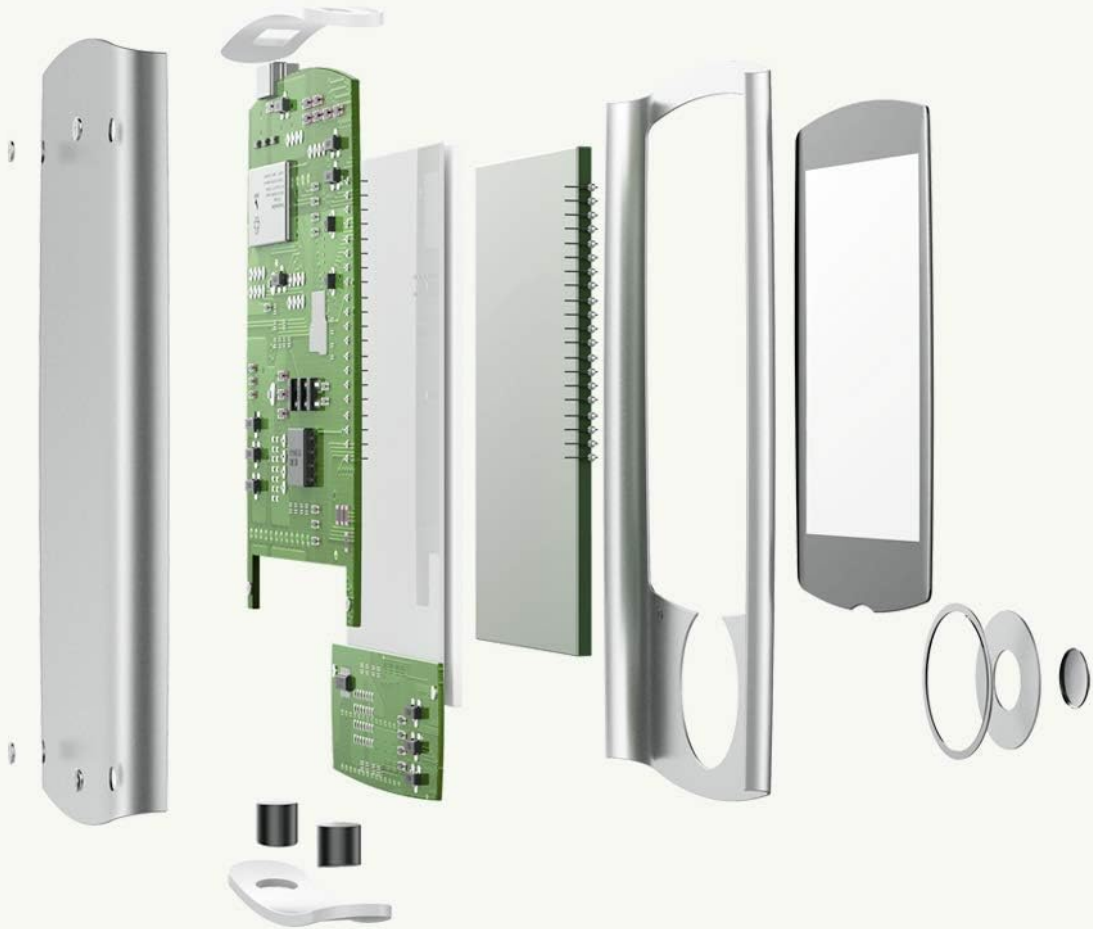


Image: Exploded diagram showing the internal components and precise electronic circuitry of the AUVON TENS unit.

SETUP

1. **Charge the Device:** Connect the TENS unit to a power source using the provided USB cable. The power indicator will show charging status. Ensure the unit is fully charged before first use.
2. **Attach Lead Wires:** Insert the lead wires into the ports on the top of the TENS unit.
3. **Connect Electrode Pads:** Attach the electrode pads to the ends of the lead wires. Ensure a secure connection.
4. **Prepare Skin:** Clean and dry the skin area where the pads will be applied. Avoid applying pads over broken skin, rashes, or irritated areas.
5. **Apply Electrode Pads:** Carefully peel the electrode pads from their protective film and firmly place them on the desired treatment area. Ensure both pads are applied to the skin for the unit to function correctly.

Set Your TENS Easily



Step 1

Main Mode: Press “M”

Sub-mode: Press
round button



Step 2

Intensity: Press +/-



Step 3

Timer: Press T
Lock: Long press
the round button

Image: The AUVON TENS unit connected to a laptop via USB cable for charging.

OPERATING INSTRUCTIONS

1. **Power On:** Press and hold the power button (usually the 'M' button) to turn on the device.
2. **Select Main Mode:** Press the 'M' button to cycle through the main massage modes. The display will show the selected mode.
3. **Select Sub-Mode:** After selecting a main mode, press the round button (often labeled 'T' or the central button) to select a sub-mode if available.
4. **Adjust Intensity:** Use the '+' and '-' buttons to increase or decrease the intensity level. Start with a low intensity and gradually increase until a comfortable sensation is felt.
5. **Set Timer:** Press the 'T' button to adjust the treatment duration. The device will automatically shut off once the timer expires.
6. **Lock Function:** Long press the round button (often the 'T' button) to lock the screen and prevent accidental changes during treatment. Repeat to unlock.
7. **Power Off:** Press and hold the power button to turn off the device after use.



About AUVON

AUVON is dedicated to helping people live better life! We have donated over 200,000 healthcare products to nursing homes and we are always there to share our love!



Image: A three-step visual guide demonstrating how to select modes, adjust intensity, and set the timer/lock function on the AUVON TENS unit.

Electrode Placement Guidelines

Proper placement of electrode pads is essential for effective therapy. Refer to the diagram below for common application areas:



Image: Diagram illustrating recommended electrode pad placements for the neck, shoulder, back, arm, and knee. For optimal results, place pads on clean, dry skin. Ensure pads are not overlapping and are at least one inch apart. Consult a healthcare professional for specific placement recommendations for your condition.

MAINTENANCE

- **Cleaning the Device:** Wipe the TENS unit with a soft, damp cloth. Do not use abrasive cleaners or immerse the device in water.
- **Electrode Pad Care:** After each use, place the electrode pads back onto the provided pads holder or the original plastic sheet. To prolong adhesion, you may lightly dab the sticky side with a drop of water.
- **Storage:** Store the device and accessories in the provided pouch in a cool, dry place, away from direct sunlight and extreme temperatures.
- **Battery Care:** Recharge the device regularly, even if not in frequent use, to maintain battery health.

TROUBLESHOOTING

Problem	Possible Cause	Solution
No sensation or weak sensation	Pads not properly adhered to skin. Intensity level too low. Lead wires not securely connected. Pads are dry or worn out. Only one pad is applied.	Ensure pads are firmly pressed onto clean, dry skin. Increase intensity gradually using the '+' button. Check all lead wire connections to the unit and pads. Replace old or dry pads with new ones. Ensure both pads are applied to the skin.
Device does not turn on	Battery is depleted. Device malfunction.	Charge the device fully. Contact customer support if charging does not resolve the issue.
Pads losing stickiness	Pads are dirty or dry. Pads have reached end of lifespan.	Clean skin before application. Lightly moisten pads with water. Replace pads with new ones.




SPECIFICATIONS




- **Model Number:** PRO18AB
- **Manufacturer:** AUVON
- **Product Dimensions:** 9.47 x 4.98 x 3.53 cm
- **Weight:** 250 g
- **Power Source:** Rechargeable Lithium Battery (included)
- **Modes:** 24 Pre-programmed modes
- **Intensity Levels:** 20
- **Country of Origin:** China

WARRANTY AND SUPPORT

AUVON provides a 12-month warranty for this product. For any questions, concerns, or technical assistance, please contact AUVON customer service. Refer to the contact information provided with your product packaging or visit the official AUVON website for support details.

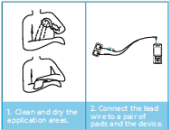Related Documents - PRO18AB

 USER MANUAL AUVON® TENS Massager Relaxation, Anytime, Anywhere Model: SM9126 Device Name: TENS & EMS	AUVON SM9126 TENS Massager User Manual: Features, Operation, and Troubleshooting Comprehensive user manual for the AUVON SM9126 TENS Massager. Learn about its features, operating instructions, safety precautions, troubleshooting tips, and maintenance for effective pain relief and muscle stimulation.
 MANUAL DE USUARIO Masajeador TENS Relajación cualquier lugar cualquier momento Modelo: AS8016 Nombre del dispositivo: TENS y Estimulador Muscular Activado	AUVON AS8016 TENS Massager User Manual - Pain Relief and Muscle Stimulation Comprehensive user manual for the AUVON AS8016 TENS massager. Learn about features, operation, safety precautions, troubleshooting, and maintenance for effective pain relief and muscle stimulation.
Quick Guide Attention The purpose of this Quick Guide is to help you understand the operation of the TENS unit much more quickly, and the contents of the quick guide are not a substitute for the user manual. Please read the user manual carefully before use to understand the contraindications and warnings. Please consult your doctor before use if you're not sure whether you can use this device. Operation Note: Do not use it while charging. 	AUVON TU2224-A TENS and Muscle Stimulator User Manual - Pain Relief Guide Comprehensive user manual for the AUVON TU2224-A TENS and Powered Muscle Stimulator. Learn how to use, troubleshoot, and care for your device for effective pain relief and muscle stimulation.

<div data-bbox="129 107 298 405"><p>Quick Guide</p><p>Attention</p><p>The purpose of this Quick Guide is to help you understand the operation of the TENS unit much more quickly, and the contents of the quick guide are not a substitute for the user manual. Please read the user manual carefully before use to understand the contraindications and warnings. Please consult your doctor before use if you're not sure whether you can use this device.</p><p>Operation</p><p>Note: Do not use it while charging.</p><div data-bbox="129 277 298 405"><p>1. Clean and dry the applicator wires.</p><p>2. Connect the lead wire to a pair of pads and the device (press the unit off).</p></div></div>	<p>AUVON TENS & EMS Unit Quick Guide and User Manual - Model AS8016</p> <p>Comprehensive quick guide and user manual for the AUVON TENS and Powered Muscle Stimulator Model AS8016. Learn how to operate the device, understand TENS and EMS therapies, proper pad placement for various body areas, battery charging, troubleshooting, and safety precautions for effective pain relief and muscle stimulation.</p>
<div data-bbox="129 750 298 1030"><p></p><p>USER MANUAL</p><p>Dual Channel TENS Massager</p><p>Relaxation, Anywhere, Anytime</p><p>Model AS8012</p><p>Device Name: TENS and Powered Muscle Stimulator</p></div>	<p>AUVON AS8012 Dual Channel TENS Massager User Manual</p> <p>Comprehensive user manual for the AUVON AS8012 Dual Channel TENS Massager, detailing features, operating instructions, safety precautions, and troubleshooting for effective pain relief and muscle stimulation.</p>
<div data-bbox="129 1124 298 1404"><p></p><p>USER MANUAL</p><p>TENS and Powered Muscle Stimulator</p><p>Please read the product literature carefully before using the product.</p><p>Model TU2124-A</p><p>Device Name: TENS and Powered Muscle Stimulator</p></div>	<p>AUVON TU2124-A TENS and Powered Muscle Stimulator User Manual</p> <p>User manual for the AUVON TU2124-A TENS and Powered Muscle Stimulator. Learn about its features, intended uses, warnings, operating instructions, placement, troubleshooting, and maintenance.</p>