

## Manuals+

[Q & A](#) | [Deep Search](#) | [Upload](#)

[manuals.plus](#) /

› [Mynt](#) /

› [Mynt Foot Massager \(Model 30906809483346\) User Manual](#)

## Mynt 30906809483346

# Mynt Foot Massager

USER MANUAL (MODEL: 30906809483346)

### Important Safety Information

Please read all instructions carefully before using the Mynt Foot Massager. Retain this manual for future reference.

- Consult a medical professional prior to using if you have any health concerns, such as a pacemaker, or are pregnant.
- If you have sensitive feet, consider wearing socks and working up from the lowest intensity setting.
- Do not use the massager in wet environments or near water.
- Keep out of reach of children.
- Unplug the device when not in use or before cleaning.
- Do not attempt to repair or disassemble the unit yourself. Contact customer support for assistance.

**Legal Disclaimer:** Statements regarding dietary supplements have not been evaluated by the FDA and are not intended to diagnose, treat, cure, or prevent any disease or health condition.

### Product Overview

The Mynt Foot Massager is designed to provide a comprehensive foot massage experience, combining Shiatsu, tapping, heat, and air compression functions. It aims to alleviate foot discomfort and promote relaxation.



Figure 1: Mynt Foot Massager (Top View)

### Key Features:

- **Multi-Functional Massage:** Combines Shiatsu, tapping, heat, and air compression for a varied massage experience.
- **Targeted Relief:** Features patented arch-tapping technology and 22 massage nodes for comprehensive foot coverage.
- **Adjustable Settings:** Customizable roller speed, massage modes, air pressure, and heat levels.
- **Hygienic Design:** Equipped with removable and washable inner sleeves for easy cleaning.

# *Immersive Experience*

*Total coverage and air compression deliver a fuller-foot massage.*



Figure 2: Internal Massage Mechanisms

# Superior Technique

*Innovative design ensures masseuse-grade, pressure point treatment.*



Figure 3: 22 Massage Nodes for Targeted Relief

# Superior Heating

*Warms quicker and more consistently.*

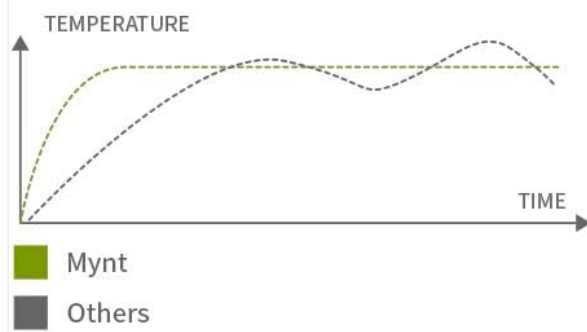


Figure 4: Heating Performance Comparison

## Product Video:

Your browser does not support the video tag.

Video 1: Demonstration of the Mynt Tapping Foot Massager's features and operation.

## Getting Started

The Mynt Foot Massager is designed for immediate use out of the box with minimal setup.

1. Unpack the foot massager from its packaging.
2. Place the massager on a flat, stable surface near an electrical outlet.
3. Ensure the power cord is free from kinks or damage.
4. Plug the power cord into a standard electrical outlet.
5. Sit comfortably and place your feet into the foot inserts.

## Operating Instructions



## Control Panel:



Figure 5: Control Panel Layout

- **Power Button:** Press to turn the unit on or off.
- **Mode Button:** Cycles through different massage programs (e.g., Shiatsu, Tapping, combination modes).
- **Roller Speed:** Adjusts the speed of the massage rollers.
- **Air Pressure:** Controls the intensity of the air compression massage.
- **Intensity:** Adjusts the overall massage intensity.
- **Heat Button:** Activates or deactivates the soothing heat function.

## Using the Massager:

1. Plug in the massager and press the **Power Button** to turn it on.
2. Select your desired massage program using the **Mode Button**.
3. Adjust the **Roller Speed**, **Air Pressure**, and **Intensity** to your preference.
4. Activate the **Heat Function** if desired.
5. Enjoy your foot massage. The device typically has an auto-shutoff feature after a set duration for

safety.

6. To turn off, press the **Power Button** again or unplug the unit.

## Care and Maintenance

Proper care ensures the longevity and hygiene of your Mynt Foot Massager.

### Cleaning the Inner Sleeves:



Figure 6: Removing Washable Inner Sleeves

- The inner sleeves are removable and washable for sanitary use.
- Unzip the sleeves from the foot inserts.
- Hand wash the sleeves with mild soap and water, or machine wash on a gentle cycle.
- Allow the sleeves to air dry completely before reinserting them into the massager.
- Ensure the sleeves are properly zipped back into place before next use.

### General Cleaning:

- Wipe the exterior of the massager with a soft, damp cloth.

- Do not use abrasive cleaners or solvents.
- Ensure the device is unplugged before cleaning.

## Troubleshooting

If you encounter issues with your Mynt Foot Massager, refer to the following common solutions:

Problem	Possible Cause	Solution
Massager does not turn on.	Not plugged in; power outlet issue; power button not pressed.	Ensure the power cord is securely plugged into a working outlet. Press the power button firmly.
Massage intensity is too weak/strong.	Incorrect intensity setting.	Adjust the "Intensity" or "Air Pressure" settings on the control panel. Consider wearing socks if too strong.
Heat function is not working.	Heat function not activated.	Press the "Heat" button on the control panel. Allow a few minutes for the heat to become noticeable.
Unusual noise during operation.	Obstruction; internal component issue.	Turn off and unplug the unit. Check for any obstructions in the footwells. If the noise persists, contact customer support.

If the problem persists after attempting these solutions, please contact Mynt customer support.

## Specifications

Feature	Detail
Brand	Mynt
Model Number	30906809483346
Product Dimensions	16.3 x 15.35 x 9.57 inches
Item Weight	14.55 Pounds
Power Source	Electric
Special Features	Heat, Shiatsu, Tapping, Air Compression
Massager Form	Foot Massager Machine
Date First Available	February 12, 2018

## Warranty and Support

For warranty information, product support, or to purchase accessories, please visit the official Mynt store or contact their customer service.



[Visit the Mynt Store on Amazon](#)

© 2023 Mynt. All rights reserved.