

## Manuals+

[Q & A](#) | [Deep Search](#) | [Upload](#)

[manuals.plus](#) /

› [Waterpik](#) /

› [Waterpik Complete Care 9.0 Electric Toothbrush with Water Flosser \(CC-01\) User Manual](#)

## Waterpik CC-01

# Waterpik Complete Care 9.0 Electric Toothbrush with Water Flosser (CC-01) User Manual

## INTRODUCTION

The Waterpik Complete Care 9.0 combines a water flosser and an electric toothbrush into a single, compact unit. This integrated oral care system is designed to save counter space and power outlets while providing comprehensive cleaning for your teeth and gums. The Triple Sonic toothbrush offers multiple cleaning modes, and the water flosser features adjustable pressure settings to effectively remove plaque and debris.





Image: The Waterpik Complete Care 9.0 unit, featuring a white base with a clear water reservoir, an electric toothbrush on the left, and a water flosser handle on the right.

## SAFETY INFORMATION

---

Keep out of reach of children. If ingested, seek medical help immediately.

## PACKAGE CONTENTS

---

Your Waterpik Complete Care 9.0 package includes the following items:

- Waterpik Complete Care 9.0 Base Unit
- Triple Sonic Electric Toothbrush
- Water Flosser Handle
- 5 Water Flossing Tips (various types for different dental needs)
- 2 Brush Heads (advanced design with high/low bristle pattern)
- Premium Toothbrush Travel Case
- Tip Storage Case

# Included in the box



Image: A top-down view of all components included in the Waterpik Complete Care 9.0 package, neatly arranged. This includes the main unit, electric toothbrush, water flosser handle, various tips, brush heads, and storage cases.

## SETUP

1. **Unpack the Product:** Carefully remove all components from the packaging.
2. **Charge the Toothbrush:** Place the electric toothbrush into its designated charging dock on the base unit. Ensure the unit is plugged into a 120VAC/60Hz outlet (for North America only). The charge indicator lights on the toothbrush will illuminate.
3. **Fill the Water Reservoir:** Lift the hinged reservoir lid and remove the clear reservoir from the base. Fill the reservoir with lukewarm water. Place the reservoir back onto the base, ensuring it is seated correctly.
4. **Attach a Water Flossing Tip:** Select a desired water flossing tip and insert it firmly into the center of the water flosser handle until it clicks into place.

# waterpik® SONIC 9.0

## completecare

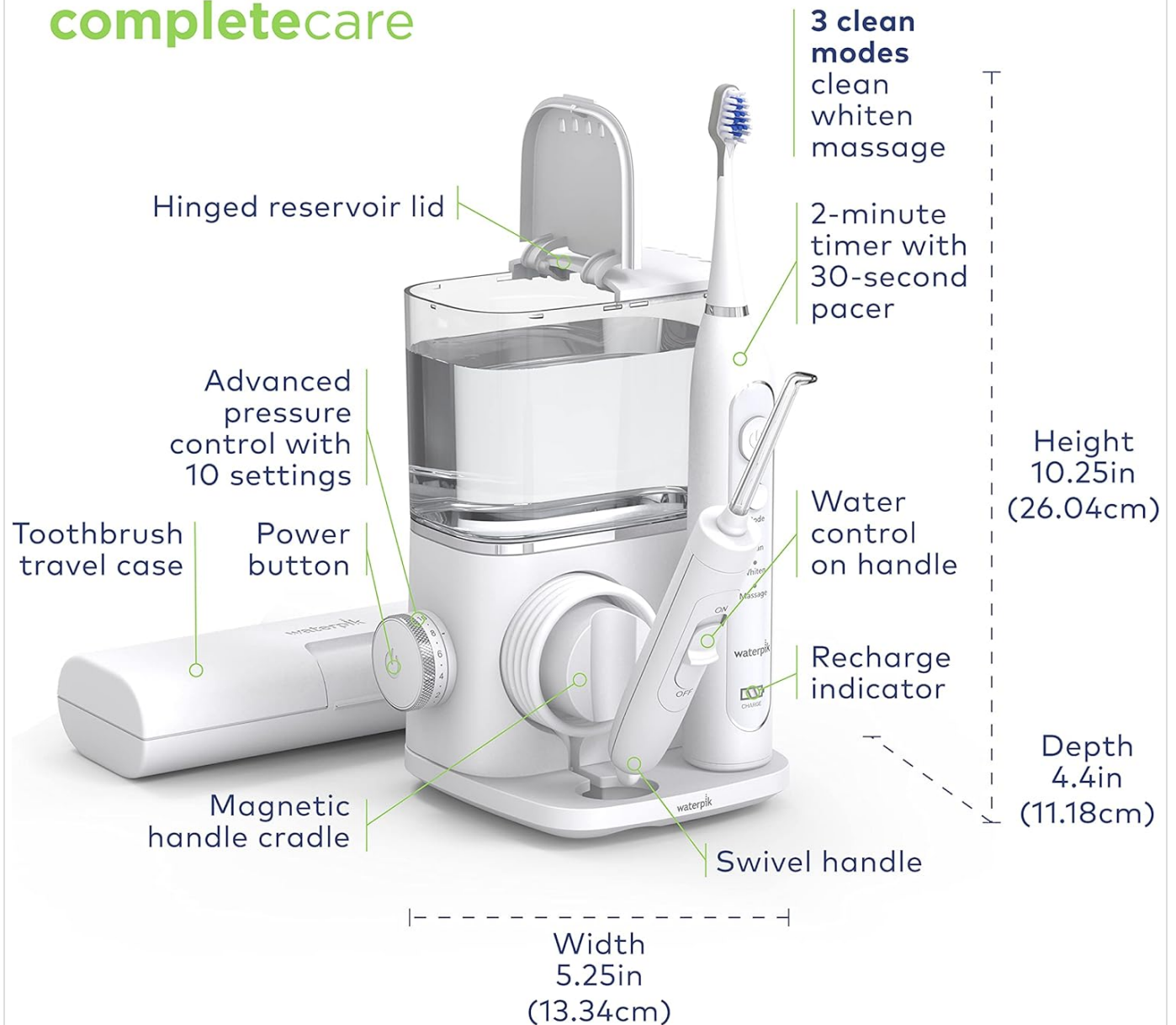


Image: A detailed diagram of the Waterpik Complete Care 9.0, highlighting key features such as the hinged reservoir lid, advanced pressure control with 10 settings, power button, magnetic handle cradle, water control on handle, recharge indicator, toothbrush travel case, and overall dimensions.

## OPERATING INSTRUCTIONS

### Electric Toothbrush Usage

1. **Attach Brush Head:** Securely attach one of the provided brush heads to the electric toothbrush handle.
2. **Apply Toothpaste:** Apply a small amount of toothpaste to the brush head.
3. **Select Mode:** Press the 'Mode' button to cycle through the three available brushing modes: 'Clean', 'Whiten', and 'Massage'. Choose the mode that suits your preference.
4. **Brush Teeth:** Place the brush head against your teeth at a slight angle towards the gumline. Press the power button to start brushing. The toothbrush features a 2-minute timer with a 30-second pacer, indicating when to move to a different quadrant of your mouth.
5. **Recharge:** After use, return the toothbrush to its charging dock on the base unit.

# On/Off water control on handle

Conveniently start  
and stop water flow



Image: A close-up view of the Waterpik electric toothbrush handle, clearly displaying the power button, mode selection button, and indicator lights for the 'Clean', 'Whiten', and 'Massage' sonic brushing modes.

## Water Flosser Usage

1. **Adjust Pressure:** Use the pressure control dial on the base unit to select your desired water pressure setting (1-10). Start with a lower setting and gradually increase as you become accustomed to the flosser.
2. **Position Tip:** Lean over the sink and place the water flossing tip in your mouth. Partially close your lips to prevent splashing, allowing water to flow from your mouth into the sink.
3. **Activate Water Flow:** Turn the unit on using the power switch on the base. Then, slide the ON/OFF switch on the flosser handle to the 'ON' position to start the water flow.
4. **Floss:** Aim the tip at the gumline at a 90-degree angle. Move along the gumline, pausing briefly between each tooth. Ensure you clean both the front and back surfaces of your teeth. The tip rotates 360 degrees for easy access to all areas.
5. **Turn Off:** When finished, slide the ON/OFF switch on the handle to 'OFF', then turn off the unit at the base.



# High volume reservoir

**90 Seconds of water capacity**



Image: A hand holding the Waterpik water flosser handle, showing the convenient ON/OFF water control switch located directly on the handle for easy operation.

Video: An official Waterpik video demonstrating the features and benefits of the Waterpik Complete Care system, including both the water flosser and electric toothbrush in use.

Video: A seller's perspective on the Waterpik Electric Toothbrush, showcasing its design and functionality as part of the Complete Care system.

Video: A demonstration of the powerful flossing action of the Waterpik water flosser, highlighting its effectiveness in cleaning between teeth.

## MAINTENANCE

### Cleaning the Unit

- Empty any remaining water from the reservoir after each use.
- Wipe down the exterior of the unit with a soft, damp cloth.

- Periodically clean the reservoir by hand with warm, soapy water or place it on the top rack of a dishwasher.

### Replacing Brush Heads and Tips

- Replace toothbrush heads every 3 months, or sooner if bristles are frayed.
- Replace water flossing tips every 3-6 months, depending on usage and tip type.

# 5 Water flossing tips included



Image: Two Waterpik electric toothbrush heads, featuring an advanced design with high and low bristle patterns, ready for attachment to the toothbrush handle.



Image: A selection of five different Waterpik water flossing tips, each designed for specific dental needs, such as plaque removal, orthodontic care, or periodontal pockets.

## TROUBLESHOOTING

### Toothbrush Not Turning On/Holding Charge

If your toothbrush is not turning on or appears to not be holding a charge, you may need to perform a battery reset procedure:

1. Unplug the flosser base from the electrical outlet.
2. Locate the mode button on the toothbrush (beneath the power button).
3. Press and hold down the mode button for a full 30 seconds.
4. Plug the flosser unit back into the wall outlet and place the toothbrush in its charging dock.
5. Allow the toothbrush to charge for a full, uninterrupted 24 hours.



If the issue persists after following these steps, please contact Waterpik customer support.

## SPECIFICATIONS

Feature	Detail
Product Dimensions	6.25 x 7.25 x 8.75 inches
Item Weight	2.47 Pounds
Model Number	CC-01
Power Source	Battery Powered (Toothbrush), Corded Electric (Base Unit)
Voltage/Frequency	120VAC/60Hz (for North America only)
Special Features	Multiple Pressure Settings (10), Multiple Tips (5), Rechargeable, Rotatable Tip, Timer (2-minute with 30-second pacer)
Product Benefits	Complete Oral Care
Unit Count	1 Count

## WARRANTY AND SUPPORT

The Waterpik Complete Care 9.0 is backed by a 2-year limited manufacturer’s warranty. For detailed warranty information, please refer to the manual included in your package. Waterpik’s Colorado-based support team is available to assist with any product questions or needs.