

Manuals+

[Q & A](#) | [Deep Search](#) | [Upload](#)

manuals.plus /

› [Silk'n](#) /

› [Silk'n Flash&Go Pro - At Home Permanent Hair Removal System User Manual](#)

Silk'n 07290016627121

Silk'n Flash&Go Pro - At Home Permanent Hair Removal System

Brand: Silk'n | Model: 07290016627121

INTRODUCTION

The Silk'n Flash&Go Pro is a hair removal device designed by physicians for professional results in the comfort of your home. Enjoy smooth skin using light pulsing technology that disables hair growth at the root.

Benefits:

- Convenient, at-home hair removal
- Silky smooth skin
- Quick, painless and safe

PRODUCT OVERVIEW



The Silk'n Flash&Go Pro device, a compact and ergonomic design for comfortable handling.



Key components of the device: Power Button, Skin Tone Sensor, and Flash Window.



Diagram illustrating the hair reduction treatment cycle over multiple sessions.



HPL™ technology targets the pigment in hair follicles to disable hair growth without damaging the skin.

SETUP AND FIRST USE

Watch this official video from Silk'n for a comprehensive guide on setting up and using your Flash&Go Pro device:

Official Silk'n Flash&Go Pro product video demonstrating setup and usage.

Quick Start Guide:

1. Remove the Silk'n Flash&Go™ HAIR REMOVAL UNIT and other components from the packaging.
2. Verify that the DISPOSABLE LAMP CARTRIDGE is inserted properly.
3. Plug the ADAPTOR cord into the Silk'n Flash&Go™ HAIR REMOVAL UNIT.
4. Plug the ADAPTOR into an electrical outlet. The STAND BY INDICATOR on your Silk'n Flash&Go™ is now ready to start.

5. Choose areas to treat that have not had recent exposure to the sun.
6. The skin of the treatment areas should be shaved, clean, dry and free of deodorants. Do not wax, pluck or tweeze the hair out. Shaving is an important step to get the results you want!
7. Press the ON/STBY SWITCH. The POWER ON INDICATOR (green) will turn on.
8. Approximately 3.5 seconds after pressing the ON/STBY SWITCH, the fan sound (similar to the sound of a hairdryer) will start. The device is then ready for you to trigger the energy level.
9. If this is your first treatment, choose the energy level. Flash&Go has 5 levels (1 the lowest to 5 the highest). To choose the right energy level for your treatment area separately.

OPERATING THE DEVICE

The Silk'n Flash&Go Pro utilizes advanced technology for effective hair removal. The device features a skin tone sensor and a movement sensor to ensure safe and efficient operation.

- **Power Levels:** The device offers 5 adjustable energy levels. Start with a lower level and gradually increase as your skin tolerates.
- **Skin Tone Sensor:** The integrated skin tone sensor automatically detects your skin complexion. If the skin is too dark for a safe treatment at the selected energy level, the device will adjust the power down or prevent flashing.
- **Touch&Glide Technology:** This feature allows for continuous flashing as you glide the device over your skin, making treatments quicker and more convenient.
- **Movement Sensor:** The device will only flash when it is in full contact with the skin and moving, ensuring precise application.




Applying the Silk'n Flash&Go Pro to the skin for treatment.

TREATMENT SCHEDULE

For optimal results, follow the recommended treatment schedule:

- **Treatments 1-4:** Plan two weeks apart.
- **Treatments 5-7:** Plan four weeks apart.
- **Treatments 8+:** Treat as needed, until desired results are achieved.

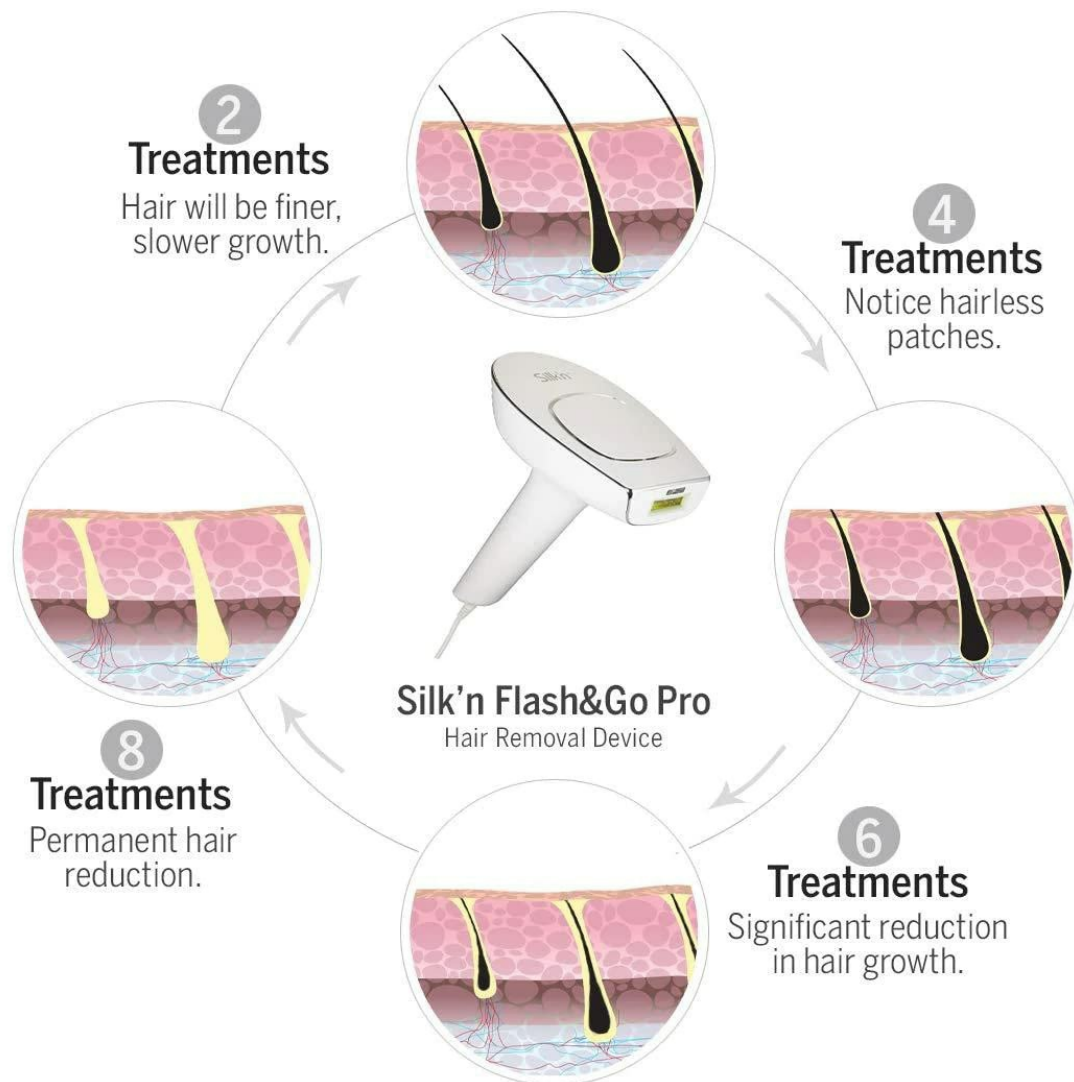

* Based on bi-weekly treatment schedule.

Consistent application according to the schedule is key for effective hair reduction.

SAFETY GUIDELINES AND CONTRAINDICATIONS

It is crucial to read and understand all safety information before using the Silk'n Flash&Go Pro. This device is not designed for everyone. Please consult a doctor or dermatologist if you are unsure.

Warnings:

- **Skin Tone:** Do not use Silk'n Flash&Go on naturally dark skin complexion. The device focuses on pigments, and darker skin tones may experience adverse effects such as redness and discomfort. Refer to the color skin chart in the inner cover of the manual for compatibility.
- **Sun Exposure:** Do not use on tanned skin or after recent sun exposure (including artificial tanning machines). Do not expose treated areas to the sun.
- **Eye Safety:** Never use Silk'n Flash&Go around or near the eyes.
- **Skin Conditions:** Do not treat on tattoos or permanent makeup, dark brown or black spots (such as large freckles, birth marks, moles or warts), nipples, genitals or lips.

- **Medical Conditions:** Do not treat on active eczema, psoriasis, lesions, open wounds or infections (cold sores), or abnormal skin conditions caused by systemic or metabolic diseases (e.g., diabetes). Light-based treatments can cause adverse reactions on already sensitive areas.



Illustration of body areas generally suitable for treatment. Always refer to safety guidelines.

MAINTENANCE

To ensure the longevity and optimal performance of your Silk'n Flash&Go Pro, regular maintenance is recommended:

- **Cleaning:** After each use, ensure the device is unplugged and cool. Gently wipe the flash window and the surrounding area with a soft, dry cloth. Do not use abrasive cleaners or harsh chemicals.
- **Storage:** Store the device in its original packaging or a protective case in a cool, dry place away from direct sunlight and moisture.
- **Lamp Cartridge:** The Flash&Go Pro is designed for lifetime flashes, meaning no refill cartridges are needed. However, ensure the lamp cartridge area is kept clean.

TROUBLESHOOTING

If you encounter any issues with your Silk'n Flash&Go Pro, refer to the following common troubleshooting tips:

- **Device Not Flashing:** Ensure the device is properly plugged into a working electrical outlet. Verify that the skin tone sensor is in full contact with the skin and that your skin tone is compatible with the device's safety parameters. The device will not flash if skin contact is not detected or if the skin is too dark.
- **Blinking Indicator Lights:** If certain indicator lights are blinking, it may signal a specific issue. Refer to the detailed manual for specific blinking patterns and their corresponding solutions. For example, a blinking light might indicate the need for a flash replacement (though this model is designed for lifetime flashes, older models or specific issues might trigger such warnings).
- **Device Getting Hot:** It is normal for the device to warm up during use. However, if it becomes excessively hot or uncomfortable, discontinue use and allow it to cool down. Ensure proper ventilation around the device during operation.




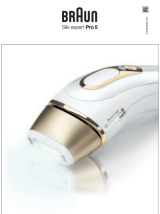


SPECIFICATIONS

Feature	Detail
Product Dimensions	4.88 x 8.13 x 8.13 inches; 11.48 ounces
Item Model Number	07290016627121
Manufacturer	Silk'n
ASIN	B07CBPHC8P
Best Sellers Rank	See Top 100 in Beauty & Personal Care

WARRANTY AND SUPPORT

Information regarding product warranty and customer support was not available in the provided data. Please refer to the product packaging or the official Silk'n website for detailed warranty terms and contact information for support.

Related Documents - 07290016627121

 <p>Silk'n Motion Premium H3220</p> <ul style="list-style-type: none"> 1. SAFETY PRECAUTIONS 2. BEFORE YOU START TREATING 3. HOW TO USE THE DEVICE 4. TREATMENT PLAN 5. TREATMENT RESULTS 6. TROUBLESHOOTING 7. CONTACT INFORMATION 8. WARRANTY 9. INDEX 	<p>Silk'n Motion Premium H3220 User Manual - Permanent Hair Reduction Device</p> <p>Comprehensive user manual for the Silk'n Motion Premium H3220, a home-use IPL device for permanent hair reduction. Covers safety, operation, treatment plans, and troubleshooting.</p>
	<p>SILK PRO IPL Hair Removal System User Manual and Guide</p> <p>Comprehensive user manual for the SILK PRO IPL Hair Removal System by LUXE. Learn about safety guidelines, operating instructions, technical parameters, and hair management for effective and safe hair removal.</p>
 <p>Silk'n Motion Premium H3220</p> <ul style="list-style-type: none"> 1. SAFETY PRECAUTIONS 2. BEFORE YOU START TREATING 3. HOW TO USE THE DEVICE 4. TREATMENT PLAN 5. TREATMENT RESULTS 6. TROUBLESHOOTING 7. CONTACT INFORMATION 8. WARRANTY 9. INDEX 	<p>Silk'n Motion Premium H3220 User Manual - Guide to Home Hair Removal</p> <p>Discover how to use the Silk'n Motion Premium (H3220) with this comprehensive user manual. Learn about safe operation, treatment plans, safety features, and troubleshooting for effective at-home IPL hair removal.</p>
 <p>BRAUN Silk-expert Pro 5</p>	<p>Braun Silk-expert Pro 5 IPL Hair Removal Device User Manual</p> <p>Comprehensive user manual for the Braun Silk-expert Pro 5 Intense Pulsed Light (IPL) hair removal device. Learn about its features, how to use it safely, contraindications, troubleshooting, and maintenance for effective and comfortable home hair removal.</p>
 <p>Silk'n Flash&Go Express</p>	<p>Silk'n Flash&Go Express Hair Removal Device FAQs</p> <p>Frequently asked questions about the Silk'n Flash&Go Express hair removal device, covering its effectiveness, treatment sessions, usage, results, and suitability for different skin types and conditions.</p>
 <p>Silk'n Infinity Fast H3104</p> <ul style="list-style-type: none"> 1. SAFETY PRECAUTIONS 2. BEFORE YOU START TREATING 3. HOW TO USE THE DEVICE 4. TREATMENT PLAN 5. TREATMENT RESULTS 6. TROUBLESHOOTING 7. CONTACT INFORMATION 8. WARRANTY 9. INDEX 	<p>Silk'n Infinity Fast H3104 User Manual: Your Guide to At-Home Hair Removal</p> <p>Discover the Silk'n Infinity Fast H3104 user manual for comprehensive instructions on safe and effective at-home hair removal using eHPL™ technology. Learn about device features, usage, safety, and troubleshooting.</p>