

VOYOR BR610

VOYOR BR610 5-in-1 Blackhead Remover Vacuum

Instruction Manual

INTRODUCTION

Thank you for choosing the VOYOR BR610 5-in-1 Blackhead Remover Vacuum. This device is designed to provide deep pore cleansing, effectively removing blackheads, oil, impurities, and acne, while also exfoliating dead skin for a more radiant complexion. Please read this manual thoroughly before use to ensure safe and optimal performance.

SAFETY INFORMATION

Important Safety Precautions:

- When using the product, slowly move the suction head in a circular or linear direction on your skin.
- Ensure the device does not stay on any single skin area for more than 2 seconds to avoid bruising or irritation.
- Actual results may vary due to different skin conditions.
- We recommend not using this device more than twice a week.
- Do not use on irritated, damaged, or broken skin.
- Keep out of reach of children.

Materials: The device is constructed from durable ABS plastic.

PRODUCT OVERVIEW

The VOYOR BR610 is a versatile skin care tool featuring innovative vacuum absorption technology and an intuitive LED display. It comes with five interchangeable functional heads to cater to various skin needs.

L'OUTIL INDISPENSABLE POUR LES POINTS NOIRS



Réduire
Les Pores



Enlève Les Points
Noirs Et La Graisse



5 Têtes
Fonctionnelles



Forte Puissance
D'aspiration



Indicateur
De Charge De
La Batterie



5 Niveaux
De Puissance



Écran LCD



Rechargeable



Image: The VOYOR BR610 device showcasing its design and key features like 5 heads, 5 suction levels, battery indicator, and LCD screen.

Key Features:

- **Deep Pore Cleansing:** Effectively removes blackheads, oil, impurities, grease, and acne from clogged pores. It also exfoliates dead skin, making your skin more radiant.
- **5 Adjustable Suction Levels:** Easily adjust the most appropriate suction levels for different skin areas. Suitable for sensitive, dry, and oily skin with stubborn blackheads.
- **5 Functional Heads Included:**
 - **Large Circular Hole Head:** For strong suction, ideal for removing stubborn blackheads and oil.
 - **Small Circular Hole Head:** With gentle suction for sensitive skin areas.
 - **Oval Hole Head:** Designed for better blackhead removal from the nose and around the mouth.
 - **Microcrystalline Head:** For dead skin exfoliation and improving skin texture.
 - **Silicone Head:** For gentle suction, suitable for delicate skin.
- **LED Display:** Elegant and easy-to-operate screen shows suction level and battery status.
- **Rechargeable and Portable:** Built-in 1200mAh USB rechargeable lithium battery provides long-lasting use. Lightweight (140 grams), compact, and ideal for travel.

5 Têtes Fonctionnelles Pratiques, Ne Souffrez Plus Des Problèmes De Peau



Image: Illustration of the five interchangeable functional heads: small circular, silicone, oval, microcrystalline, and large circular.

SETUP

1. **Charge the Device:** Before first use, fully charge the device using the provided USB cable. The LED display will indicate charging status. A full charge takes approximately 1.5 hours.
2. **Attach a Suction Head:** Select the desired functional head and gently twist it onto the device until securely attached.
3. **Perform a Test:** Before applying to your face, test the suction on your arm to familiarize yourself with the different suction levels and ensure no adverse reaction. Start with the lowest suction level (Level 1).



Image: The VOYOR BR610 connected to a USB charger, illustrating the charging process and duration.

OPERATING INSTRUCTIONS

1. Skin Preparation (Crucial Step)


Before using the product, it is essential to open your pores to maximize effectiveness and minimize skin irritation. Apply a warm cloth, hot towel, or blackhead-removing gel/cream to your skin for 5 minutes. A warm shower can also help open pores.



Image: Visual guide for skin preparation, including cleansing and steaming to open pores.

2. Device Operation

1. Press the power button to turn on the device.
2. Select your desired suction level using the '+' and '-' buttons.
 - **Levels 1-3:** Suitable for sensitive and dry skin.
 - **Levels 2-4:** Suitable for normal skin.
 - **Levels 3-5:** Suitable for combination and oily skin.
3. Place the suction head gently on your skin.
4. Move the device slowly and continuously across the skin in a circular or linear motion.
5. Do not keep the device stationary on one spot for more than 2 seconds to prevent redness or bruising.

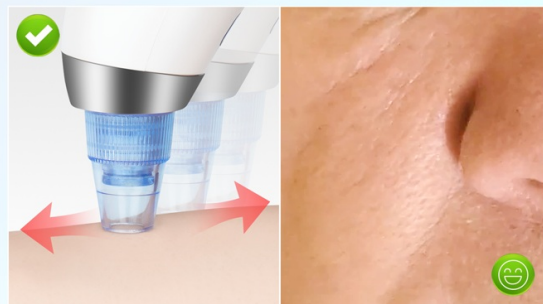

Déplacez Doucement L'astuce D'aspiration Sur La Peau

Image: Correct usage technique, emphasizing slow, continuous movement.



Image: Incorrect usage technique, showing the device held stationary, which can cause skin marks.

3. Post-Treatment Care

After using the device, rinse your face with cold water to close your pores. Apply a soothing and moisturizing cream or mask to help your skin recover and reduce any temporary redness.

MAINTENANCE

1. **Clean the Suction Heads:** After each use, remove the suction head and rinse it thoroughly under running water. You can also use a cotton swab with alcohol for deeper cleaning. Ensure the heads are completely dry before reattaching or storing.
2. **Clean the Device Body:** Wipe the main body of the device with a soft, damp cloth. Do not immerse the device in water.
3. **Storage:** Store the device and its accessories in a cool, dry place, away from direct sunlight and moisture.



Image: Cleaning a suction head under running water.

TROUBLESHOOTING

Problem	Possible Cause	Solution
---------	----------------	----------

Problem	Possible Cause	Solution
No suction or weak suction	Incorrectly attached head; clogged head; low battery; suction level too low.	Ensure the suction head is securely attached. Clean the suction head. Charge the device. Increase the suction level.
Redness or bruising after use	Device held stationary for too long; suction level too high for skin type; skin not properly prepared.	Move the device continuously. Reduce suction level. Ensure pores are open before use. Apply soothing cream after use.
Device does not turn on	Low battery; device malfunction.	Charge the device fully. If the problem persists after charging, contact customer support.

SPECIFICATIONS

- **Model:** BR610
- **Brand:** VOYOR
- **Material:** ABS Plastic
- **Color:** White
- **Item Weight:** 300 Grams
- **Battery:** 1200mAh Rechargeable Lithium Battery
- **Charging Time:** Approximately 1.5 hours
- **Suction Levels:** 5 adjustable levels
- **Included Heads:** 5 functional heads (Large Circular, Small Circular, Oval, Microcrystalline, Silicone)
- **UPC:** 792855824398
- **ASIN:** B07D4F59JD

WARRANTY AND SUPPORT

VOYOR products are designed for quality and performance. For any questions, concerns, or warranty claims, please refer to the contact information provided with your purchase or visit the official VOYOR website. Please retain your proof of purchase for warranty purposes.