

Colgate CN06919A

Colgate 360 Advanced Sonic Battery Power Toothbrush Instruction Manual

Model: CN06919A

INTRODUCTION

This manual provides essential instructions for the safe and effective use, maintenance, and troubleshooting of your Colgate 360 Advanced Sonic Battery Power Toothbrush. Please read this manual thoroughly before first use and retain it for future reference.





Figure 1: Colgate 360 Advanced Sonic Battery Power Toothbrush packaging.

SAFETY INFORMATION

WARNING:

- Adult supervision is required for persons with disabilities when using this product.
- Replace the brush head refill if it becomes loose due to damage.
- Battery choke hazard. Do not install batteries backwards, charge, put in fire, or mix with other battery types.
- Batteries may explode or leak, causing injury.
- Replace all batteries at the same time.

PACKAGE CONTENTS

- 1 x Colgate 360 Advanced Sonic Battery Power Toothbrush Handle
- 1 x Replaceable Brush Head
- 2 x AAA Batteries (included)

SETUP

Battery Installation

1. Firmly hold the bottom of the toothbrush handle.
2. Rotate the top section of the handle counter-clockwise to open the battery compartment.
3. Insert 2 new AAA batteries, ensuring correct polarity as indicated inside the brush handle.
4. Firmly rotate the top section clockwise back into position to close the battery compartment securely.



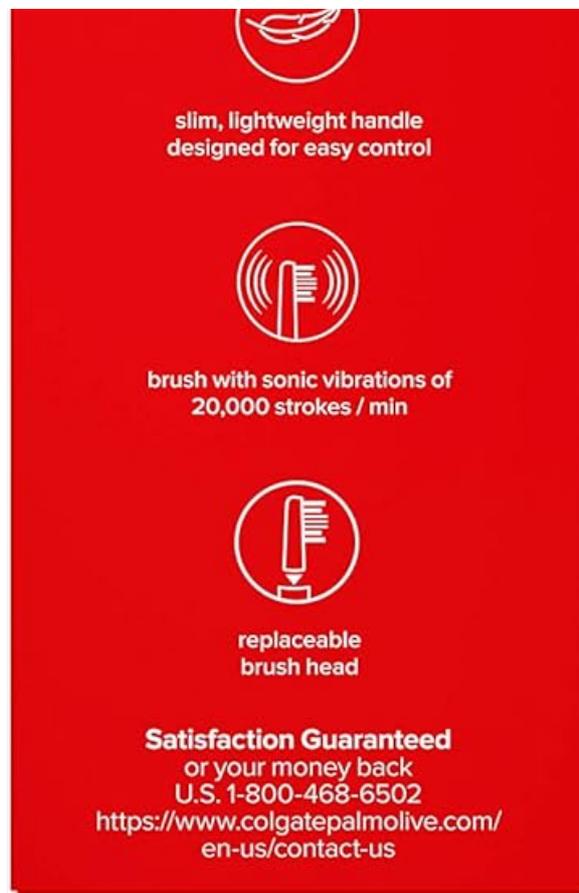


Figure 2: Battery installation instructions on packaging.

OPERATING INSTRUCTIONS

Brushing Your Teeth

1. Attach a brush head to the handle.
2. Wet the bristles and apply a small amount of toothpaste.
3. Place the brush head against your teeth at a slight angle towards the gumline.
4. Press the power button (red button) to turn on the toothbrush. The toothbrush will provide sonic vibrations at 20,000 strokes per minute.
5. Move the brush head slowly from tooth to tooth, allowing the bristles to clean each surface. Do not scrub aggressively.
6. The toothbrush includes a 2-minute timer to assist with recommended brushing duration. Continue brushing until the timer indicates completion.
7. After brushing, press the power button again to turn off the toothbrush.



Sonic Vibrations

20,000 strokes/min

2 Minute Timer

On/off button

Slim Design

Lightweight handle
for easy control

Figure 3: Toothbrush features including sonic vibrations and 2-minute timer.



Figure 4: Illustration of the 2-minute timer during brushing.

Using the Cheek & Tongue Cleaner

The back of the brush head features a cheek and tongue cleaner designed to remove odor-causing bacteria. After brushing your teeth, gently glide the textured surface over your cheeks and tongue.



Replaceable Brush Head with Cheek & Tongue Cleaner

Figure 5: Replaceable brush head with integrated cheek and tongue cleaner.

MAINTENANCE

Cleaning the Toothbrush

After each use, rinse the toothbrush head thoroughly under running water. Shake off any excess water and store the toothbrush upright to allow it to air dry.

Replacing the Brush Head

Dentists and hygienists recommend replacing your toothbrush refill every 3 months, or sooner if the brush head becomes damaged or worn.

1. Gently pull the old brush head directly away from the handle.
2. Press a new Colgate 360 Advanced brush head firmly into place on the handle until it clicks securely.

Replacing Batteries

Refer to the "Battery Installation" section under "Setup" for detailed instructions on replacing the AAA batteries. Always replace both batteries at the same time.

TROUBLESHOOTING

Problem	Possible Cause	Solution
Toothbrush does not turn on.	Batteries are depleted or incorrectly installed.	Replace batteries with new AAA batteries, ensuring correct polarity. Refer to "Battery Installation" section.
Brush head feels loose.	Brush head is not securely attached or is damaged.	Ensure the brush head is firmly pressed onto the handle. If damaged, replace the brush head.
Reduced vibration intensity.	Batteries are low.	Replace batteries with new AAA batteries.

SPECIFICATIONS

Model Number: CN06919A

Power Source: 2 x AAA Batteries (included)

Vibration Speed: 20,000 strokes per minute

Special Features: 2-minute timer, Lightweight handle, Replaceable brush head, Cheek & Tongue Cleaner

Item Firmness: Soft

Product Dimensions: 2.36 x 1.38 x 8.66 inches

Weight: 3.21 ounces

Manufacturer: Colgate-Palmolive Co

WARRANTY AND SUPPORT

Colgate offers a satisfaction guarantee. If you are not satisfied with your product, you may be eligible for a refund or replacement.

For further assistance, product support, or warranty claims, please contact Colgate-Palmolive directly:


Website: colgatepalmolive.com/us/contact-us

Phone: 1-800-468-6502

© 2023 Colgate-Palmolive Co. All rights reserved.


Related Documents - CN06919A

<p>How to Use</p> <ul style="list-style-type: none"> 1. Check your batteries <ul style="list-style-type: none"> • Check that you have 2 AAA batteries. • Check that the batteries are inserted correctly. • Check that the batteries are inserted with the correct polarity. • Make sure the batteries are inserted into the correct slots. 2. Turn on the toothbrush <ul style="list-style-type: none"> • Press the power button. • The toothbrush will turn on and the LED indicator will light up. 3. Brush your teeth <ul style="list-style-type: none"> • Brush for 2 minutes. • Use the Colgate app to track your brushing time. • The Colgate app will provide you with personalized brushing analysis and rewards. 4. Recharge the toothbrush <ul style="list-style-type: none"> • Charge the toothbrush in the charging case. • The charging case will charge the toothbrush when it is closed. • The LED indicator will light up when the toothbrush is fully charged. 5. Replace the brush head <ul style="list-style-type: none"> • Replace the brush head when it is worn. • The Colgate app will provide you with personalized brushing analysis and rewards. 	<p>How to Use Your Hum by Colgate Smart Toothbrush</p> <p>A comprehensive guide on how to use the Hum by Colgate smart toothbrush, including charging, setup, brushing modes, and cleaning instructions. Learn how to connect to the Colgate app for personalized brushing analysis and rewards.</p>
---	--

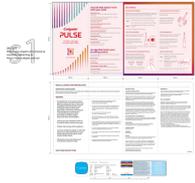


[Colgate Pulse Electric Toothbrush User Manual and Guide](#)

Comprehensive guide for the Colgate Pulse electric toothbrush, covering setup, usage, features, maintenance, warranty, and safety information.

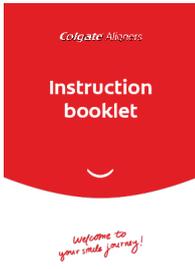

[Colgate Aligners Instruction Booklet: Your Guide to a Straighter Smile](#)

Comprehensive instruction booklet for Colgate Aligners, guiding you through the process of achieving a straighter smile with clear steps, care instructions, and important precautions.

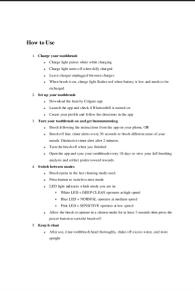


[How to Use Your Hum by Colgate Smart Toothbrush](#)

A comprehensive guide on how to use the Hum by Colgate smart toothbrush, including charging, setup, brushing modes, and cleaning instructions. Learn how to connect to the Colgate app for personalized brushing analysis and rewards.



[Colgate Pulse Electric Toothbrush User Manual and Guide](#)

Comprehensive guide for the Colgate Pulse electric toothbrush, covering setup, usage, features, maintenance, warranty, and safety information.

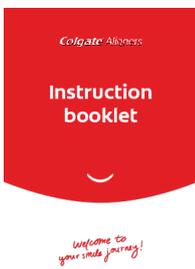

[Colgate Aligners Instruction Booklet: Your Guide to a Straighter Smile](#)

Comprehensive instruction booklet for Colgate Aligners, guiding you through the process of achieving a straighter smile with clear steps, care instructions, and important precautions.