

VOYOR MS110

VOYOR Neck Massager MS110 Instruction Manual

Model: MS110 | Brand: VOYOR

INTRODUCTION

The VOYOR Neck Massager MS110 is a manual shiatsu deep tissue massager designed for pain relief in the neck, shoulders, and other body areas. Its ergonomic handle and dual trigger points allow for controllable massage intensity, helping to relieve muscle pain and fatigue.

PRODUCT FEATURES

- **Controllable Massage Intensity:** Two soft yet firm silicone balls target pressure points, allowing users to apply desired strength for various massage needs.
- **Muscle Pain & Fatigue Relief:** Designed to press, squeeze, and knead tight, stiff, and aching muscles in the neck, shoulders, and legs, improving blood circulation and promoting relaxation.
- **Mimics Human Massage:** Dual trigger point massager balls provide a deep penetrating massage sensation, similar to human hands, for instant pain relief.
- **Lightweight & Portable:** Weighing only 310g, its portable design allows for personal deep tissue massage at home or during travel.
- **Ergonomic & Flexible Design:** Features an anti-skid ergonomically designed handle for comfortable grip. The flexible handle can be adjusted to suit different body shapes.

SAFETY INFORMATION

WARNING:

- Keep out of reach of children.
- Avoid contact with skin and eyes.
- If ingested, seek medical help immediately.
- Do not use near water or when wet.
- Do not apply excessive pressure to avoid injury.
- Always follow manufacturer instructions for use.

SETUP

The VOYOR Neck Massager MS110 comes with detachable massage balls for easy cleaning and customization. To set up your massager:

1. Unpack all components: the main massager frame and the massage balls.
2. Align the massage balls with the designated slots on the massager frame. The balls are designed to fit securely.
3. Gently push the massage balls into place until they click or feel firmly seated. Ensure they are properly installed before use.
4. The massage balls can be installed at different angles according to your needs. They are designed to be fixed and not rotate during massage, providing targeted pressure.

Your browser does not support the video tag.

Video: VOYOR Neck Massager MS110 assembly and operation. This video demonstrates how to install the massage balls and basic usage.

OPERATING INSTRUCTIONS

The VOYOR Neck Massager MS110 is a manual device that allows you to control the intensity and target specific areas. Follow these steps for effective use:

1. **Positioning:** Place the massager around your neck, shoulders, or other desired body part, ensuring the silicone massage balls are aligned with the area you wish to massage.
2. **Grip:** Hold the ergonomic handles firmly with both hands. The anti-skid design ensures a comfortable and secure grip.
3. **Apply Pressure:** Gently squeeze the handles to bring the massage balls closer together, applying pressure to your muscles. The intensity is controlled by how much pressure you apply.
4. **Movement:** Move the massager slowly up and down, or side to side, to target different points. The fixed massage balls work by pressing into the muscle, not by rolling.
5. **Target Trigger Points:** Focus on areas of tension or knots. You can hold the pressure on a specific trigger point for a few seconds to deepen the massage.
6. **Duration:** Use for approximately 5-10 minutes per session, or as needed, to relieve muscle pain and promote relaxation.

The massager is designed to mimic the hands of a masseuse, providing deep tissue relief. Adjust the pressure to your comfort level.

FOR 10 MINUTES A DAY, YOUR SHOULDERS AND NECK BECOME RELAXED

Say goodbye to neck soreness and indulge in muscle massage



Double trigger
massage



Adjustable
intensity



Home & car &
office use



Multiple body
parts use



Curved
design



Simple
operation



Image: Front view of a woman using the VOYOR Neck Massager on her neck, demonstrating the ergonomic grip and positioning.

CAN BE USED ON MULTIPLE BODY PARTS, NOT JUST THE NECK

Suitable for multiple body parts, relax your shoulder anytime, anywhere

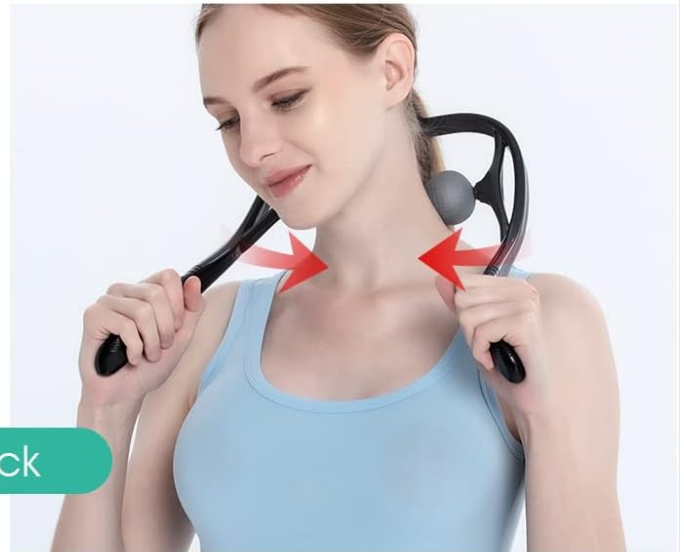


Image: Back view of a woman using the VOYOR Neck Massager, targeting the neck and upper shoulders.



Image: Illustration demonstrating the versatility of the massager for use on shoulders, thighs, and calves.

MAINTENANCE

To ensure the longevity and hygiene of your VOYOR Neck Massager MS110:

- **Cleaning:** The massage balls are removable for easy cleaning. Detach them from the frame and wash with mild soap and water. Rinse thoroughly and allow to air dry completely before reattaching.
- **Storage:** Store the massager in a clean, dry place away from direct sunlight and extreme temperatures.
- **Inspection:** Periodically inspect the massager for any signs of wear or damage. Discontinue use if any damage is observed.

Your browser does not support the video tag.

Video: "What's in the massage kit?" This video shows the components of the massager, including the removable massage balls, which aids in understanding maintenance.

TROUBLESHOOTING

The VOYOR Neck Massager MS110 is a simple, non-electronic device, minimizing potential issues. If you encounter any problems:

- **Insufficient Pressure:** Ensure you are applying enough manual force by squeezing the handles. Adjust your grip and body position to maximize pressure on the desired area.
- **Discomfort:** If you experience discomfort, reduce the pressure or reposition the massager. Avoid using on bony areas or directly on the spine.
- **Hair Entanglement:** If using on the neck, ensure long hair is tied up to prevent entanglement with the massage balls.
- **Product Damage:** If the product appears damaged or broken, discontinue use immediately and contact customer support.

SPECIFICATIONS

Feature	Detail
Brand	VOYOR
Model Number	MS110
Product Dimensions	13.78 x 7.09 x 1.3 inches
Item Weight	10.34 ounces (293 Grams)
Material	Silicone
Power Source	Manual (Battery Powered: No)
Use For	Neck, Shoulders, Back, Legs, Thighs, Calves

WARRANTY AND SUPPORT

For any questions regarding your VOYOR Neck Massager MS110, please contact VOYOR customer support. Specific warranty details may vary and are typically provided with your purchase documentation. You can visit the official VOYOR Store for more information:[VOYOR Store](#)

LEGAL DISCLAIMER

Statements regarding dietary supplements have not been evaluated by the FDA and are not intended to diagnose, treat, cure, or prevent any disease or health condition.

ASK A QUESTION ABOUT THIS MANUAL

Ask about setup, troubleshooting, compatibility, parts, safety, or missing instructions. Manuals+ will review the question and use this page's manual context to help answer it.

Name

Anonymous

Email

you@example.com

After a successful submission, we securely hash the supplied account or form email before sending it to Microsoft Advertising for conversion attribution. [Privacy details](#).

Question

Example: How do I reset this device or fix this error?

Details

Model number, symptoms, what you tried, or the section of the manual you are using.

Submit question