

Project E Beauty PE706

Project E Beauty LED Light Therapy Mask

Model: PE706

1. PRODUCT OVERVIEW

The Project E Beauty Wireless LED Photon Therapy Mask is an advanced skincare device designed for facial and neck rejuvenation. It utilizes 7 different LED light colors to address various skin concerns, including anti-aging, acne treatment, and overall skin health improvement. This lightweight and easy-to-use device offers professional-grade skincare from the comfort of your home.



Image 1.1: The Project E Beauty LED Light Therapy Mask in use, showing the facial and neck components.

2. PACKAGE CONTENTS

Please check the package carefully to ensure all items are present and undamaged:

- 1 x LED Mask for Face and Neck
- 1 x Remote Control
- 1 x Power Cord
- 1 x User Manual

3. PRODUCT FEATURES

The Project E Beauty LED Light Therapy Mask offers a range of features for effective and convenient skin treatment:

- **Professional-Grade Skincare:** Equipped with 99 UV-free LED bulbs emitting 7 colors of light for various skin

benefits.

- **Exclusive Neck LED Panel:** Includes a dedicated neck panel with 40 LED bulbs to treat signs of aging in the neck area simultaneously.
- **7 LED Colors for All Skin Types:** Customizable therapy with specific wavelengths for different skin concerns.
- **Easy Treatment Sessions:** Simple operation with adjustable intensity levels and session durations.
- **Painless and Non-Invasive:** A comfortable cosmetic procedure with no downtime, cuts, or discomfort.
- **Wireless Technology:** Designed for ease of use without being tethered by wires during treatment.
- **Lightweight Design:** Ensures comfort during use.

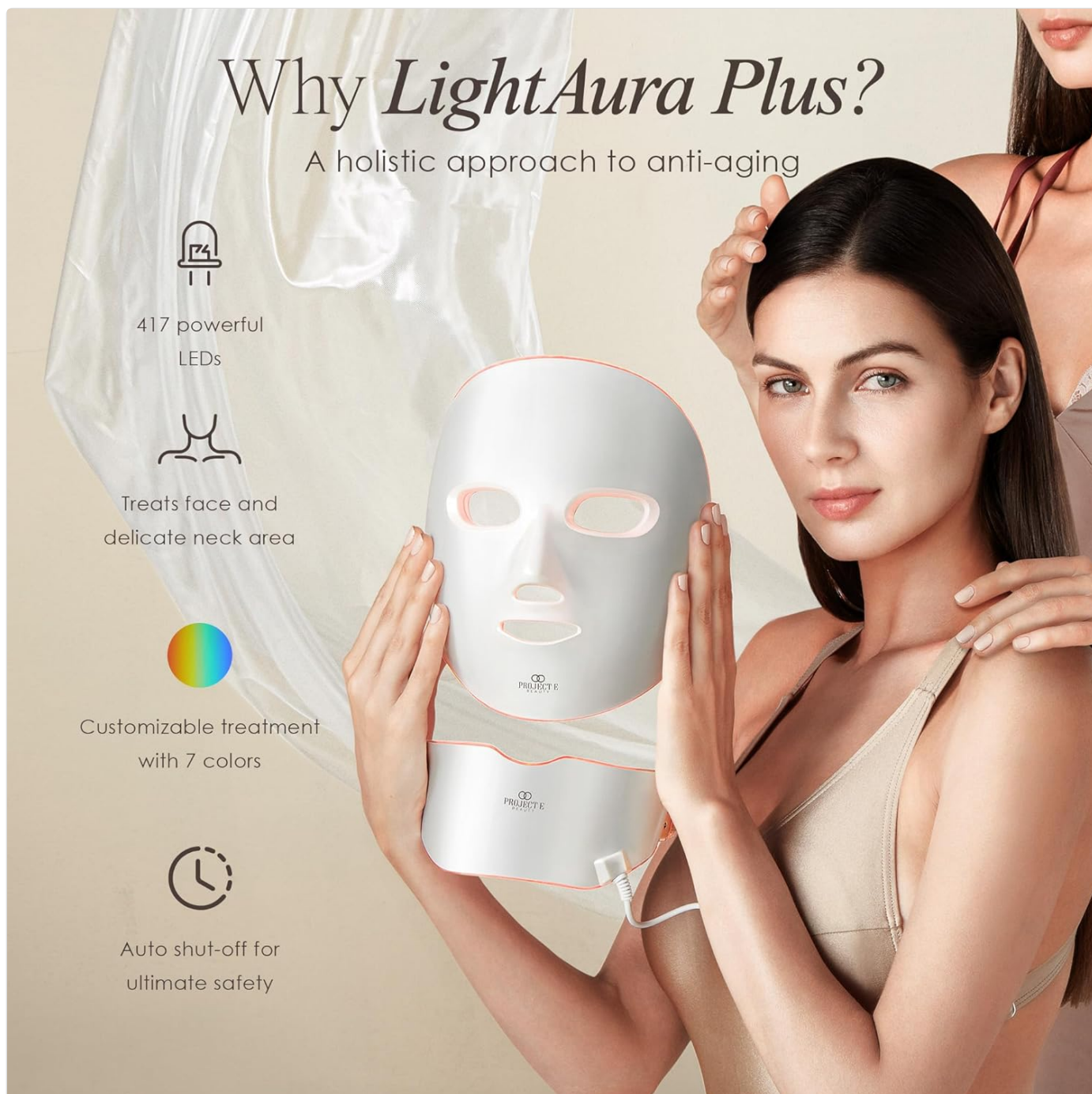


Image 3.1: Key features of the Project E Beauty LED Light Therapy Mask, highlighting its comprehensive design and functionality.

4. SETUP

Follow these steps to prepare your device for its first use:

1. **Unpack:** Carefully remove all components from the packaging.

2. **Inspect:** Check the mask, remote control, and power cord for any visible damage.
3. **Connect:** Connect the power cord to the mask and a suitable power outlet (110V-220V, 50/60Hz).

How to Use LightAura Plus?

A glow like *no other*

STEP 1

Prep the skin



STEP 2

Power on the device



STEP 3

Unwind



STEP 4

Apply AM/PM skincare



Image 4.1: Visual guide for preparing and using the LED Light Therapy Mask.

5. OPERATING INSTRUCTIONS

To operate your Project E Beauty LED Light Therapy Mask, follow these steps:

1. **Prepare Your Skin:** Cleanse your face thoroughly before each treatment session.
2. **Power On:** Press the "Power Switch" button on the remote control to activate the LED mask.
3. **Start Treatment:** Press "Start/Stop" to begin the treatment session.
4. **Adjust Time:** Use the 'Time' buttons on the remote control to set your desired session duration. For beginners, start with 10 minutes and gradually increase to 15 minutes as your skin adapts.

5. **Select Color:** Press the "Color" button to cycle through the 7 available LED light colors. Choose the color based on your specific skin concern.
6. **Adjust Intensity:** Use the "Energy" buttons to adjust the light intensity. Start with the lowest level and gradually increase as comfortable.
7. **Relax:** Wear the mask and relax during the treatment. The device is designed to be comfortable and non-invasive.
8. **Post-Treatment:** After the session, you may apply your regular AM/PM skincare products.

LED Light Color Functions:

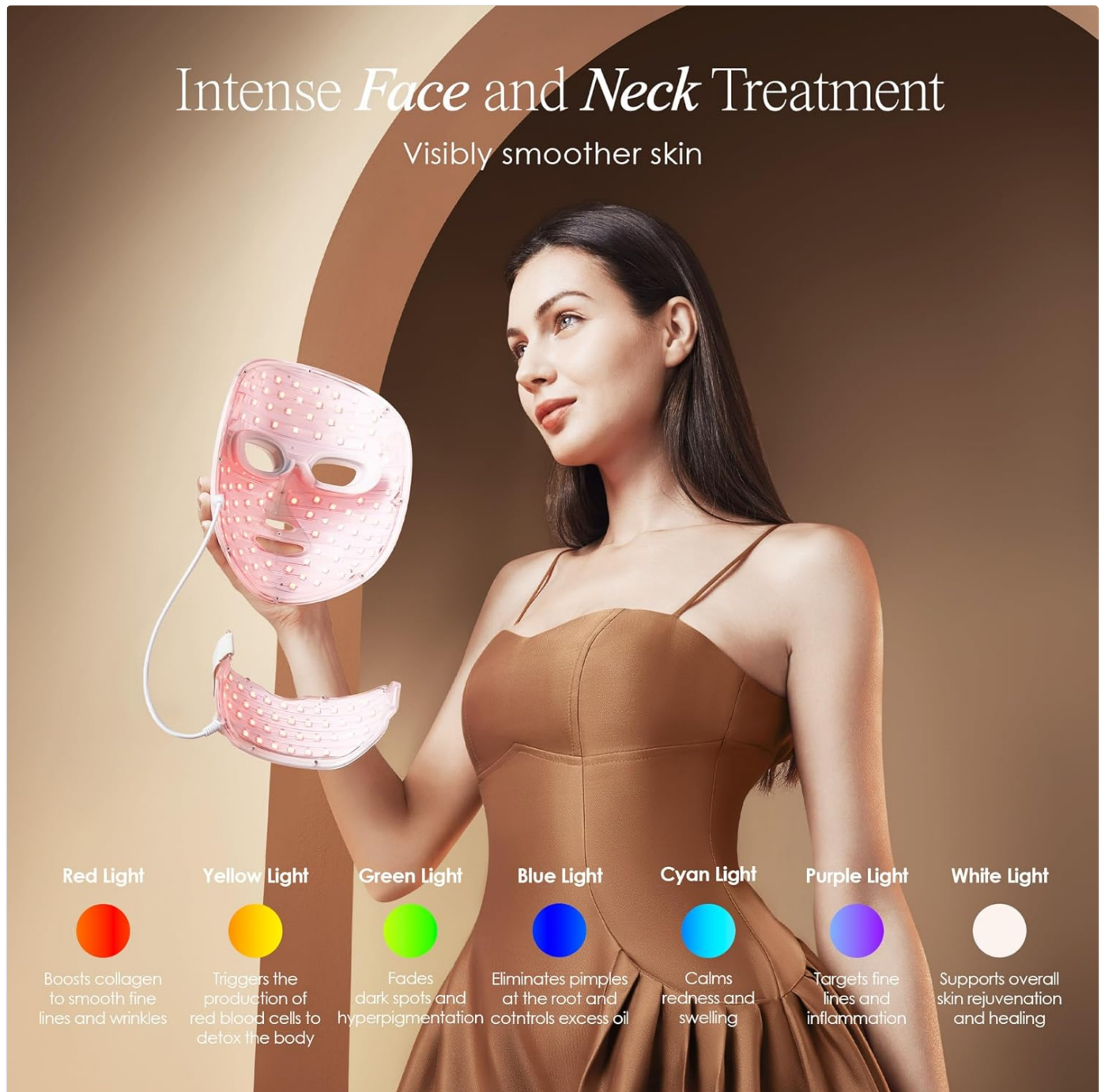


Image 5.1: Overview of the 7 LED light colors and their specific benefits for skin treatment.

- **Red Light (630nm):** Boosts collagen production, smooths fine lines and wrinkles.
- **Yellow Light (630+520nm):** Triggers red blood cell production for detoxification.
- **Green Light (520nm):** Fades dark spots and hyperpigmentation, controls excess oil.
- **Blue Light (460nm):** Eliminates pimples at the root, controls excess oil.
- **Cyan Light (460+520nm):** Calms redness and swelling.

- **Purple Light (460+630nm):** Targets fine lines and inflammation.
- **White Light (460+630+520nm):** Supports overall skin rejuvenation and healing.

6. MAINTENANCE

Proper maintenance ensures the longevity and effectiveness of your LED mask:

- **Cleaning:** After each use, gently wipe the mask surfaces with a soft, dry cloth. For deeper cleaning, use a slightly damp cloth with a mild, non-abrasive cleaner. Ensure the mask is completely dry before storing.
- **Storage:** Store the mask in a cool, dry place, away from direct sunlight and extreme temperatures. Keep it in its original packaging or a protective case to prevent dust accumulation and damage.
- **Avoid Water:** Do not immerse the mask or remote control in water or any other liquids.

7. TROUBLESHOOTING

If you encounter any issues with your Project E Beauty LED Light Therapy Mask, refer to the following common solutions:

- **Device Not Turning On:**
 - Ensure the power cord is securely connected to both the mask and a working power outlet.
 - Verify that the power outlet is functional by testing it with another device.
- **LED Lights Not Illuminating:**
 - Check if the device is powered on and the "Start/Stop" button has been pressed.
 - Ensure the intensity level is not set to zero. Increase the energy level using the "Energy" buttons.
- **Remote Control Not Responding:**
 - Ensure the remote control is pointed directly at the mask's receiver.
 - Check if the remote control's battery needs replacement (if applicable).

If the problem persists after attempting these solutions, please contact Project E Beauty customer support for further assistance.

8. SPECIFICATIONS

Specification	Detail
Brand	Project E Beauty
Model Number	PE706
Color	White
Voltage	110V-220V, 50/60Hz
Facial Mask LEDs	102 LEDs
Neck Mask LEDs	40 LEDs
Treatment Area	Face and Neck
Item Weight	1 Kilogram


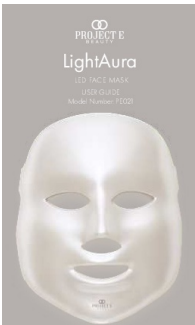
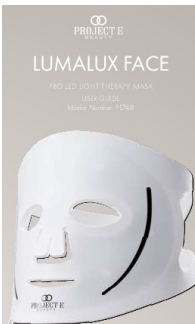
Package Dimensions	32.2 x 19 x 12.7 cm
Compatible Skin Types	Oily, Combination, Normal, Sensitive, Dry
Power Source	1 Lithium-ion battery (requires batteries)
UPC	616576442160




9. WARRANTY AND SUPPORT

Project E Beauty offers a 1-year warranty for this product. To register your product for warranty coverage, please refer to the instructions and QR code provided on the warranty card included in your product packaging. For any questions, technical assistance, or warranty claims, please visit the official Project E Beauty website or contact their customer support directly.

© 2024 Project E Beauty. All rights reserved.

Related Documents - PE706

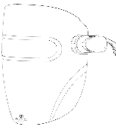	<p>Photon Skin Rejuvenation Face & Neck Mask (PE706) User Manual</p> <p>User manual for the Project E Beauty Photon Skin Rejuvenation Face & Neck Mask (Model No. PE706). Learn about LED photon therapy, its benefits, safety precautions, and how to use the device for optimal skin rejuvenation.</p>
	<p>Project E Beauty LightAura LED Face Mask User Guide</p> <p>Comprehensive user guide for the Project E Beauty LightAura LED Face Mask. Learn about its features, benefits, safety instructions, how to use, and maintenance for optimal skincare results.</p>
	<p>Project E Beauty Lumalux Face Pro LED Light Therapy Mask User Guide</p> <p>Explore the Project E Beauty Lumalux Face Pro LED Light Therapy Mask (Model PE748). This user guide details its advanced LED light therapy features for skin rejuvenation, targeting aging, acne, and more. Learn about benefits, usage, and safety.</p>

<div data-bbox="118 91 309 412"><p>Project E Beauty LED Light Therapy Mask (Model No. PE730)</p><p>User Manual</p><p>PROJECT E BEAUTY</p><p>www.projecteauty.com</p></div>	<p>Project E Beauty LED Light Therapy Mask PE730 User Manual</p> <p>User manual for the Project E Beauty LED Light Therapy Mask (Model PE730). Learn about its features, how to use it for wrinkle reduction and acne treatment, safety guidelines, and maintenance.</p>
<div data-bbox="118 660 309 981"><p>PROJECT E BEAUTY</p><p>LUMALUX FACE+</p><p>TWO LED LIGHT THERAPY FACE & NECK MASK</p><p>Model Number PE749</p></div>	<p>Project E Beauty Lumalux Face+ Pro LED Light Therapy Face & Neck Mask User Guide PE749</p> <p>User guide for the Project E Beauty Lumalux Face+ Pro LED Light Therapy Face & Neck Mask (Model PE749). Learn about its functions, benefits, safety instructions, and how to use this advanced skincare device.</p>
<div data-bbox="118 1025 309 1346"><p>Skin Rejuvenation Photon Mask (Model No. PE021)</p><p>User Manual</p><p>PROJECT E BEAUTY</p><p>www.projecteauty.com</p></div>	<p>Project E Beauty Skin Rejuvenation Photon Mask PE021 User Manual</p> <p>Discover the Project E Beauty Skin Rejuvenation Photon Mask (Model PE021). This user manual provides essential information on LED photon therapy, usage instructions, safety precautions, and benefits for skin rejuvenation.</p>