

Epilady EP-720-04

Epilady Epilaser 808 nm Absolute Permanent Laser Stylus User Manual

Model: EP-720-04

Brand: Epilady

1. INTRODUCTION

About Your Epilady Epilaser Absolute

The Epilady Epilaser 808 nm Absolute Permanent Laser Stylus is an over-the-counter device designed for the removal of unwanted facial hair. The Absolute is based on advanced infrared laser technology. It is intended for permanent facial hair reduction after initial treatment. The Absolute allows for permanent and rapid facial hair reduction using professional laser technology with unlimited treatments.

It offers lasting and long-term results with the convenience, affordability, and privacy of self-treatment, eliminating the need for costly and inconvenient appointments. Small, lightweight, and ergonomic, the Absolute is simple and comfortable to use.

Key Features:

- Laser stylus for permanent hair reduction
- Prevents undesirable effects on surrounding skin areas
- 30-minute charging time with USB plug
- Intelligent, painless, compact, and easy to use
- Ideal for small areas, treating the face, chin, neck, bikini line, stomach, fingers, and toes
- Target gender: Unisex



The Epilady Epilaser Absolute device, shown in white and silver, with its protective cap detached, highlighting its sleek design and compact size.

2. IMPORTANT SAFETY INFORMATION

WARNING: Do not use on delicate areas like the eyes, ears, and other private regions.

This device is designed for specific skin and hair types. Please consult the compatibility chart below before use to ensure safe and effective treatment.



A visual guide illustrating the Epilady Epilaser Absolute's compatibility with different skin tones. It indicates suitability for White to Light Brown skin tones (Fitzpatrick Types 1-4) and Dark Brown to Black skin tones (Fitzpatrick Types 5-6). A note indicates future availability for White skin with Blonde hair.


Wondering if **Our Laser Hair Removal** is right for you?

<div style="display: flex; align-items: center;"> <div style="writing-mode: vertical-rl; transform: rotate(180deg);">Skin Color</div> <div style="display: flex; gap: 5px;"> </div> </div>	Type 1	Type 2	Type 3	Type 4	Type 5	Type 6
Black	✓	✓	✓	✓	✗	✗
Dark Brown	✓	✓	✓	✓	✗	✗
Dark Blonde Light Brown	✓	✓	✓	✓	✗	✗
Red/Orange	✗	✗	✗	✗	✗	✗
Light Blonde	✗	✗	✗	✗	✗	✗
Grey	✗	✗	✗	✗	✗	✗
White	✗	✗	✗	✗	✗	✗

Skin Tone & Body Hair Chart | Based on the **Fitzpatrick Scale**

A comprehensive chart detailing the effectiveness of laser hair removal based on the combination of skin tone (Fitzpatrick Scale Types 1-6) and hair color. Green checkmarks indicate compatibility, while red crosses indicate non-compatibility. Generally, the device is effective on darker hair colors with lighter to medium skin tones.

The Epilady Epilaser Absolute uses true laser technology, not Intense Pulsed Light (IPL). Laser technology employs a single powerful wavelength to target a specific chromophore, resulting in a higher energy density than IPL.

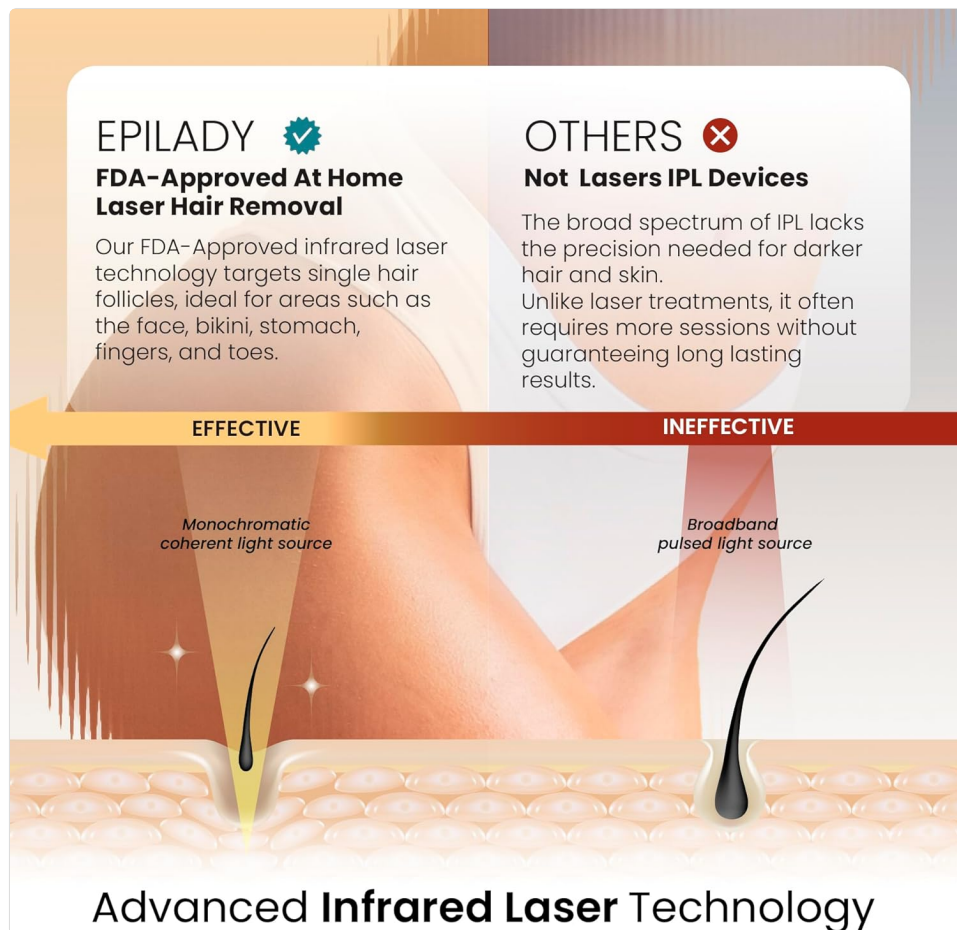

A comparative diagram illustrating the difference between Epilady's FDA-Approved At Home Laser Hair Removal, which uses a monochromatic coherent light source for effective targeting, and other IPL devices that use a broadband pulsed light source, which are often less effective for permanent hair reduction.

3. PACKAGE CONTENTS

Please ensure all items are present in your package:

- Epilady Epilaser 808 nm Absolute Permanent Laser Stylus
- USB Charging Cable
- 1 AA Battery (included)

4. SETUP

Charging the Device

Before first use, or when the battery is low, connect the device to a USB power source using the provided USB charging cable. A full charge takes approximately 30 minutes. Ensure the device is fully charged for optimal performance.

5. OPERATING INSTRUCTIONS

Preparation for Treatment

For best results, shave the targeted area before use. This helps ensure the laser effectively targets hair at the root.

Using the Device

1. Ensure the device is charged and clean.
2. Turn on the device (refer to the power button on the device).
3. Place the tip of the stylus directly over the hair you wish to treat.
4. Gently press the tip down to activate the laser. The treatment time per hair is approximately 3 seconds.
5. Move to the next hair and repeat the process.




The Epilady Epilaser Absolute device in use, demonstrating its application on facial skin. This device is FDA-approved for laser hair removal on areas including the bikini line, feet, underarms, and face.

Recommended Treatment Areas:

The Epilady Epilaser Absolute is ideal for small areas such as the face, chin, neck, bikini line, stomach, fingers, and toes.

Treatment Sessions and Expected Results:

Noticeable results can be observed within 4 weeks of consistent use. Laboratory results show that after two months of use, hair loss rates measure at around 90% reduction.

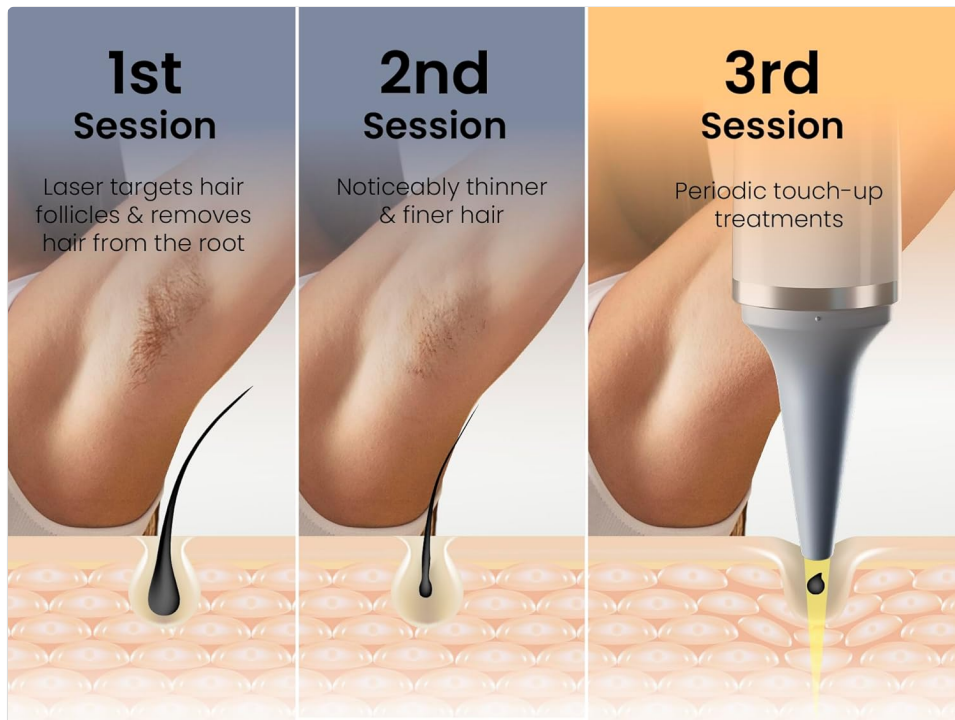

1st Session
Laser targets hair follicles & removes hair from the root

2nd Session
Noticeably thinner & finer hair

3rd Session
Periodic touch-up treatments

Safe, Effective & Quick Laser Hair Removal

Laboratory results show that after two months of use, hair loss rates measure at around 90% reduction.
(NOTE: For enhanced results, continual use and periodic touch-ups with this product are highly beneficial).

A visual representation of the expected results over multiple treatment sessions. The 1st session targets hair follicles, the 2nd session shows noticeably thinner and finer hair, and the 3rd session indicates the need for periodic touch-up treatments for sustained results. Laboratory results suggest up to 90% hair reduction after two months of use.



Approved



Patented Intelligent Skin Detector



Treatment Time



Noticeable Result



Compact & Travel-Ready

Warning:
Do not use on delicate areas like the eyes, ears, and other private regions.



For Best Results: **Shave the targeted area before use.**

This helps ensure laser targets hair at the root.

An overview of key features and usage recommendations for the Epilady Epilaser Absolute. It highlights FDA approval, a patented intelligent skin detector, a 3-second treatment time, noticeable results in 4 weeks, and its compact, travel-ready design. Important warnings include avoiding use on delicate areas like eyes, ears, and private regions. For optimal results, it is recommended to shave the targeted area before use to ensure the laser effectively targets hair at the root.

Next Generation Laser

Hair Removing Technology



- ☀ Safe (Safety Goggles Not Needed)
- ☀ Lightweight
- ☀ Painless
- ☀ Easy-To-Use
- ☀ Patented Technology

An image highlighting the key benefits of the Epilady Epilaser Absolute's next-generation laser hair removal technology. These benefits include being safe (no safety goggles needed), lightweight, painless, easy-to-use, and incorporating patented technology.

6. MAINTENANCE

Cleaning the Device

After each use, gently wipe the tip of the stylus with a clean, dry cloth. Do not use abrasive cleaners or immerse the device in water.

Storage

Store the device in a cool, dry place, away from direct sunlight and out of reach of children. Always replace the protective cap when not in use.

7. TROUBLESHOOTING

Device Not Activating

- Ensure the device is fully charged.
- Confirm the tip is pressed firmly against the skin.

Device Gets Hot During Use

- It is normal for the device to warm up after several activations. If it becomes uncomfortably hot, allow it to cool down before continuing use.

No Noticeable Hair Reduction

- Ensure you are using the device consistently as per the operating instructions.
- Verify your skin and hair type are compatible with the device using the provided charts.
- Ensure the area is shaved before treatment.

8. SPECIFICATIONS

Product Name	Epilady Epilaser 808 nm Absolute Permanent Laser Stylus
Brand	Epilady
Item Model Number	EP-720-04
Battery	1 AA (included)
Product Dimensions (L x W x H)	20.32 x 6.35 x 22.22 cm
Weight	70 grams
Manufacturer	Epilady, Inc.

9. WARRANTY AND SUPPORT

For warranty information, technical support, or any inquiries regarding your Epilady Epilaser Absolute, please contact Epilady, Inc. directly. Refer to the product packaging or the official Epilady website for the most current contact details.