

[Manuals.plus](#) /

› [TOUCHBeauty](#) /

› TOUCHBeauty AG-1767 3-in-1 EMS Facial Massager with LED Light Therapy User Manual

TOUCHBeauty AG-1767

TOUCHBeauty AG-1767 3-in-1 EMS Facial Massager with LED Light Therapy

Model: AG-1767 | Brand: TOUCHBeauty

INTRODUCTION

The TOUCHBeauty AG-1767 is a versatile 3-in-1 facial beauty device designed for comprehensive skin care. It combines Microcurrent technology for facial and eye/lip massage with LED Light Therapy to address various skin concerns, promoting skin tightening, rejuvenation, and overall skin health. This manual provides detailed instructions for safe and effective use of your device.

PRODUCT OVERVIEW

Components Included

- Main Device Body
- Microcurrent Facial Massage Head (with rollers)
- Microcurrent Eye & Lip Massage Head
- LED Light Therapy Head
- USB Charging Cable

Device Features

The AG-1767 features an intuitive LCD display for mode and intensity settings, and interchangeable heads for specialized treatments.



Image: The TOUCHBeauty AG-1767 device with its three interchangeable heads: Microcurrent Facial Massage Head, Microcurrent Eye & Lip Massage Head, and LED Light Therapy Head.



Image: Detailed diagram illustrating the dimensions of the TOUCHBeauty AG-1767 device and labeling its key components including the Light Therapy Head, Microcurrent Facial Massage Head, Microcurrent Eye & Lip Massage Head, LCD Screen, Power/Mode/Intensity Switch Button, Type-C Charging Port, and USB Charging Cable.

SAFETY INFORMATION

Please read all safety instructions before using the device to ensure safe and effective operation.

- Do not use on broken, irritated, or infected skin.
- Avoid direct contact with eyes when using the LED Light Therapy Head. Protective eyewear is recommended.
- Consult a physician before use if you are pregnant, have a pacemaker, or any other medical condition.
- Keep out of reach of children.
- Do not immerse the device in water.
- Use only the provided charging cable.

SETUP

Charging the Device

1. Connect the USB charging cable to the Type-C charging port on the device.
2. Plug the other end of the USB cable into a suitable USB power adapter (not included) or a computer USB port.
3. The battery indicator on the LCD screen will show charging status. Ensure the device is fully charged before first use.

Attaching/Changing Massage Heads

To attach a massage head, align the head with the main device body and gently push until it clicks into place. To remove, pull the head gently away from the device body.


OPERATING INSTRUCTIONS

Ensure your skin is clean and dry before use. For microcurrent treatments, apply a conductive gel or serum to the treatment area.

1. Microcurrent Facial Massage Head

This head features two rollers for facial contouring and skin tightening.

EMS V-shaped Facial Massage

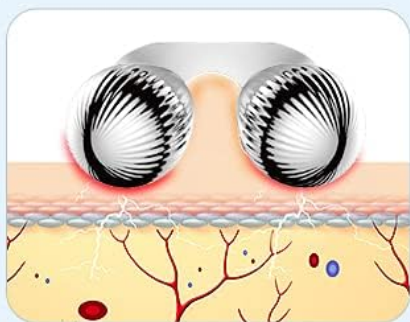

Microcurrent



Apply with gel



Image: A woman demonstrating the use of the EMS V-shaped facial massager, emphasizing the application of gel for effective microcurrent treatment.

1. Attach the Microcurrent Facial Massage Head to the device.
2. Press and hold the Power/Mode/Intensity button to turn on the device.
3. Press the Power/Mode/Intensity button repeatedly to cycle through 5 intensity levels. Select a comfortable intensity.
4. Apply a thin layer of conductive gel or serum to the desired facial area.
5. Gently glide the rollers over your skin. Recommended movements include:
 - From the jawline up to the cheekbones.
 - From the corner of the mouth towards the ear.
 - Slowly from the top of the eyebrows to the hairline (with eyes closed).
6. The treatment duration is 6 minutes. The device will emit an audible alert every 2 minutes, indicating it's time to move to a different treatment area.



Image: Overview of the TOUCHBeauty 3-in-1 Multi-Functional Facial Beauty Device, showcasing its benefits such as firming, lifting, promoting absorption, fading fine lines, relieving eye fatigue, alleviating dark circles, and offering multi-color light therapy.



Image: The LCD display of the TOUCHBeauty AG-1767, illustrating the five available intensity settings for microcurrent mode and the 5-minute smart timing feature.

2. Microcurrent Eye & Lip Massage Head

This head is designed for delicate areas around the eyes and lips to reduce wrinkles and fine lines.

EMS Eye&Lip Massage



Relieve Fine Lines





Image: A woman using the EMS eye and lip massager, demonstrating its application to relieve fine lines around the eyes and mouth.

1. Attach the Microcurrent Eye & Lip Massage Head to the device.
2. Press and hold the Power/Mode/Intensity button to turn on the device.
3. Press the Power/Mode/Intensity button repeatedly to cycle through 3 intensity levels. Select a comfortable intensity.
4. Apply a thin layer of conductive gel or serum to the desired eye or lip area.
5. Gently move the head around the contours of your eyes and lips.
6. The treatment duration is 3 minutes. The device will emit an audible alert every minute, reminding you to change the treatment area.

3. LED Light Therapy Head

This head utilizes different LED light wavelengths to target various skin concerns.

4 in 1 LED Light Therapy

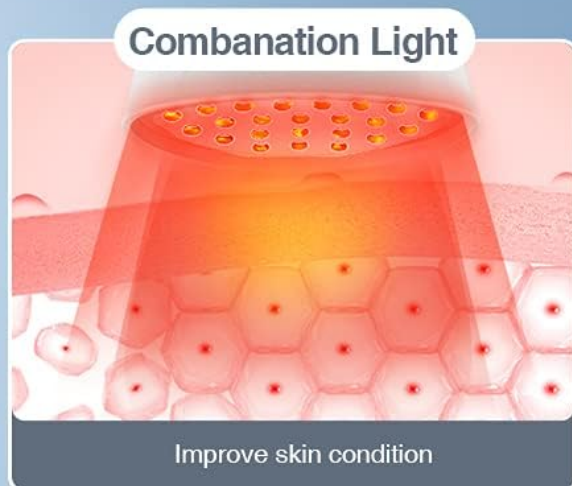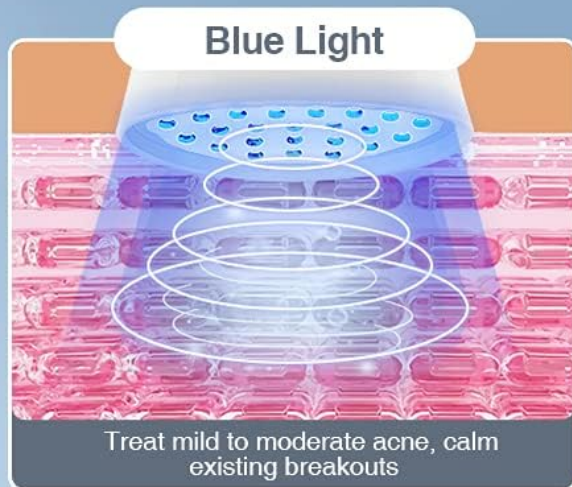

Notice: Please avoid eye area when using

Image: An illustration detailing the four types of LED light therapy offered by the device: Red Light for collagen regeneration and fading fine lines, Blue Light for treating acne and calming breakouts, Yellow Light for brightening skin and diminishing melanin, and Combination Light for overall skin improvement.

1. Attach the LED Light Therapy Head to the device.
2. Press and hold the Power/Mode/Intensity button to turn on the device.
3. Press the Power/Mode/Intensity button repeatedly to select the desired LED light mode:
 - **Red Light:** Helps regenerate collagen and fade fine lines. (6-minute smart timer)
 - **Yellow Light:** Brightens skin and diminishes melanin. (6-minute smart timer)
 - **Blue Light:** Treats mild to moderate acne and calms breakouts. (3-minute smart timer)
 - **Combination Light:** Improves overall skin condition.
4. Hold the LED head approximately 1-2 cm away from the skin, moving it slowly over the treatment area.
5. **Important:** Avoid direct eye contact with the LED light.

General Usage Tips

- For optimal results, use the device regularly as part of your skincare routine.
- Always apply a conductive gel or serum for microcurrent treatments to ensure proper conductivity and comfort.

- Clean the device and heads after each use.

MEDIA RESOURCES

Product Demonstration Video

Your browser does not support the video tag.

Video: A brief overview demonstrating the key features and usage of the TOUCHBeauty AG-1767 3-in-1 Facial Beauty Device, including its microcurrent and LED light therapy functions.

MAINTENANCE

Cleaning

- After each use, wipe the massage heads with a soft, damp cloth.
- Do not use abrasive cleaners or harsh chemicals.
- Ensure the device is completely dry before storing.

Storage

- Store the device in a cool, dry place, away from direct sunlight and extreme temperatures.
- Keep the device and its accessories in their original packaging or a protective case to prevent damage.

TROUBLESHOOTING

Problem	Possible Cause	Solution
Device does not turn on	Low battery or not charged	Charge the device fully.
No sensation during microcurrent treatment	Insufficient conductive gel or low intensity setting	Apply more conductive gel. Increase the intensity level.
Device stops unexpectedly	Treatment timer expired or low battery	Recharge the device. Restart the treatment if needed.
Skin irritation after use	High intensity setting or sensitive skin	Reduce intensity. Discontinue use if irritation persists and consult a dermatologist.

SPECIFICATIONS

- **Model Number:** AG-1767
- **Brand:** TOUCHBeauty
- **Material:** Metal
- **Power Source:** Rechargeable battery
- **Voltage:** 230 Volts
- **Product Dimensions (L x W x H):** 10 x 5.5 x 15 cm
- **Product Weight:** 650 grams
- **First Available Date:** October 10, 2018

WARRANTY AND SUPPORT

For warranty information, technical support, or customer service inquiries, please refer to the official TOUCHBeauty website or contact your retailer. Keep your purchase receipt for warranty claims.