

mornwell RST2203

Mornwell RST2203 Electric Rotary Toothbrush User Manual

Model: RST2203

1. IMPORTANT SAFETY INFORMATION

Please read all instructions carefully before using this product. Keep this manual for future reference.

- Keep out of reach of children under 8 years old.
- Do not use for purposes other than oral hygiene.
- Replace brush head every 3 months.
- Statements regarding dietary supplements have not been evaluated by the FDA and are not intended to diagnose, treat, cure, or prevent any disease or health condition.

2. PRODUCT OVERVIEW

The Mornwell RST2203 Electric Rotary Toothbrush is designed for effective oral cleaning with its high-speed rotating brush head and multiple cleaning modes. It features a smart timer and is USB rechargeable for convenience.



Image: Mornwell RST2203 Electric Rotary Toothbrush package contents, including the handle, four round replacement brush heads, a USB charging cable, and the instruction manual.

Key Features:

- **Powerful Rotary Action:** High-speed rotation of 8800 oscillations per minute for thorough cleaning.
- **3 Cleaning Modes:** Clean, Soft, and Massage modes to suit various oral care needs.
- **2-Minute Smart Timer:** Ensures dentist-recommended brushing time with 30-second interval pauses.
- **IPX7 Waterproof:** Designed for safe use in wet environments.
- **USB Rechargeable:** 3.5 hours charging provides up to 30 days of use.
- **Lockout Function:** Prevents accidental activation during travel.

3. SETUP AND CHARGING

3.1 Attaching the Brush Head

Align the brush head with the metal shaft on the toothbrush handle and push it down until it clicks into place. Ensure it is securely attached.

Your browser does not support the video tag.

Video: Demonstration of attaching the brush head to the Mornwell electric toothbrush handle and connecting the USB charging cable.

3.2 Charging the Toothbrush

Connect the USB charging cable to the charging port at the bottom of the toothbrush handle. Plug the other end into any standard USB power source (adapter not included). The battery indicator light will show charging status.



Image: The Mornwell RST2203 Electric Rotary Toothbrush highlighting its 30-day battery life and 3.5-hour fast charging via a DC charging cable.

A full charge takes approximately 3.5 hours and provides up to 30 days of regular use.

4. OPERATING INSTRUCTIONS

4.1 Power On/Off and Mode Selection

Press the power button to turn the toothbrush on. Press the power button repeatedly to cycle through the three cleaning modes: Clean, Soft, and Massage. The indicator light next to the selected mode will illuminate.



Image: Close-up of the Mornwell RST2203 Electric Rotary Toothbrush handle, illustrating the three cleaning modes: Clean (for deep cleaning), Soft (suitable for novices), and Massage (suitable for sensitive teeth).

4.2 Brushing Technique

Apply a small amount of toothpaste to the brush head. Place the brush head against your teeth at a slight angle towards the gum line. Move the brush head slowly from tooth to tooth, allowing the rotary action to clean each surface. Do not scrub vigorously.

Removes Up to

500% More Plaque

*vs. a manual toothbrush

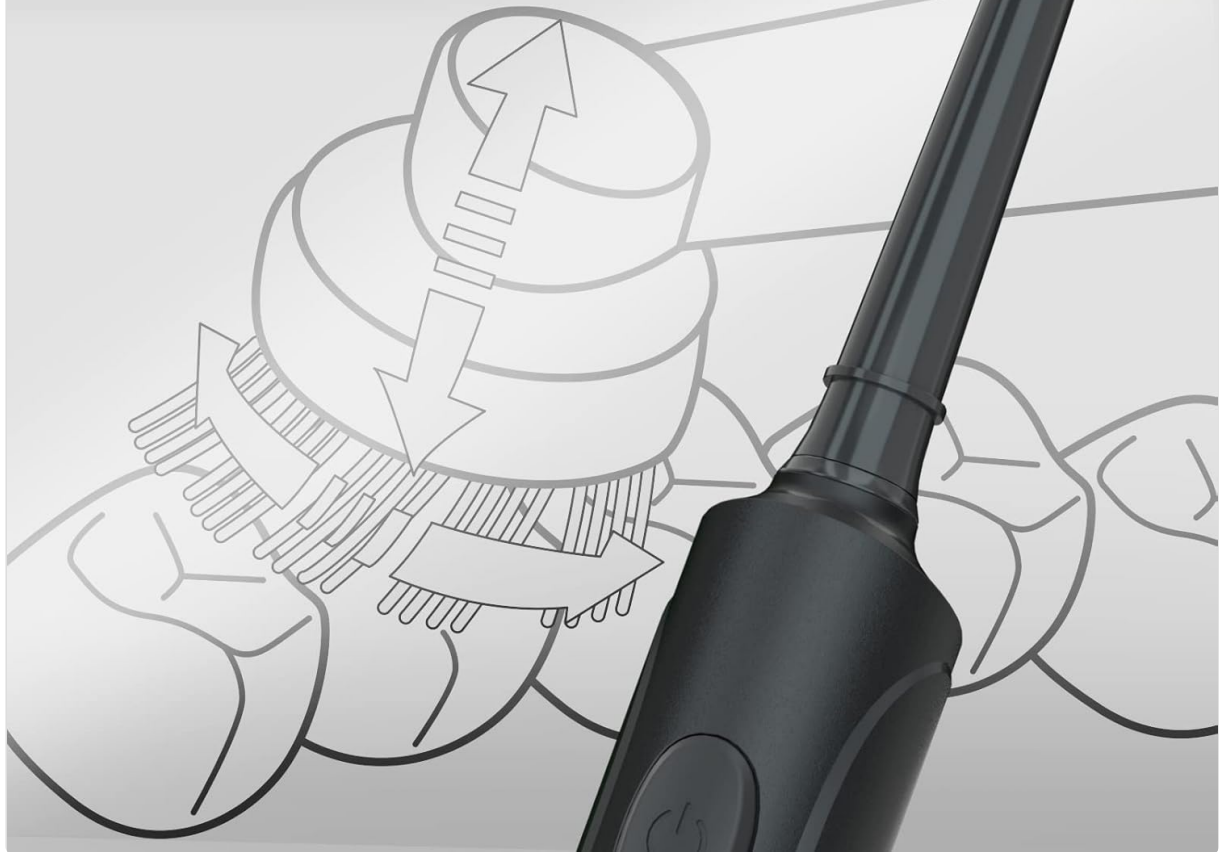


Image: Illustration showing the rotary action of the Mornwell RST2203 Electric Rotary Toothbrush effectively removing plaque from tooth surfaces.

4.3 Smart Timer Function

The toothbrush features a 2-minute smart timer, which is the dentist-recommended brushing time. It will pause briefly every 30 seconds to remind you to move to a different quadrant of your mouth, ensuring even cleaning. The toothbrush will automatically shut off after 2 minutes.



Image: The Mornwell RST2203 Electric Rotary Toothbrush in a bathroom setting, highlighting its 30-second interval pause and 2-minute automatic shut-off features.

4.4 Lockout Function

To prevent accidental activation during travel, the toothbrush has a lockout function. If the toothbrush is in lockout mode and keeps flashing, press and hold the power button for 5 seconds to release the lockout function. Ensure the battery is charged.



Lockout Function Prevents Accidental Start

Long press for 5s to turn on/off

*Tips:

If you can't turn on your toothbrush, it might be in lock mode. Just press for 5s to release.

Image: The Mornwell RST2203 Electric Rotary Toothbrush packed in a suitcase, demonstrating the lockout function to prevent accidental starts during travel. A tip indicates to long press for 5 seconds to turn on/off if in lock mode.

5. MAINTENANCE

5.1 Cleaning the Toothbrush

After each use, rinse the brush head thoroughly under running water. Wipe the handle with a damp cloth. Do not immerse the entire handle in water for extended periods, despite its IPX7 rating, to prolong its lifespan.

5.2 Replacing Brush Heads

Dentists recommend replacing your brush head every 3 months, or sooner if the bristles become frayed. To replace, pull the old brush head straight off the handle and snap a new one into place.



Image: The Mornwell RST2203 Electric Rotary Toothbrush handle displayed alongside four new round replacement brush heads, emphasizing the availability of spares.

6. TROUBLESHOOTING

6.1 Toothbrush Not Turning On or Flashing

If the toothbrush does not turn on or the indicator lights are flashing, it might be in lockout mode. Press and hold the power button for 5 seconds to release the lockout function. Ensure the battery is charged.

If the toothbrush is not turning on, ensure it is fully charged. If the issue persists, try holding the power button for 5 seconds to deactivate the lockout function. If the problem continues, contact customer support.

6.2 Reduced Performance

If the toothbrush performance seems reduced, ensure it is fully charged. Also, check the brush head for wear and tear; a worn brush head can reduce cleaning effectiveness. Replace if necessary.

7. SPECIFICATIONS

Feature	Detail
---------	--------

Model Number	RST2203
Brand	mornwell
Power Source	Battery Powered (1 A battery included)
Oscillations per Minute	8800
Charging Time	3.5 hours
Battery Life	Up to 30 days (minimum)
Special Feature	Brushing Timer, Lockout Function
Waterproof Rating	IPX7
Dimensions	9.37 x 1.24 x 1.2 inches
Weight	6.74 ounces

8. WARRANTY AND SUPPORT

For any quality problems with your Mornwell toothbrush, please contact us directly. While Amazon's return policy is 30 days, Mornwell is committed to customer satisfaction. Refer to the included instruction manual for specific warranty details or contact Mornwell customer support for assistance.

Manufacturer: Prooral Health Technology (Shenzhen) Ltd.