

iVOLCONN IN-006H

iVOLCONN Hand Massager IN-006H Instruction Manual

Model: IN-006H | Brand: iVOLCONN

1. INTRODUCTION

Thank you for choosing the iVOLCONN Hand Massager. This device is designed to provide soothing relief and relaxation for your hands, wrists, and fingers through a combination of heat therapy, air compression, and pressure point massage. It features 3 distinct modes and 3 intensity levels to customize your experience. Please read this manual thoroughly before use to ensure safe and effective operation.



Image 1: The iVOLCONN Hand Massager providing relaxation. The device is white and silver, with a hand inserted into the opening.

2. PRODUCT FEATURES

- **Cordless Design:** Enjoy portability and use anywhere.
- **Heat Therapy:** Gentle warmth to soothe muscles and improve circulation.
- **Air Compression Massage:** Applies pressure to the wrist, palm, and fingers.
- **3 Massage Modes:** Customizable massage experiences.
- **3 Intensity Levels:** Adjust pressure to your comfort.
- **Targeted Pressure Points:** Designed to stimulate key pressure points in the hand.

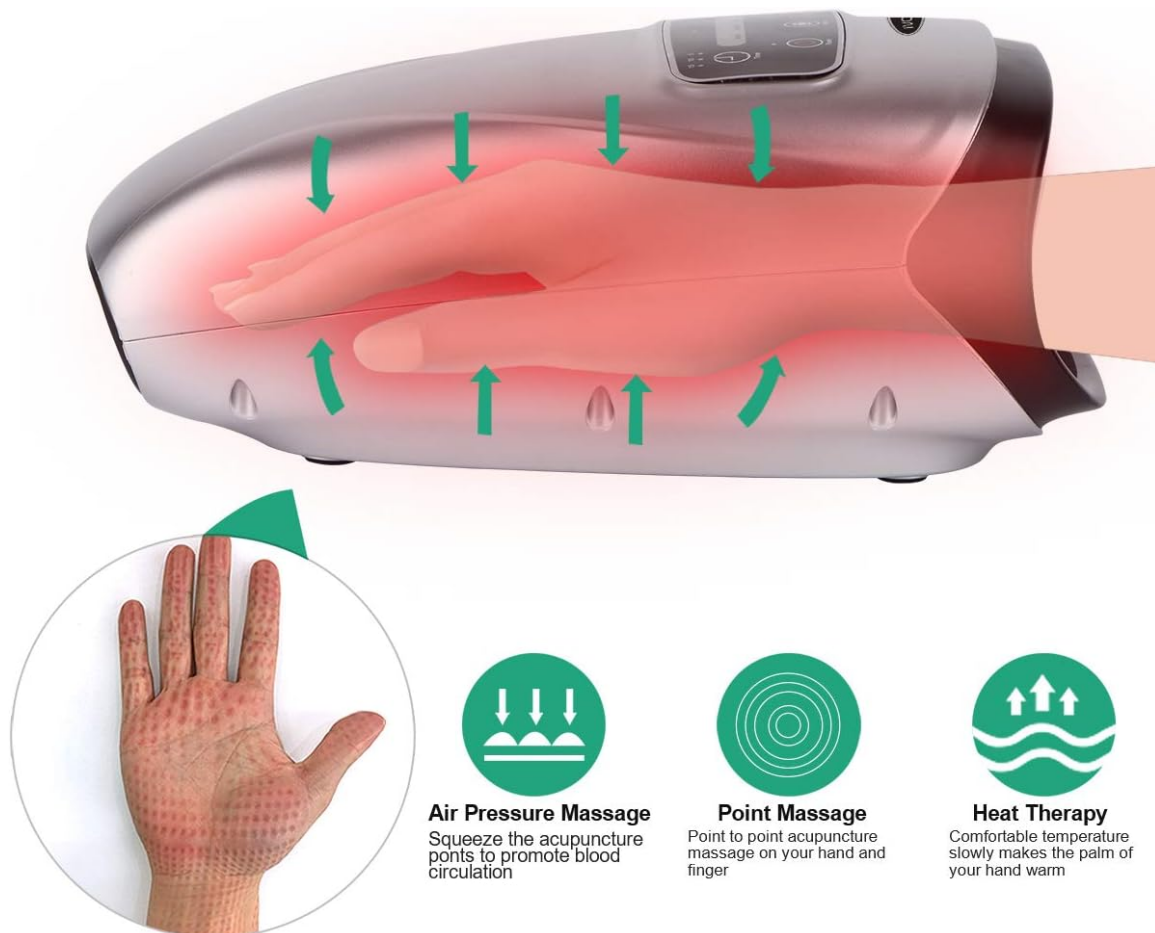


Image 2: A visual representation of the massager's functions, showing air pressure, point massage, and heat therapy applied to the hand.

3. PACKAGE CONTENTS

Ensure all items are present in the package:

- iVOLCONN Hand Massager Unit
- USB-C Charging Cable
- User Manual (this document)
- Disposable Gloves (for hygiene or use with lotions)
- Finger Massager Tool

4. SETUP

1. **Charging:** Before first use, fully charge the device using the provided USB-C cable. The charging indicator light will show charging status. A full charge typically takes a few hours.
2. **Placement:** Place the massager on a flat, stable surface.
3. **Hand Insertion:** Gently insert your hand into the massager opening. Ensure your palm is flat against the internal nodes and fingers are comfortably positioned.

5. OPERATING INSTRUCTIONS

The control panel features intuitive buttons for easy operation:

Hand Reflexology

It compresses trigger points on the palm and back of hand using a mechanical air pressure system that air pressure squeeze massage the hand. physical therapists found the massager to be effective in relieving pain and discomfort in the hands; its soothing rhythm helps to relieve the tension, stress, and fatigue

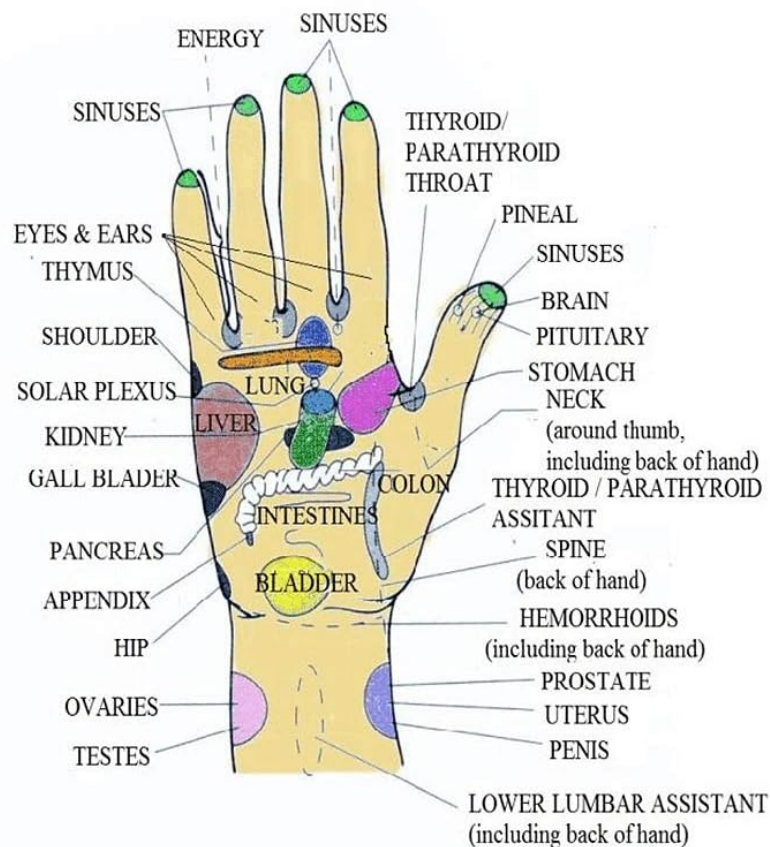


Image 3: Detailed view of the control panel, highlighting each button and indicator.

- **Power Button:** Press and hold to turn the device ON/OFF.
- **Mode Button:** Cycle through 3 different massage modes (e.g., Refresh, Recover, Relax). Each mode offers a unique massage pattern targeting different areas of the hand.
- **Intensity Button:** Adjust the air compression strength between 3 levels (Low, Medium, High).
- **Heat Button:** Toggle heat therapy ON/OFF. The heat function provides gentle warmth to enhance relaxation.
- **Timer:** The device typically operates in 15-minute cycles and will automatically shut off.

Experiment with different modes and intensity levels to find your preferred setting. It is recommended to start with a lower intensity and gradually increase as desired.

6. MAINTENANCE

- **Cleaning:** Disconnect the device from power before cleaning. Wipe the exterior with a soft, damp cloth. Do not use abrasive cleaners or immerse the device in water.
- **Storage:** Store the massager in a cool, dry place away from direct sunlight and extreme temperatures.
- **Disposable Gloves:** Use the provided disposable gloves to maintain hygiene, especially if applying lotions or creams before massage.

7. TROUBLESHOOTING

- **Device not turning on:** Ensure the device is fully charged. If the issue persists, contact customer support.
- **No or weak compression:** Check the selected intensity level. Ensure your hand is properly inserted.
- **Heat not working:** Verify the heat function is activated on the control panel. Allow a few minutes for the heat to build up.
- **Unusual noises:** If you hear any abnormal sounds, discontinue use and contact customer support.

8. SPECIFICATIONS

Attribute	Detail
Brand	iVOLCONN
Model Number	IN-006H
Power Source	Battery Powered (1 Lithium Ion battery included)
Material	Plastic
Color	Silver
Package Dimensions	12.13 x 7.64 x 5.51 inches
Item Weight	2.78 Pounds



Image 4: The physical dimensions of the hand massager, including finger window size.

9. SAFETY INFORMATION

- Do not use if you have any open wounds, skin irritations, or severe hand conditions without consulting a medical professional.
- Keep out of reach of children.
- Do not attempt to disassemble or repair the device yourself.
- Avoid using the device while driving or operating heavy machinery.
- Statements regarding this device have not been evaluated by the FDA and are not intended to diagnose, treat, cure, or prevent any disease or health condition.

10. WARRANTY AND SUPPORT

For warranty information or customer support, please refer to the contact details provided on the product packaging or visit the official iVOLCONN website.

