



Manuals.plus /

- › **Tilcare** /
- › Tilcare Male Urinals with Glow-in-the-Dark Lid (3-Pack) - Instruction Manual

Tilcare 1

Tilcare Male Urinals with Glow-in-the-Dark Lid (3-Pack)

Instruction Manual

1. PRODUCT OVERVIEW

The Tilcare Male Urinal is designed for convenience and discretion, offering a reliable solution for individuals requiring a portable or bedside urination aid. This 3-pack includes durable, 1000ml capacity plastic urinals, each equipped with a unique glow-in-the-dark, leak-proof screw lid and a sealed handle for secure handling and placement.



Image: Tilcare Male Urinal with blue screw lid. This image shows the overall design of the urinal, highlighting its compact and ergonomic shape.

Key Features:

- **Leak-Proof Screw Lid:** Ensures no spills or odors, providing peace of mind.
- **Glow-in-the-Dark Cap:** The lid glows after exposure to light, making it easy to locate and use in low-light conditions.
- **Durable Construction:** Made from thick, high-quality plastic for long-lasting use.
- **Sealed Handle:** Designed to prevent urine from entering the handle, ensuring hygiene and a secure grip. It is also shaped to fit most bed rails.
- **1000ml Capacity:** Ample volume for extended use.
- **Accurate Measurement Markings:** Translucent walls with clear markings allow for easy volume monitoring.

2. SETUP

Before first use, please follow these steps:

1. **Unpack:** Carefully remove the urinals from their packaging.
2. **Initial Cleaning:** Wash each urinal thoroughly with warm water and mild soap. Rinse completely and allow to air dry.
3. **Charge Glow-in-the-Dark Lid:** For optimal glow, expose the lid to natural sunlight or a bright artificial light source for a few minutes. This charges the phosphorescent material, allowing it to glow in the dark.

EASY TO FIND AND USE IN THE DARK



STEP 1

CHARGE THE LID
IN THE SUNLIGHT



STEP 2

GLOW IN THE
DARK CAP



Image: Instructions for charging the glow-in-the-dark lid. Step 1 shows placing the lid in sunlight, and Step 2 shows the lid glowing in the dark.

3. OPERATING INSTRUCTIONS

Using the Tilcare Male Urinal is straightforward:

1. **Preparation:** Ensure the urinal is clean and ready for use. If using in the dark, ensure the lid has been charged.
2. **Positioning:** The wide and smooth opening is designed for comfortable and secure use. Position the urinal appropriately to prevent spills.
3. **Secure the Lid:** After use, firmly screw the glow-in-the-dark lid onto the urinal. Ensure it is tightly closed to prevent any leaks or odors.

4. **Placement:** The sealed handle allows for easy carrying and can be hooked onto most bed rails for convenient access.



Image: Hand holding a Tilcare urinal with a blue liquid, showing a leak-proof icon. This illustrates the secure, leak-proof design of the urinal.

Recommended Uses:

- **Bedside Use:** Ideal for bedridden patients or individuals with limited mobility, allowing for convenient use without leaving the bed.
- **Travel:** Perfect for long car trips, camping, hiking, or any situation where immediate access to a restroom is unavailable.
- **Emergency:** A practical solution for unexpected situations.

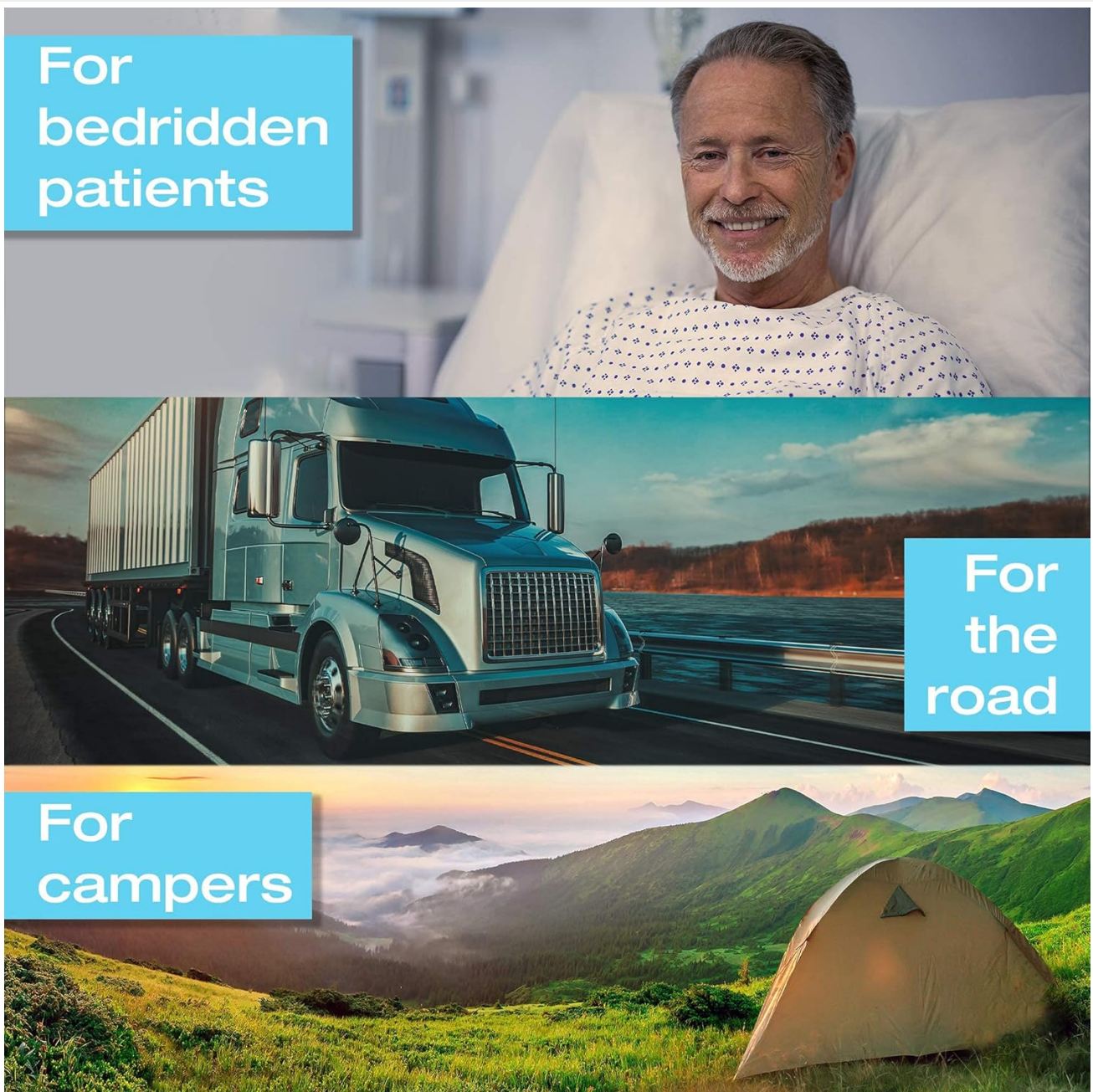


Image: Collage showing various use cases for the urinal, including for bedridden patients, during travel on the road, and for campers.

4. MAINTENANCE

Proper maintenance ensures hygiene and extends the life of your Tilcare Male Urinal:

- **Cleaning After Each Use:** Empty the urinal promptly after each use. Rinse thoroughly with cold water.
- **Sanitization:** Wash with warm water and a mild disinfectant soap. Use a brush if necessary to clean all internal surfaces.
- **Rinsing:** Rinse thoroughly with clean water to remove all soap residue.
- **Drying:** Allow the urinal to air dry completely before re-attaching the lid or storing.
- **Storage:** Store in a clean, dry place away from direct sunlight when not in use.

5. TROUBLESHOOTING

- **Leakage:** If you experience leakage, ensure the screw lid is properly aligned and tightened securely. Check for any

damage to the lid or the urinal's opening.

- **Lid Not Glowing:** If the glow-in-the-dark lid is not glowing, it may not have been sufficiently charged. Expose it to a bright light source (sunlight or artificial light) for a longer period.
- **Odor:** Persistent odor usually indicates inadequate cleaning. Ensure thorough cleaning and sanitization after each use.

6. SPECIFICATIONS

Brand	Tilcare
Model Number	1
Capacity	1000 ml (32 oz)
Material	Durable Plastic
Dimensions (approx.)	24.5 x 23.8 x 10.2 cm (9.6 x 9.4 x 4 inches)
Weight (approx.)	340 grams (0.75 lbs)
Units	3 Total



Image: Tilcare urinal with dimensions and features labeled. This image highlights the wide opening, sealed handle, and accurate measurement markings, along with specific dimensions.

7. WARRANTY AND SUPPORT

Warranty Information:

Your Tilcare Male Urinals come with a **30-day factory guarantee**. This warranty covers manufacturing defects and ensures your peace of mind regarding the quality and functionality of the product.

Customer Support:

For any questions, concerns, or assistance with your Tilcare Male Urinals, please contact our customer support. Refer to the packaging or the retailer's website for specific contact details.

