



Manuals.plus /

› Intelitopia /

› Intelitopia 5.2 Percussion Massage Gun User Manual

Intelitopia 5.2

Intelitopia 5.2 Percussion Massage Gun User Manual

INTRODUCTION

Thank you for choosing the Intelitopia 5.2 Percussion Massage Gun. This device is designed to help relieve muscle stiffness and soreness, increase blood circulation, and improve the overall health of soft tissues. Please read this manual carefully before use to ensure safe and effective operation.

SAFETY INFORMATION

- Do not immerse the device in water or use it in wet environments.
- Keep the device away from heat sources.
- Do not apply excessive pressure to the massage gun during use, as this may cause the device to stop or lead to injury.
- Consult a physician before using this product if you have any medical conditions, are pregnant, or have a pacemaker.
- Keep out of reach of children.

PACKAGE CONTENTS

- 1x Intelitopia 5.2 Percussion Massage Gun
- 4x Massage Heads
- 1x Power Charger
- 1x User Manual



Image: The Intelitopia 5.2 Percussion Massage Gun with its included accessories, neatly organized in its packaging.

PRODUCT OVERVIEW

The Intelitopia 5.2 features a durable, ergonomic design for comfortable handling. It includes a powerful brushless motor and a long-lasting battery, indicated by an LED light ring.



Image: The Intelitopia 5.2 Massage Gun displayed with its interchangeable massage heads.

SETUP AND CHARGING

1. **Initial Charge:** Before using the device for the first time, ensure it is fully charged.
2. **Power On for Charging:** Locate the power button at the bottom of the massager handle. Turn this button to the 'ON' position before connecting the charger.
3. **Connect Charger:** Plug the power charger into the device and then into a 110-220V power outlet.
4. **Charging Indicator:** The LED light ring will blink during charging. When fully charged, the blue light will illuminate completely around the ring. Charging typically takes 2-3 hours.

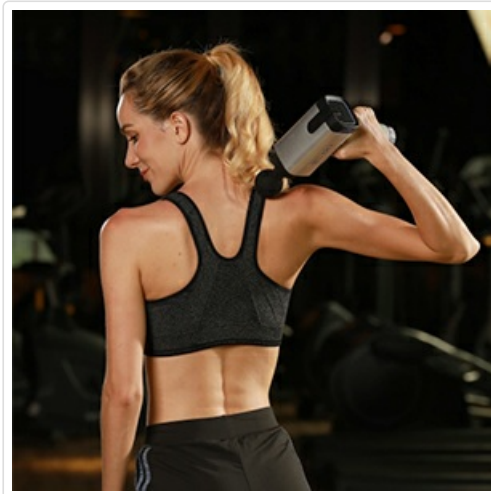


Image: Close-up of the LED battery indicator on the Intelitopia 5.2 Massage Gun, showing charge status.

OPERATING INSTRUCTIONS

1. **Attach Massage Head:** Select the desired massage head and firmly insert it into the front opening of the massage gun.
2. **Power On:** Ensure the power button at the bottom of the handle is in the 'ON' position. Press the main power button on the device to start the massage gun.
3. **Adjust Speed:** Press the main power button repeatedly to cycle through the three adjustable speed levels:

- Level 1: 1800 percussions per minute
- Level 2: 2400 percussions per minute
- Level 3: 3200 percussions per minute

4. **Apply to Body:** Gently apply the massage head to the desired muscle area. Move the device slowly over the muscles, allowing the percussion to work. Avoid applying excessive pressure.
5. **Power Off:** To turn off the device, cycle through all speed levels until the device stops, or turn the power button at the bottom of the handle to the 'OFF' position.

Your browser does not support the video tag.

Video: Demonstration of attaching a massage head and operating the Intelitopia 5.2 Percussion Massage Gun, showcasing its ease of use and adjustable speeds.

MESSAGE HEADS

The Intelitopia 5.2 comes with four specialized massage heads for targeting different muscle groups:

- **Fork Head:** Designed for the spine and neck.
- **Flat Head:** Suitable for any part of the body, offering a general massage.
- **Bullet Head:** Ideal for joints and deep tissue massage.
- **Round Head:** Best for large muscle groups like arms, back, hips, and thighs.

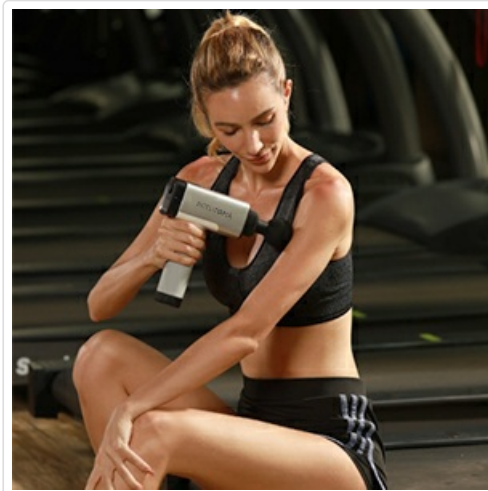


Image: The four interchangeable massage heads included with the Intelitopia 5.2 Percussion Massage Gun.

CARE AND MAINTENANCE

- **Cleaning:** Wipe the device with a slightly damp cloth after each use. Do not use abrasive cleaners or solvents.
- **Storage:** Store the massage gun in a cool, dry place, away from direct sunlight and extreme temperatures.
- **Battery Care:** To prolong battery life, avoid fully discharging the battery frequently. Charge the device regularly, even if not in constant use.

TROUBLESHOOTING

- **Device Not Turning On:** Ensure the power button at the bottom of the handle is in the 'ON' position and the device is sufficiently charged.

- **Device Stops During Use:** This may occur if excessive pressure is applied. Reduce pressure and restart the device. The high-quality motor is designed with ultra-strong heat dissipation to prevent overheating.
- **Reduced Performance:** Ensure the massage head is securely attached and the battery is adequately charged.

If you experience persistent issues, please refer to the contact information provided with your purchase for customer support.

SPECIFICATIONS

Feature	Detail
Model Number	5.2
Dimensions (L x W x H)	8.66 x 6.89 x 3.15 inches
Weight	1000g (2.2 lbs)
Speed Levels	3 (1800, 2400, 3200 RPM)
Battery Capacity	15000mAh Lithium-ion
Working Time	4-6 Hours per charge
Charging Time	2-3 Hours
Charging Voltage	110-220 V
Power Source	Battery Powered
Massage Technique	Percussion
Material	Silicone (grip)
Water Resistance	Not Water Resistant