

## Waterpik SF-02

# Waterpik Sonic-Fusion Professional Flossing Electric Toothbrush (Model SF-02)

Integrated Sonic Toothbrush and Water Flosser

## INTRODUCTION

---

The Waterpik Sonic-Fusion SF-02 combines the benefits of an advanced rechargeable sonic electric toothbrush with the proven effectiveness of a Waterpik Water Flosser. This single device is designed to provide a complete oral clean by brushing and flossing simultaneously or independently, removing harmful bacteria and debris from between teeth and below the gumline.

## WHAT'S INCLUDED

---

# Advanced pressure control

With 10 water flossing pressure levels for a custom clean



Image: Contents of the Waterpik Sonic-Fusion SF-02 box, showing the main unit, two flossing brush heads, and a deluxe travel case.

- Waterpik Sonic-Fusion Main Unit (Model SF-02)
- 2 Water Flossing Brush Heads (color-coded)
- Deluxe Toothbrush Travel Case
- Built-in Charging Base

## SETUP

1. **Charge the Unit:** Place the Sonic-Fusion handle onto the built-in charging base. Ensure the unit is plugged into a global voltage (100-240VAC, 50/60Hz) compatible outlet. The recharge indicator light will show charging status.
2. **Fill the Reservoir:** Lift the hinged lid of the reservoir. Fill with lukewarm water. You may also add your favorite mouthwash to the water.
3. **Attach Brush Head:** Firmly press a water flossing brush head onto the handle.

# Included in the box

**Waterpik Sonic-Fusion®**  
**2 flossing brush heads**  
**Deluxe toothbrush case**



Image: Filling the easy-fill reservoir with water. The reservoir has a 60-second water capacity.

## OPERATING INSTRUCTIONS

The Waterpik Sonic-Fusion offers three cleaning modes: Brush, Floss, and Brush+Floss. It also features 10 advanced water pressure settings for a customized cleaning experience.

### Modes of Operation:

- **Brush Mode:** For sonic brushing only. Apply toothpaste to the brush head.
- **Floss Mode:** For water flossing only.
- **Brush+Floss Mode:** Combines sonic brushing and water flossing simultaneously.

### Using the Device:

1. **Select Pressure Setting:** Adjust the advanced water pressure control dial to your preferred cleaning level (1-10).
2. **Position Brush Head:** Place the brush head in your mouth, aiming the bristles at the gumline.
3. **Activate Mode:** Press the desired mode button (Brush, Floss, or both).
4. **Brush and Floss:** Guide the brush head along your gumline, pausing briefly between teeth to allow

the water to clean effectively. The 2-minute timer with 30-second pacer will guide you through brushing each quadrant of your mouth.

5. **On-the-Go Brushing:** The handle can be detached from the hose for brushing without the water flosser function. Push the round button on the back of the power handle to disconnect.



Image: The control panel on the Waterpik Sonic-Fusion handle, featuring distinct buttons for Brush and Floss modes.

# Deluxe Toothbrush Case

Holds your toothbrush  
and 2 brush heads

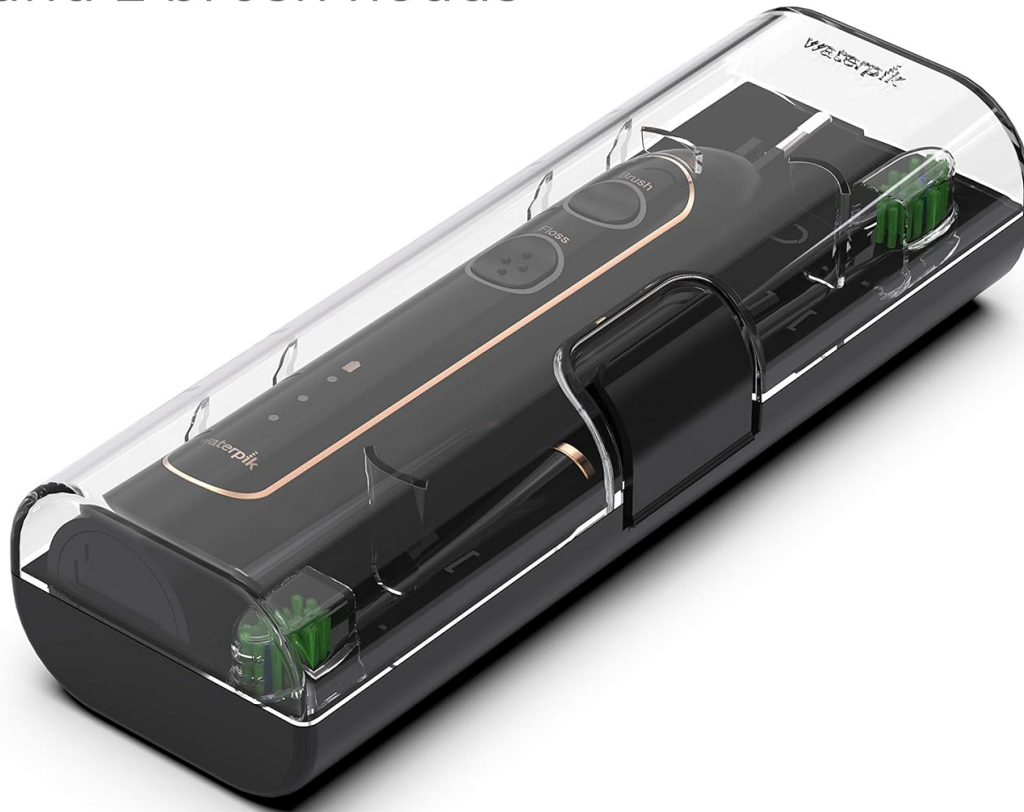


Image: The advanced pressure control dial on the Waterpik Sonic-Fusion unit, offering 10 water flossing pressure levels.

Your browser does not support the video tag.

Video: Demonstrates how the Waterpik Sonic-Fusion allows you to brush and floss simultaneously for a complete oral clean.

Your browser does not support the video tag.

Video: Highlights the dual cleaning action of the Waterpik Sonic-Fusion, combining water and sonic power for enhanced effectiveness.

## MAINTENANCE

- **Rinsing After Mouthwash:** If you use mouthwash in the reservoir, always rinse the unit thoroughly with plain water afterward to prevent buildup.
- **Charging:** For optimal performance, keep the toothbrush handle in the charging base when not in use.
- **Reactivating Handle:** If the handle is removed from the base for more than 10 minutes, return it to the base briefly to re-activate the water flossing function.

## TROUBLESHOOTING

If you encounter any issues with your Waterpik Sonic-Fusion, please refer to the detailed troubleshooting section in the full product manual. For immediate assistance, contact Waterpik Customer Service.

## SPECIFICATIONS

Feature	Detail
Model Number	SF-02
Product Dimensions	5.5 x 5.25 x 10.88 inches
Item Weight	1.75 Pounds
Power Source	Battery Powered (Rechargeable)
Voltage	100 - 240VAC, 50/60Hz (Global Voltage)
Tank Volume	473 Milliliters (14 ounces) for 60 seconds of flossing time
Color	Black
Material Type	Plastic, Stainless Steel
Product Benefits	Gum Health, Plaque Removal
UPC	073950247085

## WARRANTY & SUPPORT

The Waterpik Sonic-Fusion SF-02 is backed by a 3-year manufacturer's limited warranty. For any product questions or needs, please contact the Waterpik support team based in Fort Collins, Colorado:

**Customer Service:** 1-888-226-3042

**Hours:** Monday - Friday, 7 AM - 5 PM MST

**waterpik®**  
**sonic-fusion®**



**Waterpik offers a 3 year  
limited warranty**

**Questions?**

Please contact our  
Customer Service team for any  
product issues or questions.

**Customer Service**  
**1-888-226-3042**  
Monday-Friday  
7 am – 5pm MST



Image: Waterpik Sonic-Fusion unit with warranty information and customer support details.

## **SAFETY INFORMATION**

Please read all instructions before using this product. This device is intended for oral hygiene use only. Refer to the full instruction manual for comprehensive safety warnings and usage guidelines.

- General battery safety warnings
- Oral use safety warnings
- Electrical safety warnings
- Child safety warnings