

Manuals+

[Q & A](#) | [Deep Search](#) | [Upload](#)

Manuals.plus /

› [ANLAN](#) /

› ANLAN Ultrasonic Facial Scrubber ALCPJ07-02 Instruction Manual

ANLAN ALCPJ07-02

ANLAN Ultrasonic Facial Scrubber ALCPJ07-02 Instruction Manual

Model: **ALCPJ07-02** | Brand: **ANLAN**

1. INTRODUCTION

Thank you for purchasing the ANLAN Ultrasonic Facial Scrubber. This device is designed to enhance your skincare routine with multiple functions including ultrasonic peeling, ion cleansing, ion moisturizing, and EMS microcurrent. Please read this manual thoroughly before use to ensure safe and effective operation.

4つ 自宅でお手軽エステ級スキンケア の機能を備え

4つの機能で毛穴の汚れ・メイク残りなどを水の流れで吹き飛ばし、肌の保湿・ハリ感が更にアップ



Image 1: ANLAN Ultrasonic Facial Scrubber and its packaging.

2. SAFETY PRECAUTIONS

- Always perform a patch test before full use to ensure skin compatibility.
- Use only the dedicated charging cable provided in the package. Using other cables may cause power issues.
- Do not use on broken, irritated, or infected skin.
- Keep out of reach of children.
- Consult a doctor before use if you have any medical conditions, especially if you have a pacemaker or other electronic medical implants.

3. PACKAGE CONTENTS

- ANLAN Ultrasonic Facial Scrubber (Main Unit)
- USB Type-C Charging Cable (Dedicated)
- Charging Stand
- Protective Cap
- User Manual (This document)

4. PRODUCT OVERVIEW

Familiarize yourself with the different parts of your ANLAN Ultrasonic Facial Scrubber.

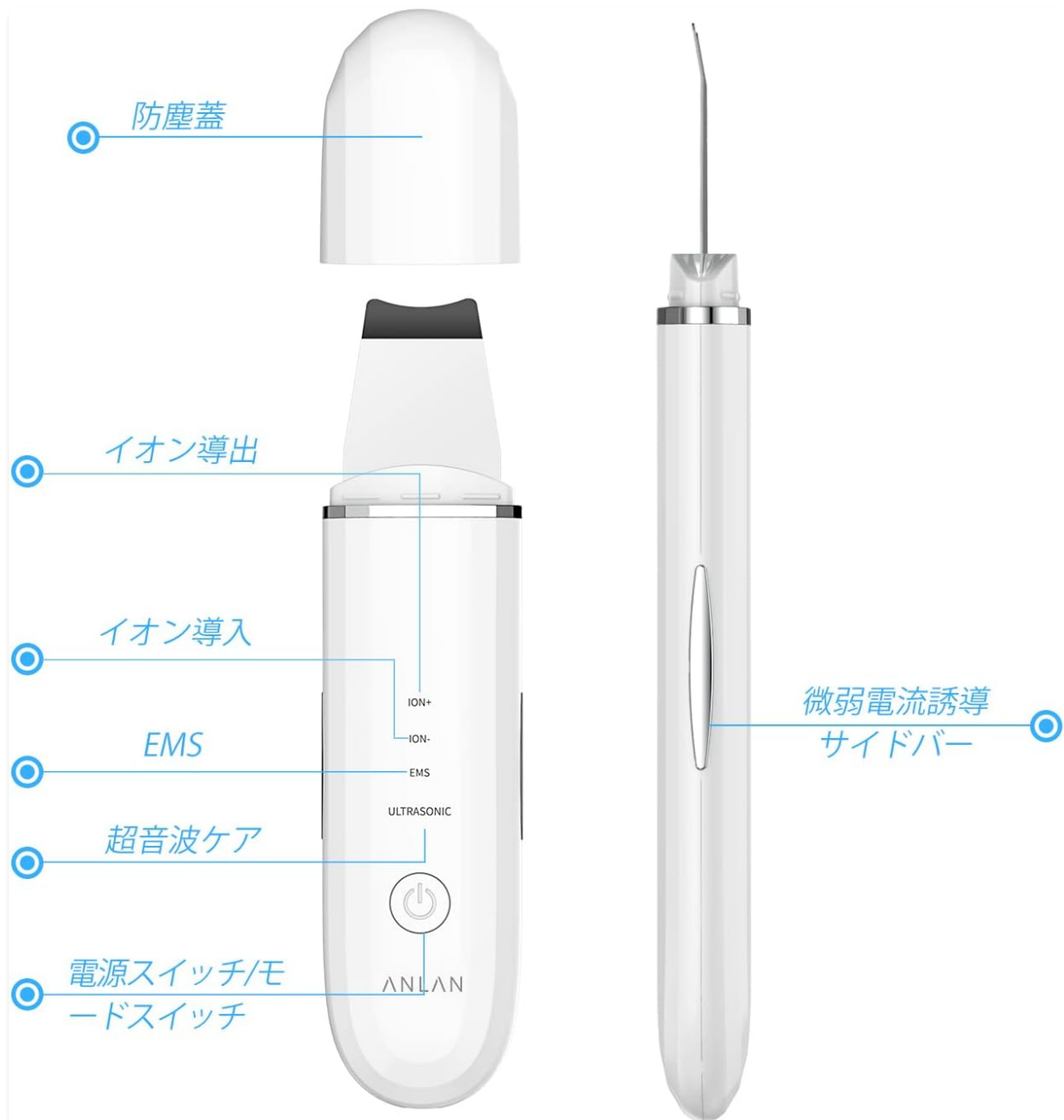


Image 2: Product Parts Diagram. Key components include the protective cap, ion export/import areas, EMS, ultrasonic care, power/mode switch, and microcurrent induction sidebar.

5. SETUP AND CHARGING

5.1 Initial Charge

Before first use, fully charge the device. Connect the dedicated USB Type-C charging cable to the charging stand and then place the device on the stand. The indicator light will show charging status.

軽量コンパクト 持ち運びは便利、 いつでも手軽にエステ可能！



Image 3: Charging the device on its stand. The magnetic charging stand allows for easy and convenient charging.

5.2 Waterproof Feature (IPX7)

The device features IPX7 waterproof rating, allowing safe use in the bath or shower. The charging port is designed to prevent water intrusion.

Your browser does not support the video tag.

Video 1: Demonstration of the ANLAN Royal Water Peeling device's IPX7 waterproof capability, showing it can be used while bathing.

Your browser does not support the video tag.

Video 2: IPX7 immersion test for the ANLAN Royal Water Peeling device, demonstrating its water resistance.

6. OPERATING INSTRUCTIONS

The device features four main skincare modes. Press the power button to turn on/off and cycle through the modes.



Image 4: The four main functions of the ANLAN Ultrasonic Facial Scrubber.

6.1 Ultrasonic Peeling Mode (Sonic)

This mode uses 24,000 vibrations per second to atomize water on the skin, effectively removing sebum, makeup residue, and dirt from pores. It helps to achieve smoother, brighter skin.

- **Atomization:** Water is turned into mist by ultrasonic vibrations.
- **Splitting:** Tiny bubbles are generated and burst, dislodging dirt.
- **Detachment:** Removes small particles of dirt from the skin.

Ultrasonic 超音波ピーリング

24000回/秒の超音波振動でお毛穴
の奥の汚れを除去します。

肌をスベスベ、光沢になる

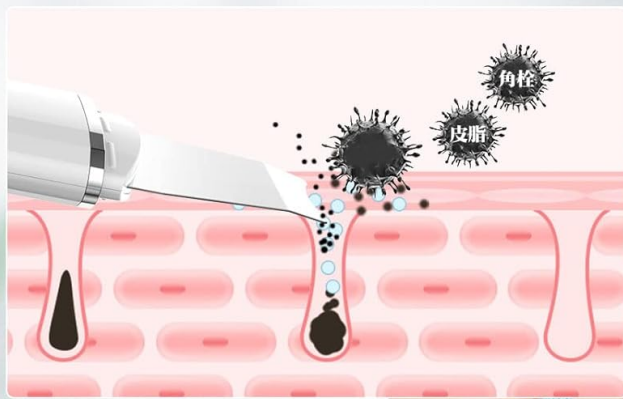


Image 5: Ultrasonic peeling in action, demonstrating the removal of dirt from pores.


6.2 Ion+ Cleansing Mode

This mode uses positive ions and ultrasonic vibrations to absorb deep-seated dirt and makeup residue from the skin. It prepares the skin for better absorption of beauty ingredients.

Ion+ : イオンクレンジング

肌深層の汚れをイオンと超音波振動のチカラで吸着します。

さらに美容成分が浸透しやすいお肌に整えます。

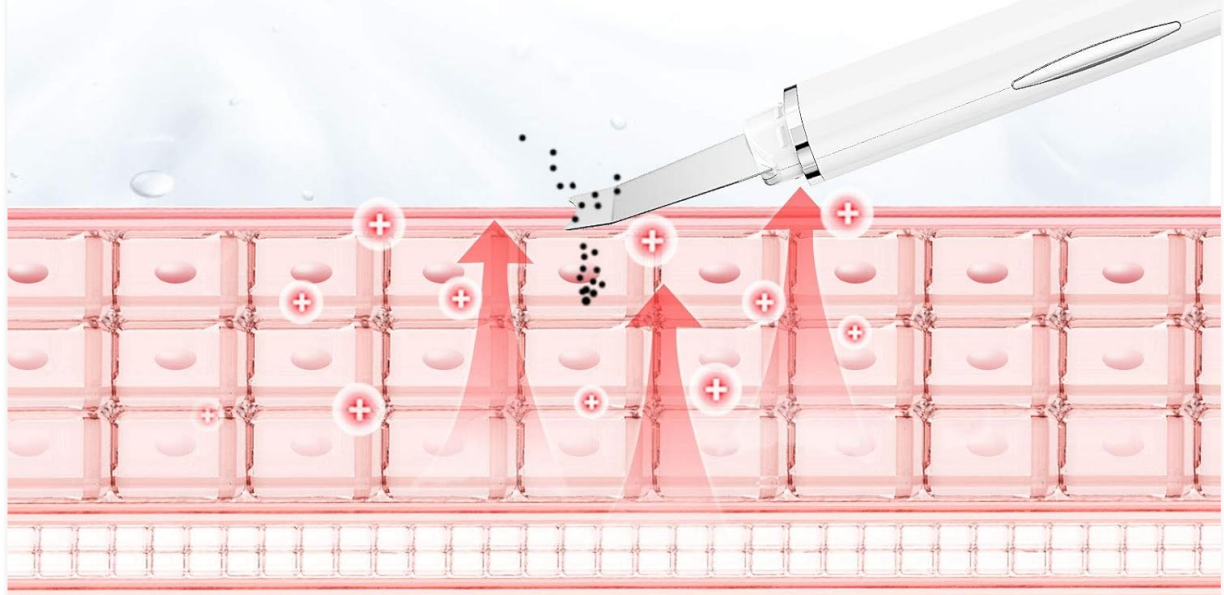

1. いつも通り洗顔をし、水気をしっかり抜き取ります。



2. 霧吹きミストなどで顔を水でたっぷり濡らします。



3. ピーリングヘッドを下から押し上げるように動かします。


Image 6: Ion+ cleansing mode, showing how positive ions attract and remove deep skin impurities.

6.3 Ion- Moisturizing Mode

This mode uses negative ions to help beauty essences and sheet mask ingredients penetrate deeper into the skin, enhancing hydration and elasticity. This can increase the penetration of active ingredients like Vitamin C by approximately 13 times compared to manual application.

Ion- : 浸透モード

美容液の有効成分をお肌の奥までに浸透し、いつもの美容液の効果が格段にアップします。保湿・ハリある肌へ。

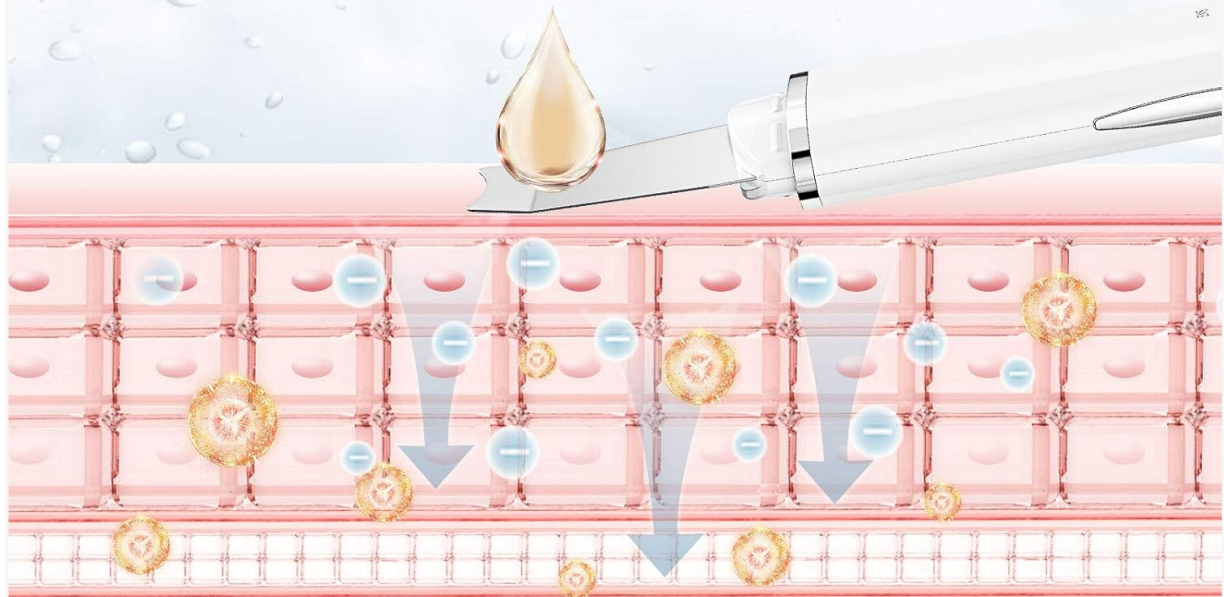

ビタミンCの浸透比較



Image 7: Ion- moisturizing mode, demonstrating enhanced penetration of skincare products.

6.4 EMS Microcurrent Mode

The EMS microcurrent mode targets facial muscles that are difficult to reach manually, helping to tighten and firm the skin, and reduce puffiness.

EMS微電流

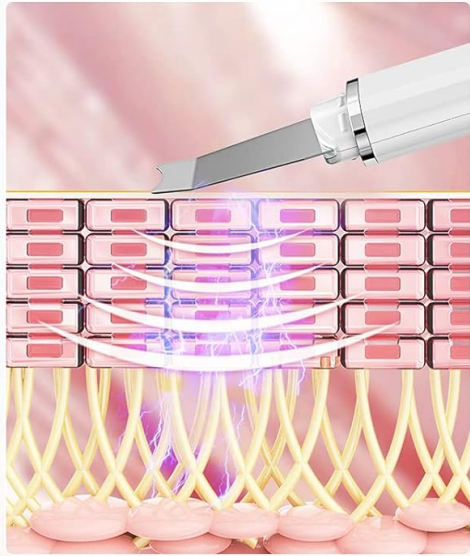


Image 8: EMS microcurrent mode, illustrating muscle stimulation for skin tightening.

7. RECOMMENDED USAGE STEPS

For optimal results, follow these steps:

1. **Daily Cleansing:** Cleanse your face as usual and gently pat dry.
2. **Ion+ Cleansing:** Apply a facial mist or toner to keep your face moist. Use the Ion+ mode to deeply cleanse.
3. **Ultrasonic Peeling:** Keep your face moist. Use the ultrasonic peeling mode, moving the spatula downwards and lifting.
4. **Ion- Moisturizing:** Apply your favorite beauty essence or sheet mask. Use the Ion- mode to help ingredients penetrate.
5. **EMS Massage:** Use the EMS mode to gently massage and firm your skin.
6. **Finish:** Complete your routine with your regular lotion or cream.



Image 9: Visual guide for the recommended usage steps.

8. MAINTENANCE

8.1 Cleaning

After each use, wipe the spatula clean with a soft, damp cloth. Do not use abrasive cleaners or immerse the entire device in water for extended periods, despite its waterproof rating, to prolong its lifespan.

8.2 Storage

Store the device in a cool, dry place, away from direct sunlight and high temperatures. Ensure the protective cap is on the spatula when not in use.

9. TROUBLESHOOTING

9.1 Device Not Turning On/Charging

- Ensure the dedicated USB Type-C charging cable and stand are used.
- Check if the charging cable is securely connected to both the power source and the charging stand.
- Allow the device to charge for at least 2 hours.

9.2 Weak Vibration/No Effect

- Ensure your skin is adequately moist during use, especially for ultrasonic peeling.
- Confirm the correct mode is selected for your desired function.
- Ensure the device is fully charged.

10. SPECIFICATIONS

| Specification | Detail |
|---------------|------------|
| Model Number | ALCPJ07-02 |
| Brand | ANLAN |

| | |
|----------------------|---|
| Dimensions | 18.6 x 12.6 x 2.4 cm |
| Weight | 160 g |
| Material | Acrylonitrile Butadiene Styrene (ABS), Stainless Steel, Polycarbonate |
| Color | White |
| Waterproof Rating | IPX7 |
| Charging Type | USB Type-C (with dedicated stand) |
| Ultrasonic Frequency | Approx. 24,000 vibrations/second |

11. WARRANTY AND SUPPORT

The product comes with a **1-year warranty** from the date of purchase.

For any questions or issues regarding the product, please contact the seller where you purchased the item or ANLAN Customer Service via email. Our customer support team aims to respond within 24 hours.

Please retain your proof of purchase for warranty claims.
