

[Manuals.plus](#) /

> [AquaFlosser](#) /

> AquaFlosser Oral Hygiene Standard Water Hose Plastic Handle for Waterpik Oral Irrigator User Manual

## AquaFlosser WP100HandleV2

# AquaFlosser Oral Hygiene Standard Water Hose Plastic Handle User Manual

Model: WP100HandleV2

## INTRODUCTION

---

This manual provides instructions for the installation, operation, and maintenance of your AquaFlosser Standard Water Hose Plastic Handle. This replacement part is designed to restore the functionality of compatible Waterpik oral irrigators. Please read this manual thoroughly before installation and use to ensure proper function and safety.


## PRODUCT OVERVIEW

---



**Figure 1:** AquaFlosser Standard Water Hose Plastic Handle Assembly. This image displays the complete replacement handle, including the hose, handle unit, and the connection point for the irrigator tip.

The AquaFlosser replacement handle assembly is a crucial component for maintaining the performance of your Waterpik oral irrigator. It includes a durable plastic handle with a pause button, a flexible water hose, and a rotating tip connection for optimal oral hygiene.

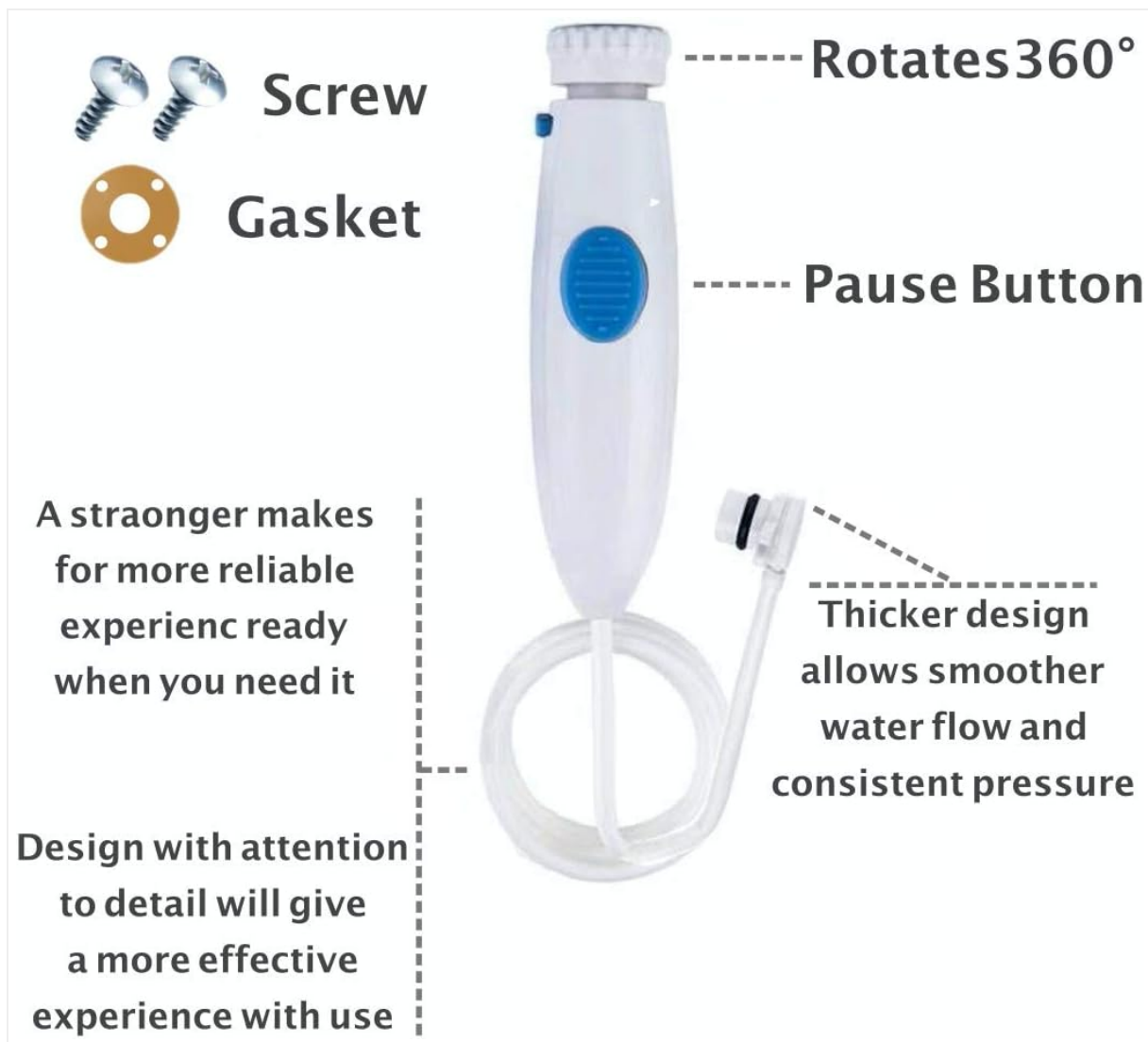

**Figure 2:** Detailed view of the AquaFlosser handle highlighting key features. The image points out the 360° rotatable tip connection, the pause button on the handle, and describes the thicker hose design for smoother water flow and consistent pressure. It also shows included screws and a gasket.

- **Slip-Resistant Design:** Enhances grip during use.
- **Pause Button:** Allows for temporary cessation of water flow directly from the handle.
- **360° Rotating Tip Connection:** Provides flexibility for directing water flow.
- **Thicker Hose:** Engineered for smoother water flow and consistent pressure.

## COMPATIBILITY

This AquaFlosser replacement handle is compatible with the following Waterpik oral irrigator models:

## *Compatible Model :*

☆ *WP-100*

☆ *WP-900*



*Gasket*



*Screw*



**Figure 3:** Visual representation of compatible Waterpik models. The image explicitly lists WP-100 and WP-900, and also shows the included gasket and screws.

- Waterpik WP-100 Series
- Waterpik WP-450 Series
- Waterpik WP-250 Series
- Waterpik WP-300 Series
- Waterpik WP-660 Series
- Waterpik WP-900 Series

Please verify your Waterpik model number before purchasing to ensure compatibility.

### **PACKAGE CONTENTS**

Your package should contain the following items:

- 1 x AquaFlosser Standard Water Hose Plastic Handle Assembly
- 1 x Gasket
- 2 x Mounting Screws

### **SAFETY INFORMATION**

Please observe the following safety precautions:

- Keep away from children.
- Consult a healthcare professional before using if you have medical conditions related to the mouth or gums.

- Do not submerge the device in water or use near a water source while the device is plugged in.
- Unplug the device when not in use.

## INSTALLATION INSTRUCTIONS

Follow these steps to replace the handle assembly on your Waterpik oral irrigator:

1. **Unplug the Unit:** Ensure your Waterpik oral irrigator is unplugged from the electrical outlet before beginning any installation or maintenance.
2. **Drain Water:** Empty any remaining water from the reservoir and the unit.
3. **Access the Base:** Carefully turn your Waterpik unit upside down to access the base where the hose connects. You may need to remove rubber feet or a cover depending on your specific model.
4. **Remove Old Handle:** Locate the screws securing the old hose/handle assembly to the unit. Use an appropriate screwdriver to remove these screws. Gently detach the old hose and handle.
5. **Replace Gasket:** It is recommended to replace the old gasket with the new one provided for optimal sealing and hygiene. Position the new gasket correctly.
6. **Install New Handle:** Align the new AquaFlosser handle assembly with the connection point on your Waterpik unit. Ensure the hose is routed correctly and does not kink.
7. **Secure with Screws:** Use the provided mounting screws to securely fasten the new handle assembly to the unit. Do not overtighten.
8. **Test Connection:** Once installed, gently tug on the hose to ensure it is firmly attached.



**Figure 4:** View of the underside of a Waterpik WP-100W unit, illustrating the hose connection point and screw locations. This image serves as a reference for where the replacement handle assembly connects to the main unit.

## OPERATING INSTRUCTIONS

After successful installation, operate your Waterpik oral irrigator as usual. The AquaFlosser handle functions identically to the original.

- **Fill Reservoir:** Fill the Waterpik reservoir with lukewarm water.
- **Insert Tip:** Place your preferred Waterpik tip into the rotating connection point on the AquaFlosser handle until it clicks into place.
- **Power On:** Plug in the Waterpik unit and turn it on.
- **Adjust Pressure:** Set the desired pressure level on your Waterpik unit.
- **Begin Flossing:** Place the tip in your mouth, lean over a sink, and press the power button on the Waterpik unit. Use the pause button on the AquaFlosser handle for temporary stops.
- **Rotate Tip:** Utilize the 360° rotating tip connection to reach all areas of your mouth effectively.

## MAINTENANCE

Regular maintenance ensures the longevity and hygienic operation of your AquaFlosser handle and Waterpik unit.

- **Gasket Replacement:** It is recommended to change the gasket every 6 months for hygiene and to maintain optimal water pressure and prevent leaks.
- **Cleaning the Handle:** Wipe the handle with a soft, damp cloth. Do not use abrasive cleaners or submerge the handle in water.
- **Draining After Use:** After each use, drain any remaining water from the reservoir and run the unit briefly to clear water from the hose and handle.

## TROUBLESHOOTING

If you encounter issues with your AquaFlosser handle, refer to the following common problems and solutions:

Problem	Possible Cause	Solution
<b>Reduced Water Pressure</b>	Clogged tip, kinked hose, worn gasket, or issue with the main Waterpik unit.	Ensure the tip is clean. Check the hose for kinks. Replace the gasket if it's old. Refer to your Waterpik unit's manual for further troubleshooting.
<b>Water Leakage from Handle/Hose</b>	Loose connection, damaged hose, or improperly installed gasket.	Verify the handle is securely attached to the Waterpik unit. Inspect the hose for visible damage. Ensure the gasket is correctly seated and not damaged.
<b>Tips Pop Off Easily</b>	Worn tip connection point or incompatible tip.	Ensure you are using genuine Waterpik tips or compatible alternatives. Push the tip firmly until it clicks. If the issue persists, the connection point might be worn.

## SPECIFICATIONS

- **Model Number:** WP100HandleV2
- **Brand:** AquaFlosser

- **Material:** Plastic
- **Compatibility:** Waterpik WP-100, WP-450, WP-250, WP-300, WP-660, WP-900 Series
- **Package Dimensions:** Approximately 7.13 x 2.64 x 1.85 inches
- **Item Weight:** Approximately 2.08 ounces

## WARRANTY AND SUPPORT

---

For any questions regarding your AquaFlosser replacement handle, please contact the manufacturer, AquaFlosser. Specific warranty details may vary; please refer to your purchase documentation or contact customer support for more information.

Manufacturer: AquaFlosser