

Safety 1st TH0940600

Safety 1st Quick Read 2-In-1 Thermometer Instruction Manual

Model: TH0940600

1. INTRODUCTION

The Safety 1st Quick Read 2-In-1 Thermometer is designed for fast and accurate temperature readings. It offers both oral and underarm measurement options and includes features like a Reset Timer and Fever Light™ Bar for ease of use. Please read this manual thoroughly before use to ensure proper operation and care.

2. PRODUCT OVERVIEW

The thermometer features a large digital display, a power button, and a flexible tip for comfort. It is designed for quick readings and includes a protective cover for hygiene and storage.



Image 1: The Safety 1st Quick Read 2-In-1 Thermometer displaying a temperature. This image shows the thermometer's main body, digital screen, and power button.

3. SETUP

3.1. Unpacking

1. Carefully remove the thermometer from its packaging.
2. Inspect the device for any damage. Do not use if damaged.

3. Retain the protective cover for storage.

3.2. Battery Installation (if required)

The thermometer comes with batteries pre-installed. If replacement is needed, refer to the battery compartment usually located at the back or side of the device. Ensure correct polarity when inserting new batteries.



Image 2: The thermometer with its clear protective probe cover removed, exposing the metal tip. This illustrates how the thermometer appears when ready for use or when accessing the battery compartment.

4. OPERATING INSTRUCTIONS

4.1. Taking a Temperature (Oral or Underarm)

1. Press the power button to turn on the thermometer. The display will show the last reading, then 'Lo' and a flashing '°F' or '°C' indicating it's ready.
2. **For Oral Use:** Place the tip of the thermometer under the tongue, ensuring good contact. Close the mouth and hold the thermometer in place.
3. **For Underarm Use:** Place the tip of the thermometer in the armpit, ensuring it touches the skin directly. Hold the arm firmly against the body to keep the thermometer in place.
4. Wait for approximately 8 seconds until the thermometer beeps, indicating the reading is complete.
5. Read the temperature displayed on the screen.
6. Press the power button to turn off the thermometer, or it will automatically shut off after a period of inactivity.

4.2. Reset Timer Feature

The thermometer includes a 60-second Reset Timer. After taking a reading, allow 60 seconds before taking another reading to ensure the thermometer has reset to room temperature for optimal accuracy.

4.3. Fever Light™ Bar

The Fever Light™ Bar will glow amber if an above-normal temperature is detected, providing a quick visual alert.

4.4. Switching Between °F and °C

To switch between Fahrenheit (°F) and Celsius (°C), ensure the thermometer is off. Press and hold the power button for approximately 3 seconds until the unit changes. Release the button, and the new unit will be set.

Your browser does not support the video tag.

Video 1: Safety 1st Lifestyle Video. This video provides a general overview of the Safety 1st brand and its commitment to baby safety and development, rather than specific instructions for this thermometer model.

5. MAINTENANCE

5.1. Cleaning

To clean the thermometer, wipe the tip and body with a soft cloth dampened with rubbing alcohol or mild soap and water. Do not immerse the entire thermometer in water. Allow it to air dry completely before storing.

5.2. Storage

Always store the thermometer in its protective cover when not in use. Keep it in a cool, dry place, away from direct sunlight and extreme temperatures.



Image 3: The Safety 1st Quick Read 2-In-1 Thermometer shown with its clear protective cover. This cover helps maintain hygiene and protects the thermometer when stored.

6. TROUBLESHOOTING



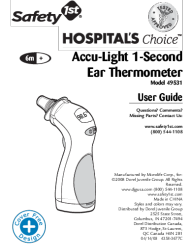

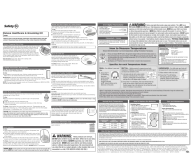

Issue	Possible Cause	Solution
Inaccurate or inconsistent readings	Not waiting for Reset Timer; improper placement; low battery.	Ensure 60-second wait between readings. Verify correct placement (oral/underarm). Replace batteries if needed.
Thermometer does not turn on	Dead batteries; device malfunction.	Replace batteries. If issue persists, contact customer support.
No beep sound when reading is complete	Device malfunction.	Observe the flashing '°F' or '°C' symbol; when it stops flashing, the reading is stable. Contact customer support if the lack of beep is a concern.

7. SPECIFICATIONS

Feature	Detail
Brand	Safety 1st
Model Name	Quick Read 2-In-1 Thermometer
Item Model Number	TH0940600
Display Type	Digital
Response Time	8 seconds
Power Source	Battery Powered
Special Feature	High Accuracy, Fever Light™ Bar, Reset Timer
Product Care Instructions	Wipe Clean
Included Components	Protective Probe Cover

8. WARRANTY AND SUPPORT

For warranty information or customer support, please refer to the contact details provided with your product packaging or visit the official Safety 1st website. You can also visit the [Safety 1st Store on Amazon](#) for more information.

	<p><u>Safety 1st 3-in-1 Nursery Thermometer User Guide</u></p> <p>User guide for the Safety 1st 3-in-1 Nursery Thermometer (TH063, TH064, TH075), providing instructions on usage, care, troubleshooting, and specifications for oral, rectal, and underarm temperature measurement.</p>
	<p><u>Safety 1st Prompt Read 3-in-1 Digital Thermometer User Guide TH097</u></p> <p>User guide for the Safety 1st Prompt Read 3-in-1 Digital Thermometer (Model TH097). Provides instructions on how to use, care for, and troubleshoot the thermometer for accurate oral, rectal, and underarm temperature readings, including normal temperature ranges and warranty information.</p>
	<p><u>Safety 1st Accu-Light 1-Second Ear Thermometer Model 49531 User Guide</u></p> <p>Comprehensive user guide for the Safety 1st Accu-Light 1-Second Ear Thermometer (Model 49531). Learn how to take accurate temperature readings, understand fever information, care for your device, and troubleshoot common issues. Features quick, gentle, and hygienic measurement.</p>
	<p><u>Safety 1st Fever Light™ Ear Thermometer 49551 & 49659 User Guide</u></p> <p>User guide for the Safety 1st Fever Light™ Ear Thermometer (models 49551 and 49659), providing instructions on operation, troubleshooting, care, and battery replacement.</p>
	<p><u>Safety 1st Deluxe Healthcare & Grooming Kit IH496 User Manual</u></p> <p>Comprehensive user manual for the Safety 1st Deluxe Healthcare & Grooming Kit (Model IH496), detailing the use, care, and features of each included item for infant and child health and hygiene. Includes instructions for nail clippers, brush and comb, toothbrush, digital thermometer, medicine dispenser, nasal aspirator, and more.</p>
	<p><u>Safety 1st Deluxe Healthcare & Grooming Kit User Guide and Instructions</u></p> <p>Comprehensive guide for the Safety 1st Deluxe Healthcare & Grooming Kit, including instructions for the digital thermometer, nasal aspirator, grooming tools, and more. Learn about usage, care, and troubleshooting.</p>