

**breo BRE1200H**

# BREO iSee M Eye Massager Instruction Manual

Model: BRE1200H | Brand: breo

## 1. INTRODUCTION

---

The breo iSee M is a portable eye massager designed to provide a soothing experience for the eyes. It incorporates features such as warmth and gentle massage to help alleviate eye strain and promote relaxation. This manual provides essential information for the safe and effective use of your device.



Figure 1.1: User wearing the breo iSee M eye massager, demonstrating its intended use.

## 2. SAFETY INFORMATION

---

Please read all safety instructions carefully before using the device. Failure to do so may result in injury or damage to the product.

- Do not use the device if you have any eye conditions, injuries, or have recently undergone eye surgery. Consult a medical professional if unsure.
- Keep out of reach of children.
- Do not immerse the device in water or use it in wet environments. The device is not waterproof.
- Use only the provided charging cable and adapter.
- Do not disassemble, repair, or modify the device yourself. Contact customer support for service.
- If you experience any discomfort or pain during use, discontinue immediately.
- Avoid prolonged use. Recommended session duration is typically 10-15 minutes.

- Do not use while driving or operating heavy machinery.

### 3. PACKAGE CONTENTS

---

Verify that all items are present in the package:

- breo iSee M Eye Massager Unit
- USB Charging Cable
- Storage Bag
- User Manual (this document)



Figure 3.1: The breo iSee M eye massager main unit.

### 4. PRODUCT FEATURES

---

The breo iSee M is equipped with several features designed for optimal eye relaxation:

- **Warm Compress Function:** Provides gentle heat to the eye area, promoting comfort and relaxation.
- **Portable Design:** The device can be folded for easy storage and transport, making it ideal for travel.
- **Ergonomic Fit:** Designed to conform to the contours of the face for a comfortable and effective massage experience.
- **Simple Controls:** Intuitive buttons for power and mode selection.

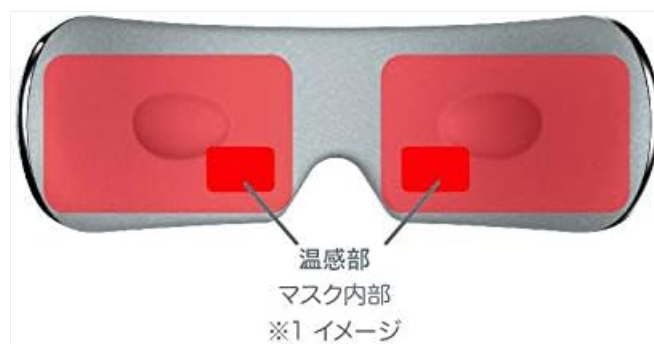


Figure 4.1: Illustration of the warm compress areas within the eye massager mask.



Figure 4.2: Diagram illustrating the slim design and how the device fits the contours of the face, including the temples.

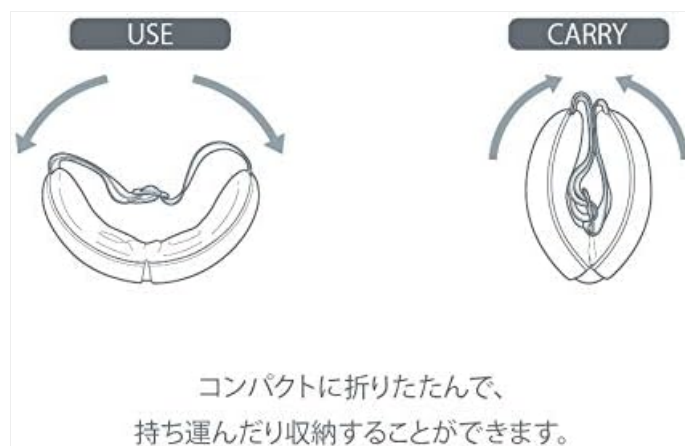


Figure 4.3: Instructions on how to fold the device compactly for portability and storage.

## 持ち運びにも便利な専用保管袋付き

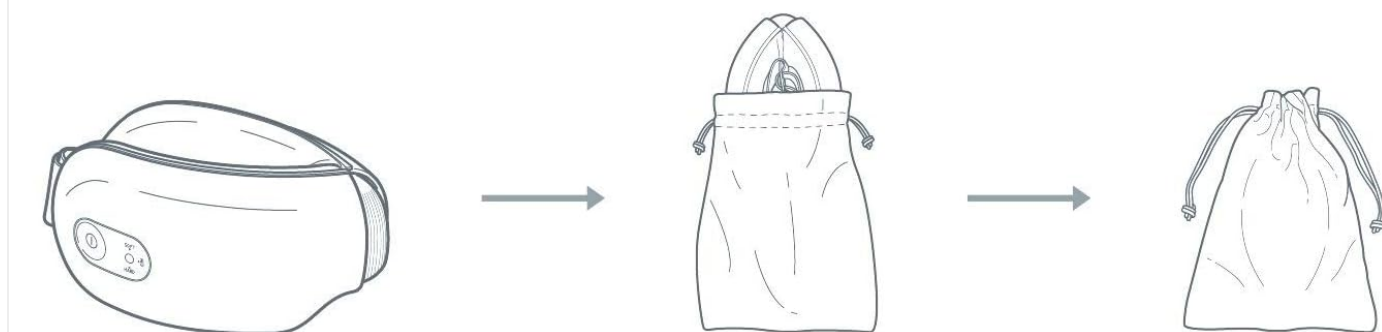


Figure 4.4: The dedicated storage bag for convenient carrying and protection of the device.

## 5. SETUP

---

### 5.1 Charging the Device

1. Connect the USB charging cable to the charging port on the device.
2. Connect the other end of the USB cable to a compatible USB power adapter (not included) or a computer's USB port.
3. The indicator light will show the charging status. Refer to the indicator light section in Operating Instructions for details.
4. Once fully charged, disconnect the charging cable.

### 5.2 Wearing the Device

1. Unfold the eye massager if it is in its compact storage state.
2. Place the massager over your eyes, ensuring the soft padding rests comfortably around your eye sockets.
3. Adjust the elastic strap to secure the device firmly but comfortably around your head. It should not feel too tight or too loose.

## 6. OPERATING INSTRUCTIONS

---

### 6.1 Controls Overview



Figure 6.1: Location of the power and mode buttons on the device.

- **Power Button ( ):** Press and hold to turn the device ON/OFF. Short press to cycle through modes (if applicable).
- **Mode Button ( ):** Short press to change massage modes or intensity levels.

## 6.2 Basic Operation

1. **Power On:** Press and hold the Power button ( ) for 2-3 seconds until the device powers on.
2. **Select Mode:** Short press the Mode button ( ) to cycle through available massage modes. The device may have different combinations of air pressure, vibration, and heat.
3. **Adjust Heat (if applicable):** Some models may have a dedicated heat button or heat integrated into modes. Refer to the specific mode descriptions.
4. **Power Off:** Press and hold the Power button ( ) for 2-3 seconds until the device powers off. The device may also automatically shut off after a set duration (e.g., 15 minutes).

## 7. MAINTENANCE

---

### 7.1 Cleaning

- Ensure the device is powered off and unplugged before cleaning.
- Wipe the surface of the device with a soft, slightly damp cloth.
- Do not use abrasive cleaners, alcohol, or harsh chemicals.
- Allow the device to air dry completely before storing or next use.

### 7.2 Storage

- Store the device in a cool, dry place, away from direct sunlight and extreme temperatures.
- Use the provided storage bag to protect the device from dust and scratches.
- Fold the device compactly as shown in Figure 4.3 for efficient storage.

## 8. TROUBLESHOOTING

---

Problem	Possible Cause	Solution
---------	----------------	----------

Problem	Possible Cause	Solution
Device does not turn on.	Battery is low or depleted.	Charge the device fully.
No heat or weak heat.	Heat function not activated or device malfunction.	Ensure the correct mode with heat is selected. If problem persists, contact support.
Device stops unexpectedly.	Automatic shut-off or low battery.	This is normal if the session duration is reached. Recharge the device if the battery is low.
Discomfort during use.	Improper fit or sensitive skin.	Adjust the strap for a comfortable fit. Discontinue use if discomfort persists.

## 9. SPECIFICATIONS

Parameter	Value
Model Number	BRE1200H
Brand	breo
Color	Gray
Power Source	Battery Powered (Lithium Ion)
Product Weight	576 g
Package Dimensions	29.7 x 23.4 x 10.1 cm
Water Resistance Level	Not Waterproof
Application Area	Face (Eyes)

## 10. WARRANTY AND SUPPORT

For warranty information and customer support, please refer to the warranty card included with your product or visit the official breo website. Keep your purchase receipt as proof of purchase.




If you encounter any issues not covered in the troubleshooting section, please contact breo customer service for assistance.

**Contact Information:**

- Official breo Website: [www.breo.com](http://www.breo.com) (Please check for regional websites)
- Customer Service Email/Phone: Refer to product packaging or official website.



**Related Documents**

 <p>Breo iSee M</p> <p>SIMPLY RELAX</p> <p>breo</p>	<p><a href="#">Breo iSee M Eye Massager User Manual</a></p> <p>User manual for the Breo iSee M eye massager, detailing product features, functions, safety precautions, troubleshooting, and maintenance.</p>
 <p>Breo WOWO S</p> <p>SIMPLY RELAX</p> <p>breo</p>	<p><a href="#">Breo WOWO S Hand Massager User Manual   Features, Functions, and Safety</a></p> <p>Comprehensive user manual for the Breo WOWO S electric hand massager. Learn about its features, functions, modes, troubleshooting, and safety precautions for optimal use.</p>
 <p>Breo</p> <p>USER MANUAL</p> <p>breo</p>	<p><a href="#">Breo iSee4 Eye Massager User Manual</a></p> <p>User manual for the Breo iSee4 Eye Massager, detailing product configuration, functions, operation, precautions, and troubleshooting.</p>
	<p><a href="#">Breo iDream 5S Head Massager User Manual - Features, Functions, and Safety</a></p> <p>Comprehensive user manual for the Breo iDream 5S head massager. Learn about its features, functions, operating instructions, troubleshooting, and safety precautions.</p>